

HW4 Instruction Manual

Version 2

1	OVERVIEW	7
1.1	HW4 Standard Edition	7
1.2	HW4 Professional Edition	7
1.3	HW4 Professional with OPC Server	8
1.4	HW4 Professional with AwQuick	8
1.5	HW4 Professional Trial Version	8
1.6	HW4 Validated	8
2	NEW IN VERSION 2.0.0	8
2.1	Microsoft .NET version 2.....	8
2.2	Water activity functions	9
2.3	Online graph – alarm band	9
2.4	User login security	9
2.5	Signature of protocols and log files	9
2.6	Protocol file format	9
2.7	Alarm function	9
2.8	Probe adjustment.....	10
2.9	Mean kinetic temperature added to statistical calculations	10
2.10	Windows notification area	10
2.11	Device Search at the beginning of each session	10
2.12	Online buffer and online graph.....	10
2.13	Product key HW4 ID number	10
2.14	Changes to the HW4 main menu bar.....	11
2.15	Multiple and concurrent HW4 sessions.....	11
2.16	Devices compatible with HW4.....	11
3	INSTALLATION.....	12
3.1	Software Licensing Agreement.....	12
3.2	HW4 requirements	13
3.2.1.	Computer/Operating System Requirements	13
3.2.2.	Operating System Compatibility	13
3.2.3.	Important information to review prior to installing HW4.....	13
3.3	Installing HW4.....	15
3.4	Updating HW4 to a newer version	19
3.4.1.	Partial update: HW4.exe.....	19
3.4.2.	Full update: Install_HW4.exe	19
3.5	Uninstalling HW4 (full uninstall)	19
4	INITIAL START-UP	20
4.1	Language selection	20
4.2	Product key and registration	21
4.3	Automatic device discovery	22
4.4	Creating the first user (HW4 Professional)	23
4.5	Access to functions and screens (HW4 Professional).....	23
5	CONNECTING THE FIRST DEVICE	24
5.1	Device connection methods.....	24
5.1.1.	Physical Serial port.....	24
5.1.2.	Bluetooth serial port (HygroLog NT with the Bluetooth option).....	24
5.1.3.	USB port	25
5.1.3.1.	Installing the ROTRONIC USB driver	25
5.1.4.	LAN (TCP/IP) - wired connection	28
5.1.4.1.	Automatic configuration (new devices only)	28
5.1.4.2.	Manual configuration	28
5.1.5.	LAN (TCP/IP) - wireless connection.....	32
5.1.5.1.	Wireless security.....	32
5.1.5.2.	Automatic configuration	32
5.1.5.3.	Manual configuration	32
6	CONNECTING AND DISCOVERING ADDITIONAL DEVICES.....	39
6.1	Definition: masters and slaves (HW4 Professional).....	39
6.2	Masters	39
6.3	RS-485 slaves (HW4 Professional)	39

6.3.1.	Slaves connected to a COM or USB master	39
6.3.2.	Slaves connected to a LAN master	40
6.3.3.	RS-485 address	40
7	HW4 MAIN SCREEN	42
7.1	Device Tree (left pane)	43
7.1.1.	Tree populated with individual devices	43
7.1.2.	Tree populated with device group(s) - HW4 Professional	43
7.1.3.	Instrument group: renaming and changing the contents	44
7.2	Right pane in Device View mode	45
7.2.1.	Current Values tab	46
7.2.1.1.	Customizing the Current Values tab	46
7.2.2.	Online graph	49
7.2.3.	Online Buffer tab	50
7.2.4.	Log to PC tab	52
7.2.5.	View tab	53
7.2.6.	OPC Server tab	53
7.2.7.	AwQuick / AwE Mode tab	54
7.3	Right pane in Group View mode	55
7.3.1.	Current Values tab	55
7.3.1.1.	Customizing the Current Values tab	55
7.3.2.	On-line graph	58
8	HW4 MAIN MENU BAR	59
8.1	File	59
8.2	Devices and Groups	59
8.3	View	61
8.4	Settings and Tools	62
8.5	Help	63
9	SHORTCUT BUTTONS	63
10	ETHERNET CONFIGURATION TOOL	64
11	USERS AND PASSWORDS (HW4 Professional)	65
11.1	Minimum user rights	65
11.2	Additional user rights	66
11.3	Creating and logging-on the first user	67
11.4	User table form	69
11.5	Menu bar	70
11.6	User login security	71
12	COLUMN HEADERS	72
13	HW4 GLOBAL SETTINGS	74
13.1	View Tab	74
13.2	General Tab	76
13.3	Language/Unit System Tab	79
13.4	File Locations Tab	81
13.5	Graph Settings Tab	83
13.6	Alarm Settings Tab (HW4 Professional)	86
13.7	Events Tab (HW4 Professional)	88
14	USER EVENTS (HW4 Professional)	89
14.1	Menu bar	90
15	ALARM TABLE (HW4 Professional)	91
15.1	Alarm acknowledgement	92
15.2	Alarm table menu bar and buttons	93
16	PSYCHROMETRIC CONVERSIONS	95
17	DEVICE MANAGER - HygroLog NT	96
17.1	Device Manager Menu Bar	97
17.2	Device Owner	98
17.3	Device Information	99
17.4	Language / Unit System	100
17.5	Memory Card	101
17.6	Standard Inputs	102
17.7	Display	103

17.8	LEDs.....	104
17.9	Keypad	105
17.10	Sounds.....	106
17.11	Power Source	107
17.12	Digital Interface	108
17.13	Docking Station.....	108
17.13.1.	Probe type: HygroClip.....	109
17.13.2.	Probe type: Analog probe (single channel):.....	110
17.13.3.	Probe type: Pressure probe	111
17.13.4.	External Contacts.....	112
17.13.5.	Relay 1 / Relay 2 (docking station DS-R-1)	112
18	DEVICE MANAGER - Transmitters and indicators	116
18.1	Device Manager Menu Bar	117
18.2	Device Information	118
18.3	Language / Unit System.....	119
18.4	Probe Inputs	120
18.4.1.	HygroClip probe.....	120
18.4.2.	Analog Probe (single channel):	121
18.4.3.	Pressure probe	122
18.5	Analog Outputs.....	123
18.6	Display.....	124
18.7	Keypad	125
18.8	Digital Interface	126
18.9	Custom Calculation	127
19	DEVICE MANAGER – HygroClip Digital Interface.....	128
19.1	Device Manager Menu Bar	129
19.2	Device Information	129
19.3	Language / Unit System.....	130
19.4	Probe Inputs	131
19.5	Digital Interface	132
20	DEVICE MANAGER – HygroClip Alarm.....	133
20.1	Device Manager Menu Bar	134
20.2	Device Information	134
20.3	Logical Inputs	135
20.4	Relays	136
20.5	Programming.....	137
20.6	Digital Interface	138
20.7	Buttons and Switches.....	139
21	DEVICE MANAGER – HygroStat MB	140
21.1	Device Manager Menu Bar	141
21.2	Device Information	141
21.3	Language / Unit System.....	142
21.4	Probe Inputs	143
21.5	Relays	144
21.6	Display.....	145
21.7	Digital Interface	146
22	DEVICE PROTECTION	147
23	DATA LOGGING – with the HygroLog NT.....	149
23.1	Menu Bar.....	150
23.2	Programming the HygroLog NT for data logging	150
23.2.1.	Recording the status of external contacts (logical inputs).....	152
24	AUTOMATIC LOG FILE DOWNLOAD - HygroLog NT	153
25	DATA LOGGING - directly on the PC	156
25.1.1.1.	Start logging data to the PC.....	156
25.1.1.2.	Stop logging data to the PC.....	158
26	ACCESSING LOG FILES, PROTOCOLS AND EVENT FILES.....	159
26.1	HW4 main menu bar	159
26.2	Access Data (HygroLog NT)	162
26.2.1.	HW4 Explorer	162

26.2.2.	Menu Bar	164
26.2.3.	Files located on the PC	165
26.2.4.	Files located on the logger	166
27	VIEW / SIGN A PROTOCOL	168
27.1	Protocol viewer menu bar	169
28	VIEW / SIGN A LOG FILE	171
28.1	Opening a log file with View Data	171
28.1.1.	From HW4 Explorer	171
28.1.2.	From the HW4 Main menu bar	173
28.2	View Data menu bar	176
28.2.1.	File	176
28.2.2.	View	179
28.2.3.	Tools	180
28.2.4.	Help	184
28.3	View Data Toolbar	185
28.4	Graph Tab	186
28.4.1.	Using the Data Cursor	186
28.4.2.	Displaying hidden data labels	187
28.4.3.	Selecting which traces are displayed by the graph	187
28.4.4.	Attaching a text / data box to a trace	188
28.4.5.	Adding a text note directly in the graph	189
28.4.6.	Deleting a text note or text/data box	189
28.5	Data Table / File Name Tab	190
28.6	Working with log files and graphs (How To)	190
29	PROBE ADJUSTMENT	191
29.1	Overview	192
29.2	Menu Bar	192
29.3	Temperature adjustment	193
29.4	Adjustment against a reference relative humidity value	195
29.5	Adjustment against a reference chilled mirror	197
30	ALARM INDICATION & REPORTING - Overview	199
30.1	Standard alarm notification (all HW4 versions)	199
30.1.1.	Device tree	199
30.1.2.	Right Pane - Current Values Tab	199
30.2	Optional alarm notification (HW4 Professional)	200
30.3	Alarm Reporting (HW4 Professional)	201
30.3.1.	Viewing the Alarm Table	201
30.3.2.	Alarm summary report	201
30.4	Additional sources of alarm data	201
30.4.1.	Log files	201
30.4.2.	HygroLog NT Logger event file	202
30.4.3.	HygroClip Alarm event memory and event viewer	203
31	FDA / GAMP COMPLIANCE (HW4 Professional)	205
31.1	Required settings and selections	205
31.2	Electronic records	205
31.3	Log File Format	206
32	RECORD KEEPING - Overview (HW4 Professional)	206
32.1	Event Files	206
32.1.1.	HW4 events	206
32.1.2.	User events	207
32.1.3.	Logger events	209
32.2	Protocols	211
32.2.1.	Protocol file names	212
32.2.1.1.	HW4 configuration protocols	212
32.2.1.2.	Instrument configuration protocols	212
32.2.1.3.	Log function programming protocols	212
32.2.1.4.	Probe adjustment protocols	212
32.3	Cross referencing protocols and event files	212
33	FILE PROTECTION	214

33.1	Authentication stamp	214
33.2	Overview	215
34	RELOCATING THE HW4 USER FOLDER	215
34.1	Method 1 (using Windows Explorer)	216
34.2	Method 2 (after the initial HW4 start-up)	219
34.3	Retrieving your previous settings and other data	222
35	CHANGING THE BAUD RATE OF AN ETHERNET DEVICE	224
36	CONCURRENT HW4 SESSIONS ON DIFFERENT WORKSTATIONS	227
37	MULTIPLE HW4 SESSIONS ON THE SAME PC	228
38	STARTING HW4 AUTOMATICALLY with Windows	229
39	PROBE ADJUSTMENT BASICS	230
39.1	Temperature	230
39.2	Humidity	230
39.2.1.	Calibration device	230
39.2.2.	ROTRONIC certified humidity standards	231
40	BASIC TCP/IP CONCEPTS	232
40.1	Compatibility requirements	232
40.2	DHCP	233
40.3	MAC Address	233
41	WATER ACTIVITY MEASUREMENT WITH HW4	234
41.1	Water Activity: definition and applications	234
41.1.1.	Definitions	234
41.1.2.	Aw and Temperature	235
41.1.3.	Applications	235
41.2	Instruments and probes for measuring water activity	236
41.3	Water activity measurement modes in HW4	236
41.4	Mode selection and settings	237
41.5	Using the AwE mode	238
41.6	Using the AwQuick mode	239
41.7	Water activity measurement report	240
41.7.1.	Signing the report	241

Document revisions:

Document Version	Software Version	Date	Notes
-10	2.0.0	Dec 7, 2006	Original release See: New in version 2.0.0

1 OVERVIEW

The HW4 software was developed by ROTRONIC AG for use with the ROTRONIC line of digital instruments and devices. HW4 is available in the following versions:

1.1 *HW4 Standard Edition*

- Unlimited number of instruments (depends on capabilities of PC and available ports)
- On-line display of the measured and calculated values (dew point or other), limited to one instrument at a time.
- Automatic retention of the most recent data in a temporary on-line buffer (multiple instruments).
- Logging of measured and calculated values to the PC (multiple instruments).
- Easy-to-read graphs and data tables.
- Statistical data tools: mean, standard deviation, minimum and maximum of the recorded data (off line only).
- Data print-out (table or graph)
- Automatic device/instrument recognition and identification
- Device/instrument configuration
- Access to the data recorded by a data logger and data transfer to the PC
- Adjustment (calibration) of the HygroClip digital probes
- Psychrometric Conversion Tool
- Built-in security to protect against data manipulation

1.2 *HW4 Professional Edition*

HW4 Professional edition meets the FDA / GAMP requirements. This version of HW4 allows multiple users, with either administrator or standard rights and password protection.

- All the functionality of HW4 Standard Edition
- RS-485 multi-dropped sub-networks of up to 64 instruments per sub-network
- On-line display of the data from multiple instruments and probes
- Possibility to overlay and synchronize data from several log files into a single graph
- User Event Logging: automatically records user main operations.
- Multiple users distributed into two groups: administrator and standard, each with different rights
- Self Event Logging: automatically records any software problem to facilitate troubleshooting
- Logger Event Logging: automatically records the data logger internal events and configuration changes.
- Automatic creation of protocols detailing instrument configuration and programming changes as well as and probe adjustments.
- Optical or acoustical indication of an alarm condition, tracking of alarm conditions in a table, and possibility of printing a report, sending an e-mail, etc.
- Visual Indication an alarm conditions when viewing log file data
- Password protected log-in
- Meets the requirements of FDA 21CFR Part II for electronic records and electronic signatures
- Meets the EU GMP requirements regarding pharmaceuticals.

1.3 HW4 Professional with OPC Server

HW4 Professional with OPC server is a special version of HW4 Professional which includes all features of HW4 Professional while also providing OPC tags that can be enabled for each individual device communicating with HW4. This allows transferring data with practically any OPC client (requires configuration / programming of the OPC client by the user).

1.4 HW4 Professional with AwQuick

HW4 Professional with AwQuick is a special version of HW4 Professional designed to facilitate measurement of the water activity (Aw) of foods, pharmaceuticals, etc. This version includes all the features of the regular HW4 Professional Edition. Two modes are available for measuring water activity: AwE and AwQuick. Both modes can be used with any of our instruments. In the AwE mode, HW4 monitors the natural equilibration of the product being measured and automatically stops the measurement process when equilibrium is reached. With most products, natural equilibrium requires from 45 to 90 minutes. The AwQuick mode reduces the time required to measure water activity to a few minutes, usually with almost the same accuracy as the AwE mode.

1.5 HW4 Professional Trial Version

HW4 Trial is a fully functional version of HW4 Professional, including the OPC server and water activity measurement functionality. A compressed installation file (zip) can be downloaded free of charge from the ROTRONIC website at:

http://www.rotronic.ch/downloads/feuchte_software.php?language=ch-en

Note: this web page is subject to change

Activation requires an HW4 product key which can be requested by filling a form on the ROTRONIC web site. A trial product key will be sent to your e-mail address. This must be entered in the registration form which appears when starting HW4 for the first time. After 30 days, the trial product key expires and HW4 can no longer be started.

1.6 HW4 Validated

HW4 Validated offers the same functionality as HW4 Professional with OPC server. Additionally, the "HW4 e-compliance Package" is available. This extensive collection of documents (including template of qualification documents) is designed to support the "regulated user" by qualifying/validating HW4-based solutions.

2 NEW IN VERSION 2.0.0

HW4 version 2.0.0 includes all the functions and features of HW4 version 1.2.2 with the following modifications and additions:

2.1 Microsoft .NET version 2

Previous versions of HW4 required the Microsoft .NET framework version 1 to be installed on the PC. HW4 version 2.0.0 requires the Microsoft .NET framework version 2. Both versions of the .NET framework can coexist on the same PC. A copy of .NET framework version 2 is provided on the HW4 installation CD. The change to .NET version 2 is required for HW4 to be compatible with the coming replacement for the Windows XP operating system (Vista).

2.2 Water activity functions

See Overview - HW4 Professional for Aw Measurement

2.3 Online graph – alarm band

HW4 version 2.0.0 allows displaying an alarm band on the online graph

2.4 User login security

- Attempts to login with the wrong password are now being recorded.
- HW4 startup will be aborted after 3 attempts with the wrong password

2.5 Signature of protocols and log files

To provide additional authentication, protocols and log files can be electronically signed by a user having the necessary rights. The user can also add text / comments to the file. The user must enter his or her password before signing the document.

Previous versions of HW4 used an “electronic signature” to verify that a document is an original, unaltered HW4 document. To distinguish between the user signature and the so called “electronic signature”, we have replaced “electronic signature” with the words “authentication stamp”.

2.6 Protocol file format

The Microsoft Word format (file extension DOC) is no longer supported by HW4. All protocol files are now regular text files (extension TXT) and can be opened with Notepad, independently from HW4. Integrity of both the file content and signature area are protected by a checksum generated by HW4.

2.7 Alarm function

The entire HW4 alarm function has been completely redesigned:

- A priority level (low, medium or high) can be attached to each type of alarm
- A time delay can be specified for the following alarm types: no communication with device, measured values out-of-limits or not available and HygroLog NT error.
- An alarm test can be performed manually at any time or automatically at regular intervals of time
- The alarm functions can be enabled during regular workdays only and / or during some hours of the day.
- Sending an alarm by email, sending an alarm to the Windows default printer or generating an acoustic signal on the PC can be restricted depending on the alarm priority level.

To reflect these changes, both the alarm table and the summary report that can be printed from the table have also been redesigned. Both the configuration of the alarm table and alarm summary report can be customized.

Other new features include:

- A separate path can be specified for the HW4 file that retains a record of all alarm events (HW4 alarm file)
- The contents of the HW4 alarm table (or part of it) can be saved to a file in XLS format.

2.8 Probe adjustment

The probe adjustment function (HygroClip probes only) has been completely redesigned to make it more similar to the equivalent function in each instrument firmware.

For both humidity and temperature adjustment the following choices are now provided:

- 1-point adjustment: this causes the instrument to which the probe is attached to apply a global offset across the entire measuring range of the probe
- Multi-point adjustment (from 2 to 4 points): the value of the reference is sent to the instrument and the instrument applies a correction to the probe depending on brackets that are predefined in the instrument firmware.

In addition, HW4 now allows the adjustment of humidity directly against a chilled mirror instrument. The dew or frost point is entered in HW4 as the reference value and HW4 automatically calculated the corresponding humidity based on the temperature measured by each probe.

2.9 Mean kinetic temperature added to statistical calculations

View Data (the HW4 log file viewer.) offers the possibility of displaying statistical data such as the average, the minimum, maximum and standard deviation of the data. In addition to the average, the temperature data now also includes the mean kinetic temperature (MKT).

The MKT is calculated based on the value of the activation energy [kJ / mole] entered in the HW4 Global Settings – Language / Unit System Tab. The activation energy should be entered using the international unit system, even when the English unit system is otherwise being used.

2.10 Windows notification area

The Windows notification area is located on the taskbar, immediately to the left of the clock. HW4 can now be sent (minimized) as a shortcut icon to the Windows notification area while it is still keeps running.

Sending HW4 to the Windows notification area can be done as follows:

- Automatically after each startup (HW4 Global Setting)
- Manually from the main menu bar

The alarm table will still appear on the PC monitor when HW4 is running in this mode.

2.11 Device Search at the beginning of each session

Upon starting a new session, HW4 now automatically does a new search for the devices that were present in the device tree at the end of the previous HW4 session. This ensures that users will not attempt to work with devices before HW4 has established a communication.

2.12 Online buffer and online graph

The General tab of HW4 Global Settings now offers the possibility of clearing the online buffer and online graph. When adjusting a probe, clearing the online buffer when the probe appears to have reached equilibrium forces the automatic scale of the graph to the highest resolution.

2.13 Product key HW4 ID number

HW4 version 2.0.0 uses a new set of product keys / ID number. The product keys used with previous versions will not be accepted by this version.

2.14 Changes to the HW4 main menu bar

Access to some functions from the HW4 main menu bar has been changed: the alarm table, user events and psychrometric conversions are now located under View. Column headers is now located under Settings and tools.

2.15 Multiple and concurrent HW4 sessions

HW4 version 2.0.0 handles both multiple sessions on the same PC and concurrent sessions on several workstations differently than prior versions and no longer relies on the attributes of the HW4.ini file to keep track of sessions.

Multiple sessions on the same PC are generally not recommended as they can produce unexpected results. HW4 version 2.0.0 allows opening multiple sessions on the same PC and simply gives a warning that HW4 is already open.

Concurrent sessions of HW4 can be run on separate workstations connected to a file server. The workstations can use a common folder located on the file server so as to have access to the same recorded data, protocols, event files, etc. The first session to be opened on a workstation is the Master session. When a concurrent session is opened on another workstation, HW4 issues a warning. This session has only limited functionality and the user has only minimum rights (essentially limited to viewing the measurement data), regardless of the rights given to this user by the HW4 administrator.

HW4 version 2.0.0 uses a file located in the HW4 User Folder (HW4RunningStatus.dat) to keep track of the first session (Master session) and of which PC or workstation has opened this session.

2.16 Devices compatible with HW4

HW4 version 2.0.0 can be used with the following devices:

- HygroLog NT data logger
- HygroFlex 2, HygroFlex 3 and M3 transmitters (same icon in device tree)
- HygroLab 2 and HygroLab 3 bench indicators
- HygroPalm 2 and HygroPalm 3 portable indicators
- HygroClip DI digital interface
- HygroClip Alarm programmable logic
- HygroStat MB

Recommended minimum firmware version:

HygroLog NT	: v 1.2
Docking station for HygroLog NT	: v 1.4
HygroPalm, HygroLab, HygroFlex and M3	: v 4.1
HygroClip DI	: v 1.0
HygroClip Alarm	: v 2.1
HygroStat MB	: v 1.0

Note: HW4 version 2.0.0 can also be used to configure the HygroFlex 1, HygroLab 1, HygroPalm 1 and HygroPalm 0. These instruments are not designed to be used with HW4 for any other purpose.

3 INSTALLATION

3.1 *Software Licensing Agreement*

1) Rights to software, know-how and process

Without any specific agreement, the customer may make use of the available software HW4, the know-how, the data carriers and the documentation to the extent provided for, but must not pass it on to third parties.

Any improvement, amendment or copying of the software by the customer requires permission in writing from ROTRONIC. The customer must include the same property rights notations on all modifications and copies as there are on the original.

In all cases, ROTRONIC or its licensors retain the ownership of the software and the know-how and also retain the right to continued utilization, even if ROTRONIC supplies the source codes or if the customer subsequently makes amendments to the software programs or know-how recordings.

2) Warranty against defects

The customer's warranty rights presume that he or she has fulfilled his or her legal obligation to examine and make a complaint in respect of a defect immediately on receipt of the goods as is required.

ROTRONIC shall be released from all obligations under all warranties either expressed or implied, if the software covered hereby is modified by persons other than its own authorized personnel, unless such modification by others is made with the written consent of ROTRONIC or unless such modification in the sole opinion of ROTRONIC is minor or unless such modification is merely the installation of a new ROTRONIC update or patch for the software. ROTRONIC makes no warranties which extend beyond the description of the software covered hereby other than as expressly stated herein, ROTRONIC expressly and specifically disclaims the implied warranty of merchantability and makes no warranty with respect to the fitness of the software covered hereby for any particular purpose or use unless such a warranty is expressly set forth.

The buyer or anyone claiming under any warranty relating to the software sold hereunder agrees that if ROTRONIC breached any such warranty, or any warranty implied either in fact or by operation of law, or if the software warranted hereunder proves defective in any manner whatsoever, ROTRONIC sole liability hereunder is limited to either replacement of the defective software or at the option of ROTRONIC, refunding to the buyer the purchase price paid for such defective software. The buyer and anyone else claiming under any warranty relating to the software sold hereunder expressly and specifically agree that ROTRONIC is not responsible for, and the buyer or such other claimant or claimants shall assume, any liability for property damage, prospective profits, special, indirect, or consequential damage, or other commercial or economic loss arising out of use or possession of the software sold hereunder. ROTRONIC shall not be liable for, and a buyer or anyone else claiming under any warranty relating to the software sold hereunder further agrees, and shall assume, any liability for personal injury arising out of use or possession of the software sold hereunder.

Representations and warranties made by any person, including dealers and representatives of ROTRONIC which are inconsistent or in conflict with the terms of this warranty (including but not limited to the limitations of the liability of ROTRONIC as set forth above), shall not be binding upon ROTRONIC unless given in writing and approved by an expressly authorized representative of ROTRONIC.

3.2 HW4 requirements

3.2.1. Computer/Operating System Requirements

The following are the minimum values required to install and run HW4 on a computer. It is highly recommended to exceed these values.

- Processor: Pentium II, 450 MHz
- RAM: 128 MB
- Available hard disk space: 50 MB
- Monitor: SVGA, 1024 x 768, 256 colors
- Ports: one free serial (COM) port or one free USB port or Network Interface Card / Ethernet LAN with one free port (RJ45 connector)

3.2.2. Operating System Compatibility

- Windows XP, NT4 with SP 6a or higher, Vista
- Windows 2000 with SP 2 or higher
- Windows Server 2003

HW4 was written for the Microsoft .NET framework (version 2.0) and requires this framework to be installed on the computer.

The .NET framework offers significant improvements in the areas of networking and user security. When new software is being installed, the .NET framework also eliminates the potential problem of conflicting dynamic library files (DLL). According to Microsoft, the .NET framework will be used by all future Microsoft operating systems.


3.2.3. Important information to review prior to installing HW4

Microsoft .NET framework

Prior to installing HW4 you should check if the Microsoft .NET framework (version 2.0) is installed on the PC. To do this, open Control Panel in Windows and select Add or Remove Programs. Windows displays an alphabetical list of installed programs. If Microsoft .NET framework is not listed it is not installed on the PC. Microsoft .NET framework v. 2.0 is included in the HW4 CD-ROM and can be installed as part of the HW4 installation procedure.

Location of the HW4 User Folder

During the initial startup, HW4 creates a User Folder named ROTRONIC_HW4. By default, this folder is created in C:\Documents and Settings\Windows User\Application Data, where Windows User is the name that was used to log into the current Windows session.

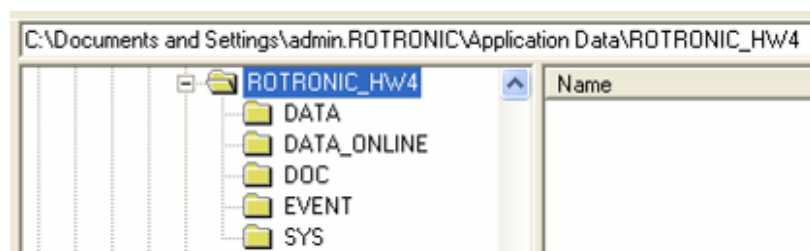

This folder is used to hold the HW4 configuration file as well as the different data, event and protocol files created by HW4.

Windows automatically gives *Windows User* **exclusive access** to the directory C:\Documents and Settings*Windows User* and subfolders. For this reason, the default path for the HW4 user folder may not be suitable for the intended use of HW4.

For instructions about changing the path of the HW4 User Folder, see [Relocating the HW4 User Folder](#) before installing HW4.

Running HW4 from several workstations with a Windows 2003 server

Generally, we do not recommend installing HW4 directly on a Windows 2003 file server. When you wish to run HW4 on multiple workstations viewing and sharing the same data, you should proceed as follows:

- Log in on one of the workstations as an administrator of the Windows 2003 server.
- Install HW4 on the workstation following the procedure described under **Installing HW4**
- After installing HW4 and prior to starting HW4, go to the HW4 installation directory on the workstation and change the path used by HW4 to locate the **HW4 User Folder** to the root directory of a file server drive that is mapped on each workstation (see [Relocating the HW4 User Folder](#)).
- Start HW4 following the procedures described in this manual under **Initial Startup**
- Close HW4 before installing HW4 on the next workstation
- Repeat this procedure for each workstation

Note: It is important to distinguish between a Windows workstation user and a HW4 user. It is equally important to distinguish between the Windows permissions to a folder and HW4 rights. When the HW4 User Folder is located on the File server and is shared by all workstations, you should give all Windows users sufficient permissions to the HW4 User Folder and subfolders (Folder Properties - Security).

3.3 Installing HW4

Insert the HW4 CD-ROM into the CD drive or your PC. The installation program should start automatically. If the installation program does not start, use My Computer in Windows to open the CD drive and double click on the file start.exe located in the root directory of the CD.



IMPORTANT: HW4 requires a PC with the Microsoft .NET Framework version 2.0 installed. If the Microsoft .NET Framework is not already installed on your PC, please click with the mouse on step 1 to install the framework.


Note: save your work and close all open files because the computer will have to restart to complete the installation process.

After installing the Microsoft .NET Framework, click on step 2 with the mouse to install HW4.


Installation begins with the following form:



Click on the NEXT button:

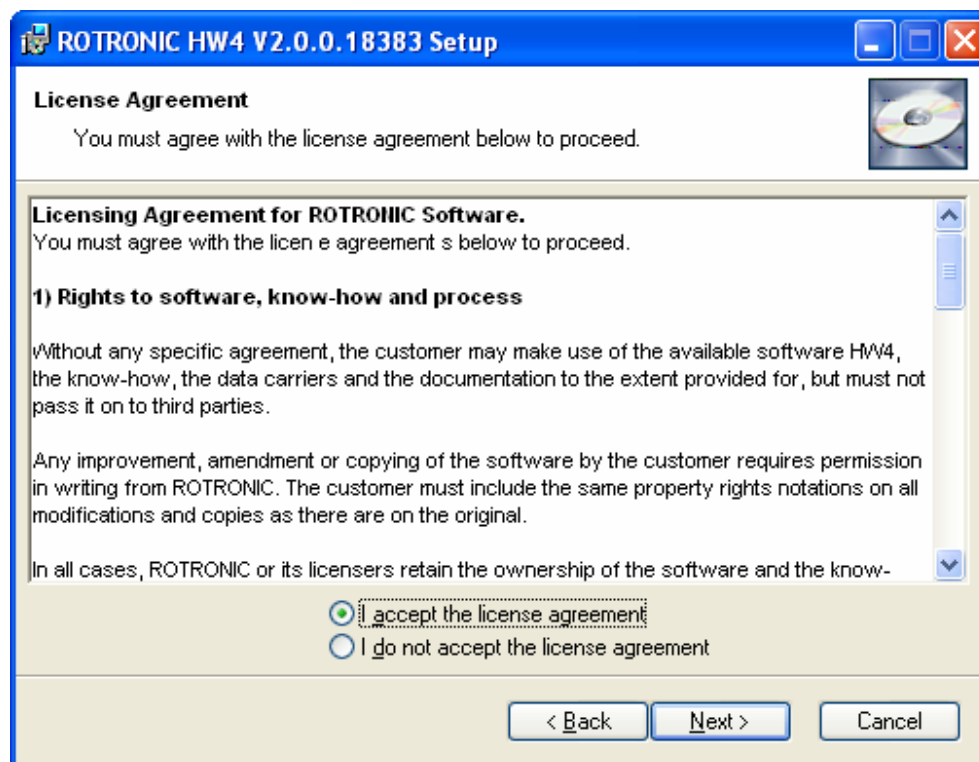

Click on the NEXT button to move through the different installation steps and follow the instructions provided on the screen.

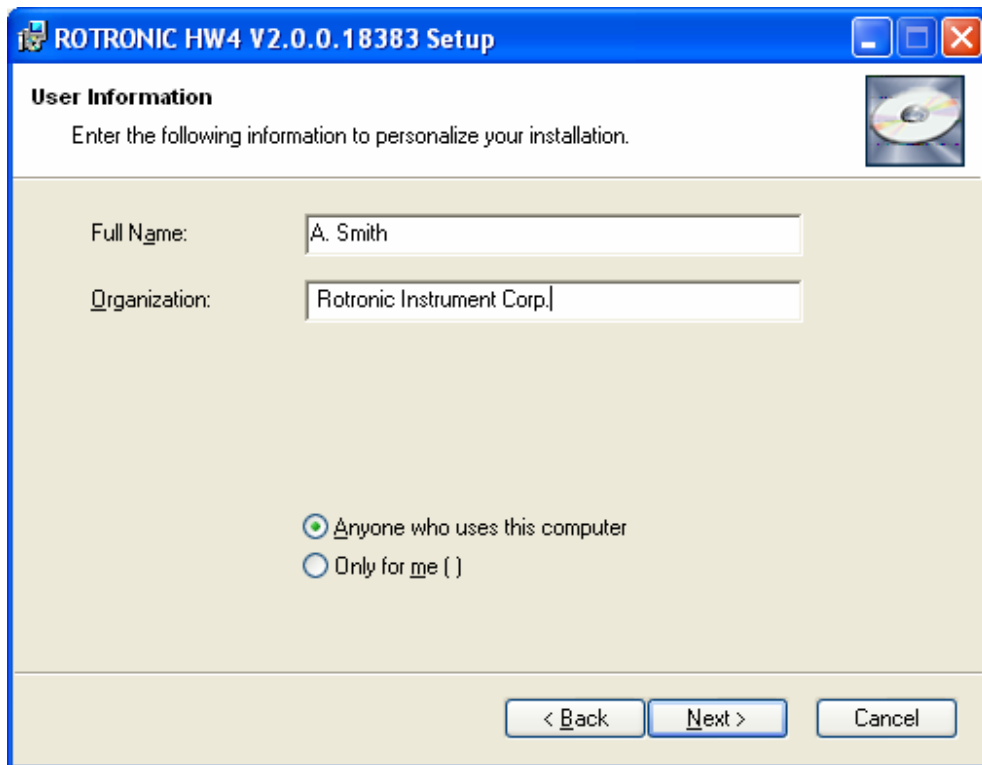

ROTRONIC HW4 V2.0.0.18383 Setup

User Information

Enter the following information to personalize your installation.


Full Name:

Organization:

☒ Anyone who uses this computer
☐ Only for me ()

< Back Next > Cancel

The choice made above determines how the HW4 install program will make entries in the Windows registry (entries for all users or for the current user only).

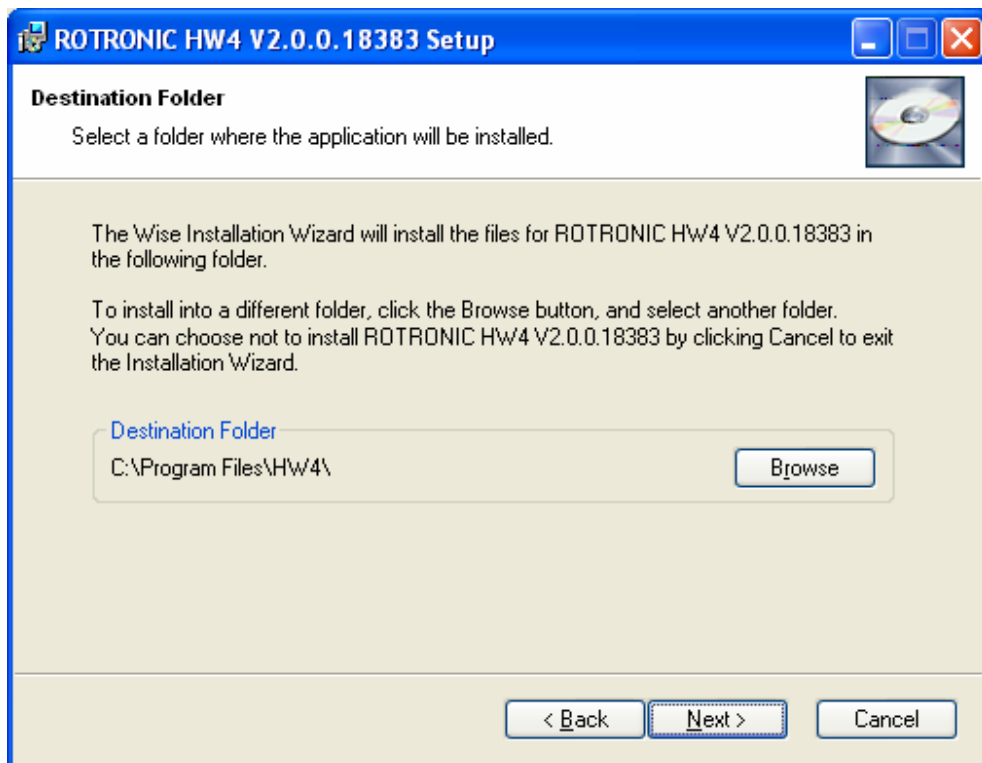

ROTRONIC HW4 V2.0.0.18383 Setup

Destination Folder

Select a folder where the application will be installed.

The 'Wise Installation Wizard will install the files for ROTRONIC HW4 V2.0.0.18383 in the following folder.

To install into a different folder, click the Browse button, and select another folder. You can choose not to install ROTRONIC HW4 V2.0.0.18383 by clicking Cancel to exit the Installation Wizard.

Destination Folder


< Back Next > Cancel

The default installation directory is C:\Program Files\HW4 and can be changed during installation.

The shortcut HW4 is automatically created on the desktop during installation:

Note: Upon starting HW4 for the first time, the folder ROTRONIC_HW4 will be automatically created in C:\Documents and Settings\Windows User\Application Data, where *Windows User* is the name that was used to log into the current Windows session.

When installation is complete, the following screen should appear:

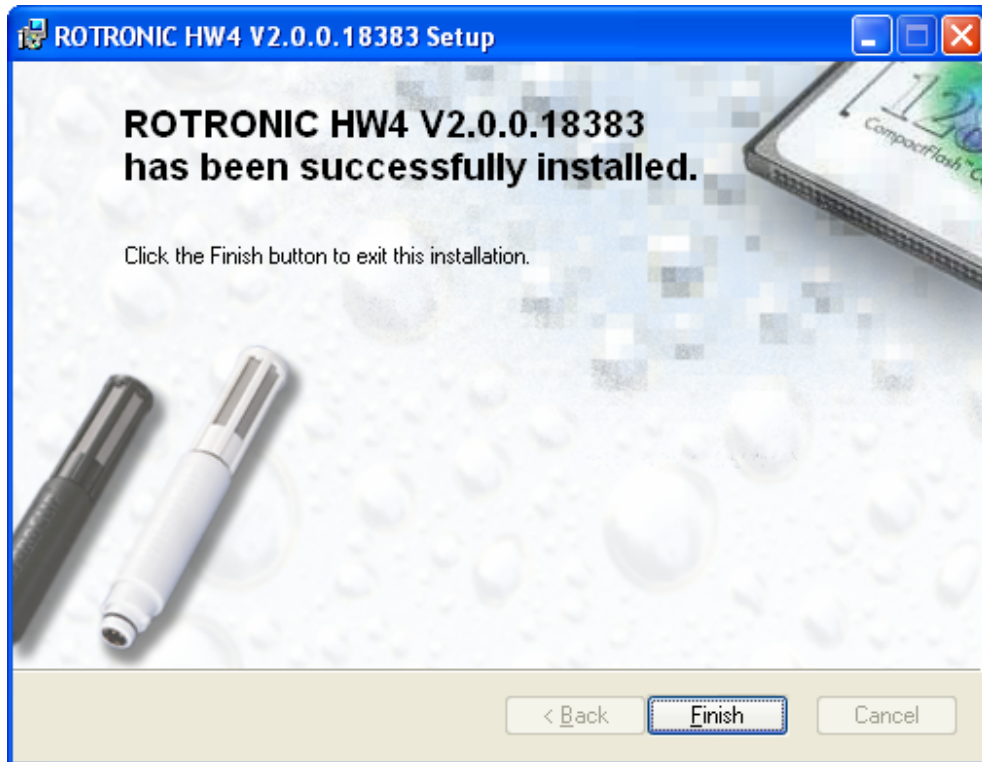

STOP HERE!

Do not start using HW4 without reading the following:

- [Relocating the HW4 User Folder](#)
- [Connecting the first device](#)

3.4 *Updating HW4 to a newer version*

When updating HW4 to a newer version, two distinct situations may present themselves. In both cases, it is advisable to make a note of your old HW4 product key prior to updating, even if in principle re-entering this information and registering the new version of HW4 are not required.

3.4.1. Partial update: HW4.exe

In the case of a partial update, simply copy the new HW4.exe file over the old file, normally located in the directory C:\Program Files\HW4.

3.4.2. Full update: Install_HW4.exe

In the case of a full update, you should first uninstall the old version of HW4. To do this properly, click on Start in Windows and open the Control Panel. Select Add or Remove Programs. ROTRONIC HW4 should be listed as one of the installed programs. Select ROTRONIC HW4 and click on the Remove button.

Note: uninstalling HW4 from the Windows Control Panel does not delete the HW4 User Folder.

IMPORTANT: do not delete the HW4 User Folder (ROTRONIC_HW4) created by the old HW4 version. Unless otherwise mentioned, this folder can be used with the updated HW4.

The default path of the HW4 User Folder is C:\Documents and Settings\XXXX\Application Data, where XXXX is the Windows user who initially installed HW4.

Deleting the HW4 User Folder would erase all users, all configuration and measurement data as well as all event tracking and protocols. It is advisable to make a temporary back-up copy of this folder.

3.5 *Uninstalling HW4 (full uninstall)*

IMPORTANT: Not following this procedure may result in problems if HW4 is installed again.

To properly uninstall HW4 and remove its main components, you should click on Start in Windows and open the Control Panel. Select Add or Remove Programs. ROTRONIC HW4 should be listed as one of the installed programs. Select ROTRONIC HW4 and click on the Remove button.


To complete a full uninstall you should manually delete the HW4 User Folder (ROTRONIC_HW4). The default path of the HW4 User Folder is C:\Documents and Settings\XXXX\Application Data, where XXXX is the Windows user who initially installed HW4.

4 INITIAL START-UP

Click on the HW4 shortcut created on your desktop by the installation program.

4.1 *Language selection*


The first form to open lists the different language files available within HW4:

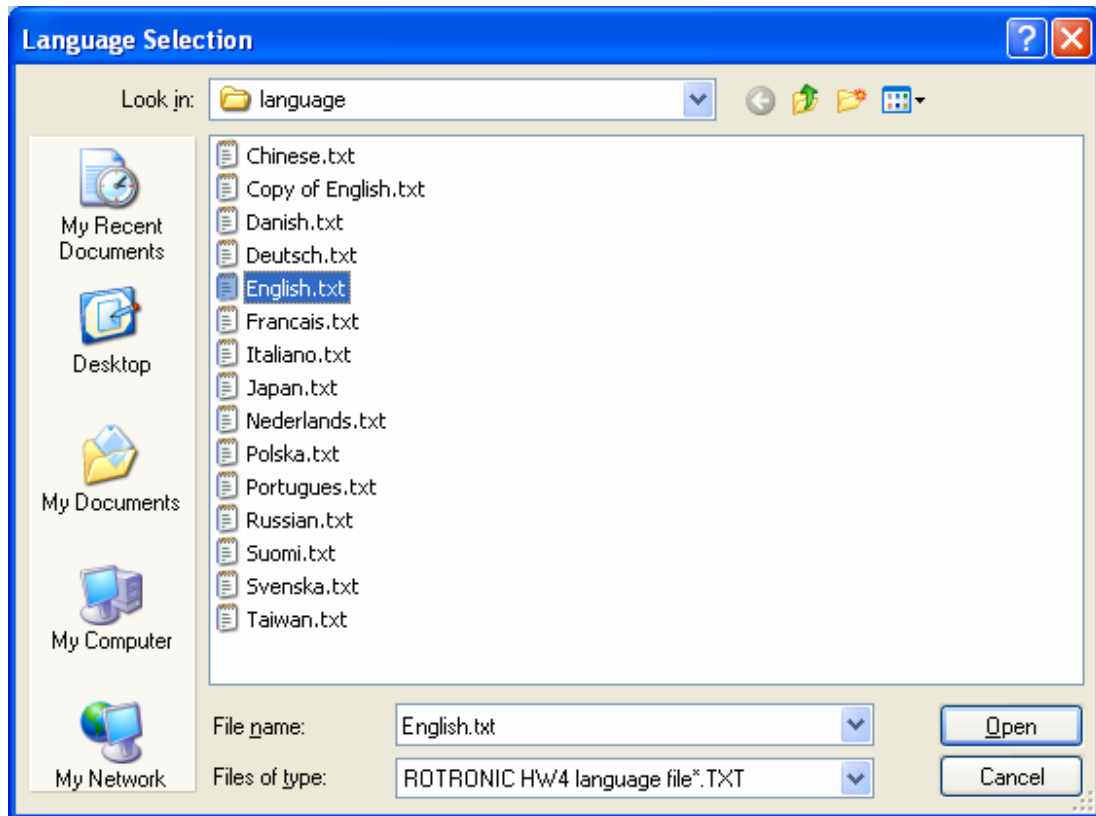

Use the mouse to highlight the desired language and click on Open.

4.2 Product key and registration

The next form is the HW4 start-up form:

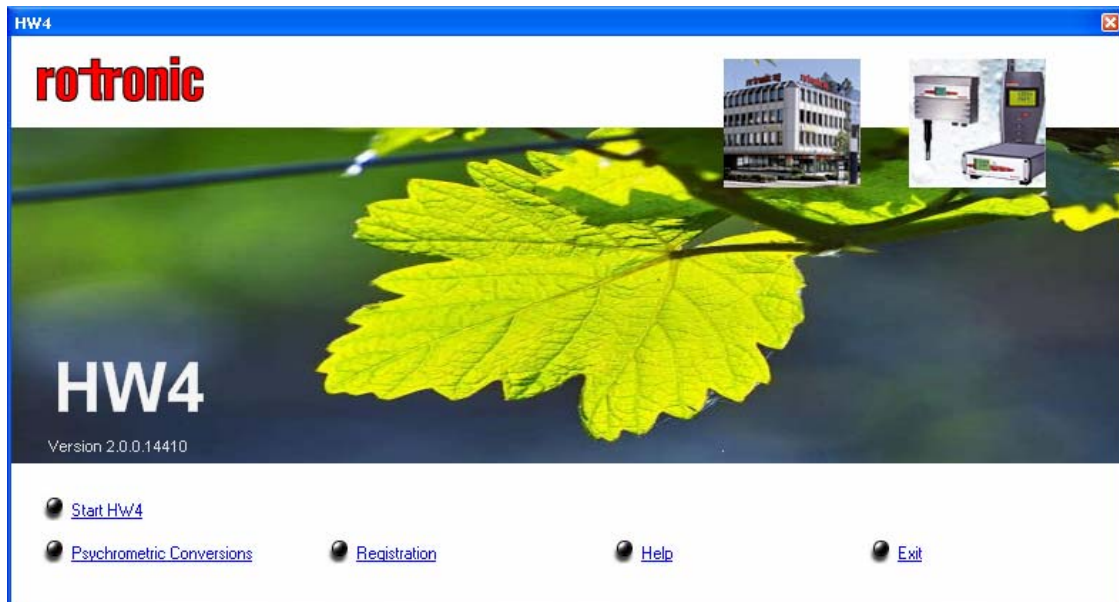


First Time Users: HW4 requires you to register and you should now click with the mouse on the underlined Registration link. This opens the following form:

Enter the HW4 product key, which is printed on a sticker affixed to the protective cover of the HW4 CD. Fill in the required information. If the box Register on-line is checked, clicking on Continue will connect you to the ROTRONIC web site, where your registration data will be automatically entered. If you do not have an internet connection or if you have already registered your copy of HW4, uncheck the Register on-line box and click on Continue.

Registering on-line offers benefits such as free updates / patches and product information.

Upon completing the registration, HW4 returns to the Start-up Screen. Depending on the HW4 product key entered in the previous form, the edition of HW4 that you are using is shown on the screen (in the example below: the validated edition of HW4 Professional).

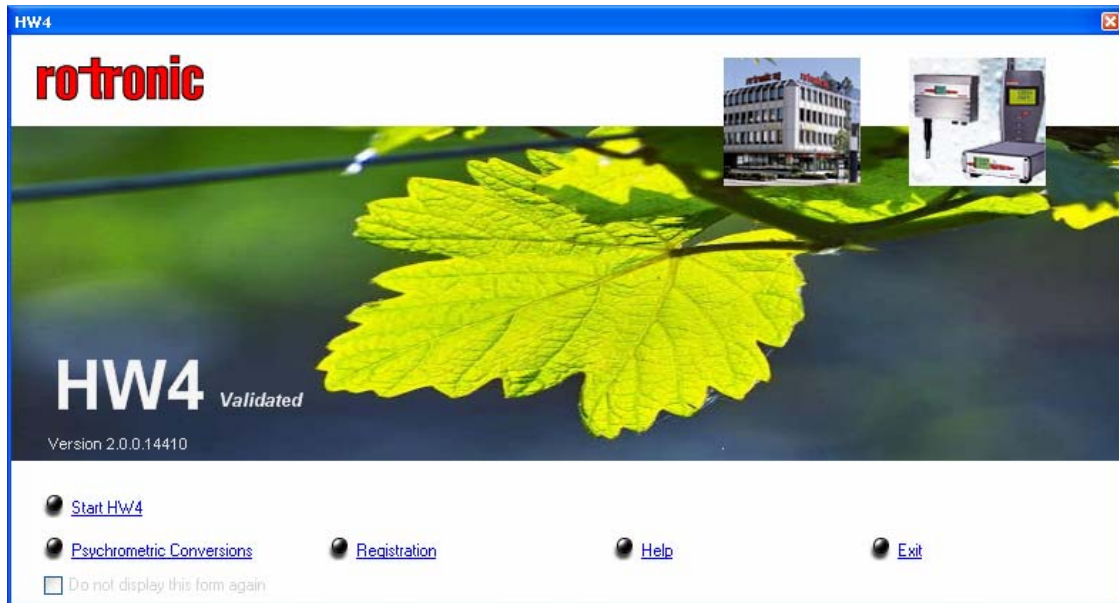

To open the HW4 main screen, click with the mouse on the blue link labeled Start HW4.


4.3 Automatic device discovery

During the initial start-up, HW4 automatically searches for any device that is connected either to a **COM port** or to a **USB port** of the PC. HW4 stops searching after finding one such device.

Devices connected by a LAN (Ethernet) must be discovered manually after HW4 completes its initial automatic search. Use **Search for Ethernet masters** or **Search for Specific IP address** under [Devices and Groups](#) in the HW4 main menu bar.

Instruments and devices connected by a Bluetooth connection must also be discovered manually, after declaring the Bluetooth serial port(s) in HW4 Global Settings.

When there is no device connected to the PC either by a COM port or by a USB port, HW4 stops searching after a single pass. In that situation, HW4 creates a fictitious HygroLog NT in the device tree (see "**HW4 main screen overview**"). This instrument is named "Offline" and should be deleted after manually discovering the actual devices connected to the PC. The red cross on top of the device icon indicates that there is no communication with the device.

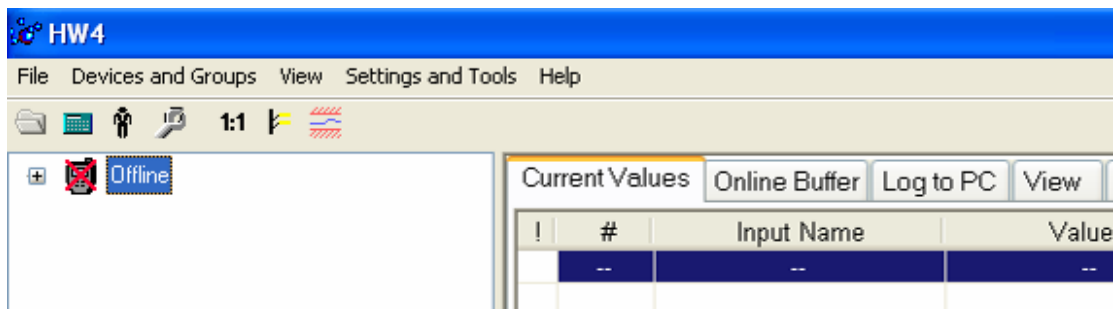


Note: whenever the device tree is empty HW4 automatically does the same search as during the initial start-up.

4.4 *Creating the first user (HW4 Professional)*

If you are using HW4 Professional, you should create at least one user with administrator rights and a password in order to comply with the FDA / GAMP requirements.

See: [Users and passwords](#).

Once the first user has been created, you will need to log in each time that you open HW4.

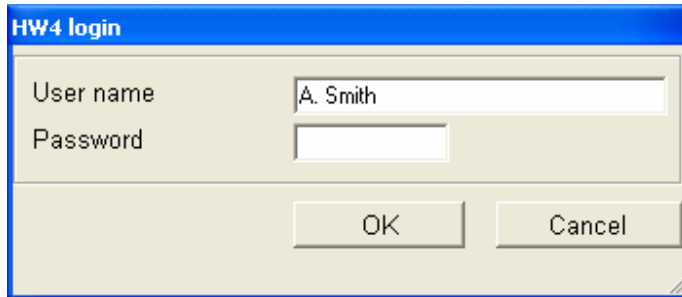The image shows a screenshot of a Windows-style dialog box titled "HW4 login". The dialog has a blue title bar. Inside, there are two text input fields. The first field is labeled "User name" and contains the text "A. Smith". The second field is labeled "Password" and is currently empty. Below the input fields, there are two buttons: "OK" and "Cancel". The dialog box has a standard Windows XP aesthetic with a light beige background and a blue border.

4.5 *Access to functions and screens (HW4 Professional)*

In all versions of HW4 Professional, access to most functions and screen views requires the current user to have sufficient rights. Whenever a function is not accessible or a screen view is not available, please verify that you have been granted the necessary rights.

For an overview of user rights, including details for each individual right, please refer to the chapter [Users and Passwords](#).

5 CONNECTING THE FIRST DEVICE

IMPORTANT:

- The HygroLog NT requires a docking station to provide an interface with the PC.
- If you want to run HW4 on multiple workstations connected to a Windows 2003 file server, use only devices connected by Ethernet (wired or wireless) because this is the only type of connection that can be shared by all workstations. Devices connected by RS-485 to an Ethernet device can also be shared.

5.1 *Device connection methods*

All versions of HW4 are compatible with the following connection methods:

- **Physical serial port**
- **Bluetooth serial port (virtual COM port)**
- **USB port**
- **LAN (TCP/IP) - cable or wireless connection**

5.1.1. **Physical Serial port**

Connecting an instrument or device to a physical COM ports does not require any driver or any special configuration. Simply connect the device to any available serial (COM) port on the PC.

5.1.2. **Bluetooth serial port (HygroLog NT with the Bluetooth option)**

Note: check the local menu of the HygroLog NT - Settings, and make sure that Bluetooth is on.

Communication via Bluetooth requires a PC equipped with an internal or external Bluetooth transceiver.

Typically, the Bluetooth driver has provisions for defining one or several Bluetooth serial ports (look for Local Services and / or Bluetooth Client Applications). Bluetooth serial ports are virtual COM ports and should be given a number that is distinct from the physical COM ports already present on the PC. Both the transceiver attached to the PC and the HygroLog NT should be given a unique COM port number.

In addition, it is usually necessary to “pair” the transceiver with the Bluetooth enabled HygroLog NT. For security reasons, the process of “pairing” requires a password to be entered. By default, all ROTRONIC Bluetooth devices are identified as HygroBlue and use 1234 as the pairing password.

In order to allow detection by HW4, any Bluetooth serial port must be manually declared in **HW4 Global Settings – [General Tab](#)** in the text box labeled **Bluetooth serial ports**, and then discovered manually using **Search for RS-232 masters** under [Devices and Groups](#) in the HW4 main menu bar.

HW4 detects and displays a Bluetooth device in the same manner as any device connected to the PC by way of a physical serial port.

5.1.3. USB port

Use of any device connected to a USB port requires installing the ROTRONIC USB driver prior to using. When this device is the only one that will be connected to the PC you should install the USB driver prior to starting HW4.

5.1.3.1. Installing the ROTRONIC USB driver

The ROTRONIC USB driver (ROTRONIC USB Option) is required only if you plan on using an instrument or docking station with USB port. This driver can be installed before or after installing HW4.

IMPORTANT: do not run HW4 while installing the USB driver on the PC.


The following example shows how to install the USB driver with Windows XP.

1. Prior to connecting the device to the PC, insert the HW4 CD-ROM in any available PC drive. The HW4 installation screen starts automatically. Exit this screen.
2. Connect the instrument or docking station to any available USB port. Upon detecting the presence of a device connected to the USB port, Windows XP automatically starts the following wizard:

If you are using Windows XP Professional SP2, the following screen appears. Select "No, not this time" and click on Next. This screen does not appear if you have not installed SP2.




3. The following screen appears:

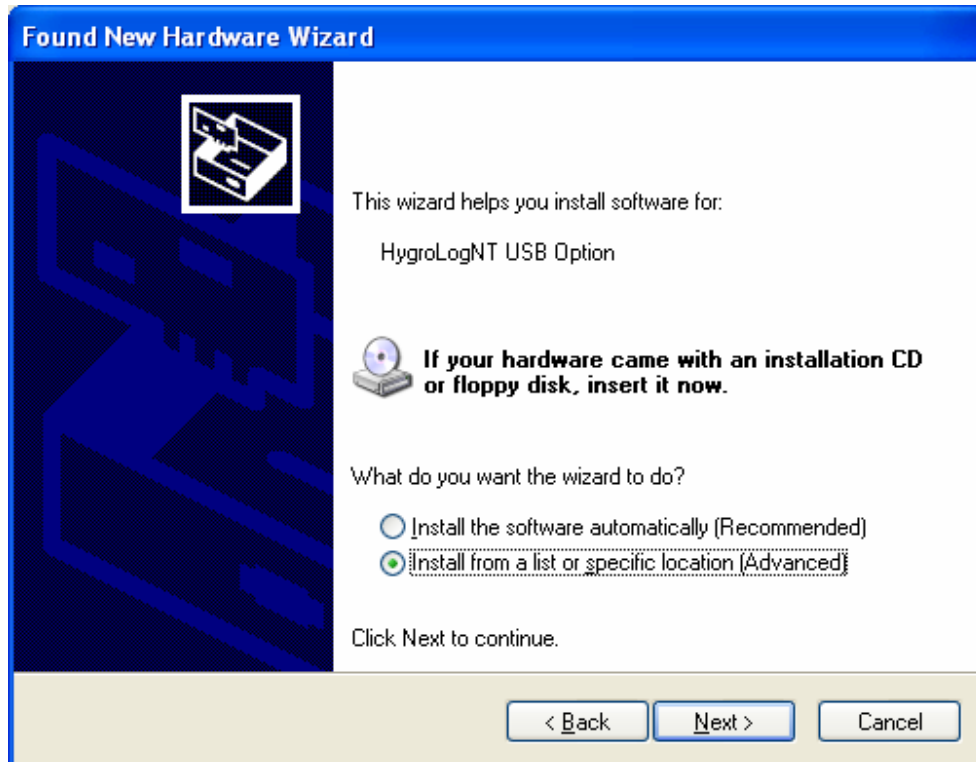

Select "Install from a list or specific location" and click on Next.

4. Select "Search for the best driver in these locations". Note that at this time the HW4 CD_ROM should have been already inserted in the PC.

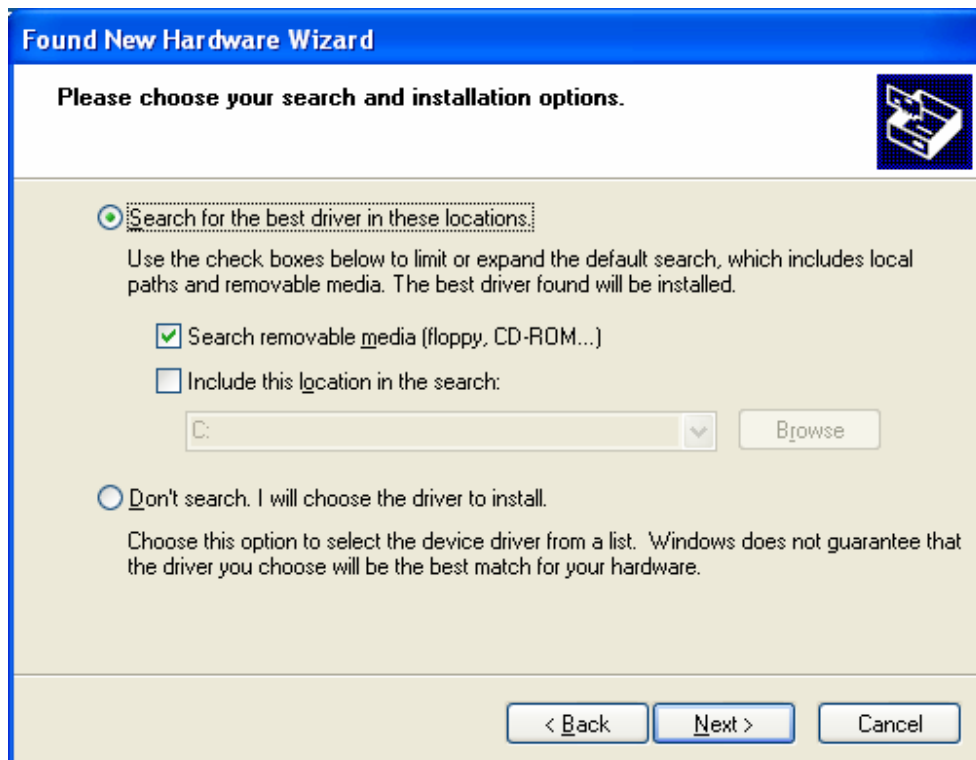


Click on Next

5. You will receive the following warning. Disregard this warning and click on “Continue anyway”.



6. Windows copies the USB driver located on the HW4 CD-ROM and displays the following message when done:

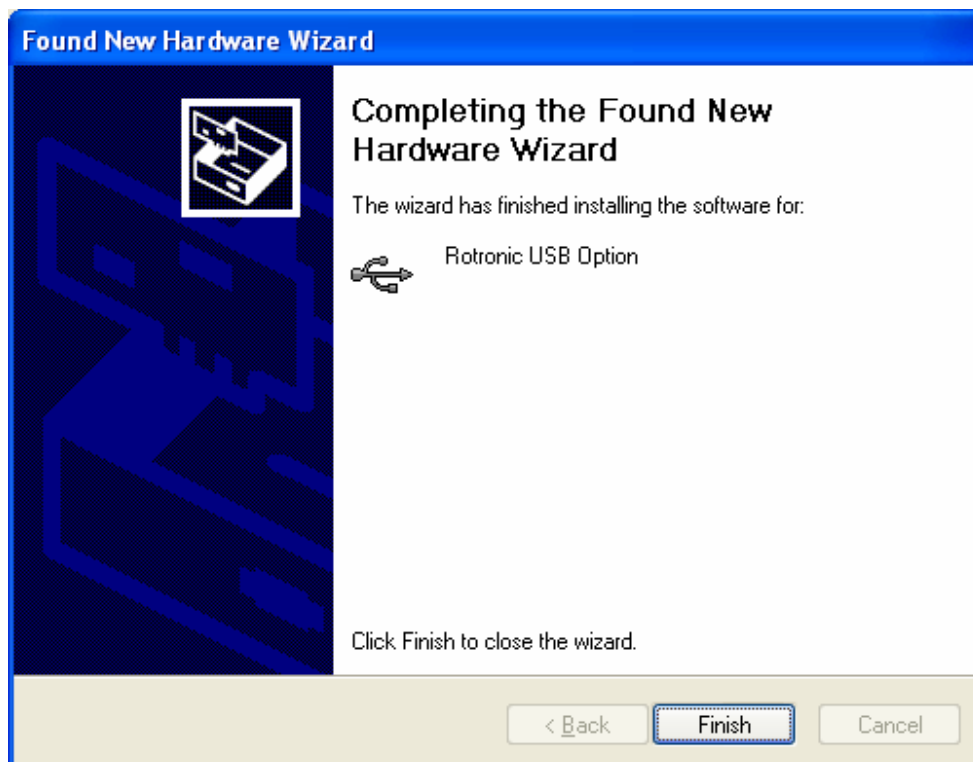

Click on “Finish” to complete the process.

Note: If your device has only a serial port and the PC has a USB port but no COM port, you may install a USB Serial Adapter on one of the PC USB ports. In this situation, install the driver supplied with the USB serial adapter (not the driver supplied with HW4) and connect the device to this adapter. HW4 will assume that the device is connected to a serial (COM) port.

5.1.4. LAN (TCP/IP) - wired connection

Usually this type of device requires some configuration to ensure compatibility with the LAN to which the PC is connected (IP address, sub-net mask, gateway, etc.). You may want to read “**Basic TCP/IP concepts**” in this manual or consult with your network administrator prior to connecting any device to your local area network. In particular, you should verify that the method used for the automatic configuration of TCP/IP settings by HW4 (see below) is compatible with your LAN.

5.1.4.1. Automatic configuration (new devices only)

Connect the ROTRONIC device to your LAN and open HW4 on a PC that is connected to your Ethernet LAN. Each device discovered by HW4 is configured based on the TCP/IP settings of the network interface card (NIC) present on the PC and is automatically given the same subnet mask and gateway address as the NIC. Each device is also given a static IP address equal to the IP address of the NIC, incremented by 1 to n (depending on the number of devices discovered by HW4). For example if the NIC has the IP address 192.168.57.100, two devices discovered by HW4 would get the IP addresses 192.168.57.101 and 192.168.57.102.

Note: the method used by HW4 for generating IP addresses may result in an address conflict on the LAN, especially if the NIC present on the PC is configured to automatically get an IP address from a DHCP server (dynamic IP address). In this situation, or if you are in doubt, you should manually configure each device as explained below. See also: [Device Protection](#)

5.1.4.2. Manual configuration

Before you begin, you should also make sure that the following is available to you:

- An unused static IP address that is compatible with your network
- The subnet mask of your network
- The local IP address of your wireless router

All ROTRONIC devices for LAN connection via cable use an internal server module from Digi International (Digi Connect ME).


A device discovery utility (**Digi Device Discovery**) is provided in the HW4 installation directory (usually, C:\Program Files\HW4). The name of the corresponding file is **dgdiscvr.exe**. The Digi Device Discovery utility can be started by double clicking with the mouse on the file.

Connect the device to be configured to your LAN and double click with the mouse on dgdiscvr.exe to start Digi Device Discovery.


When activated, Digi Device Discovery automatically detects any ROTRONIC device present on the LAN and provides a list of all such devices¹.

¹ With the possible exception of products from Digi International, third party TCP servers (or device servers) are generally not detected.


The following example shows the initial screen for one device.

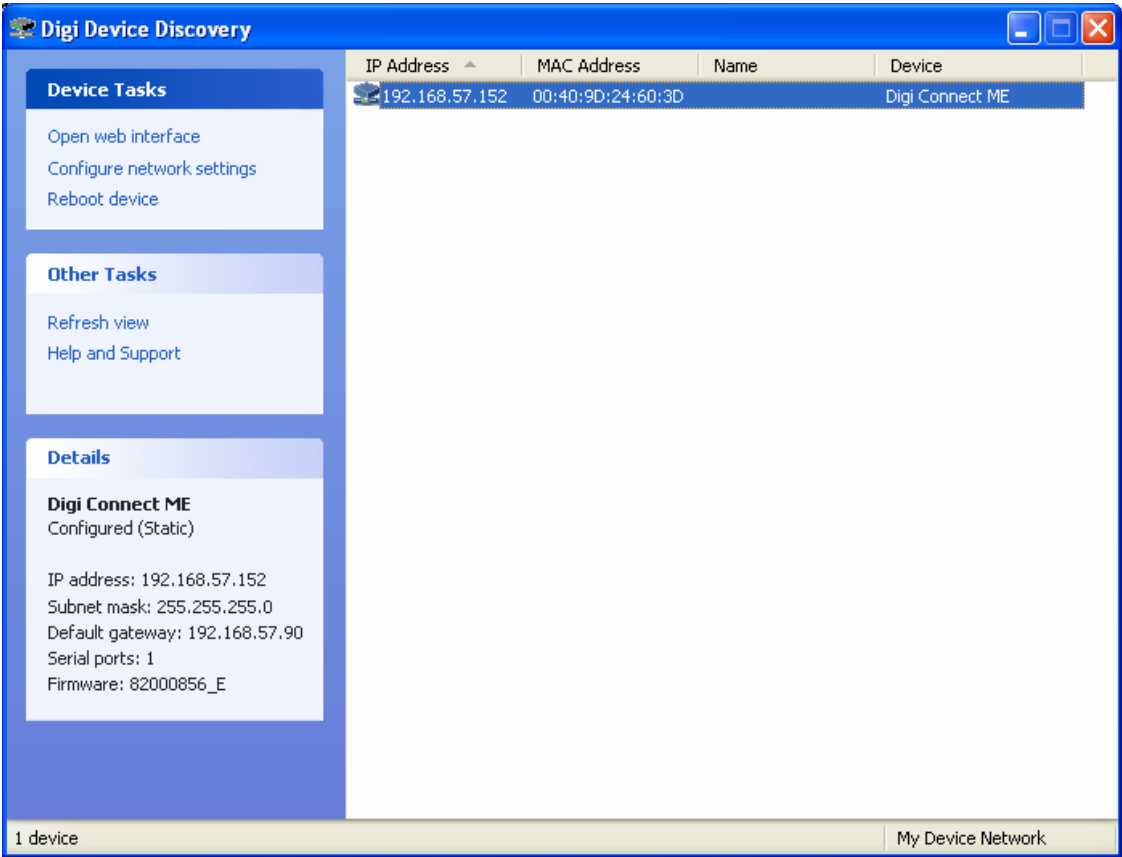

Any device with a default configuration that is not compatible with the LAN is listed with a warning in red (on the left side of the screen).

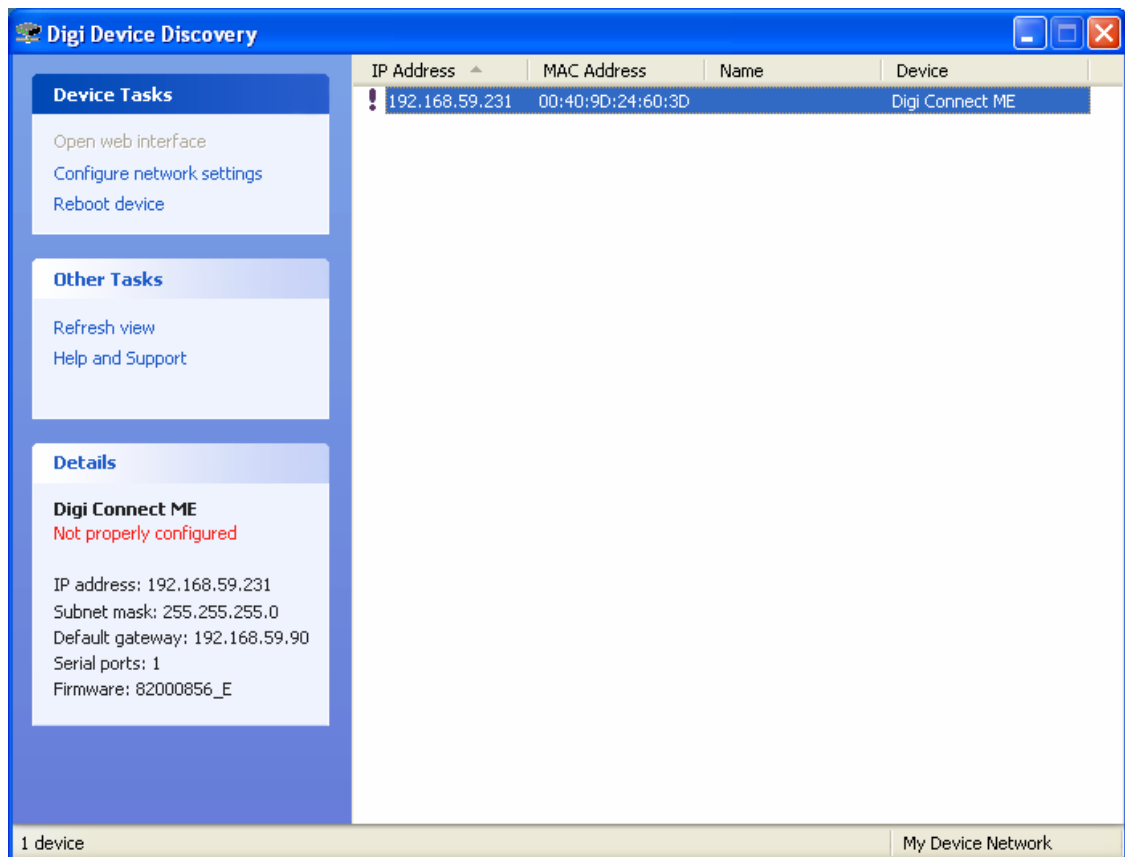


With the device highlighted, click on “Configure network settings” to open the dialog box shown below.

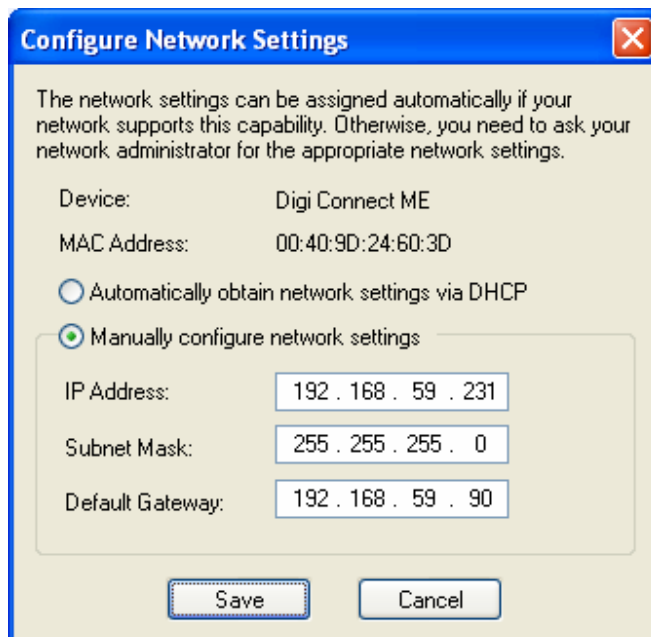

If in doubt as to which setting to enter manually, select “automatically obtain network settings via DHCP” (requires a LAN with DHCP server such as a router). Using DHCP ensures that the device settings that are compatible with the LAN (subnet mask and default gateway) and that the IP address does not generate any conflict.

Click on Save, Reboot the device and click on Refresh View.

A dynamic IP address is subject to change whenever the device is powered down and then powered up. For this reason, you should now change the device configuration from dynamic IP address (DHCP) to static IP address. To do this, click again on “Configure network settings”. The current settings of the device (IP address, subnet mask, etc.) appear in light gray.

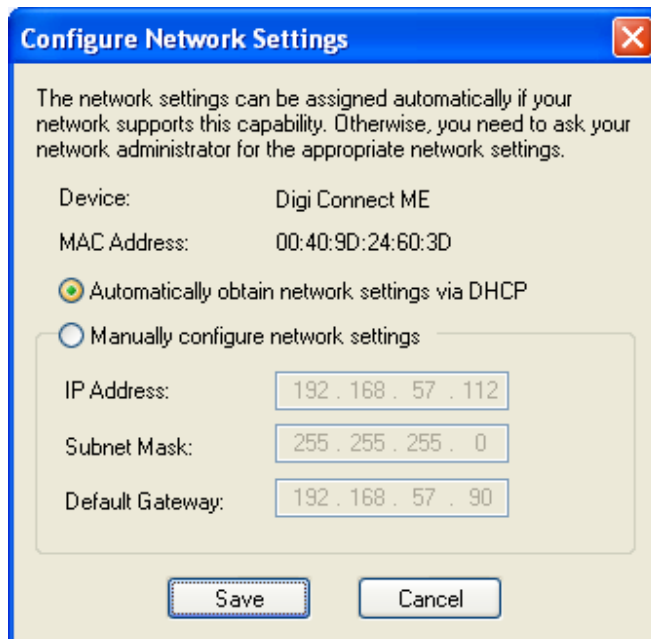

Configure Network Settings

The network settings can be assigned automatically if your network supports this capability. Otherwise, you need to ask your network administrator for the appropriate network settings.


Device: Digi Connect ME
MAC Address: 00:40:9D:24:60:3D

☒ Automatically obtain network settings via DHCP
☐ Manually configure network settings

IP Address: 192 . 168 . 57 . 112
Subnet Mask: 255 . 255 . 255 . 0
Default Gateway: 192 . 168 . 57 . 90

Save Cancel

Click on “Manually configure network settings”.

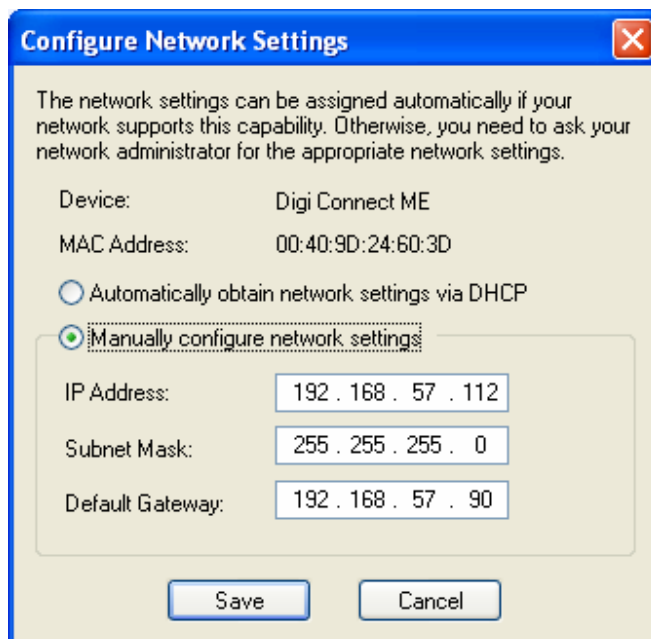

Configure Network Settings

The network settings can be assigned automatically if your network supports this capability. Otherwise, you need to ask your network administrator for the appropriate network settings.

Device: Digi Connect ME
MAC Address: 00:40:9D:24:60:3D

☐ Automatically obtain network settings via DHCP
☒ Manually configure network settings

IP Address: 192 . 168 . 57 . 112
Subnet Mask: 255 . 255 . 255 . 0
Default Gateway: 192 . 168 . 57 . 90

Save Cancel

Change the device IP address to an address that is both currently unused on the LAN and outside of the range of dynamic IP addresses used by the DHCP server. To select a proper address, you may want to consult your network administrator.

The Digi International modules feature a **web server** which can be accessed either directly from Internet Explorer (type in the IP address of the module and click on Go) or from Digi Device Discovery (right click on the device with the mouse and select "**Open web interface**").

When the device is used as the master device in a RS-485 multi-drop network, it may be necessary to use the web server to change the factory default baud rate because all devices in the RS-485 multi-drop must use the same baud rate.

See: [Changing the baud rate of an Ethernet device](#)

5.1.5. LAN (TCP/IP) - wireless connection

Note: the following instructions address the situation of a wireless network with a router. Instructions for establishing a peer to peer wireless connection (or Ad-Hoc connection) are not part of this manual.

5.1.5.1. Wireless security

It is highly advisable to protect a wireless network by requiring that an encoded key be provided by any device that tries to access the network. Deciding to activate network security is part of the configuration of a wireless router. Most wireless routers offer a choice of several types of network security. One such method consists in using WEP keys. WEP keys provide a high level of encryption and are easy to use with a Windows wireless connection.

Typically a wireless router can automatically generate up to four **WEP keys**, based on a single plain text word entered by the user during the router configuration. A WEP key consists of a sequence of characters. Only those characters used in the hexadecimal numerical notation are permitted: numbers from 0 to 9 and letters from A to F.

The WEP keys that are generated and used by the router should be either printed or written down for future configuration of whichever ROTRONIC devices will access the wireless network.

Additionally, a wireless router can be made to broadcast a network name (SSID) that is defined during configuration of the router. During configuration of a ROTRONIC wireless device the router SSID can also be used to point the device to the proper network.

5.1.5.2. Automatic configuration

HW4 does not automatically configure wireless network devices.

5.1.5.3. Manual configuration

Before you begin, you should also make sure that the following is available to you:

- An unused static IP address that is compatible with your network
- The subnet mask of your network
- The local IP address of your wireless router
- The wireless network name (SSID)
- The WEP key normally selected on your wireless router (one key out of four)


The ROTRONIC devices for wireless network connection use an internal server module from Digi International (Digi Connect WI-ME).

A device discovery utility (**Digi Device Discovery**) is provided in the HW4 installation directory (usually, C:\Program Files\HW4). The name of the corresponding file is **dgdiscvr.exe**. The Digi Device Discovery utility can be started by double clicking with the mouse on the file. When HW4 has already discovered at least one device, Digi Device Discovery can also be started directly from the **HW4 main menu bar, Settings and Tools, Ethernet Configuration Tool**.

IMPORTANT: prior to starting the configuration process, access your wireless router configuration and disable the wireless network security (WEP Keys). Failing to do so will prevent any communication with the ROTRONIC wireless device (network security will be re-established later).

To begin the configuration process, power up the ROTRONIC wireless network device. In the case of a HygroLog NT docking station, be sure that a logger is installed on the docking station.

Start Digi Device Discovery. Digi Device Discovery should automatically detect any ROTRONIC device (cable or wireless connection) present on the LAN and provides a list of all such devices. Note that devices with a cable connection use the device name **Digi Connect** while wireless devices use the name **Digi Connect Wi-ME**.

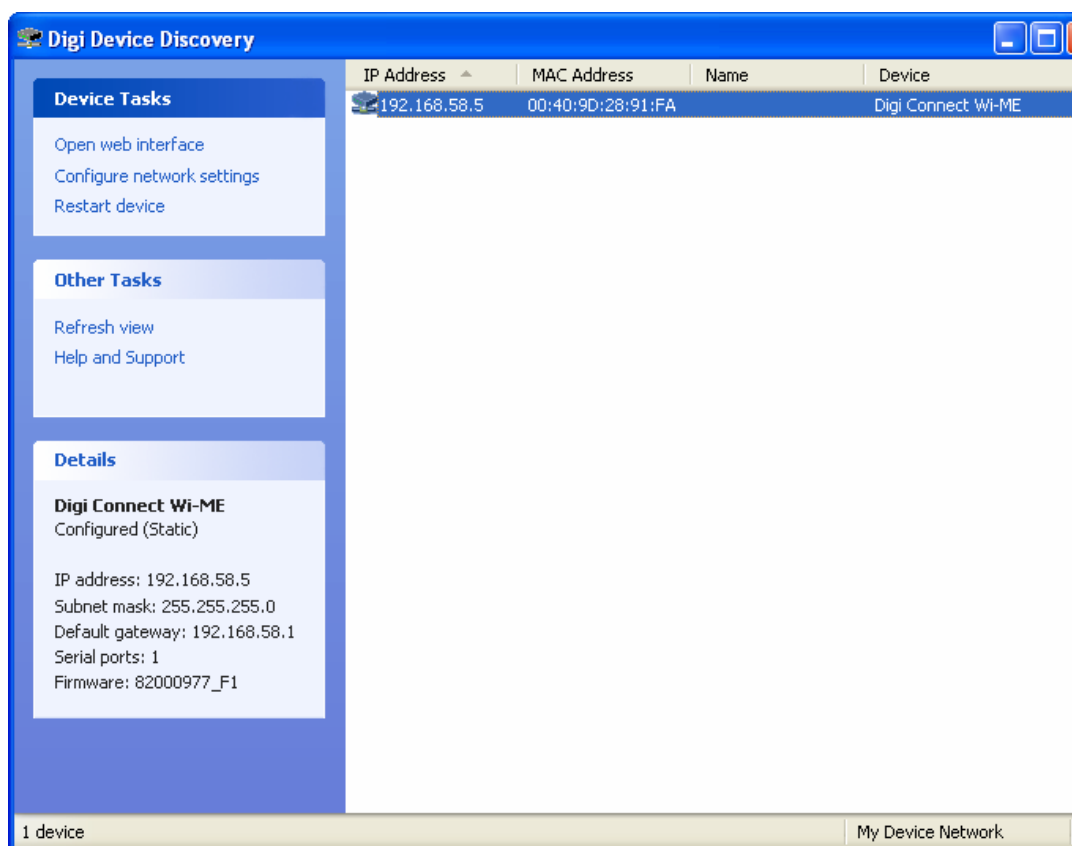


NOTE: Should Digi Device Discovery fail to detect your wireless device(s), use the HW4 main menu bar, Devices and Groups, Search for Ethernet Masters.

Because your device is not yet configured to be compatible with your wireless network, Digi Device Discovery should display the following warning message: Not properly configured.

Right click with the mouse on the device to be configured. A small menu appears. Select "Configure network settings".

The following dialog box opens:

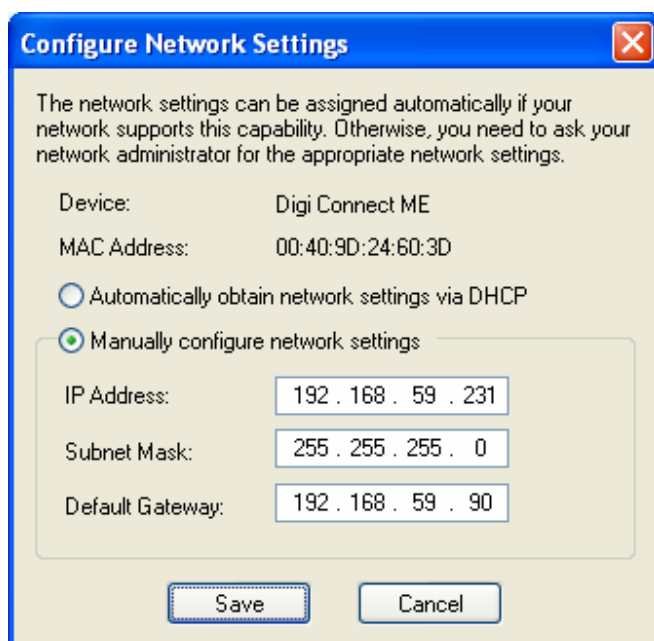

Select “Manually configure network settings” and enter an unused static IP address compatible with your wireless network, the network subnet mask and the IP address of the wireless router (Default Gateway).

Click on Save, Reboot the device (power down / power up) and click on Refresh View in Digi Device Discovery.

At this time, communication between the wireless device and HW4 should be possible after searching again for Ethernet masters.

The next steps are required to re-establish security of the wireless network.

The Digi International module features a web server which can be accessed either directly from Internet Explorer (type in the IP address of the module and click on Go) or from Digi Device Discovery (right click on the device with the mouse and select “Open web interface”). The web server can be used for configuring the module.

In the case of a wireless Digi module, a user name and password are required in order to access the home page of the module web server. The factory defaults from ROTRONIC are as follows:

- User name: **rotronic**
- Password : **wlan**

Note: should you have to reset to WI-ME module to the original manufacturer defaults, the user name reverts to **root** and the password **dbps**.

In the Digi module web server home page select **Network** under Configuration and **Wireless LAN Settings** in the Network page. The following page opens:

Digi Connect Wi-ME Configuration and Management - Microsoft Internet Explorer

File Edit View Favorites Tools Help

Back Forward Stop Home Search Favorites Reload Print Mail News RSS Feeds

Address http://192.168.58.5/config/network/wireless_config.htm

msn Search Highlight Options Pop-ups Blocked (1) My MSN

Digi Connectware™

Digi Connect Wi-ME Configuration and Management

Home

Configuration

- Network**
- Serial Ports
- GPIO
- Alarms
- System
- Remote Management
- Users

Management

- Serial Ports
- Connections

Administration

- File Management
- Backup/Restore
- Update Firmware
- Factory Default Settings
- System Information
- Reboot

Logout

Network Configuration

- IP Settings
- Wireless LAN Settings**
- Wireless Security Settings
- Wireless 802.1x Authentication Settings
- Network Services Settings
- Advanced Network Settings

Network name: (SSID)

☐ Connect to any available wireless network
☒ Connect to access point (infrastructure) networks only
☐ Connect to peer-to-peer (ad-hoc) networks only

Country:

Channel:

☐ Enable Short Preamble

Enter your wireless network name (SSID) to point the wireless module to that network. Select “Connect to access point (infrastructure) networks only”. Click on Apply.

Click on "Wireless Security Settings" (bottom of the page) to open the following page:

Digi Connect Wi-ME Configuration and Management - Microsoft Internet Explorer

File Edit View Favorites Tools Help

Back Forward Stop Reload Home Search Favorites Refresh Print Mail News RSS Feeds

Address http://192.168.58.5/config/network/wireless_security_config.htm

msn Search Highlight Options Pop-ups Blocked (1) My MSN

Digi Connectware™

Digi Connect Wi-ME Configuration and Management

Home

Configuration

- Network
- Serial Ports
- GPIO
- Alarms
- System
- Remote Management
- Users

Management

- Serial Ports
- Connections

Administration

- File Management
- Backup/Restore
- Update Firmware
- Factory Default Settings
- System Information
- Reboot

Logout

Network Configuration

- IP Settings
- Wireless LAN Settings
- Wireless Security Settings**

Network Authentication

☐ Use any available authentication method

☒ Use the following selected method(s):

- ☒ Open System
- ☐ Shared Key
- ☐ WEP with 802.1x authentication
- ☐ WPA with pre-shared key (WPA-PSK)
- ☐ WPA with 802.1x authentication
- ☐ Cisco LEAP

Data Encryption

☐ Use any available encryption method

☒ Use the following selected method(s):

- ☐ Open System (no encryption)
- ☒ WEP
- ☐ TKIP
- ☐ CCMP

WEP Keys

Transmit key: ☒ 1 ☐ 2 ☐ 3 ☐ 4

Encryption Keys:

1:	534F200B20
2:	919E080881
3:	1BEC8AADBC
4:	DA916859D9

Under Network Authentication, select “Use the following selected methods” and “Open Network”.

Under Data Encryption, select “Use the following selected method(s)” and WEP (assuming that your wireless network is set to use WEP. Under WEP keys, enter the value of each WEP key (for 64-bit encryption, each key has 20 hexadecimal characters). The WEP keys must be identical to those used by the wireless router.

Click on Apply and exit the module web server.

Access your wireless router configuration and re-establish wireless security and the WEP keys.

To change the User Name and Password used to open the web server:

Click on Configuration – Users

Digi Connect Wi-ME Configuration and Management

Home

Configuration

- Network
- Serial Ports
- GPIO
- Alarms
- System
- Remote Management
- Users**

Users Configuration

User Name	Action
rotronic	

New...

Click on the user name text box:

Digi Connect Wi-ME Configuration and Management

Home

Configuration

- Network
- Serial Ports
- GPIO
- Alarms
- System
- Remote Management
- Users

Management

- Serial Ports
- Connections

Administration

- File Management
- Backup/Restore
- Update Firmware

User Configuration - rotonic

User Configuration

User Name:

New Password:

Confirm Password:

[User Access](#)

[User Permissions](#)

Change the user name and password as desired. Click on Apply.

When the device is used as the master device in a RS-485 multi-drop network, it may be necessary to use the web server to change the factory default baud rate because all devices in the RS-485 multi-drop must use the same baud rate.

See: [Changing the baud rate of an Ethernet device](#)

6 CONNECTING AND DISCOVERING ADDITIONAL DEVICES

6.1 *Definition: masters and slaves (HW4 Professional)*

HW4 Professional makes use of the following definitions:

Master: any device / docking station that is directly connected either to a physical or to a virtual port of the PC or that is directly connected to a local area network (LAN) either by cable or by wireless.

Slave: any device / docking station that is directly or indirectly connected to a master by means of an RS-485 multi-drop. Slaves cannot be used in conjunction with HW4 Standard Edition.

6.2 *Masters*


To connect additional masters to the HW4 PC use the procedures described under “Connecting the first device”.

Additional masters must be discovered manually. Use the various search commands located in the HW4 main menu bar (see “HW4 Main Menu Bar, [Devices and Groups](#)”).

6.3 *RS-485 slaves (HW4 Professional)*


HW4 Professional allows the use of one RS-485 multi-drop with each master. Any RS-485 multi-drop is limited to a maximum of 64 devices (1 master and up to 63 slaves). All devices in the same multi-drop must be set to the same baud rate. **Slaves must be discovered manually.** Use the **Search for RS-485 slaves** command located in the HW4 main menu bar (see “HW4 Main Menu Bar, [Devices and Groups](#)”).

6.3.1. Slaves connected to a COM or USB master

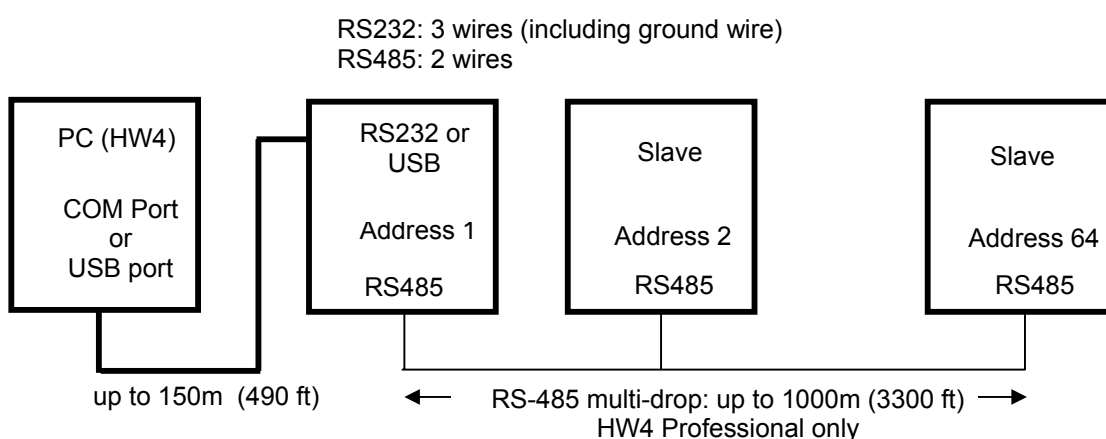

Note: the master can also be a Bluetooth device

6.3.2. Slaves connected to a LAN master

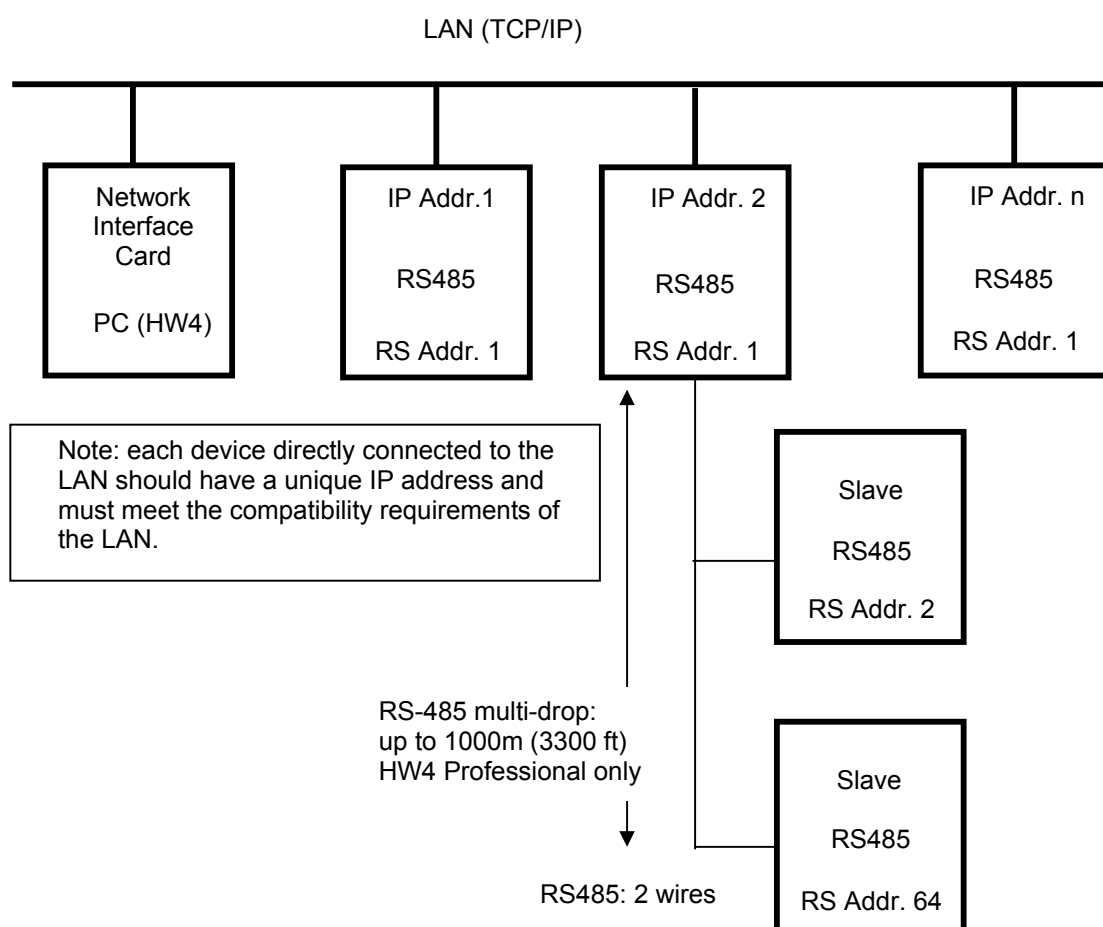

6.3.3. RS-485 address

All devices that are compatible with HW4 have an internal RS-485 address ranging from 00 to 64. The factory default for the RS-485 address is 00.

The RS-485 address is part of the communications protocol used by HW4 and is always included in the command sent to a device and in the device response. In the case of an RS-485 multi-drop, this address is used to identify each device and should be unique within each multi-drop.

When searching for RS-485 slaves, HW4 automatically changes the factory default RS-address of both masters and slaves as explained below:

- **Masters:** HW4 automatically changes the RS-485 address from 00 to 01 of any new device that is directly connected to the PC (master). As a result, all masters end up having the same RS-485 network address (01).
- **Slaves:** in the situation where an RS-485 multi-drop is detected during a search, HW4 automatically changes the address of each slave with address 00 to a unique address ranging from 02 to 64. The same address range (02 to 64) is used again when there is more than one RS-485 multi-drop.

IMPORTANT:

- With the exception of the two situations described below, do not change manually the factory default RS address (00). Allow HW4 to change this address automatically.

Change manually the RS address to 00 in the following situations:

- A master will now be used as a slave
- A slave is being moved from one multi-drop to another multi-drop.

The change of address should be done prior to changing the connection (use Device Manager). If need be, connect the devices as masters to do the address change.


Using an instrument with an address other than 00 in an RS-485 multi-drop may prevent HW4 from detecting the instrument. When using the command "Search for RS-485 slaves", HW4 will automatically change the RS address of the master to 01 and the address of each slave to 02 or higher.

- Within an RS-485 multi-drop, do not use duplicate RS addresses (with the exception of temporary address 00). Duplicate addresses will give unpredictable results in the data transmission. It is OK to use the same RS address in different multi-drops.
- All devices within a multi-drop should use the same baud rate

7 HW4 MAIN SCREEN

The HW4 main screen varies slightly depending on which edition of HW4 is running. The following description applies to the Professional edition which has more menu items and screens than the standard edition.


The HW4 main screen is subdivided into a left pane and a right pane. The width of each pane can be adjusted with the mouse. The height of the two sub-panes located in the right pane (table and graph) can also be adjusted with the mouse. To change the size of a pane, go over the separation line with the mouse cursor. When the cursor changes, click and hold with the mouse. Drag the separation line where you want it to be.

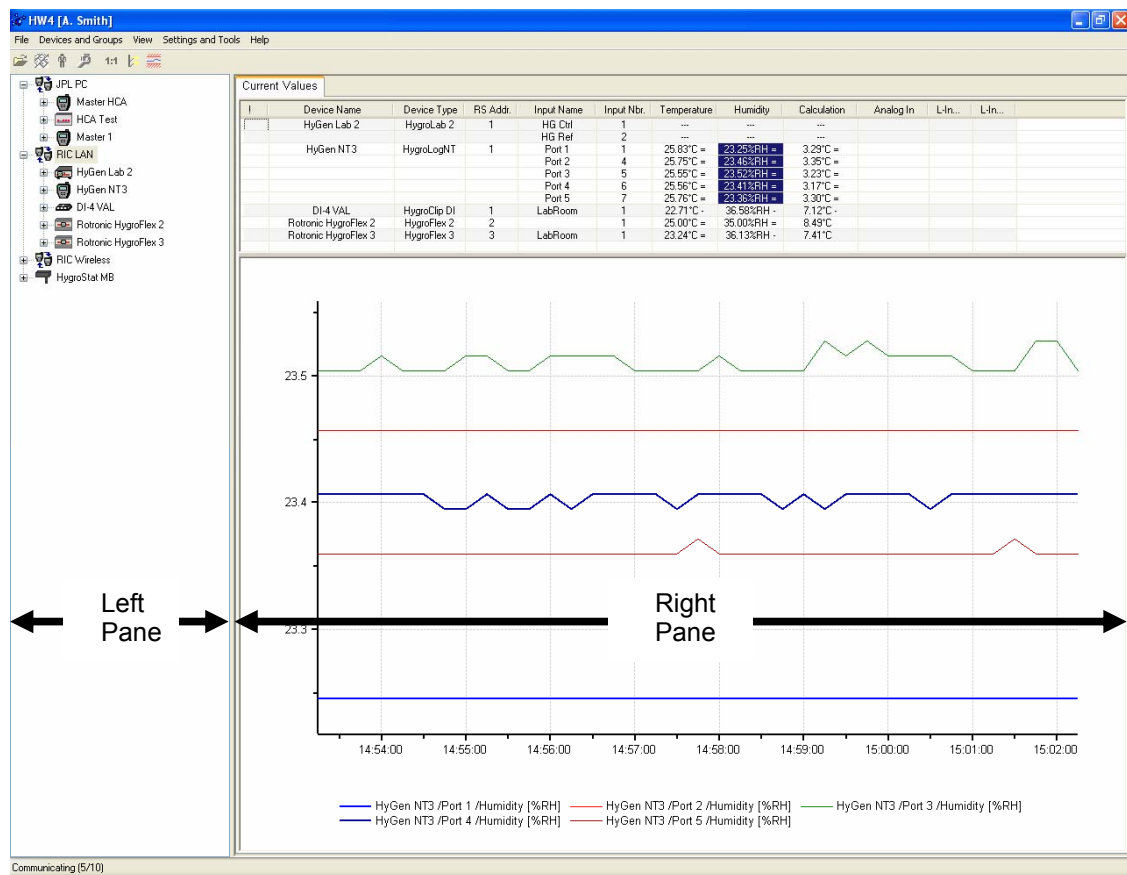

7.1 Device Tree (left pane)

7.1.1. Tree populated with individual devices

All versions of HW4 display on the left pane of the screen a tree listing the devices that have been discovered by HW4.

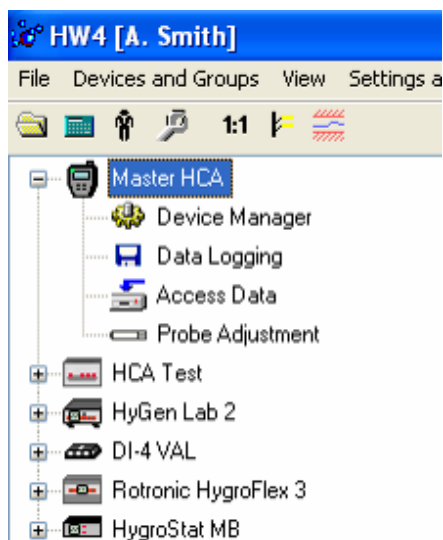

Each individual device is displayed in the left pane (for example: Master HCA, HCA Test, etc.). This is the default view for all versions of HW4 and the only view available with HW4 Standard edition

Functions available for each device can be seen when the tree is expanded. These functions are specific of each type of device.


To select a device or a function, click on it with the left mouse button.

7.1.2. Tree populated with device group(s) - HW4 Professional

A device group is a logical grouping of devices (as opposed to a physical grouping) created and defined by the user. Both masters and slaves can be part of a group, with no restriction other than not duplicating devices. Individual devices can be moved to any existing group or removed from the group (see **"Devices and Groups"**).

In addition to facilitating network management within HW4, using device groups makes it possible to simultaneously display selected data from several individual devices within a group. The data can be viewed either as a table or as a live graph.


By default HW4 displays a device group as opposed to displaying each individual device in the group. Individual devices can be displayed by expanding the group.

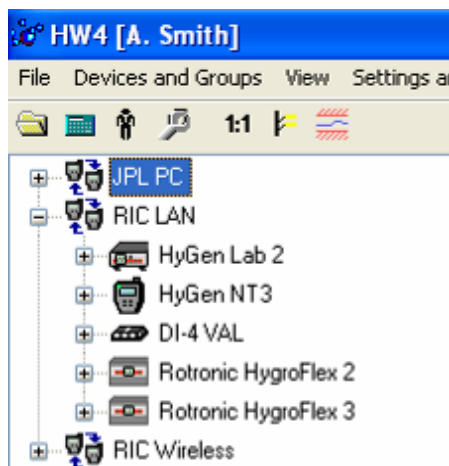

HW4 Professional edition allows the creation of device groups. When one or more groups have been created (example JPL PC, RIC LAN and RIC Wireless), and all devices have been moved to a group, the left pane displays a tree listing only groups.

A group can be expanded to show the devices within the group and each device can be expanded to show its functions.

To select a group, a device or a function, click on it with the left mouse button.

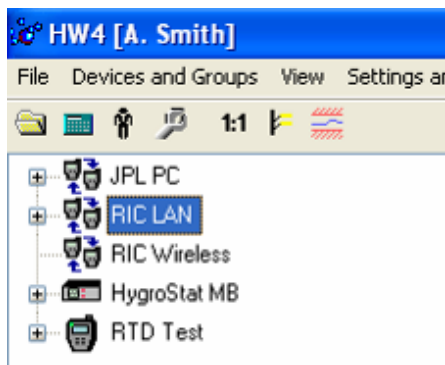


The left pane can also display one or several device groups (JPL PC, RIC LAN and RIC Wireless) as well as individual devices that are not part of any group (HygroStat MB, RTD Test).

Note: after being detected, devices appear individually in the tree regardless of the manner in which they are connected to the PC (master or slave).

7.1.3. Instrument group: renaming and changing the contents

The following menu can be opened by right clicking on a group:

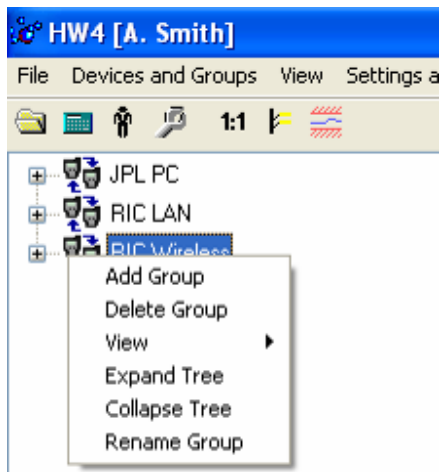

Rename Group: after creating a group, the group can be renamed.


Adding a device to a group: left click on the device and drag it on top of the group. Instruments within a group are displayed in the reverse order in which they were added to the group.

Removing a device from a group: to remove a device from a group, left click on the device and drag it outside of the group.

7.2 Right pane in Device View mode

The right pane is in the **Device View mode** whenever an individual device is selected in the device tree (left pane).

In the Device View mode, the right pane is as illustrated below. To select any of the available tabs, click on the corresponding label with the mouse left button.

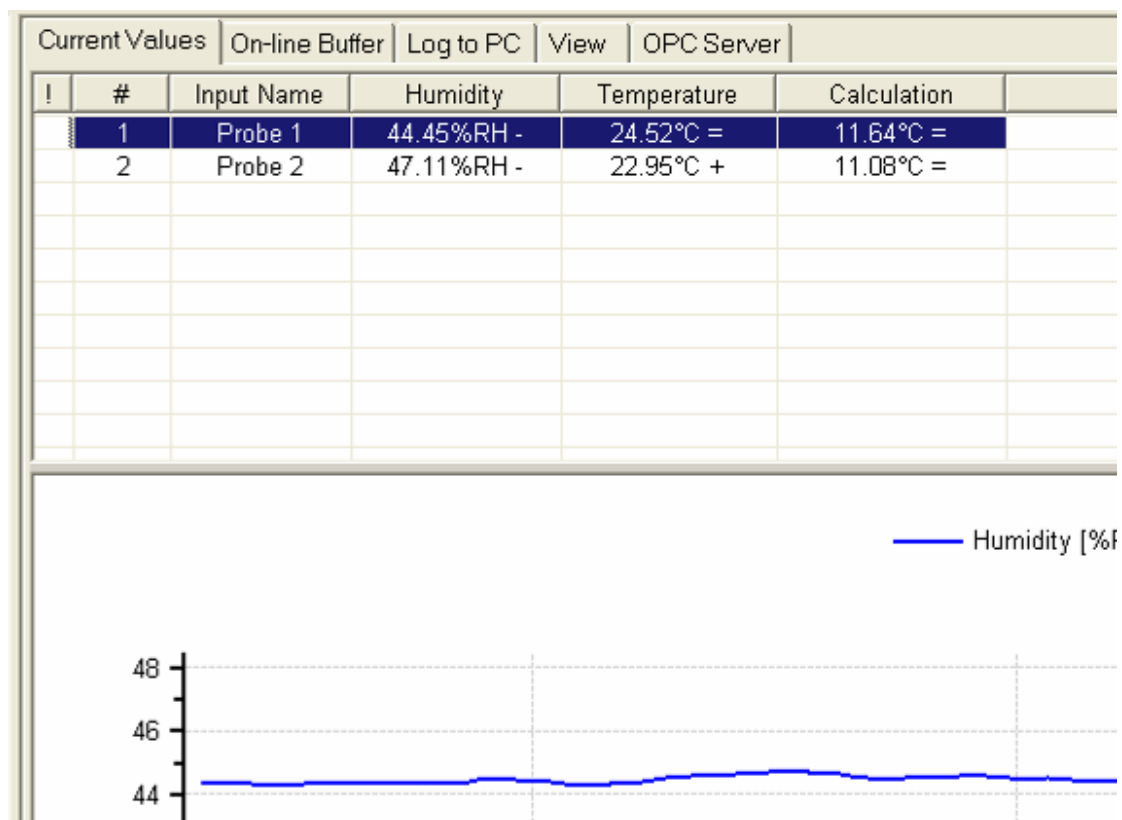

7.2.1. Current Values tab

With this tab selected, the top of the right pane displays a table of the most recent value for the different parameters associated with the device currently selected in the device tree. The values are constantly updated.

The table shows data from the different probe inputs of the device, as well as other information. Click and drag a column header with the mouse to change in which order fields are being displayed.

Current Values On-line Buffer Log to PC View OPC Server						
I	#	Input Name	Humidity	Temperature	Calculation	
	1	Probe 1	44.45%RH -	24.52°C =	11.64°C =	
	2	Probe 2	47.11%RH -	22.95°C +	11.08°C =	

The table displayed in the Current Values tab can be customized.

7.2.1.1. Customizing the Current Values tab

Data columns in the Current Values tab


Depending both on the type of device and on the configuration of the device, HW4 can display a number of data elements in the Current Values tab. Different data appear in different columns of the Current Values tab as explained in the following table:

Data	Column Headers	Notes
Humidity digital	Humidity Value 1	HygroClip probe digital output
Humidity analog	Value 1	HygroClip probe analog output (channel 1)
Temperature digital	Temperature Value 2	HygroClip probe digital output
Temperature analog	N/A	HW4 reads only single channel analog probes
Calculation	Calculation	Dew point, humidity ratio, etc.
Analog probe (single channel)	Analog Input Value 1	3d party analog probe
Pressure probe	Analog Input Value 1	3d party analog pressure probe
Custom Calculation	Custom Calc. Value 1	Requires an instrument that can be configured to perform a user defined calculation such as the difference between two temperatures.
Contact 1	Contact 1 Value 1	Logical input used to monitor an external contact
Contact 2	Contact 2 Value 2	Logical input used to monitor an external contact
Relay 1	Relay Value 1	Status of relay output 1 (energized or de-energized)
Relay 2	Relay Value 1	Status of relay output 2 (energized or de-energized)


Note: the columns "Value 1" and "Value 2" are provided for those users who wish to limit the number of columns displayed by HW4.

Column selection in the Current Values tab

To make HW4 use only specific columns in the Current Values tab, select any individual device present in the device tree (for example HyGen Lab 2 in the left pane). In the HW4 main menu bar, click on **Settings and Tools** and on **Column Headers**.

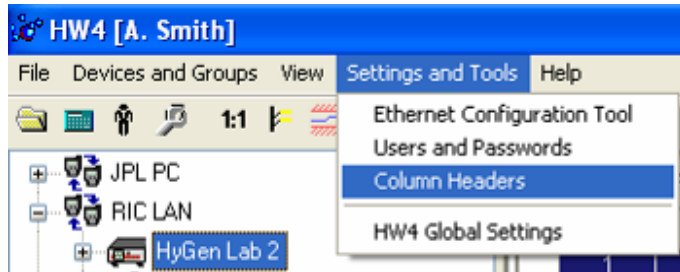

HW4 opens the following form. Make your selection using the mouse. **The selection applies globally to all devices displayed in the Current Values tab.**

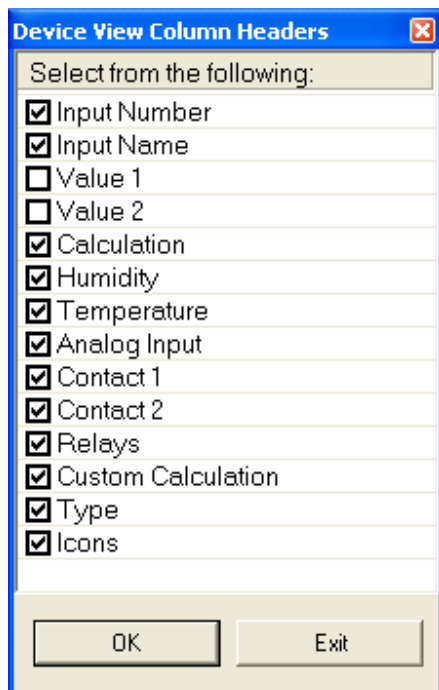

Order of the columns in the Current Values tab

The order of the columns in the Current Values tab can be changed. To move a column to the right or to the left, left click with the mouse on the column header and drag it to its new position.

Row selection in the Current Values tab

Select a device present in the device tree and use **View Tab** to select the rows to be displayed in the Current Values tab. The row selection applies only to the device currently selected in the device tree.

Current Values | On-line Buffer | Log to PC | **View** | OPC Server

Select here which items (inputs or values) to display for this device

- ☒ Input 1
- ☒ Input 2
- ☒ Input 3
- ☒ Input 4
- ☒ Input 5
- ☒ Input 6
- ☒ Input 7
- ☒ External Contact
- ☐ Custom Calculation
- ☐ Relay 1
- ☐ Relay 2
- ☐ Relay 3
- ☐ Relay 4


NOTES:

- The selections made in the View Tab affect the Current Values Tab in the Group View mode.
- The selections made in the View Tab also affect the following tabs and forms: On-line buffer tab, Log to PC tab, Device Manager form, Data Logging form and Probe Adjustment form.

If you do not see an input, a contact, etc. in any of these tabs and forms, remember to check the selections made in the View Tab.

7.2.2. Online graph

By default, the bottom of the right pane displays data in graphic form. The graph is constantly updated as new data is being received and retained in the on-line buffer. The scale of the graph is automatically selected by HW4.

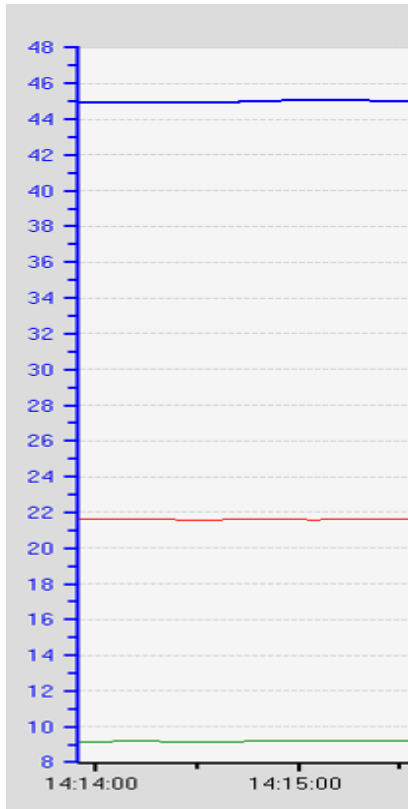


The following selections can be made with the mouse:

- To display all the parameters associated with a probe, left click on the desired probe number in the table located above the graph.
- To display only one probe parameter at a time, left click on the desired parameter in the table.
- To display only one parameter for all probes, left click on the header of the parameter.

If so desired, the graph can be hidden so as to provide more room for the probe data table (see HW4 Global Settings - View Tab).

The appearance of the graph can be customized. For details see “**Settings and Tools - HW4 Global Settings – Graph Settings tab**”.

To adjust the height of the data table and the height of the graph, bring the mouse cursor over the horizontal edge separating the graph and the table. When the mouse cursor changes, left click, hold and drag to the desired location.

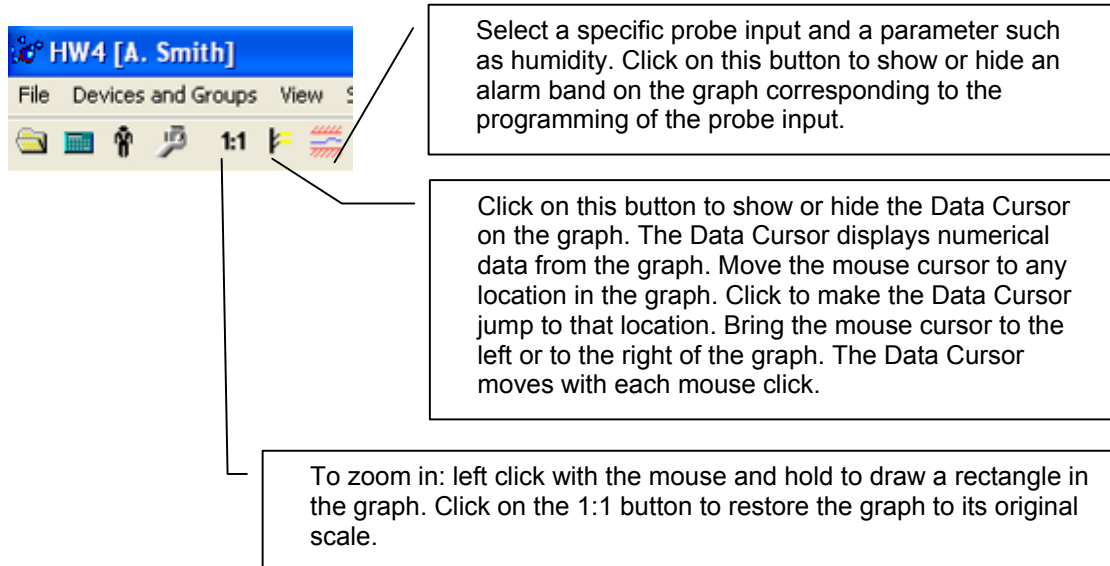

7.2.3. Online Buffer tab

With this tab selected, the right pane displays a table that shows the most recent historical data (limited to humidity, temperature and the calculated parameter). The time interval used to capture data to the on-line buffer is the value specified in **HW4 Global Settings – General tab**. This tab is also used to specify the maximum number of lines in the buffer. When the table is full, the oldest data are dumped from the bottom of the table as the most recent data is added to the top of the table.

The probe input to be displayed is selected from the text box located on top of the data table (left click on the arrow to the right of the box to display a list of the inputs that can be viewed – this list depends on the selections made in the View tab).

The on-line graph is based on the contents of the on-line buffer.

Current Values On-line Buffer Log to PC View OPC Server				
Input 1				
Date	Time	Humidity %RH	Temperature °C	Calculation [Dp] °C
10/6/2005	12:25:45	60.70	24.80	16.70
10/6/2005	12:25:30	60.94	24.80	16.76
10/6/2005	12:25:15	60.48	24.77	16.62
10/6/2005	12:25:00	60.39	24.80	16.62
10/6/2005	12:24:45	60.50	24.75	16.59
10/6/2005	12:24:30	60.86	24.71	16.66
10/6/2005	12:24:15	60.55	24.75	16.61
10/6/2005	12:24:00	60.54	24.75	16.60
10/6/2005	12:23:45	60.26	24.69	16.48
10/6/2005	12:23:30	60.19	24.70	16.47

Note: the contents of the on-line buffer are lost upon exiting HW4.

7.2.4. Log to PC tab

Current Values On-line Buffer Log to PC View OPC Server					
Input	Description	Log Interval	Log File Management	Log File Name	File Location
☐ 1	Probe 1	00:00:15	No file size limit		
☐ 2	Probe 2	00:00:15	No file size limit		
☐ 3	Probe 3	00:00:15	No file size limit		
☐ 4	Probe 4	00:00:15	No file size limit		
☐ 5	Probe 5	00:00:15	No file size limit		
☐ 6	Probe 6	00:00:15	No file size limit		
☐ 7	Probe 7	00:00:15	No file size limit		
☐ External Contacts	Contacts 1/2	00:00:15	No file size limit		
☐ Custom Calc.	Custom Calc.	00:00:15	No file size limit		
☐ Relay					

The Log to PC tab is used to record data directly on the PC from any type of instrument that is selected in the Device Tree. All input types recognized by HW4 are listed in this tab, whether they exist or not on the actual device. When the list of inputs is either empty or incomplete, please go to the **View Tab** and check the selections (see Right Pane: Device View Mode). Only those items selected in the View Tab are visible in the Log to PC tab.

The log interval and the logging mode (for example one file per day) are entered from within this tab. The data file name, its path, and the type of the file (protected or unprotected) are also entered here. Data from more than one input can be recorded, using a different file for each individual input.

Instructions for using this tab are provided in this manual under [Data Logging – directly on the PC](#).

HW4 Professional: when a HygroLog NT is selected in the device tree is a HygroLog NT, the bottom part of the Log to PC tab is used to enable the automatic downloading of log files to the PC:

Automatic Log File Download

Status Disabled

Download Interval

Download Path

Downloaded Files

Last Download

Next Download

Instructions for using this function are provided in this manual under:

[Automatic Log File Download – HygroLog NT](#)

7.2.5. View tab

The View tab is used to select which items (inputs, custom calculation, and status of optional relays) are displayed by HW4 in the Current View, On-line buffer and Log to PC tabs.

Use the mouse to select the fields that correspond to your actual use of the various device inputs.

IMPORTANT: the selections made in this tab also affect what can be seen in the following forms: Device Manager, Data Logging and Probe Adjustment.

7.2.6. OPC Server tab

This tab is available only with HW4 Professional with OPC server. Unless your HW4 product enables this edition of HW4, you cannot use the OPC feature of HW4. With OPC server enabled, HW4 can communicate with any OPC client application to the purpose of bringing together data from HW4 and data from other sources. The OPC client application needs to be configured / programmed to read or write to the OPC tags provided by HW4.

Current Values	On-line Buffer	Log to PC	View	OPC Server
Enabled Tag ...	Tag Group Name	Tag Group Description		
☐ 1	DeviceInfo	Device information tags		
☐ 2	MeasuringInfos	Data tags		
☐ 3	Alert	Alarm tags (device condition, out-of-limits values and s		
☐ 4	File	Log file tags		

A full description of the HW4 OPC server is available in the following document, available from ROTRONIC: **E_OPC-Tags_v13.doc**

7.2.7. AwQuick / AwE Mode tab

This tab is available only with HW4 Professional with AwQuick. Unless your HW4 product key enables this edition of HW4, you cannot use the water activity functions of HW4.

HW4 features two modes for measuring water activity:

- **AwE mode:**

In this mode HW4 monitors the stability of both temperature and humidity. The measurement is automatically ended as soon as both humidity and temperature reach equilibrium. The natural (or static) equilibration of most products typically requires from 45 to 60 minutes and can take as long as a couple of hours.


The criteria used to determine that temperature and humidity have reached equilibrium can be modified by the user.

Current Values		Online Buffer	Log to PC	View	AwE Mode			
Input	Input Name	Probe S/N	Humidity [Aw]	Temperature [°C]	Elapsed Time	Water Activity [Aw]	Status	
☒ 1	Port 1		0.170 <=>	24.65 <=>	00:00:00		Running	
☐ 4	Port 2		0.171	24.58				
☐ 5	Port 3		0.172	24.47				
☐ 6	Port 4		0.171	24.41				
☐ 7	Port 5		0.170	24.61				
	</							

7.3 Right pane in Group View mode

Note: device groups are available only HW4 Professional edition.

The right pane is in the **Group View mode** whenever a group of devices is selected in the device tree (left pane). In the Group View mode, the right pane is as illustrated below.

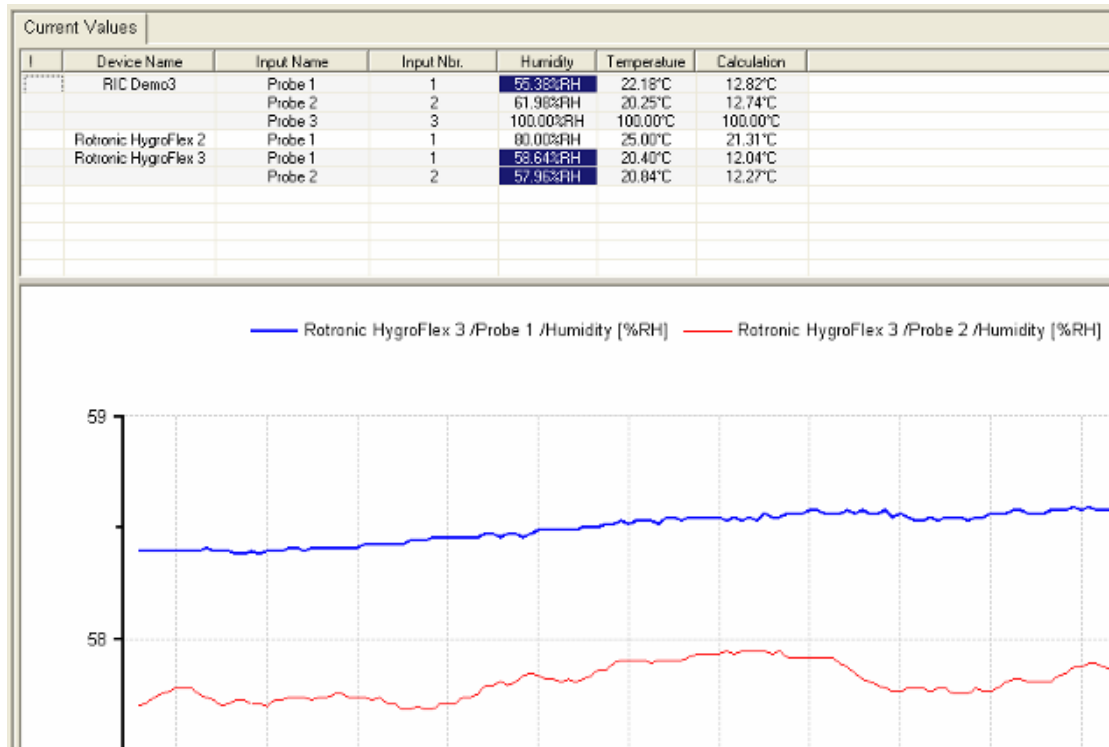

7.3.1. Current Values tab

The top of the right pane displays a table of the most recent value for the different parameters associated with the devices in the group selected in the device tree. The values are constantly updated.

7.3.1.1. Customizing the Current Values tab

Data columns in the Current Values tab

Depending both on the type of device and on the configuration of the device, HW4 can display a number of data elements in the Current Values tab. Different data appear in different columns of the Current Values tab as explained in the following table:

Data	Column Headers	Notes
Humidity digital	Humidity Value 1	HygroClip probe digital output
Humidity analog	Value 1	HygroClip probe analog output (channel 1)
Temperature digital	Temperature Value 2	HygroClip probe digital output
Temperature analog	N/A	HW4 reads only single channel analog probes
Calculation	Calculation	Dew point, humidity ratio, etc.
Analog probe (single channel)	Analog In Value 1	3d party analog probe
Pressure probe	Analog In Value 1	3d party analog pressure probe
Custom Calculation	Custom Calc. Value 1	Requires an instrument that can be configured to perform a user defined calculation such as the difference between two temperatures.
Contact 1	Contact 1 Value 1	Logical input used to monitor an external contact
Contact 2	Contact 2 Value 2	Logical input used to monitor an external contact


Note: the columns “Value 1” and “Value 2” are provided for those users who wish to limit the number of columns used by HW4.

Column selection in the Current Values tab

To make HW4 use only specific columns in the Current Values tab, select a device group that is present in the device tree (for example RIC LAN in the left pane). In the HW4 main menu bar, click on **Settings and Tools** and on **Column Headers**.



HW4 opens the following form. Make your selection using the mouse. **The selection applies globally to all groups displayed in the Current Values tab.**

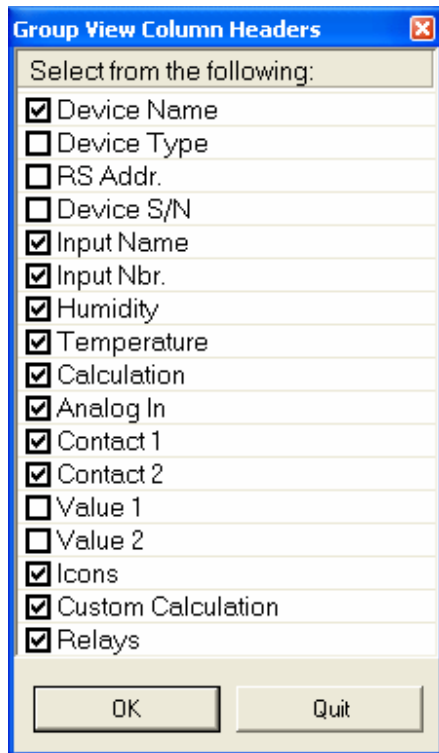

The image shows a dialog box titled "Group View Column Headers" with a close button (X) in the top right corner. Inside the dialog, there is a label "Select from the following:" followed by a list of 18 items, each with a checkbox. The items are: Device Name (checked), Device Type (unchecked), RS Addr. (unchecked), Device S/N (unchecked), Input Name (checked), Input Nbr. (checked), Humidity (checked), Temperature (checked), Calculation (checked), Analog In (checked), Contact 1 (checked), Contact 2 (checked), Value 1 (unchecked), Value 2 (unchecked), Icons (checked), Custom Calculation (checked), and Relays (checked). At the bottom of the dialog, there are two buttons: "OK" and "Quit".

Item	Selected
Device Name	Yes
Device Type	No
RS Addr.	No
Device S/N	No
Input Name	Yes
Input Nbr.	Yes
Humidity	Yes
Temperature	Yes
Calculation	Yes
Analog In	Yes
Contact 1	Yes
Contact 2	Yes
Value 1	No
Value 2	No
Icons	Yes
Custom Calculation	Yes
Relays	Yes

Order of the columns in the Current Values tab

The order of the columns in the Current Values tab can be changed. To move a column to the right or to the left, left click with the mouse on the column header and drag it to its new position.

Row selection in the Current Values tab


The rows displayed in the current values table depend on the selections made for each individual device. See **"Right pane in Device View mode - View Tab"**.

7.3.2. On-line graph

One of the reasons for creating an instrument group is the ability to display a graph showing data from different instruments within the group.

For example, to display relative humidity for all instruments and probes within a group, click on the Humidity header in the data table. To select other data for the graph, press the CONTROL key and click on any desired parameter(s) in the data table.

The graph is constantly updated as new data is being received and retained in the on-line buffer. The scale of the graph is automatically selected by HW4.

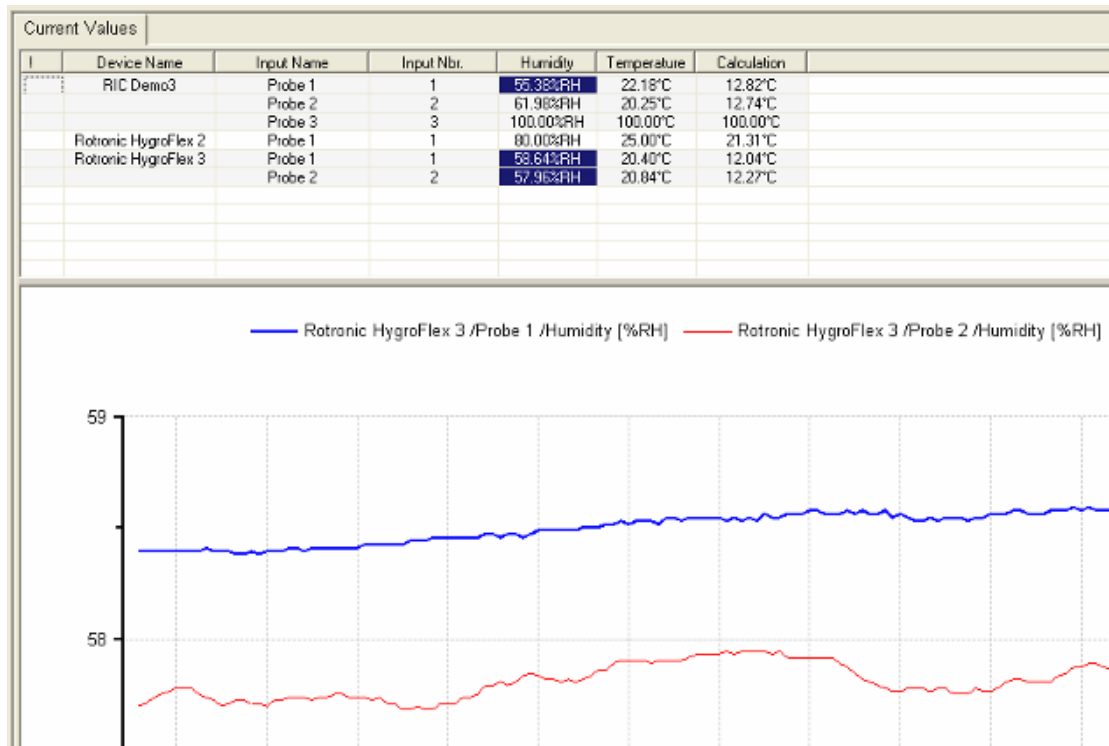

The appearance of the graph can be customized. For details see “**Settings and Tools - HW4 Global Settings – Graph Settings tab**”.

8 HW4 MAIN MENU BAR

The Main Menu Bar is located at the top of the HW4 main screen.

8.1 *File*

- **Open**

Used to open any file present on the PC such as a log file or an event file without having to open any individual instrument in the device tree (see also HygroLog NT Functions - Access Data). By default File - Open points to the default directory used by HW4 when creating files during the initial start-up (first time use): C:\Documents and Settings\Windows User\Application Data\ROTRONIC_HW4.

- **Send HW4 to the Windows notification area**

Sends / minimizes HW4 as a shortcut icon to the Windows notification area. This area is located on the taskbar, immediately to the left of the clock. HW4 keeps running.

- **Exit** exit HW4

Note: if you attempt to exit HW4 while Log to PC is active, you will get a warning from HW4.

8.2 *Devices and Groups*


- **Search for USB Masters** - HW4 searches for devices directly connected to a USB port
- **Search for RS-232 Masters** - HW4 searches for devices directly connected to either a physical COM port or to a Bluetooth virtual serial port.
- **Search for Ethernet Masters** - HW4 searches for devices directly connected to the same local area network as the PC (cable and / or wireless connection).
- **Search for Specific IP Address** - HW4 searches for a master device at a specific IP address specified by the user.



- **Search for RS-485 Slaves**¹ - HW4 searches for RS-485 slaves attached to any master that has been already discovered
- **Add Group**¹ - creates a new empty group of instruments in the left pane tree (see below)
- **Delete Group**¹ - deletes the instrument group that is currently selected in the left pane. The instruments from the group are not deleted.

¹ available only with HW4 Professional edition


- **Add Device** – allows the user to manually add a device to the device tree. This saves time in some situations. All information required to fill-in the form shown below must be known beforehand.

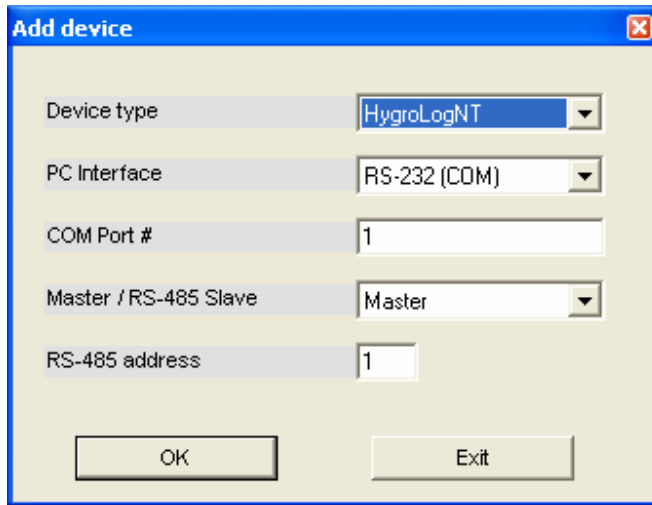

The 'Add device' dialog box has a blue title bar with a close button. The main area is light beige. It contains five rows of controls: 'Device type' with a dropdown menu showing 'HygroLogNT'; 'PC Interface' with a dropdown menu showing 'RS-232 (COM)'; 'COM Port #' with a text box containing '1'; 'Master / RS-485 Slave' with a dropdown menu showing 'Master'; and 'RS-485 address' with a text box containing '1'. At the bottom are two buttons: 'OK' and 'Exit'.

Click with the mouse on the arrow located to the right of a box to display a drop down menu of available choices. Note that if the device is not connected or if the wrong information is entered, HW4 will put a red cross over the device icon in the device tree.

- **Delete Device** - deletes the selected device from the left pane tree (this command is grayed out when no device is selected). Use this menu item after device has been disconnected
- **Delete all Devices from the Tree** – globally deletes all devices from the left pane tree. This function should be used with caution...
- **Connection over Modem** - HW4 allows one connection over a modem connected to one of the PC COM ports. In this manner, a device connected by RS-232 to the COM port of a remote modem can communicate with HW4. Additional remote devices can be connected as well, using an RS-485 multi-drop.

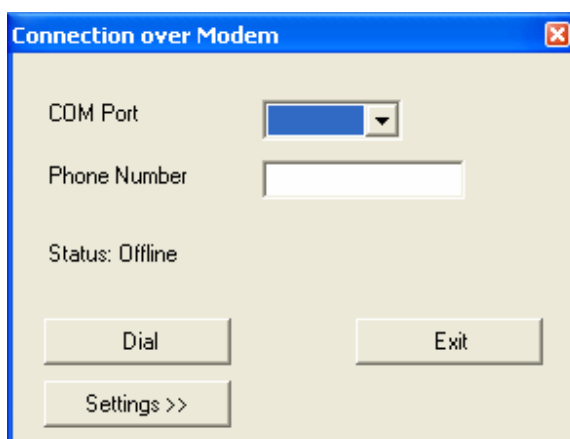

The 'Connection over Modem' dialog box has a blue title bar with a close button. The main area is light beige. It contains three rows of controls: 'COM Port' with a dropdown menu; 'Phone Number' with a text box; and 'Status: Offline' with a text label. At the bottom are three buttons: 'Dial', 'Exit', and 'Settings >>'.

Select the COM port to which the modem is connected and enter the remote phone number. The local modem can be configured after clicking on the Settings button. These settings should be:

Baud Rate : 57,600
Data bits : 8
Stop bit : 1
Parity : none

Be sure to configure the remote modem in the same manner.

8.3 *View*

- **Alarm Table**¹ – displays the alarm table
- **User Events**¹ – displays the contents of the user event file for the current session
- **Psychrometric Conversions** – opens the Psychrometric Conversions Tool

For details, see the following

[Alarm Table](#)

[User Events](#)


[Psychrometric Conversions](#)

- **Expand Tree** - expands all items in the tree (left pane)
- **Collapse Tree** - collapse all items in the tree (left pane)
- **Device Name** - show the device name in the tree (left pane)
- **Device Type** - show the device type in the tree (left pane)
- **Device Address** - show the device port / IP address and RS-485 address in the tree (left pane)
- **Device Serial Number** - show the device serial number in the tree (left pane)


¹ available only with HW4 Professional edition

Example:

View - Device Name

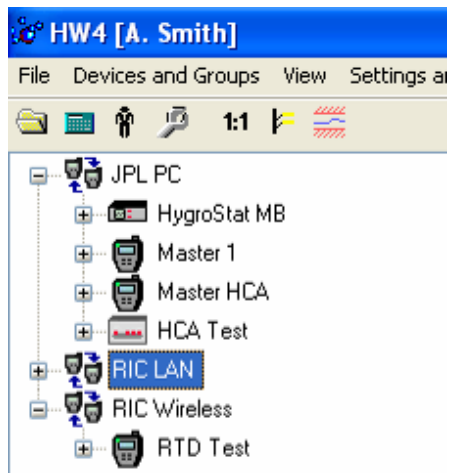

View - Device Address

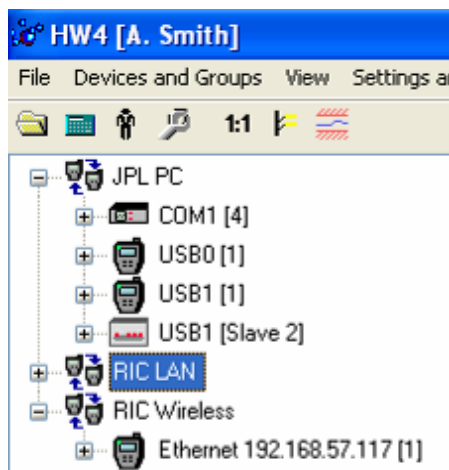

8.4 Settings and Tools

For details click on the following links:


- [Ethernet Configuration Tool](#)
- [Users and Passwords](#)¹
- [Column Headers](#) - selects the columns to be displayed in the current values table (right pane)
- [HW4 Global Settings](#)

¹ available only with HW4 Professional

8.5 Help

- **HW4 Help:** Opens HW4 Help
- **On-line Services:** When an internet connection is available, this menu item accesses the ROTRONIC web site where technical support, updates and product information are available.
- **About HW4:** Displays the version number and ID number of HW4

9 SHORTCUT BUTTONS

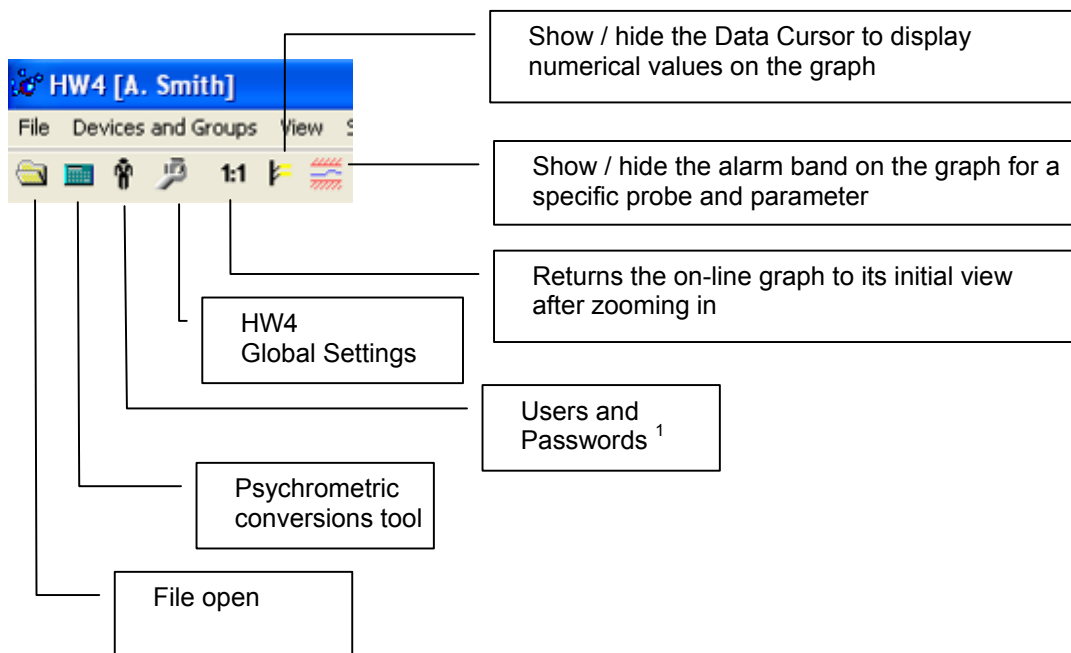


¹ Available only with HW4 Professional

10 ETHERNET CONFIGURATION TOOL

Warning: use of this tool requires basic knowledge of your LAN. If necessary, consult your network administrator.

HW4 comes with an Ethernet device discovery and configuration tool (Digi Device Discovery). When activated, the configuration tool automatically detects any ROTRONIC device present on the LAN (both cable and wireless connections) and provides a list of all such devices. The following example shows the initial screen with just one device.

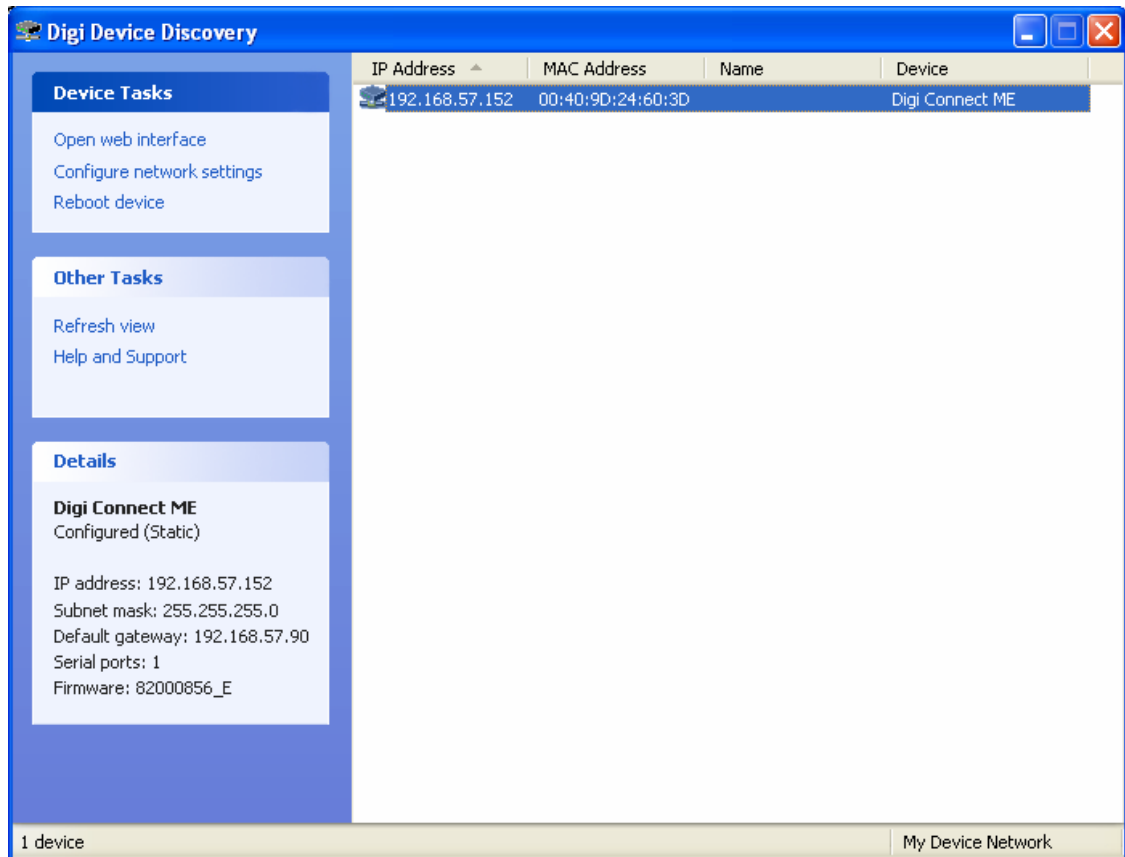

For detailed instructions see [LAN \(TCP/IP\) wired](#) or [LAN \(TCP/IP\) wireless](#).

11 USERS AND PASSWORDS (HW4 Professional)

HW4 Professional allows the creation of several user accounts each with a unique user name. Each user must have a password but passwords do not have to be unique. At any time only one user can be logged on to an HW4 session.

11.1 *Minimum user rights*

All HW4 users have automatically the following minimum rights, which essentially allow the user to view information and data without being able to change anything important:

- Start / Exit HW4
- Log off from HW4 without exiting
- HW4 main menu bar - Send (minimize) HW4 to the Windows notification area
- Select which information is displayed next to the devices in the tree: name, serial number, etc.
- Expand / collapse the device tree.
- Access to Device Manager, Data Logging and Probe Adjustment - functionality limited to viewing data and other information.
- View data in the Current Values tab (both in Device View or Group View mode)
- View data in the On-line Buffer tab.
- Access Data (HygroLog NT function) : limited to opening files located on the PC and viewing the file contents
- HW4 Settings and Tools – HW4 Global Settings - General tab: clear both the online buffer and online graph
- Users and Passwords: functionality is limited to changing own password
- HW4 Settings and Tools – Users and Passwords: functionality is limited to changing own password
- HW4 Settings and Tools – View user events
- HW4 Settings and Tools – Psychrometric conversions
- HW4 Settings and Tools – View Alarm table – the table can be viewed but alarms cannot be acknowledged

The following is not available to users with minimum rights

- Device tree: modification of the content of device groups
- Device tree - Device Manager: any change to the configuration of a device.
- Device tree – Data Logging: any kind of programming of a device log function.
- Device tree – Access Data: opening, copying, moving or deleting files present on the logger. Copying, moving or deleting files present on the PC.
- Device tree – Probe Adjustment – ADJUST button
- HW4 main menu bar - File – Open
- HW4 main menu bar - Devices and Groups menu: the entire menu is disabled
- HW4 main menu bar - View menu: the Column Headers menu item is disabled both in Device View or Group View mode.
- HW4 main menu bar - HW4 Settings and Tools menu - HW4 Global Settings: the functionality of all tabs is disabled
- HW4 main menu bar: Settings and Tools – Ethernet configuration tool
- Right pane of the HW4 main screen in Device View mode: any change to the Log to PC tab, View tab and OPC Server tab. These tabs can be selected for viewing but have no other functionality

Note: when no user has been created or when no user is logged in, HW4 gives only minimum rights to whoever is using the PC.

11.2 Additional user rights

Any user can be given a combination of the following additional rights:


User Right	Definition
Manage user accounts	Create, modify and delete user accounts
Add / Remove devices	Access to HW4 main menu - Devices and Groups - (add / remove devices and groups to or from the device tree). Access to HW4 Settings and Tools - Ethernet Configuration Tool
HW4 global settings and view settings	Access to HW4 Global Settings . Ability to change the following: Device Tree , View Tab and Column Headers (Group or Device View).
Log-to-PC function	Start and stop Log to PC operations
Change OPC settings	Enable, disable the OPC tags for individual devices (OPC version of HW4 required)
Device configuration / Log function	Access to Device Manager (configure individual devices). Access to Data Logging
Adjust probes	Access to the Probe Adjustment function.
Acknowledge alarms	Acknowledgement of alarms . The following two rights are required for defining the conditions that trigger an alarm: HW4 Global Settings and Device Configuration / Log function.
Read / Copy / Sign files	HW4 Main menu – File – Open (enables all operations on files located on the PC) HygroLog NT - Access Data: Open and copy / or just copy files from the HygroLog NT to the PC. Sign log files and protocols See also (1) and (2) below.
Delete files	Delete files (HygroLog NT or PC) See (1) and (2) below.

(1) HW4 cannot fully protect files located on the PC, since it is always possible to use Windows to open, delete or move any file unless adequate security is organized within Windows, based on the Windows user login and Windows directory access rights.

(2) **Moving files** from a device to the PC disk requires both the Read / Copy and Delete rights.

11.3 *Creating and logging-on the first user*

From the HW4 main menu, click on the following menu item:

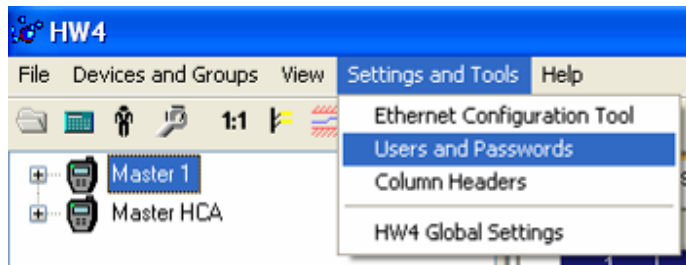

Because no user has been created yet, HW4 automatically opens an empty user form:

The screenshot shows the 'Add User' form in the HW4 application. The form has a title bar 'HW4'. It contains several input fields: 'User name', 'User description', 'Password valid until:' (with a masked value '#####'), 'Password', and 'Confirm password:'. Below these fields is a section titled 'User Rights' containing a list of permissions, each with an unchecked checkbox: 'Manage user accounts', 'Add / remove devices', 'HW4 global settings and view settings', 'Log-to-PC function', 'Change OPC settings', 'Device configuration / log function', 'Adjust probes', 'Acknowledge alarms', 'Read / Copy / Sign files', and 'Delete files'. At the bottom of the form are two buttons: 'Add' and 'Exit'.

Enter the name of the user, the user password and click on the box next to each additional right that the user will have. Be sure to make a note of the password. The user form may now look as follows:

HW4

User name: A. Smith

User description: Administrator

Password valid until: #####

Password: xxxxxxxx

Confirm password: xxxxxxxx

User Rights


- ☒ Manage user accounts
- ☒ Add / remove devices
- ☒ HW4 global settings and view settings
- ☒ Log-to-PC function
- ☒ Change OPC settings
- ☒ Device configuration / log function
- ☒ Adjust probes
- ☒ Acknowledge alarms
- ☒ Read / Copy / Sign files
- ☒ Delete files

Add Exit


Click on the Add button to complete the process.

IMPORTANT: the first user created in HW4 should as a minimum have the right to manage user accounts.

Creating a new user does not result in HW4 automatically logging on this user. To log on the first user, select again Users and Passwords:

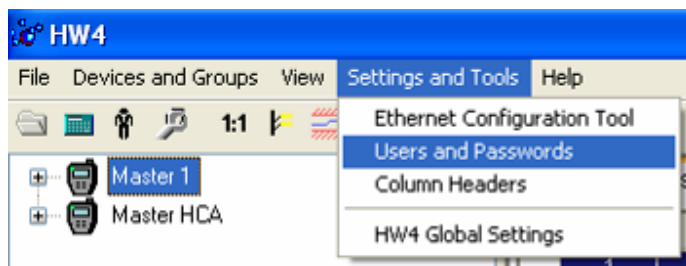


HW4 automatically opens the following box:

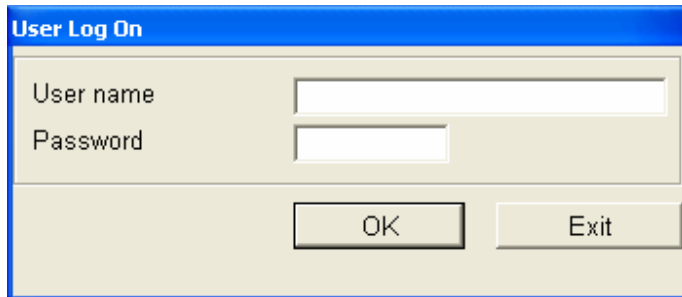

The 'User Log On' dialog box has a blue title bar. It contains two input fields: 'User name' and 'Password'. Below these fields are two buttons: 'OK' and 'Exit'.

Enter the name of the newly created user and the password for this user. Click on OK to log on the user.

HW4 acknowledges this by displaying the name of the current user above the main menu bar:

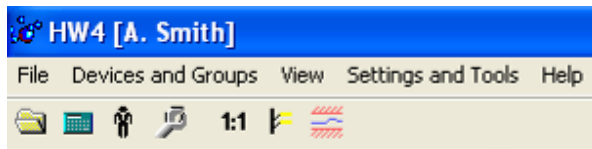


Note: when no user is logged in, there is no user name above the main menu bar. HW4 gives only minimum rights to whoever is using the PC.



11.4 User table form

The user table form is available only when at least one user exists in HW4. To open the table, click on User and Passwords:

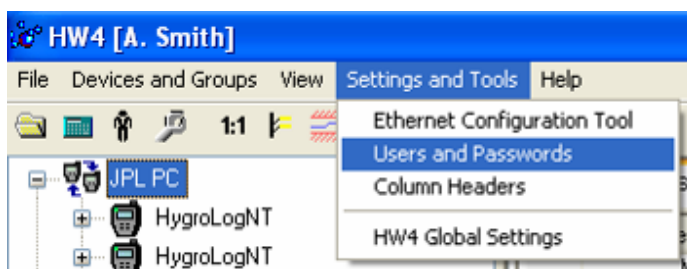

HW4 opens the user table form. When the current user has the right to manage user accounts, the table shows a list of all users currently known to HW4. When the current user does not have this right, only the current user is shown:

Status	User	Description
Logged off	J. Doe	Minimum Rights
Logged on	A. Smith	Administrator

11.5 Menu bar

View: use the View menu on the form to change the contents of the form by selecting one of the options shown below:


User Operations: this menu item gives access to the following:

View / change user data and rights, Create a new user and Delete user: Any of these selections will require the current user to re-enter his / her password:

View / change user data and rights: full access to this function requires the current user to have the right to manage user accounts. For other users, the function is limited to being able to change the current user password.

Create a new user and Delete user: any of these selections require the current user to have the “manage user accounts” right.

Creating a new user is done in the same manner as creating the first user. When changing or deleting a user, use the mouse to select the user from the list provided by the user table form:

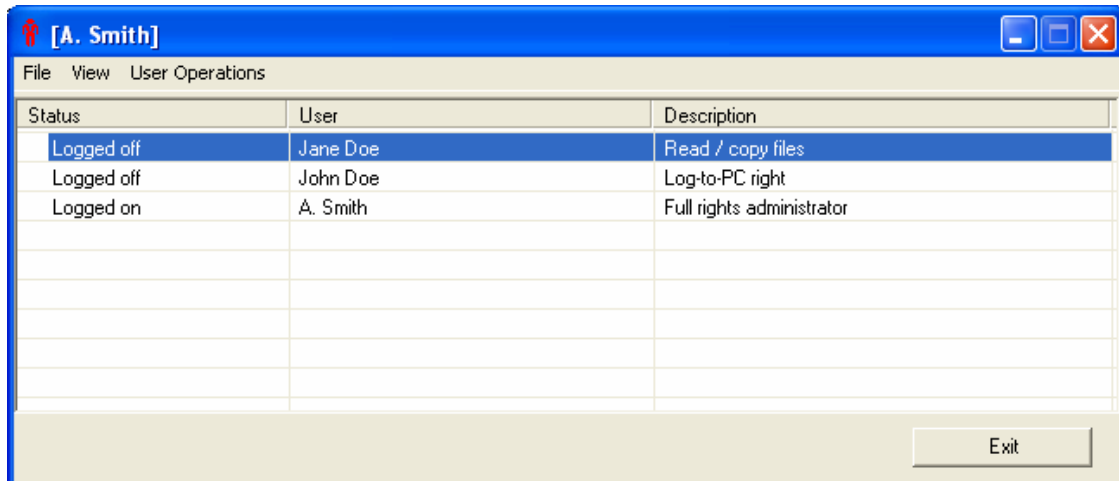


The screenshot shows a window titled "[A. Smith]" with a menu bar containing "File", "View", and "User Operations". Below the menu is a table with three columns: "Status", "User", and "Description". The table contains three rows of data:

Status	User	Description
Logged off	Jane Doe	Read / copy files
Logged off	John Doe	Log-to-PC right
Logged on	A. Smith	Full rights administrator

At the bottom right of the window is an "Exit" button.

IMPORTANT: HW4 keeps track of all users that have been created, even users that have been subsequently deleted. **Please note that the name of a user that has been deleted cannot be used again.**

Log-on a different user: selecting this item open the user log-on form. By default, the form shows the name of the user that is currently logged-on. To change the user, use the mouse to highlight the user name and type in the name of the user to be logged-on. Enter the password for that user. Click on the OK button. The current user is logged off and the new user is logged on.

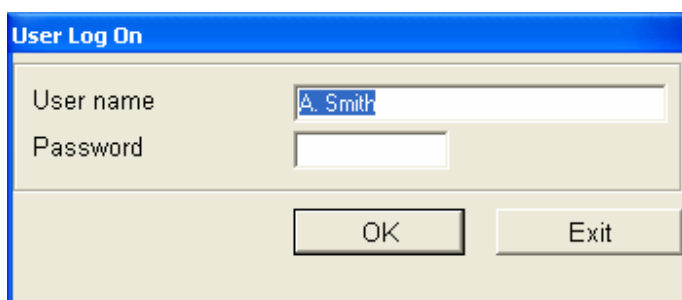

The screenshot shows a dialog box titled "User Log On". It has two input fields: "User name" and "Password". The "User name" field contains the text "A. Smith". Below the input fields are two buttons: "OK" and "Exit".

Note: clicking on the Exit button logs off the current user but does not log on the new user.

Log off current user: selecting this item logs off the current user, without logging on any specific user. When using HW4 Professional, selecting this menu item without logging-on a user may prevent access to most of the HW4 functions and screen views. To return HW4 to the desired level of functionality, log-on as a user with sufficient rights.

11.6 User login security


- Attempts to login with the wrong password are now being recorded.
- HW4 startup will be aborted after 3 attempts with the wrong password

12 COLUMN HEADERS


Column Headers is used to customize the Current Values Tab both in the Device View Mode and in the Group View Mode.

To make HW4 display only specific columns in the Current Values tab:

- **In Device View Mode:** select any individual device present in the device tree (for example HyGen Lab 2 in the left pane). In the HW4 main menu bar, click on **Settings and Tools** and on **Column Headers**.

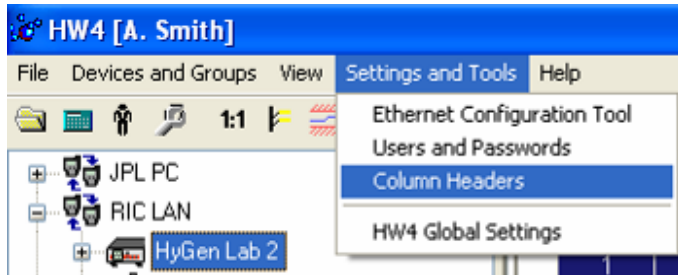


HW4 opens the following form. Make your selection using the mouse. **The selection applies globally to all devices displayed in the Current Values tab.**

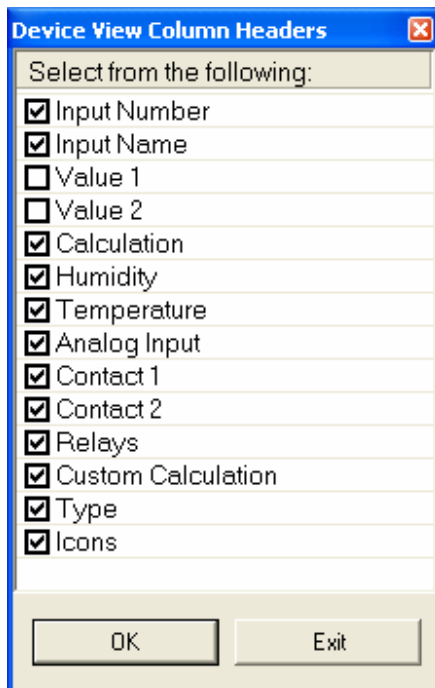

- **In Group View Mode:** select a device group that is present in the device tree (for example RIC LAN in the left pane). In the HW4 main menu bar, click on **Settings and Tools** and on **Column Headers**.



HW4 opens the following form. Make your selection using the mouse. **The selection applies globally to all groups displayed in the Current Values tab.**

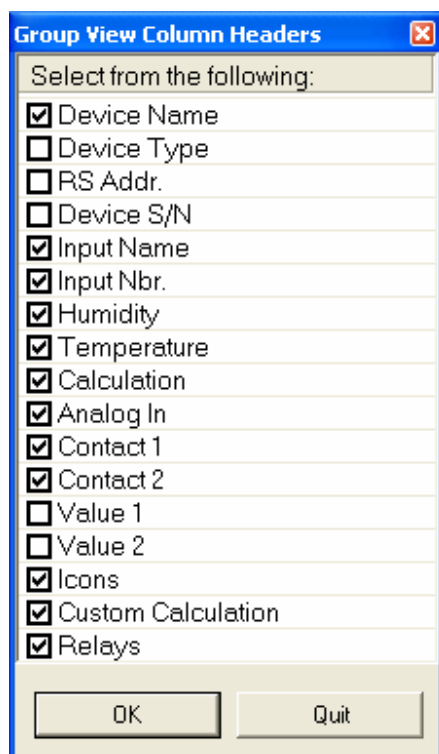

13 HW4 GLOBAL SETTINGS

This menu item opens the HW4 Global Settings Form.

13.1 View Tab

HW4 Global Settings Form

View | General | Language/Unit System | File Locations | Graph Settings | Alarm Settings | Events

Start-up Form

☒ Bypass

General

☒ Online graph

☒ HW4 background pictures

Start-up Tab

☒ Current Values ☐ Online Buffer

☒ Enable trend indicators

10 Font size in the Current Values data table [10...50]

☐ After starting HW4, send as a shortcut icon to the Windows notification area

OK Cancel Help

Start-up Form:

- Bypass: display or bypass the HW4 start-up form (log-in with password is still required)

General:

- On-line graph: shows or hide the graph in the Current Values tab
- HW4 background pictures: select or deselect to modify the appearance of several forms used in Device Manager, HW4 Explorer, etc.

Start-up Tab:

- Current Values / On-line Buffer: select the desired default tab for the right pane of the HW4 main screen
- Trend Indicators: HW4 can display to the right of each current value an indicator that shows if the value is stable (equal symbol), increasing (+ symbol) or decreasing (- symbol). The status of this indicator is read by HW4 directly from the device/instrument being monitored.

Note: this function works only for instrument models that send the status of the trend indicator ¹ to HW4 and is not available in Group View mode.

¹ Instruments such as the HygroFlex and HygroLab send the status of the trend indicator only for the two parameters that are being shown on the local LC display.

- Font Size in the Current Values data table: the default value is value 10. Use a value of 20 or larger to make it easier to read the PC monitor from a distance. Use the mouse to resize each column of the data table as required.

Probe	Description	Measurement 1
1	Probe 1	45.40%RH
2	Probe 2	45.30%RH

Probe	Description	Measurement 1
1	Probe 1	45.60%RH
2	Probe 2	46.29%RH

- After starting HW4, send as a shortcut to the Windows notification area: The Windows notification area is located on the taskbar, immediately to the left of the clock. HW4 can now be sent (minimized) as a shortcut icon to the Windows notification area while it is still keeps running. The alarm table will still appear on the PC monitor when HW4 is running in this mode.

13.2 General Tab

HW4 Global Settings Form

View **General** Language/Unit System File Locations Graph Settings Alarm Settings Events

Communication

☒ Enable Device Polling

00:00:15 Polling interval

00:00:15 Online buffer time interval

3000 Online buffer size [Clear both the online buffer and online graph](#)

Bluetooth

Bluetooth serial ports

Automatic Functions

☒ Enable system monitoring

☐ Force password change after 30 Days

☐ Enter password again after 30 minutes

☐ HW4 updates notification

☐ ROTRONIC NEWS

☐ Automatic synchronization with PC clock

☐ Auto-start mode

OK Cancel Help

Communication:

- **Enable Device Polling:** Put a check mark in this box to enable HW4 to establish the initial (start-up) communication with the devices present in the device tree and to read / update data from these devices at regular intervals of time (polling interval). The data read by HW4 is used to populate the right pane of the HW4 main screen (current values table, on-line buffer and on-line graph).

If you are going to use HW4 only to the purpose of configuring or programming instruments and are not interested in any measurement data, you may want to temporarily uncheck this box by clicking on it with the mouse. HW4 will no longer read and update the data from the devices in the tree and the right pane of the HW4 main screen will be blank. We recommend restoring the check mark prior to closing HW4 so as to allow HW4 to communicate at least once with the devices in the device tree during the next start-up.

- **Polling interval:** enter here the time interval to be used by HW4 when polling the devices connected to the HW4 PC. The minimum value of 5 seconds is also the HW4 default.

Note: the other adjustable time intervals used by HW4 (on-line buffer / graph update and log-to-PC) should be a whole multiple of the polling interval and should be at least equal to the polling interval. When enabled, OPC tags are updated as per the polling interval.

- **On-line buffer time interval:** enter here the time interval to be used by HW4 when writing data to the on-line buffer. The minimum value is 5 seconds. The on-line graph is updated using this time interval.

- **On-line buffer size:** specify here the number of lines for the on-line buffer. The size of the on-line buffer also determines the amount of data shown in the on-line graph.
- **Clear both the online buffer and online graph:** click on this link to clear both the buffer and the graph. After clicking on the link, you can exit HW4 Global Settings by clicking on the Cancel button (HW4 will not generate a protocol).


When adjusting a probe, clear the online buffer when the probe appears to have reached equilibrium. This forces the automatic scale of the graph to the highest resolution.

Bluetooth serial ports:

- The Bluetooth serial ports text box is used to declare the virtual serial port(s) used by Bluetooth enabled devices. Only the number of the serial port should be entered in the box (example: enter the number 8 for COM8). When several Bluetooth devices are used, enter each Bluetooth serial port number, separating each number with a semi-column. Do not use a space in between numbers (example: 8;10).

Automatic Functions:

- **Enable system monitoring:** check this box to have HW4 keep a record of any software problem that may occur. Keeping this type of record is an FDA / GAMP requirement. HW4 event files have the extension ERR and are located in the subfolder EVENT of the HW4 user folder. When support is required from ROTRONIC, troubleshooting the problem is greatly facilitated when the corresponding HW4 event file is available.
- **Force password change:** this feature is only available with HW4 Professional. Enable this feature to force users to change their password on a regular basis.
- **Enter password again after 30 minutes:** this feature is only available with HW4 Professional. Enable this feature to force users to enter their password again after the time indicated. This is an FDA / GAMP requirement which provides additional protection in the event that HW4 is running and left unattended for a period of time.
- **HW4 updates notification:** check this box to be automatically notified of new updates. This feature works only if the PC has access to the Internet.

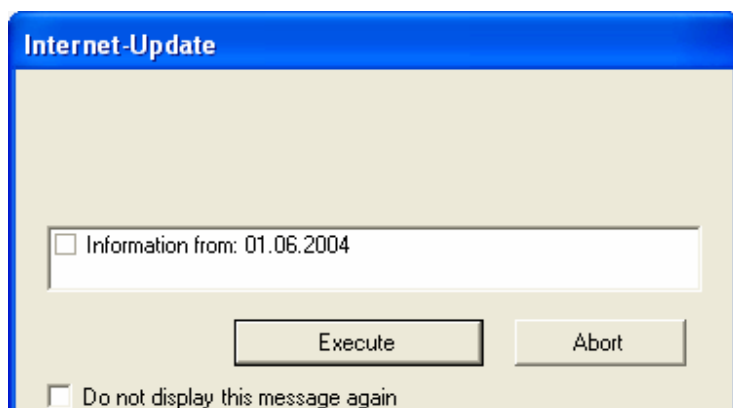

- **ROTRONIC NEWS:** check this box to be automatically notified of new issues of ROTRONIC NEWS. This feature works only if the PC has access to the Internet.

- **Automatic synchronization with PC clock:** Enable this feature to allow HW4 to automatically synchronize any device internal clock with the HW4 PC clock. Synchronization takes place daily at 2:00 am and each time that HW4 is started. In order to be synchronized, the clock of the device must differ from the PC clock by at least 1 minute.

Do not use this function when the PC and some of the polled devices are located in different time zones. See also: [Device Protection](#)

- **HW4 auto-start mode:** the auto-start mode allows HW4 to start without requiring a user to log in. HW4 logs-in the fictitious user "autostart". This user does not have to be created in HW4 and is given only the minimum HW4 rights. The auto-start mode can be used to allow a client application that connects to the HW4 OPC-server (OPC tags) to automatically start HW4. For more details, see document **E-OPC-Tags-v14.pdf**, available from ROTRONIC.

13.3 Language/Unit System Tab

HW4 Global Settings Form

View General **Language/Unit System** File Locations Graph Settings Alarm Settings Events

HW4 Settings

Language Selection [English](#)

Unit System (applies only to Psychrometric Conversion Tool)

[Use metric system](#)
[Use English system](#)

Humidity	%	RH
Temperature	Celsius	
Pressure	hPa	
Enthalpy	kJ/kg	
Density per volume	g/m ³	
Density per weight	g/kg	
Elevation (ref. sea level)	m	

Constants

Fixed pressure value	1013.25
Activation Energy [kJ/mole]	83

OK Cancel Help

The selections and values entered in this tab affect only the language used in HW4 menus and the unit system and fixed pressure value used by the Psychrometric Conversion Tool. This tab has no effect on the instruments and devices in the tree. When HW4 displays or records data from any instrument, both the unit system and any fixed pressure value are determined by the configuration of the instrument.

- **Language:** click on the underlined link to change the language used in the HW4 menus and forms to one of the available languages. The link opens the folder where the language files are located (these files have the extension .txt as in English.txt). To change the language, simply click on the appropriate file. The selection made here has no effect on the devices in the tree.

Note: After changing the language selection, you should close and re-open HW4 to ensure a complete updating of all HW4 texts.

- **Unit System:** use the underlined links to the right of the form to globally change the unit system used within the HW4 Psychrometric Conversion Tool. The English the unit system offers two choices: In Hg or PSI for pressure. Choose a unit by clicking on the arrow to the right of each text box. Relative humidity: the symbol to be used after the % symbol (RH) should be typed in the text box.

- **Constants:**

Fixed pressure value: enter here the default fixed pressure value that will be used within the HW4 Psychrometric Conversion Tool to compute calculated parameters such as wet-bulb temperature, mixing ratio, etc. The barometric pressure value can be temporarily changed from within the Tool but it reverts to the default value each time the Tool is opened.

Note: the numerical value entered here should be consistent with the pressure unit that was selected under Unit System.

Activation energy: enter here the value to be used in the calculation of the mean kinetic temperature (MKT). MKT is part of the statistical data that can be seen when using the HW4 log file viewer.

Note: the activation energy should always be entered in kJ/mole, regardless of which unit system (metric or English) is being used. Typically, a value between 60 and 100 kJ/mole may be used for liquids and solids.

13.4 File Locations Tab

The screenshot shows the 'HW4 Global Settings Form' with the 'File Locations' tab selected. The form is organized into three main sections: Measurement Data, Configuration Data, and Protocols and Events. Each section contains a list of file paths that can be edited.

Section	File Path
Measurement Data	Logger Data Files: C:\ROTRONIC_HW4\DATA
	Log to PC Files: C:\ROTRONIC_HW4\DATA_ONLINE
Configuration Data	User Folder Path: C:
	Device Configuration Files: C:\ROTRONIC_HW4\SYS
Protocols and Events	Protocols: C:\ROTRONIC_HW4\DOC
	Device Event Files: C:\ROTRONIC_HW4\EVENT
	Alarm Records: C:\ROTRONIC_HW4\ALARM
	HW4 Event Files: C:\ROTRONIC_HW4\EVENT
	User Event Files: C:\ROTRONIC_HW4\EVENT


At the bottom of the form are three buttons: OK, Cancel, and Help.

Use this tab to review or change the location of the files created by HW4. Bring the mouse cursor on top of any file path to display the full path. Left click on any file path to open HW4 Explorer and change the path.

- **Logger data files:** log files copied from the logger to the PC. These files have either the extension XLS (unprotected) or LOG (protected binary format)
- **Log to PC files:** log files directly recorded on the PC. These files have either the extension XLS (unprotected) or LOG (protected binary format)
- **User Folder Path:** This folder is used to hold the HW4 configuration file (HW4.ini) as well as the data, event and protocol files created by HW4. If you wish to change the path of the HW4 User Folder, see: [Relocating the HW4 User Folder](#).
- **Device configuration files:** these files have the extension DAT and are used to retain typical instrument configuration data so as to be able to quickly configure a number of instruments.
- **Protocols:** these files have the extension TXT and keep track of configuration changes to an instrument, instrument programming (such as data logging) and probe adjustments.
- **Device event files:** these files have the extension EVT and are generated by the data logger to keep track of the main internal logger events.

- **Alarm records:** HW4 maintains a file named HW4Alarm.alr to hold all alarm record.
- **HW4 event files:** these files have the extension ERR and keep track of any software problem encountered by HW4.
- **User event files:** these files have the extension EVT and keep track of the main operations performed by each HW4 user (HW4 Professional only). See also [Record keeping by HW4](#)

13.5 Graph Settings Tab

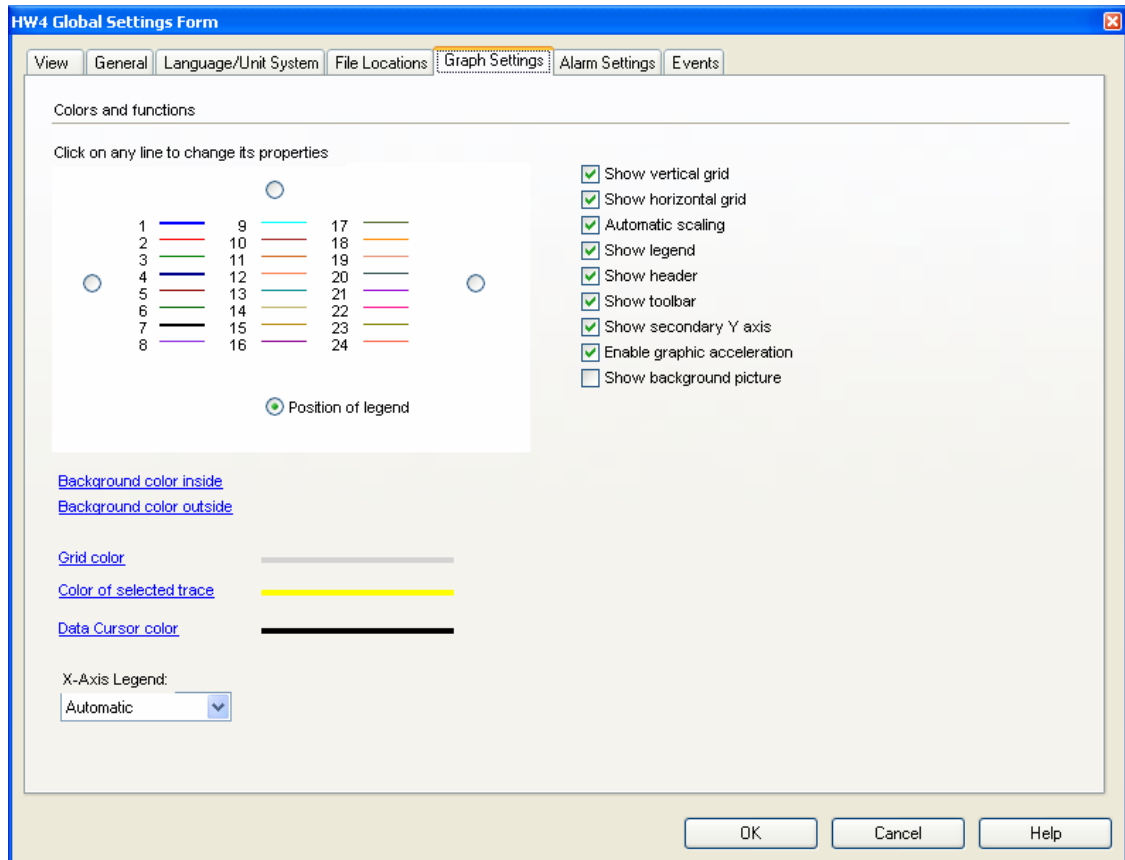

Use this tab to customize the default appearance of the graphs generated from log files. Graphs settings can also be customized for each individual log file when it is being displayed in HW4 View Data. Most of the settings in this tab do not apply to the on-line graph displayed in the Current Values tab.

- **Show Background Picture:** enables the background picture.
- **Line box:** when generating a graph, HW4 uses lines 1 to 24 as follows: line 1 = humidity probe 1, line 2 = temperature probe 1, line 3 = calculated parameter probe 1, line 4 = humidity probe 2, etc. If a probe is not connected or if a parameter is not present on the graph, HW4 does not make use of the corresponding line. The properties of each line (type, color, weight, etc.) can be customized after right clicking with the mouse on the line.
- When the box Show legend is checked, use one of the four radio buttons around the line box to set the position of the legend on the graph.
- The color of the grid, selected trace and data cursor can be changed after clicking with the mouse on the corresponding line.

- Automatic scaling: leave this box unchecked if you wish to manually define the scale for the primary and secondary Y axis of the graphs. Each axis can have its own scale. The secondary axis is used for temperature and the calculated parameter (see example below) and is not available with the on-line graph.

Example of manual scaling:

When the box “Automatic scaling is unchecked”, the form changes and two sets of two test boxes appear to let the user specify the minimum and maximum values for the scale of the primary Y axis (left) and for the scale of the secondary Y axis (right).

HW4 Global Settings Form

View | General | Language/Unit System | File Locations | **Graph Settings** | Alarm Settings | Events

Colors and functions

Click on any line to change its properties

51	1	2	3	4	5	6	7	8	9	10	11	12	13	14	15	16	17	18	19	20	21	22	23	24	25	26	27	28	29	30	31	32	33	34	35	36	37	38	39	40	41	42	43	44	45	46	47	48	49	50	51	52	53	54	55	56	57	58	59	60	61	62	63	64	65	66	67	68	69	70	71	72	73	74	75	76	77	78	79	80	81	82	83	84	85	86	87	88	89	90	91	92	93	94	95	96	97	98	99	100
----	---	---	---	---	---	---	---	---	---	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	----	-----

Position of legend: ☒ Bottom

Background color inside:

Background color outside:

Grid color:

Color of selected trace:


Data Cursor color:

X-Axis Legend:

☒ Show vertical grid
☒ Show horizontal grid
☐ Automatic scaling
☒ Show legend
☒ Show header
☒ Show toolbar
☒ Show secondary Y axis
☒ Enable graphic acceleration
☐ Show background picture

OK Cancel Help

The results may be as illustrated in the graph below:

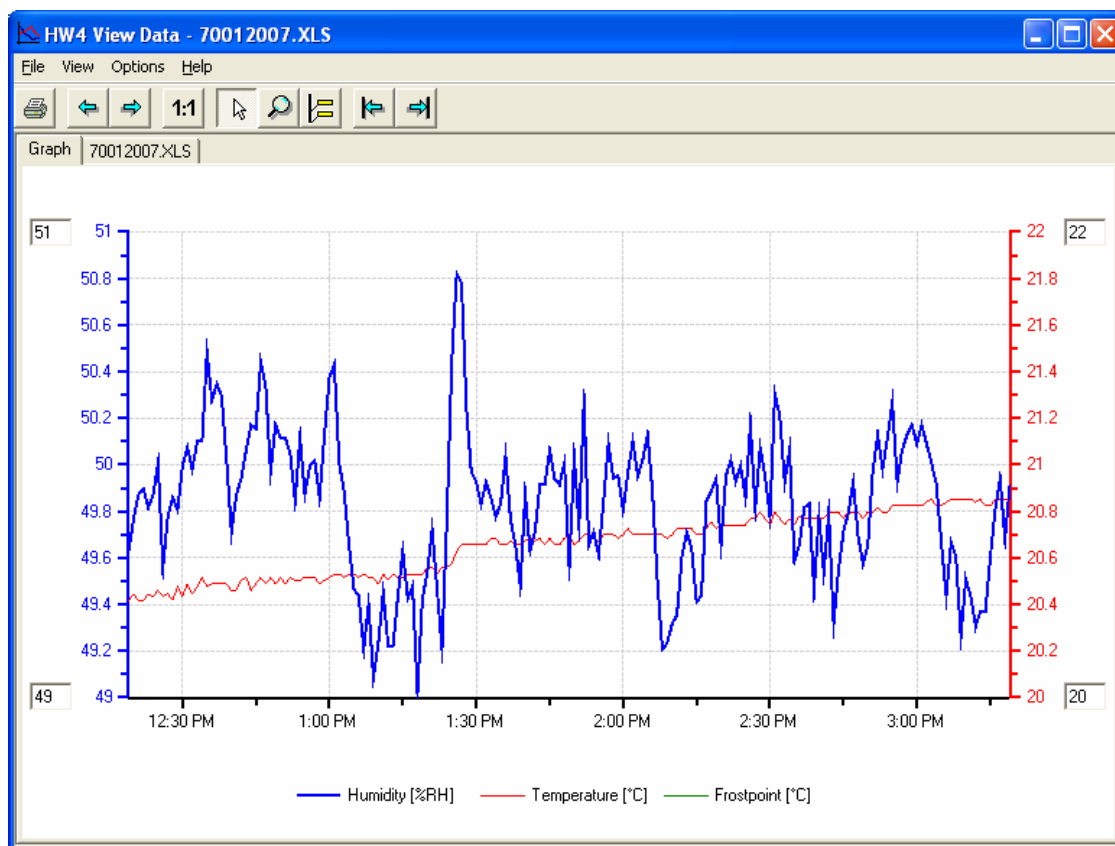

13.6 Alarm Settings Tab (HW4 Professional)

For a complete overview, see “Alarm indication and reporting by HW4” in this manual.

HW4 Global Settings Form

View General Language/Unit System File Locations Graph Settings **Alarm Settings** Events

Trigger an alarm for the following:

☐ No communication with device	Priority	Low	Time delay [hh:mm:ss]	00:00:00
☐ Measured values out of limits or N/A	Priority	High	Time delay [hh:mm:ss]	00:00:00
☐ HygroLog NT errors	Priority	High	Time delay [hh:mm:ss]	00:00:00
☐ Unavailable log-to-PC or HW4 path	Priority	High		
☐ HW4 software errors	Priority	High		
☐ Alarm Test	Priority	Low	Interval	Daily

[Run an alarm test now](#)

Time limits for the alarm function

☐ Enable workdays only [Mo-Fr]

☐ Enable during limited hours Enable time [HH:mm] 00:00 Disable time [HH:mm] 00:00

When an alarm is triggered:

☐ Send e-mail	For priority	☒ Low	☒ Medium	☒ High
☐ Print report	For priority	☒ Low	☒ Medium	☒ High
☐ Give acoustic signal on PC (once)	For priority	☒ Low	☒ Medium	☒ High
☐ Display the alarm table				

E-mail settings

Recipient: asmith@aol.com

☐ MAPI

☒ SMTP: mail.optonline.net

OK Cancel Help

This tab is used to set HW4 to track and report conditions and events that are deemed to be abnormal, such as:

- Loss of communication with any device present in the device tree
- Measured values that are out limits or are not available (see Device Manager - Input and Optional Input)
- HygroLog NT errors (or alarm condition). This includes monitoring the status of the power supply and memory card (see Device Manager)
- Unavailability of the file location (path) used by HW4 while recording data to the PC (example: disconnected network drive)
- Internal HW4 software errors / malfunctions

Click on the box next to each type of event to have HW4 trigger an alarm.

- **Alarm Test:** this function allows you to verify the reporting of alarms

Click on the box labeled Alarm Test to have HW4 simulate an alarm at regular intervals of time (daily, weekly or monthly).

Click on the link “Run an alarm test now” to have HW4 simulate an alarm right away.

- **Priority:** select from low, medium or high to fine tune the way HW4 displays alarms in the alarm table as well as the different actions that HW4 may take.
- **Time delay:** specify here the amount of time during which an event must be on before HW4 triggers an alarm.

HW4 can report an alarm condition in several ways:

- **Send an e-mail:** HW4 sends an email to the address entered in the **Recipient** box.


MAPI: HW4 will use Microsoft Outlook to send an email. Recent versions of Microsoft Outlook require a manual authorization each time that a third party software attempts to use Outlook to send a message. This will prevent HW4 from sending mail when the PC is unattended.

SMTP: HW4 will send an email directly to an outgoing-mail (SMTP) server. This can be done automatically, even when the PC is unattended.

- **Print report:** HW4 sends a short report to the default printer selected in Windows.
- **Give acoustic signal:** HW4 causes the PC to emit a one time bell sound
- **For priority:** specify here for which priority level(s) each of the above actions should be carried out.
- **Display the alarm table:** HW4 automatically records the alarm into the alarm table (see HW4 settings and tools - View alarm table) and displays the table on the PC monitor.

Even when none of the above actions is enabled, HW4 either displays a brief message or an alarm indication on the PC monitor (such as change to the color red in the device tree) in the event of an alarm condition. Abnormal conditions and events which are selected in this tab are always recorded to the HW4 alarm table (see **HW4 settings and tools - View alarm table**).

13.7 Events Tab (HW4 Professional)

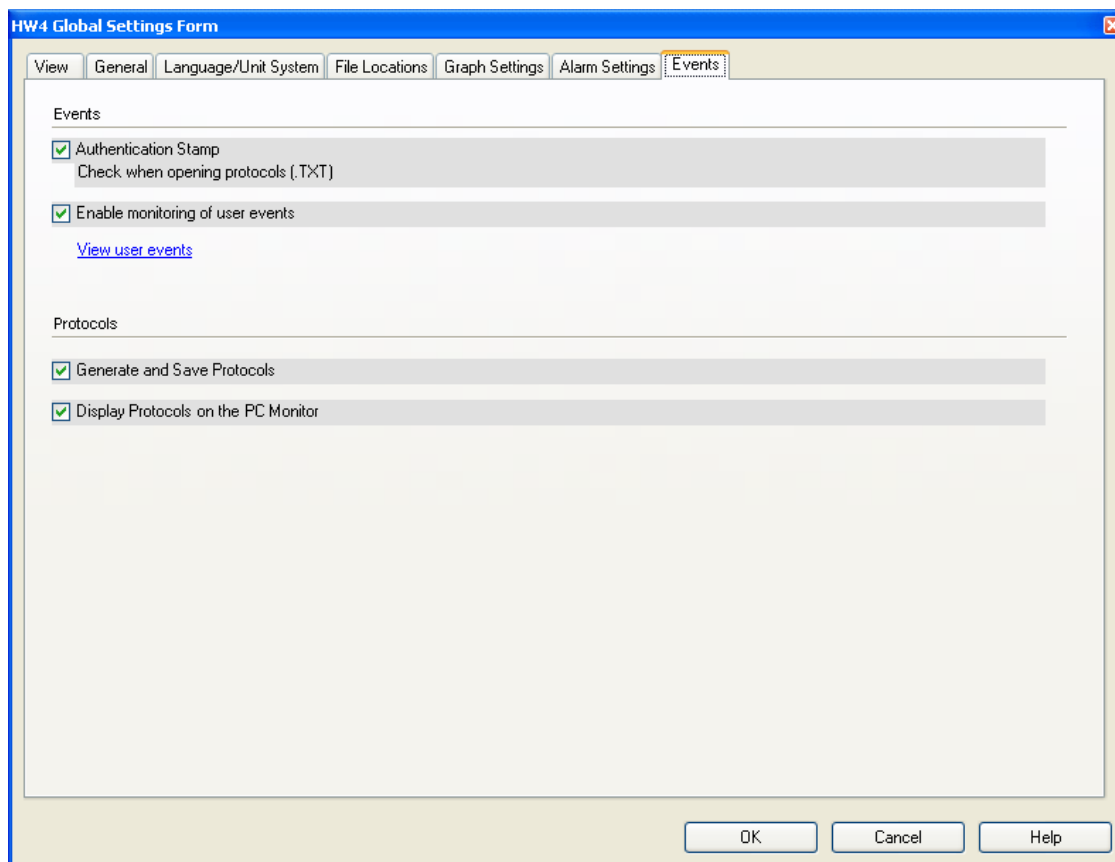

- Authentication stamp:** HW4 Professional automatically generates an authentication stamp at the end of protocols and this stamp is cross referenced in the user event files. Check this box to have HW4 validate the authentication stamp when opening a protocol (file with the extension TXT): instrument configuration, instrument programming and probe adjustment. HW4 uses the authentication stamp to verify that the protocol file has not been altered.
- Enable monitoring of user events:** check this box to have HW4 generate and maintain a user event file keeping track of the main operations performed by any user. This is a requirement of both FDA and GAMP.
- View user events:** click on this link to view the contents of the user events file.
- Generate and Save Protocols:** check this box to have HW4 automatically **generate and save to disk** a protocol file recording the details of the following operations: device configuration, log function programming (except log-to-PC) and probe adjustment.

 HW4 generates protocols in simple text format (Editor / Microsoft Notepad - file extension TXT).
- Display Protocols on the PC Monitor:** This selection causes HW4 to automatically display the protocol on the PC monitor. This option cannot be selected without selecting **Generate and Save Protocols**.

Note: for a complete overview, see "[Record keeping](#)" in this manual.

14 USER EVENTS (HW4 Professional)

When HW4 is enabled to do so (HW4 Global Settings – Events Tab), a user event file is created and maintained. For more information, see “**Record keeping by HW4**” in this manual.

The following are examples of user events:

- User login / logout
- Changes to the HW4 Global Settings
- Instrument configuration
- Log function operations, etc.

Proceed as follows to view the user event file for the current session of HW4:



HW4 opens the user event form. By default, HW4 displays only the events of the current session:


File			
Date	Time	Event	Description
4/25/2006	13:37:50	Start HW4 ID#	12 00000002 V1.2.2.28258
4/25/2006	13:37:55	User Log on	A. Smith
4/25/2006	13:37:55	Electronic signature of language file not valid	C:\Program Files\HW4\language\English.txt
4/25/2006	14:13:20	User log off	A. Smith
4/25/2006	14:15:30	User Log on	A. Smith

14.1 Menu bar


Load previous user events:

To view every user event recorded by HW4 since a specific date, proceed as follows:

In the menu bar of the user events form, click on File and then on Load previous user events.

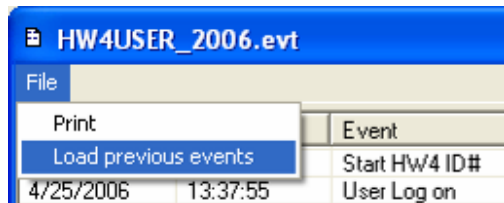

HW4 displays the following dialog box:

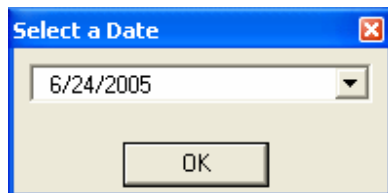

Click on the arrow to the right of the date field to open a calendar.

Select the oldest date from which you want to see all user events on record and click on the OK button. HW4 loads in the user event form all user events on record since the selected date.

Note: the next time that the user event form is opened, HW4 will display again only the events of the current session.


Print

This command prints the current contents of the user event form (current session only or every event starting from a specific date).

15 ALARM TABLE (HW4 Professional)

For an overview of the alarm function in HW4 see [Alarm indication and reporting](#).

The HW4 alarm table can be displayed on the PC monitor at any time from the HW4 main menu bar:

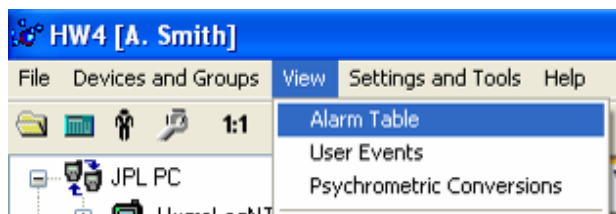

Notes:

- The HW4 alarm table remains empty until the tracking of alarm conditions is enabled in HW4 Global Settings – Alarm Settings Tab.
- The Alarm Settings Tab in HW4 Global Settings offers the option of having HW4 automatically display the alarm table on the PC monitor whenever a new alarm condition occurs.
- A general overview of alarm indication and reporting by HW4 is provided under “Alarm indication and reporting by HW4”.

By default, HW4 initially displays only alarms that have occurred during the current session. Alarms from previous sessions can be displayed by clicking on **View** in the alarm table menu bar.

Alarm conditions that are presently in effect are shown with a red exclamation mark. The red exclamation mark is automatically removed when the alarm condition no longer exists.


	Start Time	End Time	Device	Device name	Input name	Alarm Description
●	11/24/2006 4:11 PM		HygroLogNT	Upstairs		No communication with device
✓ ●	11/24/2006 3:54 PM	11/24/2006 3:54 PM	HygroLogNT	Upstairs		No communication with device
! ●	11/24/2006 3:51 PM		HygroLab 2	HyGen Lab 2	HG Ctrl	No communication with probe
●	11/24/2006 3:49 PM	11/24/2006 3:49 PM	HW4 Software			Alarm Test

Acknowledgement Text	Value	Threshold	Acknowledged by	Acknowledged
Wireless problem	--- %RH	---	A. Smith	11/24/2006 3:58 PM

Note: the alarm table is quite wide and is split in two parts in the above representation.

15.1 Alarm acknowledgement

Acknowledging an alarm requires the current user to have to the right to do so. To acknowledge an alarm, double click on the alarm line in the table. The following dialog box opens:

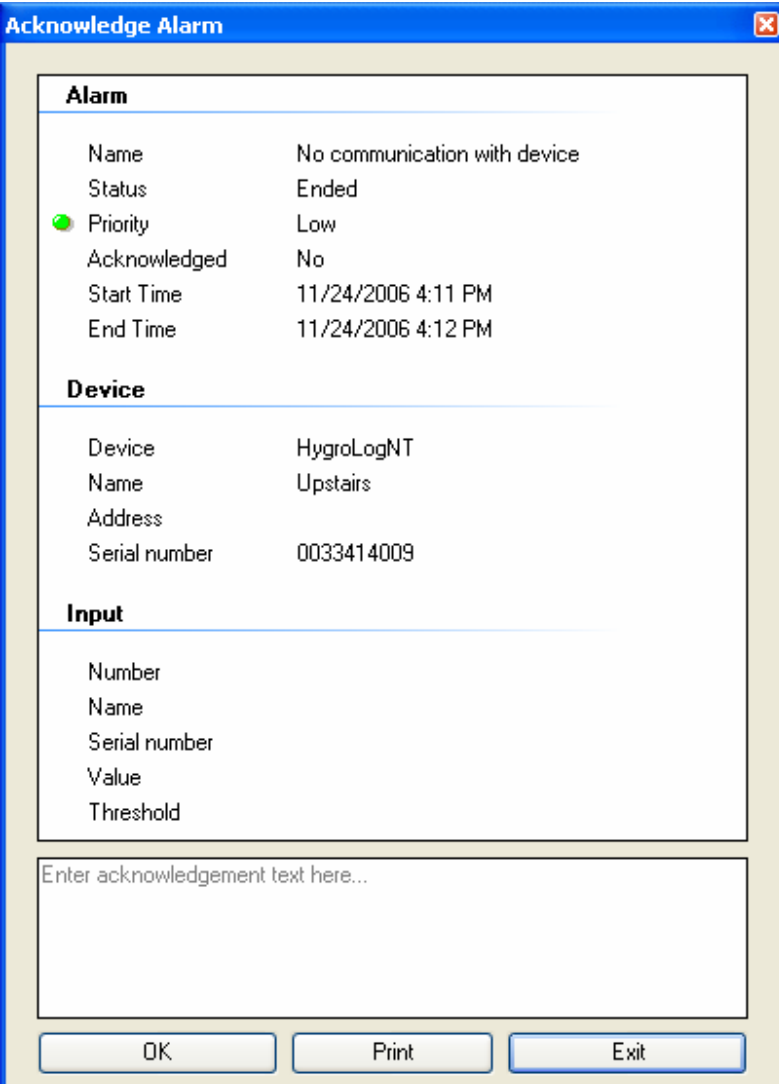

The dialog box is titled "Acknowledge Alarm" and contains three sections: Alarm, Device, and Input. The Alarm section displays details for an alarm that ended. The Device section shows information for a HygroLogNT device. The Input section lists fields for acknowledgment. At the bottom, there is a text box for entering acknowledgment text and three buttons: OK, Print, and Exit.

Alarm	
Name	No communication with device
Status	Ended
Priority	Low
Acknowledged	No
Start Time	11/24/2006 4:11 PM
End Time	11/24/2006 4:12 PM

Device	
Device	HygroLogNT
Name	Upstairs
Address	
Serial number	0033414009


Input	
Number	
Name	
Serial number	
Value	
Threshold	

Enter acknowledgement text here...

OK Print Exit

Prior to clicking on the OK button, you may want to enter a comment in the acknowledgement text box.


When the alarm has been acknowledged, this is shown in the table with a check mark as illustrated below:

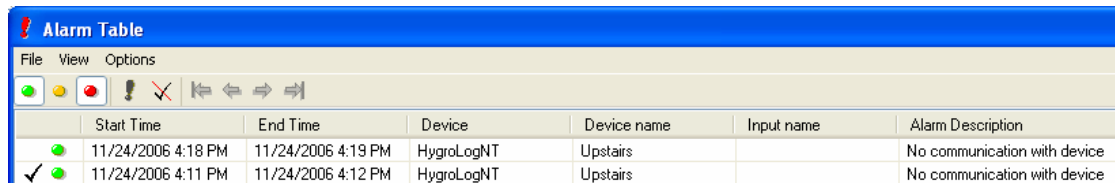

	Start Time	End Time	Device	Device name	Input name	Alarm Description
✓	11/24/2006 4:18 PM	11/24/2006 4:19 PM	HygroLogNT	Upstairs		No communication with device
✓	11/24/2006 4:11 PM	11/24/2006 4:12 PM	HygroLogNT	Upstairs		No communication with device

Acknowledging an alarm or clicking on the EXIT button does not reset the alarm. The acknowledged check mark simply means that a qualified user knows about the alarm and may have entered a comment. An alarm that is current (active) remains visualized in the main HW4 screen (using the color red in the group / device tree and in the current values table).

15.2 Alarm table menu bar and buttons

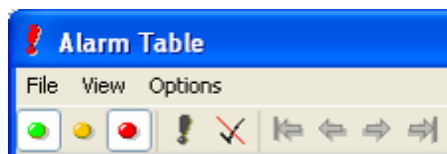

File

- **Print Summary Report:** prints the current contents of the alarm table as a report
- **Save As:** saves the current contents of the alarm table in the Microsoft Excel format
- **Exit:** closes the alarm table

View

Use **View** to change the contents of the alarm table and to view alarms selectively:

- **Alarms from the current HW4 session**
- **Alarms for one day**¹⁾
- **Alarms for one week**¹⁾
- **Alarms for one month**¹⁾
- **Low priority alarms**
- **Medium priority alarms**
- **High priority alarms**
- **Current alarms only** (shows only the alarms that are in effect)
- **Non acknowledged alarms only**

¹ Use the arrow buttons to change the day, week or month:




The information line at the bottom left corner of the table shows the number of alarm records and the date range presently covered in the table:

5 Entries on 11/24/2006	76 Entries on 11/23/2006
81 Entries from 11/20/2006 to 11/26/2006	

Options

- **Customize contents:** the information appearing both in the alarm table and in the alarm report can be customized. Clicking on "Customize contents" opens the following form:

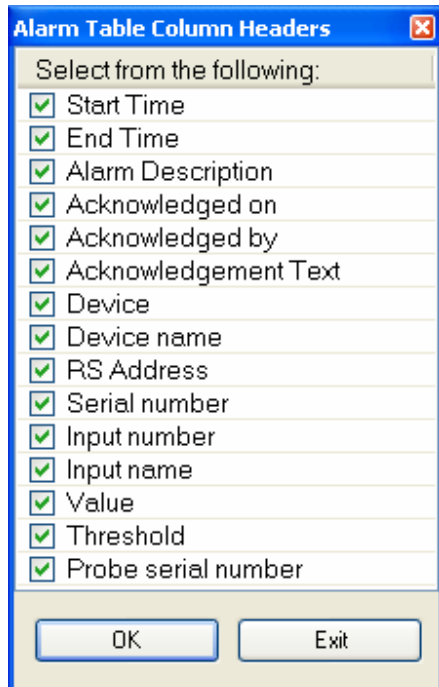

Alarm Table Column Headers

Select from the following:

- ☒ Start Time
- ☒ End Time
- ☒ Alarm Description
- ☒ Acknowledged on
- ☒ Acknowledged by
- ☒ Acknowledgement Text
- ☒ Device
- ☒ Device name
- ☒ RS Address
- ☒ Serial number
- ☒ Input number
- ☒ Input name
- ☒ Value
- ☒ Threshold
- ☒ Probe serial number

OK Exit

The selected items appear as columns in the alarm table and as lines in the alarm report.

The order of the columns within the alarm table can be changed by clicking and dragging them with the mouse.

When the value that triggers an alarm is a calculated parameter such as dew point, the symbol for that parameter appears in front of the numerical value in both the alarm table and alarm report.

Example: Fp: -16.58°C

Alarm Description	Input name	Threshold	Value
Relay Output Alarm	Relay D	1.00	1
Relay Output Alarm	Relay C	1.00	1
Value out-of-limits	Probe 2	-25.00	Fp -16.58°C

Alarm table buttons



Change the date or date range of the alarm events to be displayed

Display only non-acknowledged alarms

Display only current alarms

Select which alarm priority levels are to be displayed in the alarm table
(white background = selected)

Green : low priority
Yellow : Medium priority
Red : High priority

16 PSYCHROMETRIC CONVERSIONS

This menu item opens a form that can be used to convert humidity parameters. To change the parameter to be converted, click on Input Parameter and choose from the list.


Parameter	Symbol	Value	Unit
Humidity	H	35.00	%RH
Temperature	T	23.00	°C
Pressure	P	1013.25	hPa
Elevation (ref. sea level)		0.00	m
Wet-bulb temperature	Tw	13.98	°C
Psychr. difference	T-Tw	9.02	°C
Dew point	Dp	6.73	°C
Frost point	Fp	N/A	°C
Vapor concentration (moist)	Dv	7.19	g/m³
Specific humidity (moist)	Q	6.06	g/kg
Mixing ratio by Wt. (dry)	R	6.09	g/kg
Vapor partial pressure	E	9.83	hPa
Vapor saturation pressure	Ew	28.08	hPa
Saturation vapor concentration	Dvs	20.54	g/m³
Enthalpy	H	38.60	kJ/kg
Volume Mixing Ratio (dry)		9795.44	ppm

Calculation : WMO World Meteorological Organization
Status: OK

- Click on the blue arrow to compute
- Click on the diskette symbol to transfer the results to the default Windows text editor (such as Notepad) and eventually save the data to a text file.

Note: to change the unit system, go to the HW4 main screen. In the menu bar, click on Settings and Tools and on HW4 Global Settings. In the HW4 Global Settings form, select the tab "Language / Unit System".

17 DEVICE MANAGER - HygroLog NT

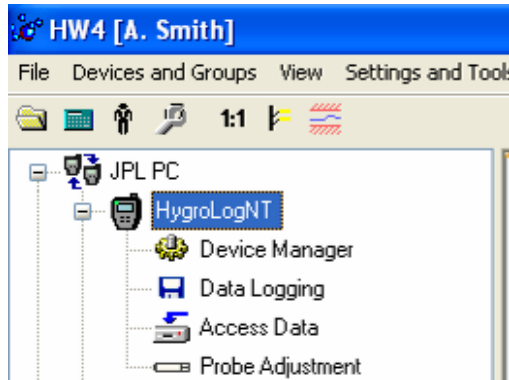

When HW4 has detected a HygroLog NT, the device appears in the left pane of the HW4 main screen. Expanding the device displays a list of the available device functions.

To select the Device Manager function, click on it with the left mouse button. HW4 opens the Device Manager form.

The Device Manager form is used to configure the HygroLog NT and to read instrument specific information (but not the data recorded by the logger). When Device Manager is started, it automatically interrogates the instrument and downloads its current configuration.

Device Manager / HygroLog NT3 / 0034277001 / HyGen NT3

File Device Info Help

Device Owner

Name: A. Smith

Company: Rotronic Instrument

Address: 160 E, Main Street

ZIP/City: 11743 Huntington, NY

Country: USA

Telephone: 631 427-3898

[Custom logo file](#)
[Display custom logo](#)

☒ Include in log file

OK Cancel

The different forms that are available within the Device Manager form are listed in a tree located on the left pane of the form. To select a form, click on it with the left mouse button.


17.1 Device Manager Menu Bar

The Device manager menu bar is located at the top of the form.

File

- **Open:** opens the device configuration directory specified in the HW4 Settings Form - File Locations Tab - and displays all available device configuration files (extension .DAT). Any device configuration file that was previously saved can be opened to quickly configure an instrument. If so desired, any directory and any file type may be opened.
- **Save As:** saves the current configuration to a file (extension .DAT) in the device configuration directory specified in the HW4 Settings Form - File Locations Tab. If so desired, any directory and any file type may be specified.
- **Exit:** exits Device Manager

Device Info

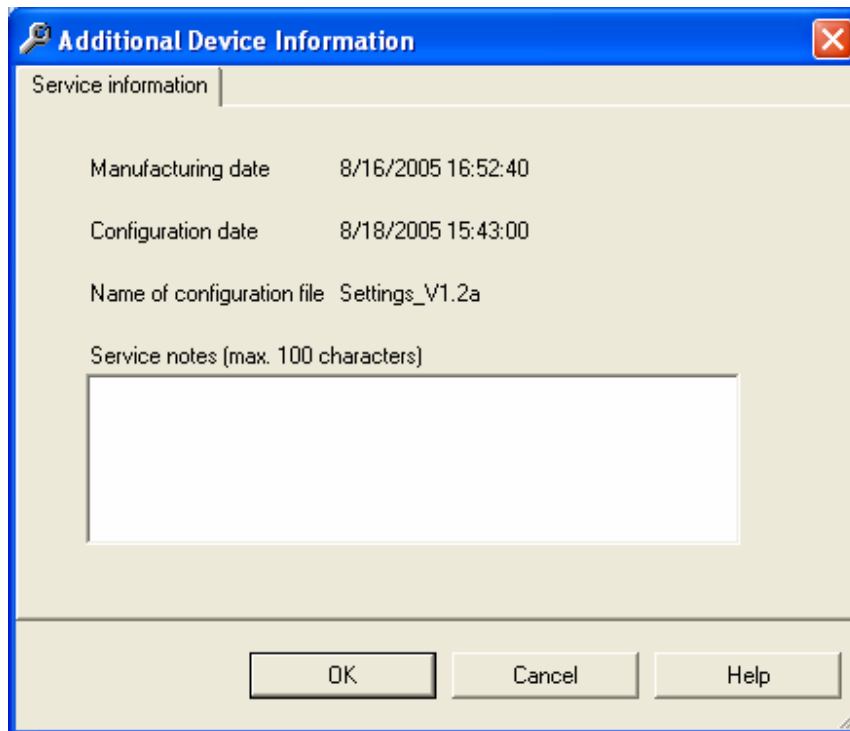

Additional Device Information

Service information

Manufacturing date 8/16/2005 16:52:40

Configuration date 8/18/2005 15:43:00

Name of configuration file Settings_V1.2a

Service notes (max. 100 characters)

OK Cancel Help

This form provides additional product information that may have been entered by the factory.

Help

The Help menu consists of:

- **HW4 Help:** Opens HW4 Help
- **About HW4:** Displays the version number and ID number of HW4

17.2 Device Owner

The screenshot shows a software window titled "Device Manager / HygroLog NT3 / 0034277001 / HyGen NT3". The window has a menu bar with "File", "Device Info", and "Help". On the left is a tree view with the following items: "Device Owner" (selected), "Device Information", "Language/Unit System", "Memory Card", "Standard Inputs" (expanded), "Display", "LEDs", "Keypad", "Sounds", "Power Source", "Digital Interface", and "Docking Station" (expanded). The main area is titled "Device Owner" and contains the following fields:

Name	A. Smith		Custom logo file
Company	Rotronic Instrument		Display custom logo
Address	160 E, Main Street		
ZIP/City	11743	Huntington, NY	
Country	USA		
Telephone	631 427-3898		

At the bottom left of the main area is a checkbox labeled "Include in log file" which is checked. At the bottom right are "OK" and "Cancel" buttons.

User information will be included in the log file when the corresponding box is checked. This is a requirement of both FDA and GAMP.

17.3 Device Information




- **Device name:** this text will be displayed next to the instrument model information. As far as possible use a unique device name.
- **RS485 address:** click on the underlined link to change the instrument address to be used in conjunction with an RS-485 network (multi-drop). Each network address should be unique and within the values of 0 to 63. Note: the default factory RS-485 address is 0. Unless necessary, do not manually modify this address. HW4 will automatically change the RS-485 address of the device, if so required.
- **Synchronize with PC clock:** click on the underlined link to set the instrument date and time to match the PC clock.

Note: the device date and time are automatically synchronized with the PC each time that the log function is being programmed.

- **Device Protection:** this function is available from firmware version 1.2b. For a description, see [Device Protection](#)

17.4 Language / Unit System

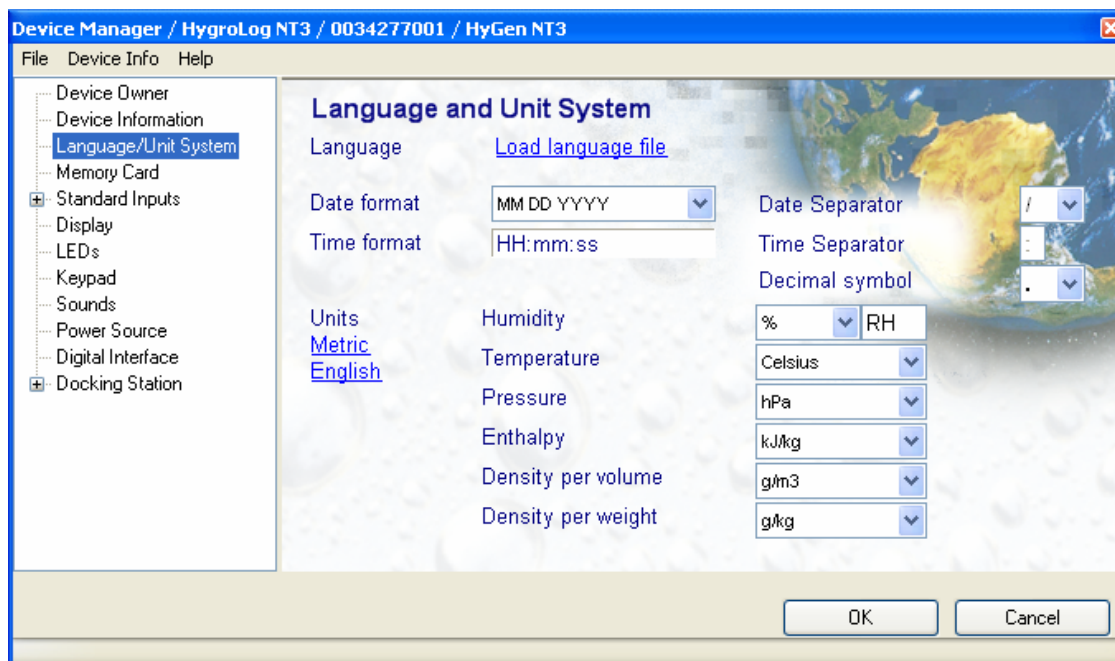

- **Language:** click on the underlined link to change the language used in the device local menus and files to one of the available languages. The link opens the folder where the language files are located (these files have the extension .LNG as in English.LNG). To change the language, simply click on the appropriate file.
- **Units:** use the underlined links to the right of the form to globally change the unit system used by the instrument. If the unit system offers several choices such as In Hg or PSI for pressure, choose a unit by left clicking on the arrow to the right of each text box.


Humidity: the symbol to be used after the % symbol (RH) should be typed in the text box.

Note: the units selected for density per volume and density per weight also apply to vapor concentration and specific humidity.

IMPORTANT:

- Do not change the language / units settings while either the logger or the PC is recording data.
- Make sure that the date format and separator symbols are identical to the format and separators used on the PC (Windows Control Panel, Regional and Language Options, Regional Options tab, Customize button, Time and Date tabs). Not doing so may prevent you from reading log files.

17.5 Memory Card

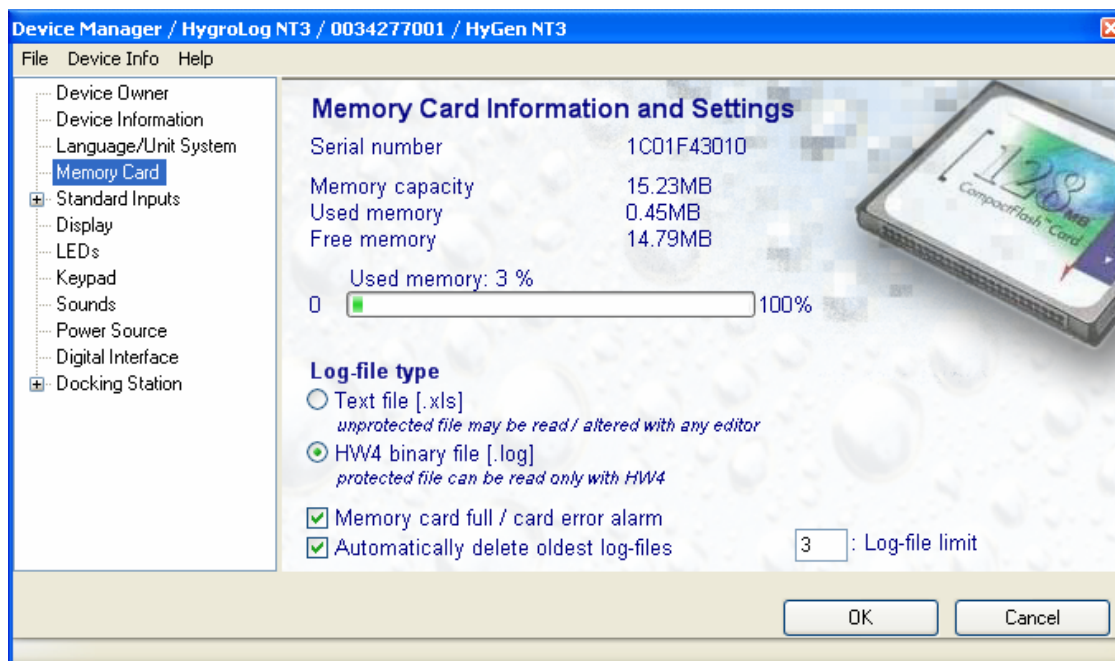

- Log-file type:

Text file (XLS): HW4 does not verify this type of file when opening the file.

HW4 binary file (LOG): files of this type can be opened with HW4 but not with an editor or with Microsoft Excel. When opening a protected file, HW4 verifies that the file has not been tampered with. Note that protected binary files use less disk space than unprotected files.

Select HW4 binary file to meet the FDA / GAMP requirements.

- Memory card full / card error alarm:

Check this box to have the HygroLog NT issue an error warning when the memory card is full or has a problem. All versions of HW4 will show a memory card alarm condition in red on the monitor screen. HW4 Professional can also be configured (HW4 global settings - Alarm settings tab) to display an alarm table and generate a report whenever any alarm or error condition issued by the HygroLog NT.

- Automatically delete oldest log files:


This option works only when the logger automatically closes the current log file and starts a new file because the current log file has reached a size limit before the logging stop time. The size of log files can be limited when programming the log function of the logger. For example, the logger can be programmed to close the current log file and to start a new one at the top of each hour, until the logging stop time has been reached.

Enter the maximum number of full size files to be retained in the text box labeled Log-file limit. In practice, this number will be exceeded by one, unless all files created during the logging process happen to have exactly the maximum file size limit.

17.6 Standard Inputs

Clicking on the plus sign to the left of “Standard Inputs” expands the tree and displays probe inputs 1 to 3. Inputs 1 to 3 are listed whether or not they actually exist on the logger.

Only those inputs selected in the View Tab are visible in Device Manager. When the list is either empty or incomplete, please exit Device manager and check the selections made in the View Tab (see Right Pane: Device View Mode).

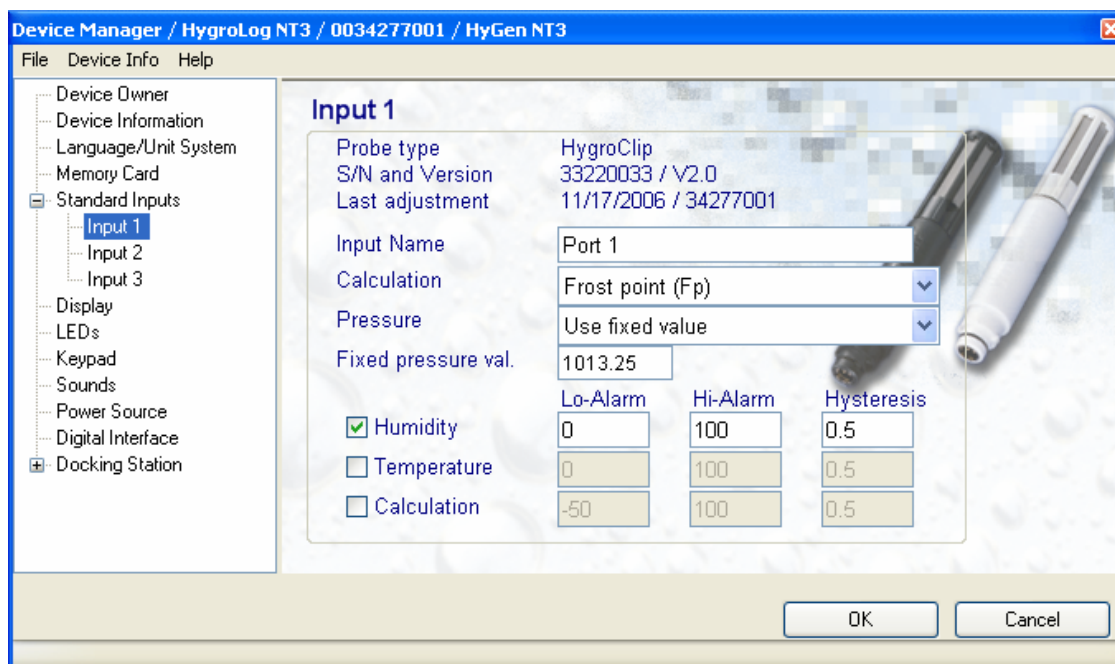


- **Input Name:** use a maximum of 12 characters
- **Calculation:** left click on the arrow to the right of the list box and select the parameter to be calculated by the HygroLog NT for this particular probe. This is also the calculated parameter that HW4 will display for this probe

Dew point (Dp) or Frost point (Fp): depending on the selection made here, both the HygroLog NT and HW4 will display either the symbol Dp or the symbol Fp. The symbol Fp indicates that any value below freezing is a frost point as opposed to a dew point. When selected, the symbol Fp is also be displayed for values above freezing. This of course is to be understood as being the same as dew point.

- **Pressure:** left click on the arrow to the right of the list box and select which barometric pressure will be used by the HygroLog NT to compute parameters such as wet-bulb, mixing ratio, etc. The HygroLog NT can use a fixed pressure value or use a measurement of the actual barometric pressure when an analog pressure probe is connected to the optional docking station.
- **Fixed pressure value:** enter here the fixed pressure value that will be used by the HygroLog NT. Note that this numerical value should be consistent with the pressure unit that was selected under Language / Unit System.

- **Alarm:** alarm conditions can be defined for humidity, temperature and the calculated parameter. Values that are below the low alarm value or above the high alarm value will trigger an alarm. A value can be specified for the alarm function hysteresis. This value is used for both the low and the high alarm. To trigger an alarm only in the event that a probe is missing, disconnected or not transmitting data, simply enable the alarm for one of the parameters and set the Hi and Lo values to the range of the probe (for example 0 and 100 for humidity). All versions of HW4 will show an out-of-limits value alarm in red on the monitor screen. HW4 Professional can also be configured (HW4 global settings - Alarm settings tab) to display an alarm table and generate a report whenever an of-of-limits condition occurs.


17.7 Display

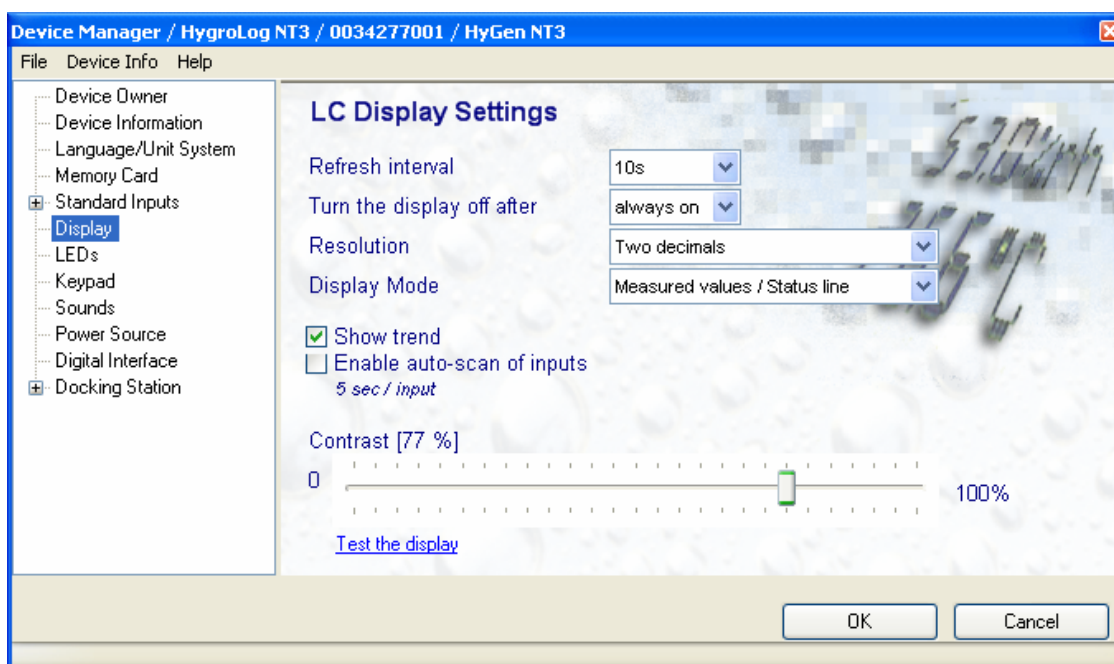

- **Refresh interval:** left click on the arrow to the right of the list box and select the refresh interval for the display of the HygroLog NT. Values range from 5 seconds to 60 minutes. To conserve battery power, use a high setting.
- **Turn the display off after:** use this item enable the Display Sleep function so as to extend the battery lifetime. Left click with the mouse on the arrow to the right of the list box and select one of the available choices. The display of the HygroLog NT goes blank after the specified time period, when no key is being pressed or when there is no communication with the PC.
- **Resolution:** left click on the arrow to the right of the list box and select the number of decimals to be used on the display of the HygroLog NT.
- **Display mode:** left click on the arrow to the right of the list box and select the type of contents for the display of the HygroLog NT. Status line causes the HygroLog NT to display the date and time as well as which input source is providing the displayed data.
- **Show trend:** when this box is checked a trend indicator is displayed to the left of the measured values and calculated parameter. This indicator shows if the values are increasing, decreasing or stable. The status of each trend indicator is read directly from the HygroLog NT.

- **Enable auto-scan of probe inputs:** check this box to make the display automatically switch between probe inputs every 5 seconds.
- **Contrast:** drag the slider to adjust the contrast of the display
- **Test the display:** click on the underline link to turn on all pixels. If the LC display is good, this should cause the display to become uniformly black (no white areas) for a few seconds.

17.8 LEDs

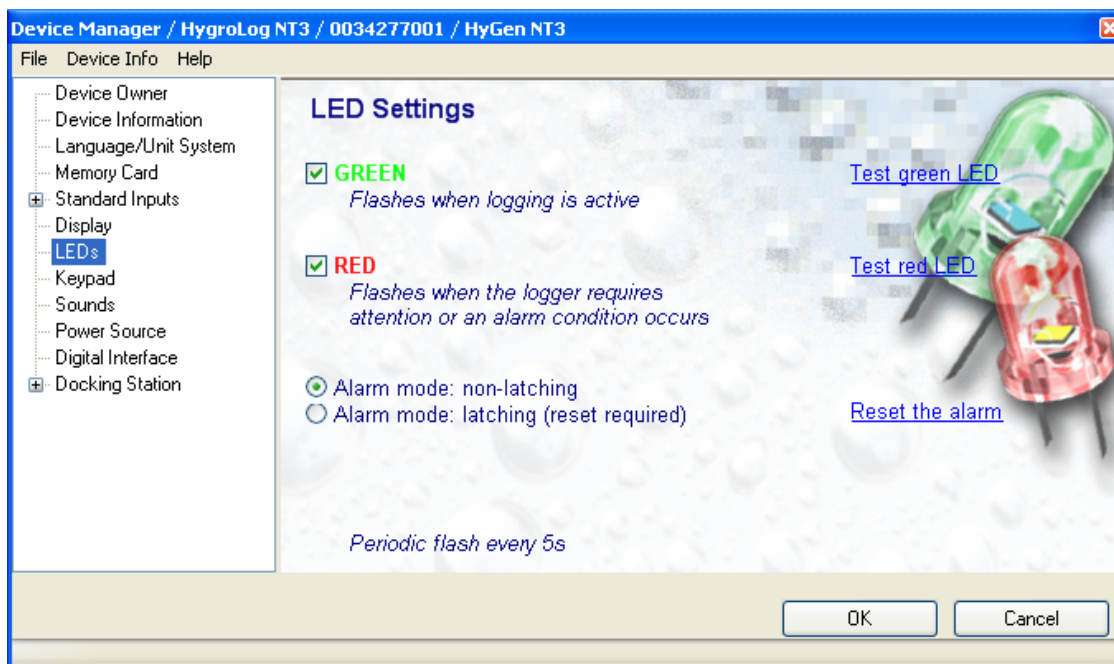


Each of the two LEDs of the HygroLog NT can be enabled or disabled. Each LED can be tested by clicking on the corresponding underlined link. During the test, the LED stays on for a few seconds.

The red LED can be programmed to flash when there is an out-of-limits value and / or when the logger requires attention (low battery, memory card full, etc.).

- **Alarm mode:**

In the non-latching mode, the red LED flashes for only as long as the data from a probe corresponds to an alarm condition. In the latching alarm mode, the red LED stays on even after the alarm condition has disappeared. When the LED is set to provide a latching alarm, it can be reset by clicking on the underlined link [Reset the alarm](#).

17.9 Keypad

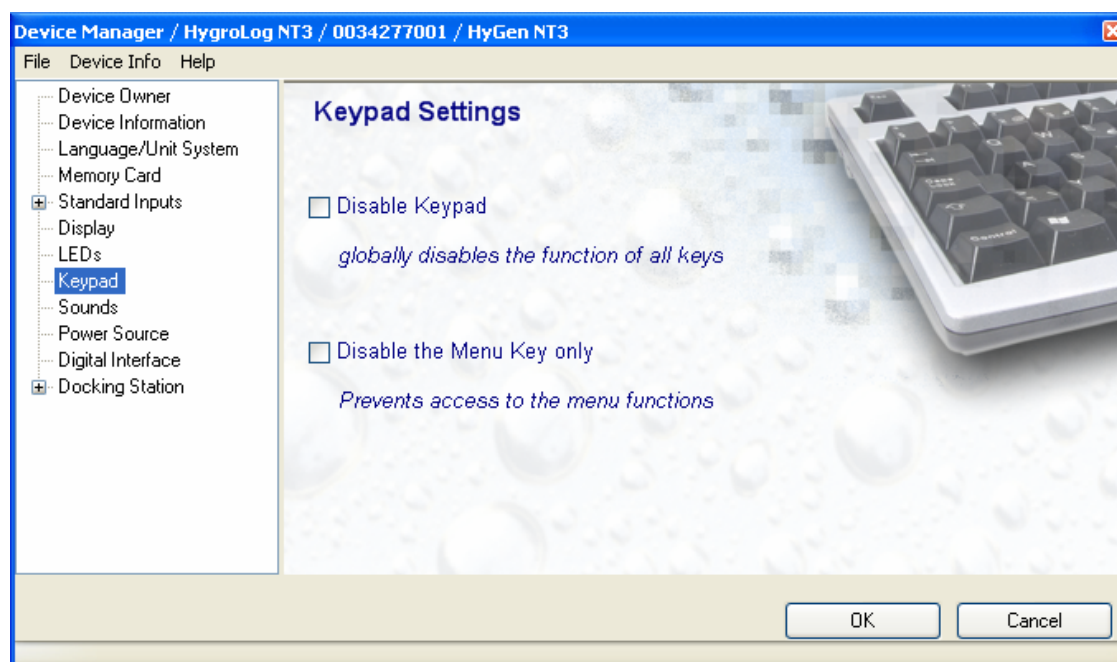

To prevent tampering with the logger or to prevent changes resulting from the accidental pressing of keys, the keypad can be completely disabled.

As an alternative, only the menu key may be disabled to protect both the logger and its probes.

IMPORTANT: you should at least disable the MENU key to comply with FDA CFR21 Part II.


17.10 Sounds



- **Keystroke:** check this box to enable a clicking feedback sound when pressing any key of the keypad.
- **System Error:** check this box to have the HygroLog NT emit a non-repeating 5 seconds sound in the event of an internal error in the HygroLog NT.
- **Alarm:** check this box to have the HygroLog NT emit a non-repeating 5 seconds sound when an alarm condition occurs.

The pitch of the sound can be set by left clicking on the arrow located to the right of the list box and by selecting the desired pitch. To test the sound function, click on the underlined link.

17.11 Power Source

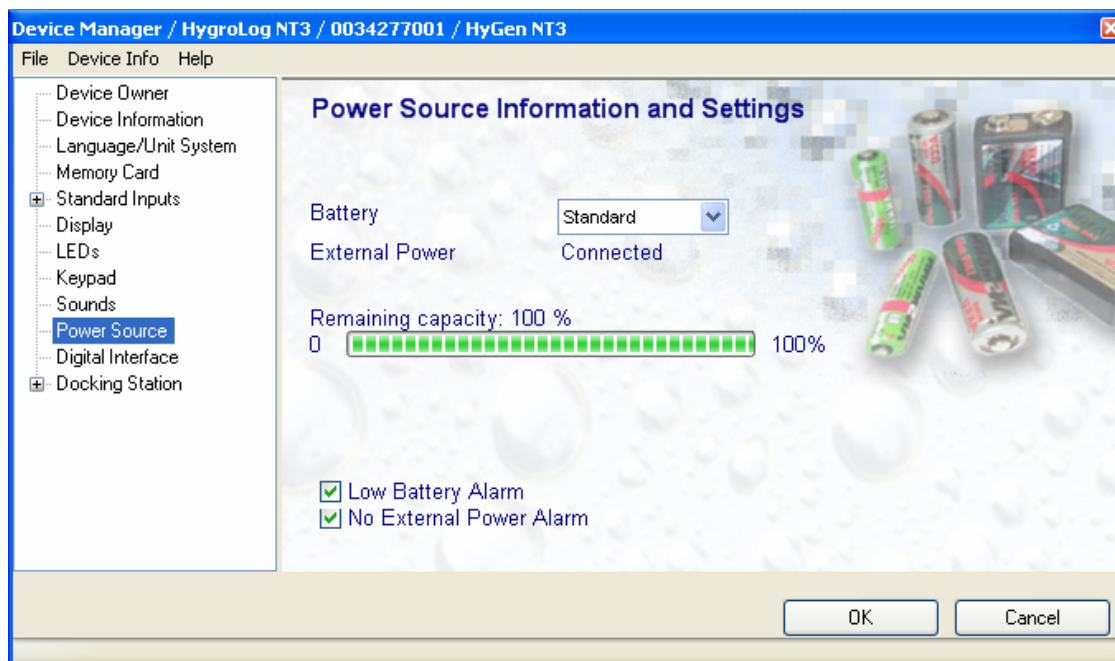

- Type of internal power source: click on the arrow located to the right of the list box and choose between battery and accumulator (rechargeable battery). Note that trying to recharge a regular battery is potentially dangerous and may damage the instrument.


Do not select accumulator if using a regular battery since this would allow the HygroLog NT to attempt recharging the battery when an external power source is connected.

Note: the only way to totally power down the HygroLog NT is to disconnect the battery.

- External Power Source: an external power source is automatically detected by the HygroLog NT
- Trigger an alarm when the battery is low: see note 1 below
- Trigger an alarm when there is no external power (AC power failure or disconnected / malfunctioning AC adapter): see note 1 below

¹ All versions of HW4 will show a power source alarm in red on the monitor screen. HW4 Professional can also be configured (HW4 global settings - Alarm settings tab) to display an alarm table and generate a report whenever any alarm or error condition issued by the HygroLog NT.

17.12 Digital Interface

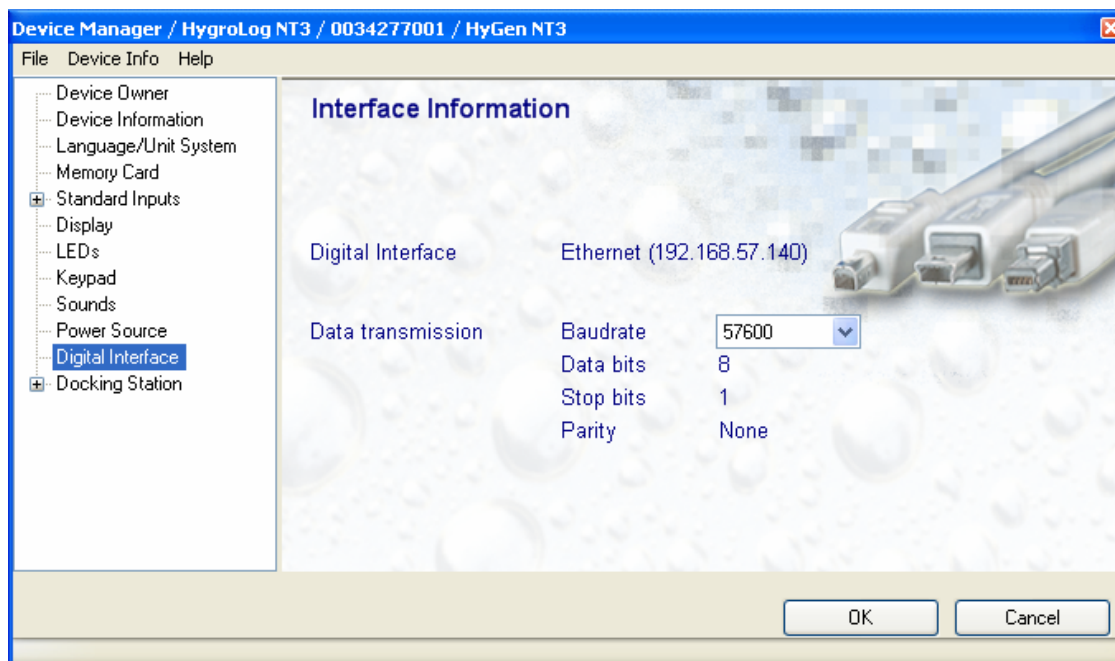

This form is used to select the RS-232 / RS-485 baud rate between 57600 (factory default) and 19200. Other data cannot be changed from within the form. Generally, leave the baud rate at 57600. Select 19200 when you plan on connecting a mix of HygroLog NT and other devices such as the HygroFlex transmitter to an RS-485 multi-drop.

IMPORTANT:

All devices within an RS-485 multi-drop should use the same baud rate. Communication within the RS-485 multi-drop will not work when the devices are configured with different baud rates.

Docking station with wired or wireless Ethernet connection: after changing the baud rate in Device Manager, you should also reflect the change in the configuration of the internal Digi International module used by the device to connect to the LAN. See: [Changing the baud rate of an Ethernet device](#)

17.13 Docking Station


Clicking on the plus sign to the left of “Docking Station” expands the tree and displays a list of optional docking station inputs (inputs 4 to 7), logical inputs and relays. These items are listed whether or not they actually exist on the docking station. Only those items selected in the View Tab are visible in Device Manager. When the list is either empty or incomplete, please exit Device manager and check the selections made in the View Tab (see Right Pane: Device View Mode). Depending on the model of docking station the following probe types may be used:

- HygroClip probe
- RTD probe (direct 4-wire connection)
- Analog probe ¹
- Pressure probe ¹

¹ Must be compatible with the requirements specified in the HygroLog NT instruction manual.

17.13.1. Probe type: HygroClip

Note: this probe type selection is also used in the case of docking station DS-PT-2 (Pt100 RTD probes) - for more information, see Pt100 RTD probes below.

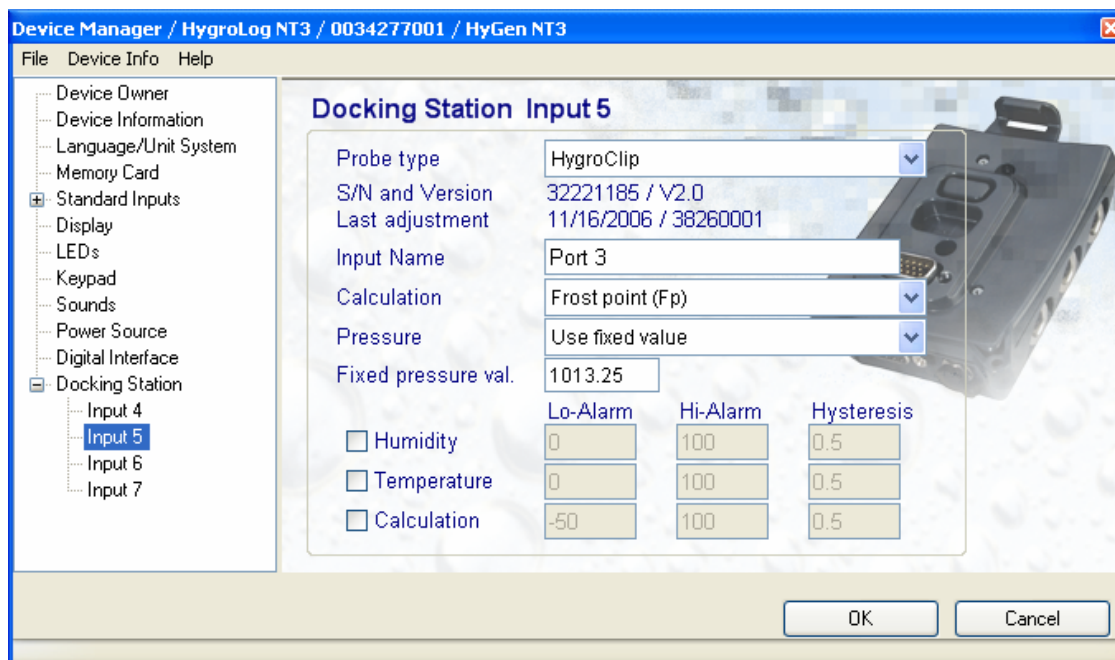

- **Input Name:** use a maximum of 12 characters
- **Calculation:** left click on the arrow to the right of the list box and select the parameter to be calculated by the HygroLog NT for this particular probe. This is also the calculated parameter that HW4 will display for this probe

Dew point (Dp) or Frost point (Fp): depending on the selection made here, both the HygroLog NT and HW4 will display either the symbol Dp or the symbol Fp. The symbol Fp indicates that any value below freezing is a frost point as opposed to a dew point. When selected, the symbol Fp is also be displayed for values above freezing. This of course is to be understood as being the same as dew point.

- **Pressure:** left click on the arrow to the right of the list box and select which barometric pressure will be used by the HygroLog NT to compute parameters such as wet-bulb, mixing ratio, etc. The HygroLog NT can use a fixed pressure value or use a measurement of the actual barometric pressure when an analog pressure probe is connected to the docking station.
- **Fixed pressure value:** enter here the fixed pressure value that will be used by the HygroLog NT.

IMPORTANT: Be sure to verify that the numerical value of the fixed pressure is consistent with the unit system of the instrument (see Language / Unit System). Neglecting to do so may result in an error on the value of any parameter that uses barometric pressure as a calculation input.

- Alarm:** alarm conditions can be defined for humidity, temperature and the calculated parameter. Values that are below the low alarm value or above the high alarm value will trigger an alarm. A value can be specified for the alarm function hysteresis. This value is used for both the low and the high alarm. To trigger an alarm only in the event that a probe is missing, disconnected or not transmitting data, simply enable the alarm for one of the parameters and set the Hi and Lo values to the range of the probe (for example 0 and 100 for humidity). All versions of HW4 will show an out-of-limits value alarm in red on the monitor screen. HW4 Professional can also be configured (HW4 global settings - Alarm settings tab) to display an alarm table and generate a report whenever an out-of-limits condition occurs.

Pt100 RTD probes: When using a docking station compatible with RTD probes, select HygroClip as the probe type. By default, temperature conditions within the range of -100 to 600°C can be measured. Enter a description for each Pt100 probe and set alarm values for temperature in the same way as for a HygroClip probe. Entries made in other fields of the form have no effect.

17.13.2. Probe type: Analog probe (single channel):

The screenshot shows the 'Device Manager' window for 'HygroLog NT3 / 0034277001 / HyGen NT3'. The left sidebar contains a tree view with the following items: Device Owner, Device Information, Language/Unit System, Memory Card, Standard Inputs (expanded), Display, LEDs, Keypad, Sounds, Power Source, Digital Interface, and Docking Station (expanded). Under 'Docking Station', 'Input 5' is selected. The main configuration area is titled 'Docking Station Input 5' and contains the following fields:


- Probe type:** Analog Probe (dropdown menu)
- Input Name:** Port 3 (text field)
- Input Voltage:** 0 - 1000 [mV] (text fields)
- Measuring Range:** 0 100 (text fields)
- Unit:** %RH (dropdown menu)
- Warm-up Time:** >3s (dropdown menu)
- Alarm:** (checkbox, currently unchecked)
- Lo-Alarm:** 0 (text field)
- Hi-Alarm:** 100 (text field)
- Hysteresis:** 0.5 (text field)

At the bottom right of the window are 'OK' and 'Cancel' buttons. An image of the docking station hardware is shown on the right side of the configuration area.

Enter a description for the probe. Specify the range of the input voltage, the measuring range of the probe, the engineering unit of the parameter measured by the probe, the time required by the probe to warm-up, and any alarm values that you may want to use.

17.13.3. Probe type: Pressure probe


The analog pressure probe is a particular case of the single channel analog probe. The difference is that the signal from the analog pressure probe can be used as an input by another probe input when calculating a psychrometric parameter that depends on the value of barometric pressure. See example below:

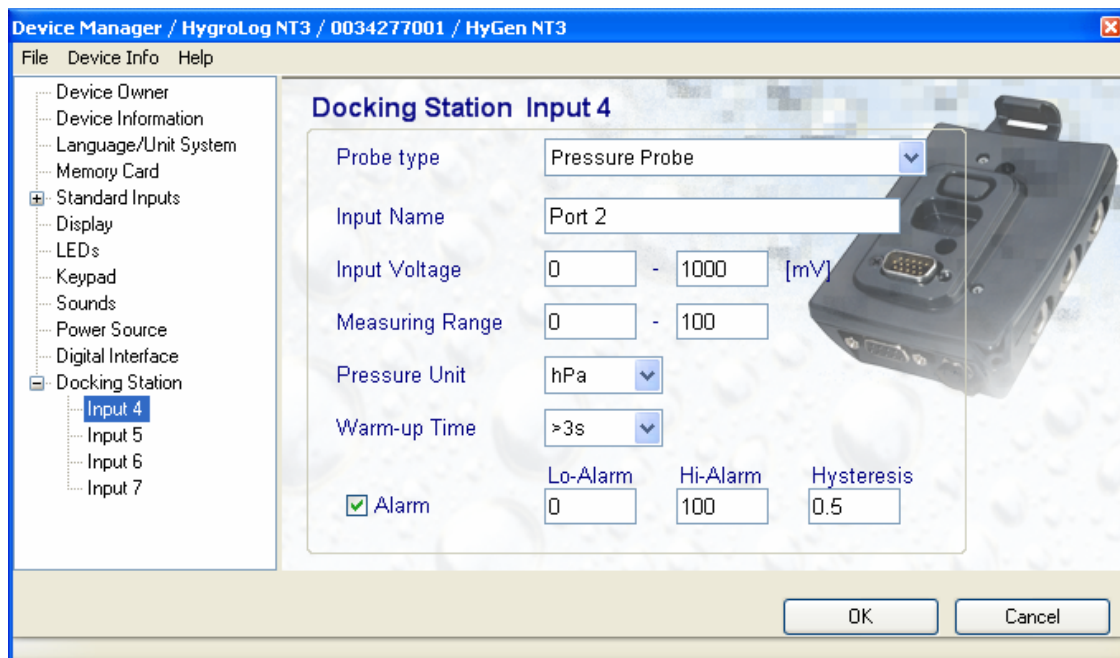

Enter a description for the probe. Specify the range of the input voltage, the measuring range of the probe, the engineering unit used for pressure, the time required by the probe to warm-up, and any alarm values that you may want to use.

17.13.4. External Contacts

Some models of docking station feature two logical inputs that can be used to monitor external contacts such as may be attached to a door (see separate HygroLog NT manual). When an external contact is open, the corresponding logical output takes the value 0.

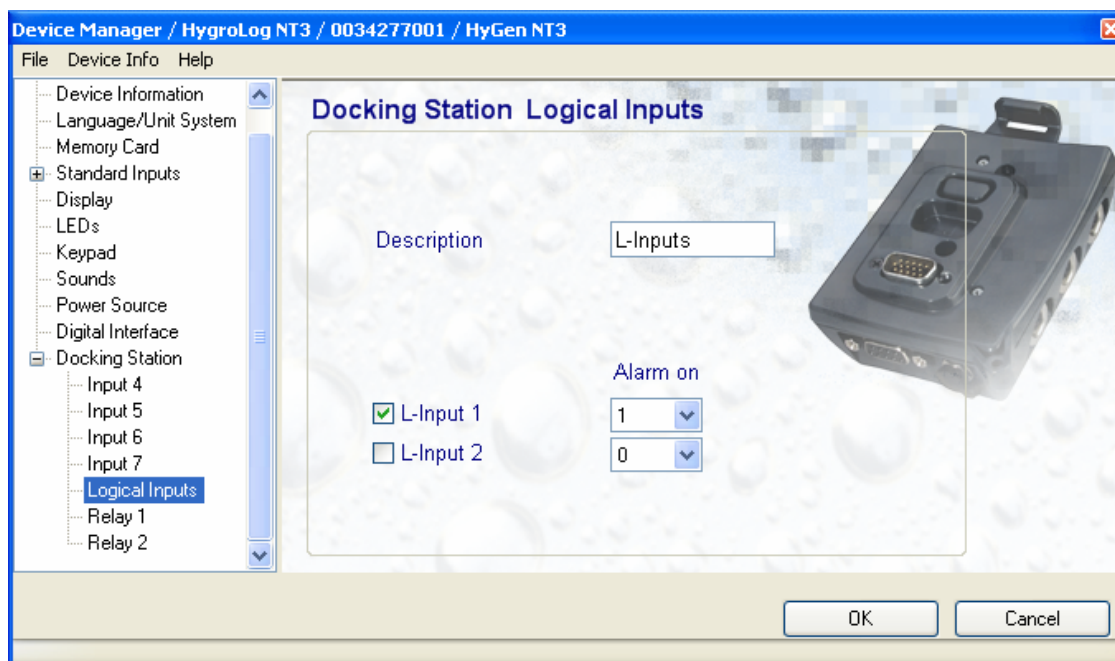


Enter a global description for the two logical inputs and define any condition that you may want to monitor for the purpose of providing an alarm.

17.13.5. Relay 1 / Relay 2 (docking station DS-R-1)

Docking station DS-R-1 is available for use with the HygroLog NT and features two independent relay outputs: Relay 1 and Relay 2 (see separate HygroLog NT manual). The following form is used to program the function of each relay and has 3 tabs: Energize Relay, De-energize Relay and Time Schedule. To select any tab, click with the mouse on the tab name with the mouse.


Note: when a relay is energized, the normally closed relay contact changes from closed to open and the normally open relay contact changes from open to closed.

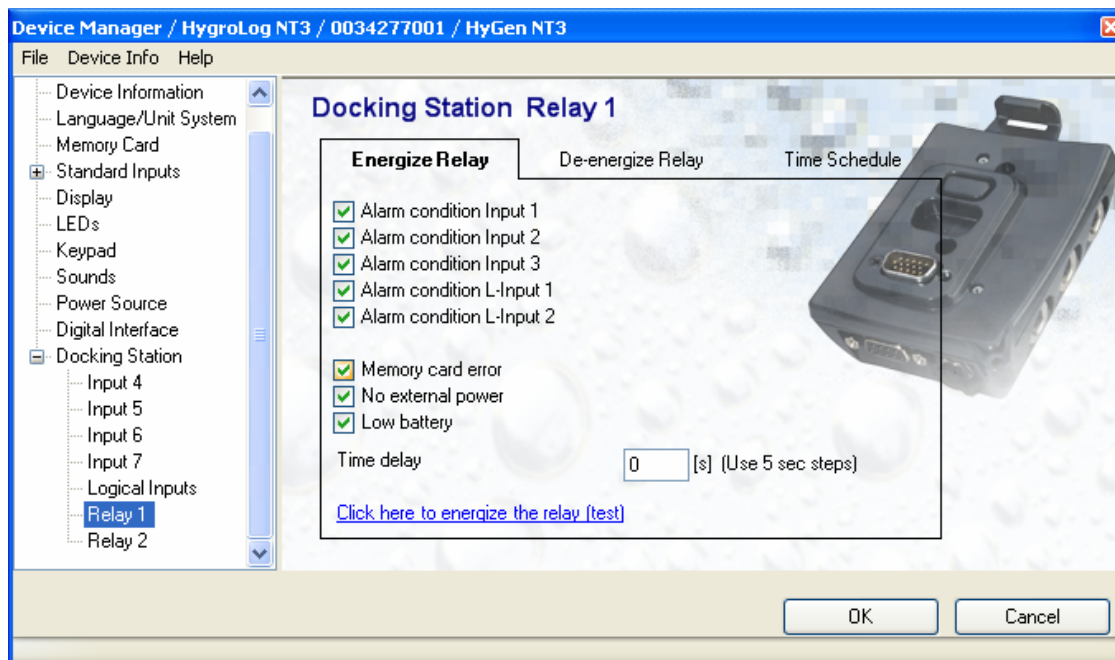Energize Relay Tab:

Use this tab to select from the list one or several alarm conditions that will cause the relay to be energized. Alarm conditions for probe inputs 1 to 3 and logical inputs 1 and 2 are defined in Device Manager under “Input” and “Logical Inputs”. To associate an alarm condition with the relay, click on the corresponding check box. A check mark symbol appears. Click again on the box to remove the check mark and deselect the condition.

If so desired, specify a time delay in seconds by typing a number in the time delay text box (use multiples of 5 seconds). When any of the selected conditions occurs, the relay will be energized after the time delay specified here, provided that the alarm condition still exists.

When the relay is not energized, click with the mouse on “Click here to energize the relay (test)” to energize and test the relay. Note: this link is not visible when the user does not have the “Device configuration / log function” right.

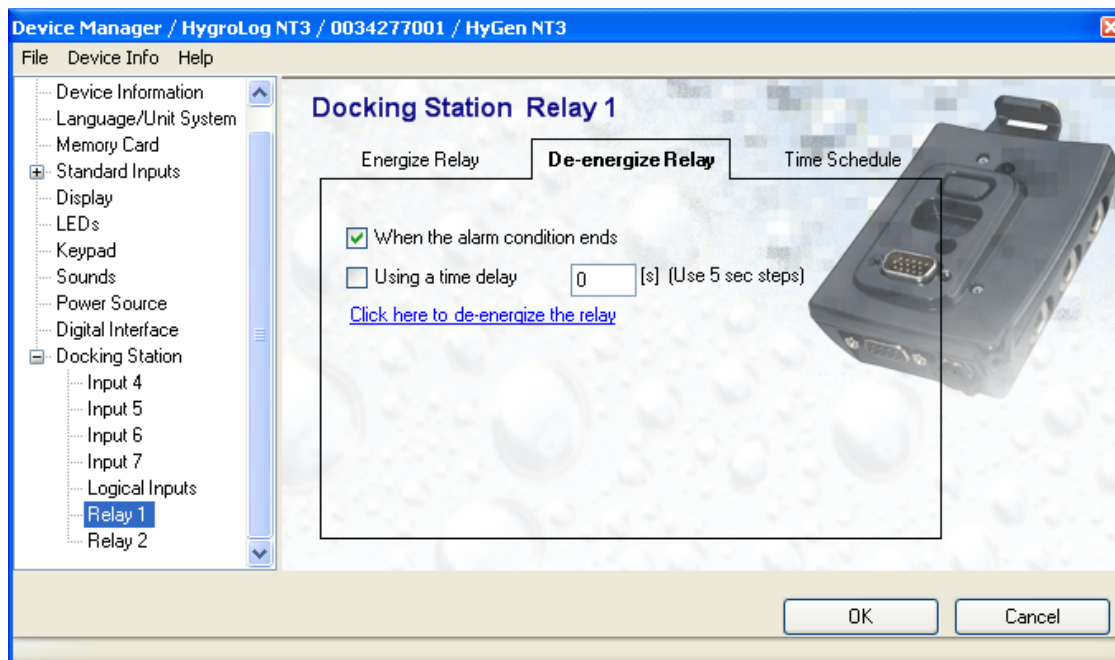De-energize Relay Tab:

Use this tab to define how the relay behaves when the alarm condition(s) ends.


For a latching relay action (relay remains energized after the alarm condition has ended), do not put a check mark in any of the boxes on the form.

For a non-latching relay action, place a check mark in one of the following boxes by left clicking with the mouse while the mouse cursor is on top of the box:

- When the alarm condition ends: the relay is automatically de-energized as soon as the alarm condition ends
- Using a time delay: the relay is automatically de-energized after the specified time delay regardless of the status of the alarm condition that caused the relay to be energized (use multiples of 5 seconds).

Note: when both boxes are check-marked, the relay is automatically de-energized depending on whichever condition occurs first.

To de-energize the relay, click with the mouse on “click here to de-energize the relay”. The relay will not be energized unless a new alarm condition occurs. Note: this link is not visible when the user does not have the “Device configuration / log function” right.

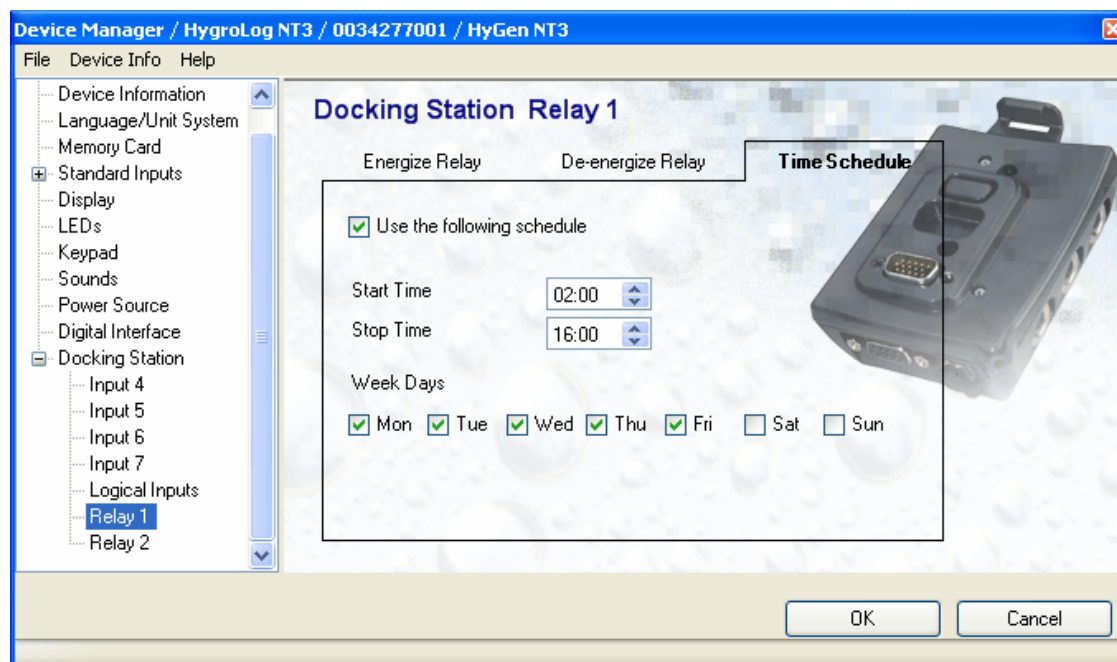Time Schedule Tab:

Use this tab to define when the relay is used to monitor the alarm conditions selected in the Energize Relay Tab.

If you do not wish to use the relay to monitor alarm conditions 7 days a week and 24 hours a day, place a check mark next to "Use the following schedule" and define a time schedule. Both the start and stop time are common to all days of the week.

18 DEVICE MANAGER - Transmitters and indicators


The following devices can be used together with HW4, provided that the firmware version of the device is 4.0a or higher:

- Indicators : HygroPalm (with docking station only) and HygroLab
- Transmitters : M3 (only as RS-485 slave) and HygroFlex

The functions available for these devices are illustrated here, using the example of the HygroFlex 3.

NOTE:

- HygroFlex 1, HygroLab 1, HygroPalm 1 and HygroPalm 0: these devices can be connected to the COM port of a PC only for configuration purposes with HW4. Connecting the HygroPalm 1 or HygroPalm 0 to a PC requires a PalmDock docking station.
- Device Manager always shows the maximum number of configuration options that are available with a specific type of device. All configuration options are not necessarily applicable to all versions of a device type. When configuring a device, the configuration options that do not apply can be safely ignored.

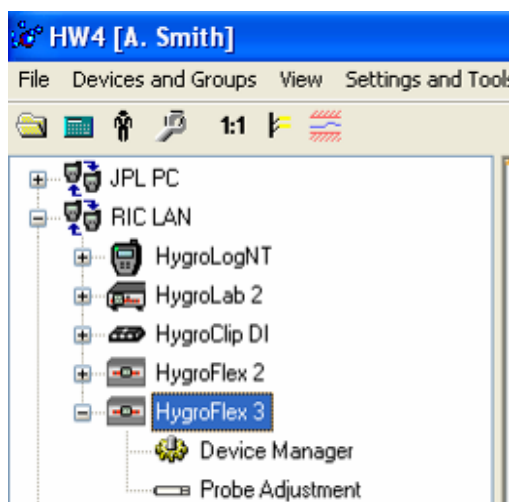

To open the Device manager form, click with the mouse on Device Manager. The Device Manager form is used to configure the HygroFlex 3 and to read instrument specific information. When Device Manager is started, it automatically interrogates the instrument and downloads its current configuration.

The different forms that are available within the Device Manager form are listed in a tree located on the left pane of the form. To select a form, click on it with the left mouse button.

18.1 Device Manager Menu Bar

The Device manager menu bar is located at the top of the form.

File


- **Open:** opens the device configuration directory specified in the HW4 Settings Form - File Locations Tab - and displays all available device configuration files (extension .DAT). Any device configuration file that was previously saved can be opened to quickly configure an instrument. If so desired, any directory and any file type may be opened.
- **Save As:** saves the current configuration to a file (extension .DAT) in the device configuration directory specified in the HW4 Settings Form - File Locations Tab. If so desired, any directory and any file type may be specified.
- **Exit:** exits Device Manager

18.2 Device Information



- **Device name:** this text will be displayed next to the instrument model information. As far as possible use a unique device name.
- **RS485 address:** click on the underlined link to change the instrument address to be used in conjunction with an RS-485 network (multi-drop). Each network address should be unique and within the values of 0 to 63. Note: the default factory RS-485 address is 0. Unless necessary, do not manually modify this address. HW4 will automatically change the RS-485 address of the device, if so required.
- **Device Protection:** this function is available from firmware version 4.1b. For a description, see [Device Protection](#)
- **Fixed pressure value:** enter here the value of the fixed pressure value that will be used by the HygroFlex when calculating parameters such as the mixing ratio (the HygroFlex can be also configured to read pressure from an analog probe). Note that this numerical value should be consistent with the pressure unit that was selected under **Language / Unit System**.
- **Reference below freezing:** select here whether the HygroFlex should provide the dew point or the frost point for values below freezing.

18.3 Language / Unit System

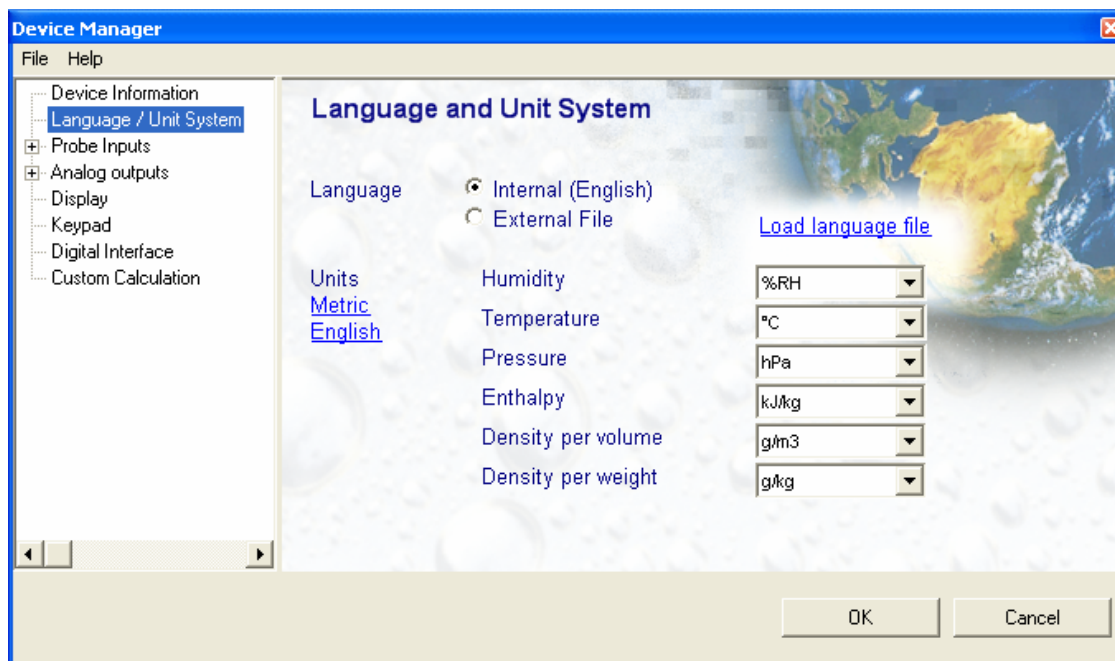

- **Language:** click on the underlined link to change the language used in the device local menus and files to one of the available languages. The link opens the folder where the language files are located (these files have the extension .LNG as in English.LNG). To change the language, simply click on the appropriate file.
- **Units:** use the underlined links to the right of the form to globally change the unit system used by the instrument. If the unit system offers several choices such as In Hg or PSI for pressure, choose a unit by left clicking on the arrow to the right of each text box.

Humidity: the symbol to be used after the % symbol (RH) should be typed in the text box.

Note: the units selected for density per volume and density per weight also apply to vapor concentration and specific humidity.

IMPORTANT: when changing the unit system used by the instrument, be sure to also change the numerical value of the fixed pressure to be consistent with the new unit system (see Device Information). Neglecting to do so will result in an error on the value of any parameter that uses barometric pressure as a calculation input.

18.4 Probe Inputs


Clicking on the plus sign to the left of “Probe Inputs” expands the tree and displays input 1 and input 2 (optional input). Both inputs are listed regardless of whether or not they do exist on the actual transmitter. Selections made in the View tab (HW4 main screen) have no effect.

The following probe types may be used:

- HygroClip probe
- Analog humidity-temperature probe ¹
- Analog probe ¹
- Pressure probe ¹

¹ Must be compatible with the requirements specified in the HygroFlex instruction manual.

18.4.1. HygroClip probe

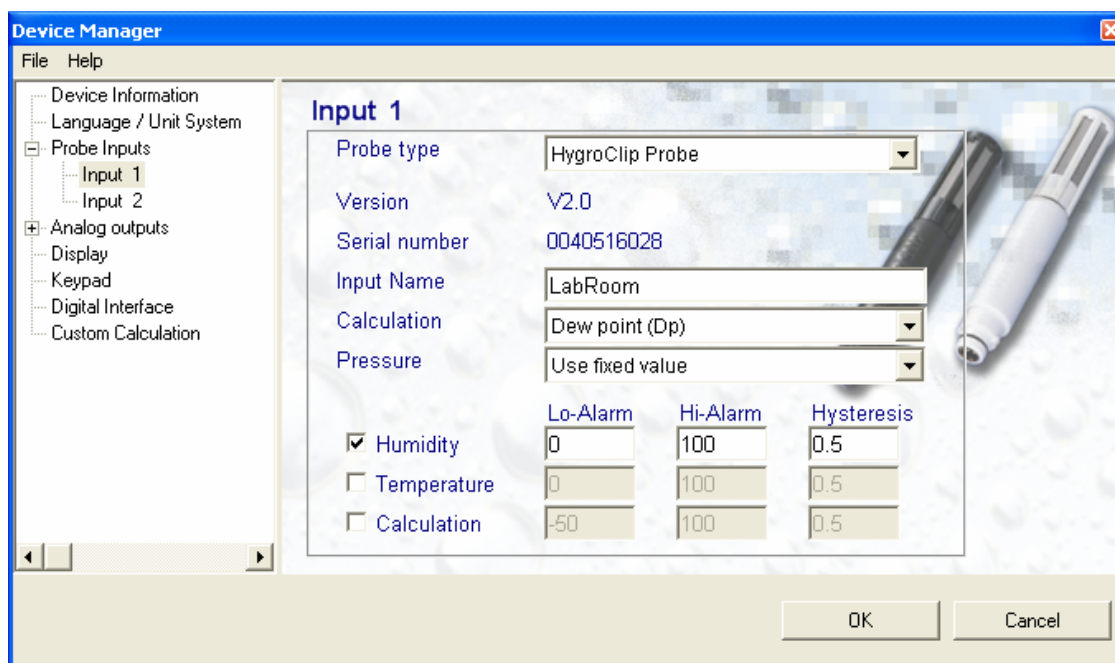


- **Probe type:** select from the pull down menu (click on the arrow)
- **Input Name:** use a maximum of 12 characters
- **Calculation:** left click on the arrow to the right of the list box and select the parameter to be calculated by the HygroFlex for this particular probe. This is also the calculated parameter that HW4 will display for this probe

Dew point (Dp): for values below freezing, the HygroFlex will calculate either the dew point or the frost point, depending on the selection made under Device Information. Both the HygroFlex and HW4 will display either the symbol Dp or the symbol Fp. The symbol Fp indicates that any value below freezing is a frost point as opposed to a dew point. When selected, the symbol Fp is also be displayed for values above freezing. This of course is to be understood as being the same as dew point.

- **Pressure:** left click on the arrow to the right of the list box and select which barometric pressure will be used by the HygroFlex to compute parameters such as wet-bulb, mixing ratio, etc. The HygroFlex can be ordered with a second probe input configured to read the analog signal from a pressure probe.
- **Alarm:** alarm conditions can be defined for humidity, temperature and the calculated parameter. Values that are below the low alarm value or above the high alarm value will trigger an alarm. A value can be specified for the alarm function hysteresis. This value is used for both the low and the high alarm. To trigger an alarm only in the event that a probe is missing, disconnected or not transmitting data, simply enable the alarm for one of the parameters and set the Hi and Lo values to the range of the probe (for example 0 and 100 for humidity). All versions of HW4 will show an out-of-limits value alarm in red on the monitor screen. HW4 Professional can also be configured (HW4 global settings - Alarm settings tab) to display an alarm table and generate a report whenever an out-of-limits condition occurs.

18.4.2. Analog Probe (single channel):

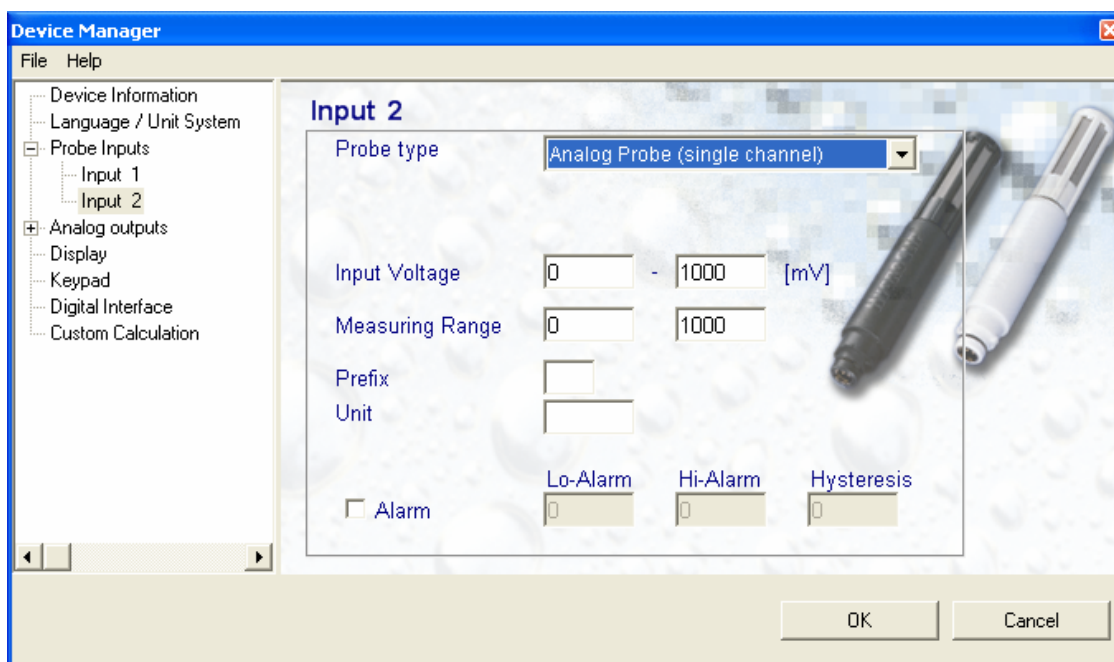


Specify the range of the input voltage (in **mV**), the measuring range of the probe, the symbol to be used on the HygroFlex display (Prefix) and the engineering unit of the parameter measured by the probe. Alarm values can be entered in the same manner as for a HygroClip probe.

Note: In the HW4 main menu bar, use **View - Column Headers** and enable the item “**Value 1**” in order to see data from the probe in the Current Values tab (do this both in Device View mode and in Device View mode).

18.4.3. Pressure probe

The pressure probe is a particular case of the single channel analog probe. The difference is that the choice of unit symbols is limited to commonly used pressure units. The signal from the pressure probe can be used as an input by another probe input (input 1) when calculating a psychrometric parameter that depends on the value of barometric pressure. See example below:

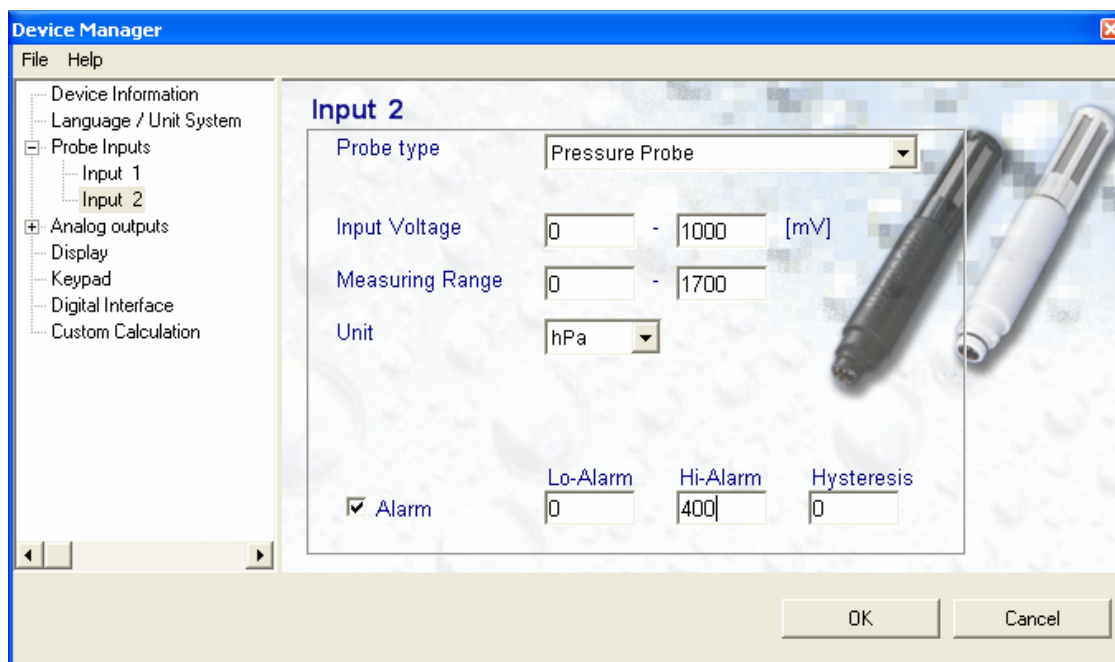


Specify the range of the input voltage (**in mV**), the measuring range of the probe and the engineering unit used by the probe (this may be different from the selection made under Language / Unit System). Alarm values can be entered in the same manner as for a HygroClip probe.

Note: In the HW4 main menu bar, use **View - Column Headers** and enable the item “**Analog In**” in order to see data from the probe in the Current Values tab (do this both in Device View mode and in Group View mode).

18.5 Analog Outputs

The HygroFlex 3 features 3 analog outputs. Each output can be configured individually.

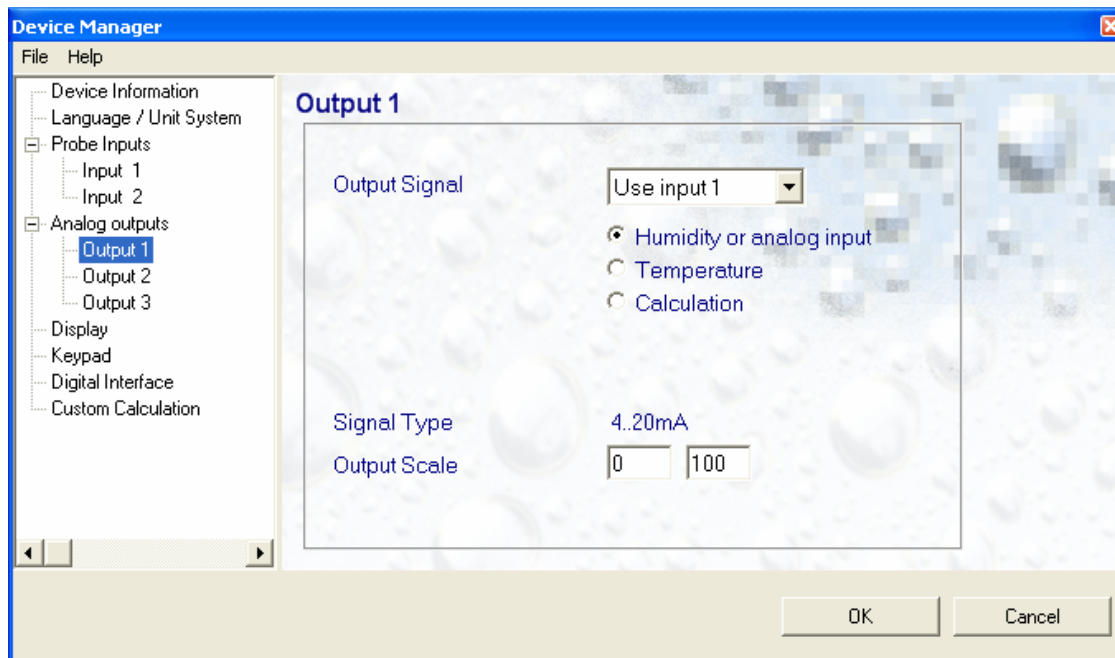

- **Output Signal:** select from the pull down menu (click on the arrow) to associate a probe input - or a custom calculation - with the analog output.

When selecting a probe input, specify which parameter should be used to generate the analog output signal: humidity, temperature or calculated parameter (such as dew point, as defined under Probe Inputs).


Note: if so desired, the same parameter may be used for more than one analog output signal

- **Signal Type:** this is purely informative and cannot be changed here. To change the signal type (4...20mA, 0...5V, etc.) the position of a number of jumpers must be changed on the PCB of the HygroFlex (see separate HygroFlex manual).
- **Output Scale:** enter the numerical values that should correspond to the lower and upper limits of the analog output signal

Example:


- Parameter for the analog output signal = humidity (%RH)
- Output scale: enter 0 and 100 for a full scale output (4...20mA = 0...100%RH)
or 50 and 100 for a partial scale output (4...20mA = 50...100%RH)

18.6 Display

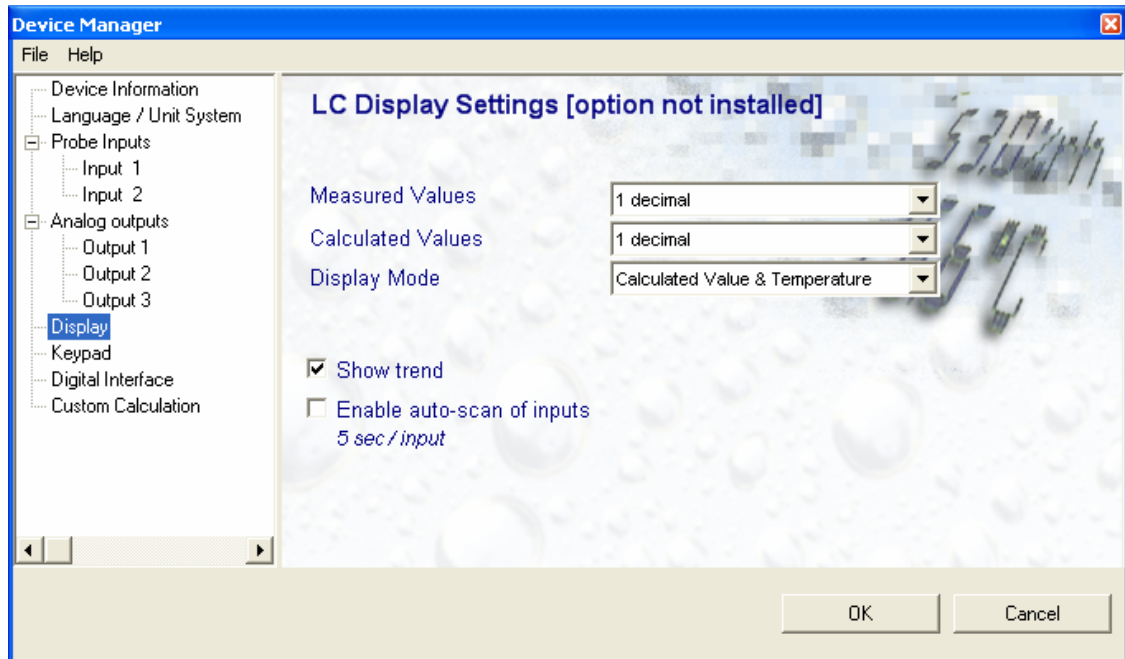


- **Measured values:** left click on the arrow to the right of the list box and select the number of decimals to be used by the optional display for directly measured values such as humidity, temperature, pressure, etc.
- **Calculated Values:** left click on the arrow to the right of the list box and select the number of decimals to be used by the optional display for calculated values such as dew point, custom calculation, etc.
- **Display mode:** left click on the arrow to the right of the list box and select the type of contents for the optional display of the HygroFlex.
- **Show trend:** when this box is checked a trend indicator appears on the optional display to the left of the measured values and calculated parameter. This indicator shows if the values are increasing, decreasing or stable. When the trend indicator is enabled, symbols such as Dp (for dew point) are no longer shown on the display).
- **Enable auto-scan of probe inputs:** check this box to make the display automatically switch between probe inputs every 5 seconds.

18.7 Keypad

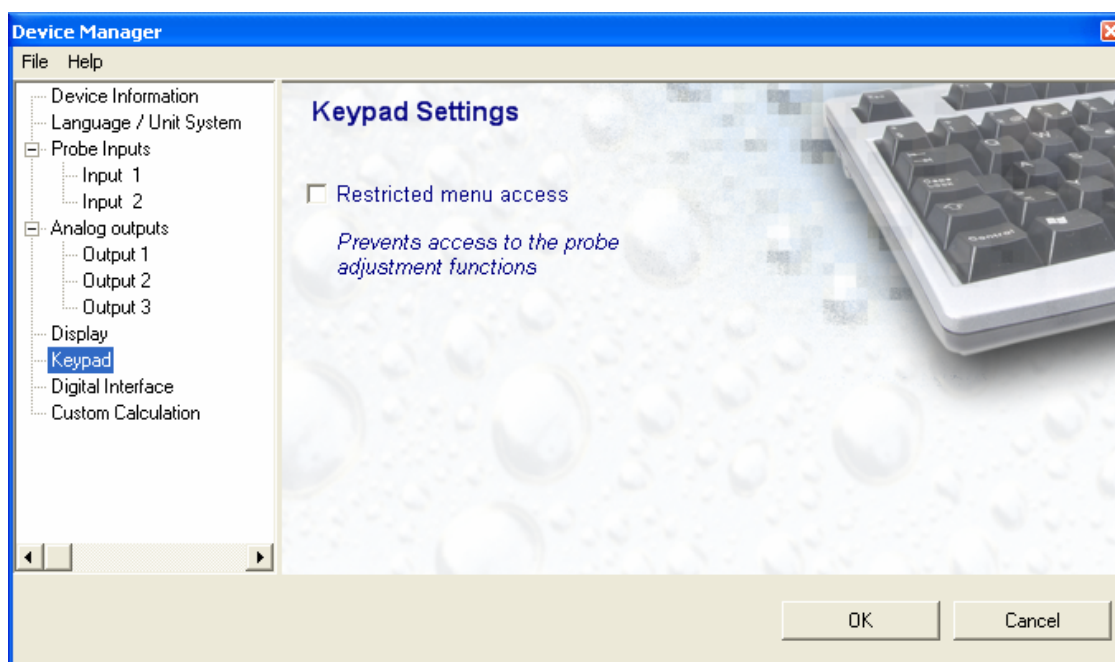

Restricted menu access: check this box to prevent tampering with the probes. Access to all probe adjustment functions from the local keypad of the HygroFlex is blocked out.

18.8 Digital Interface

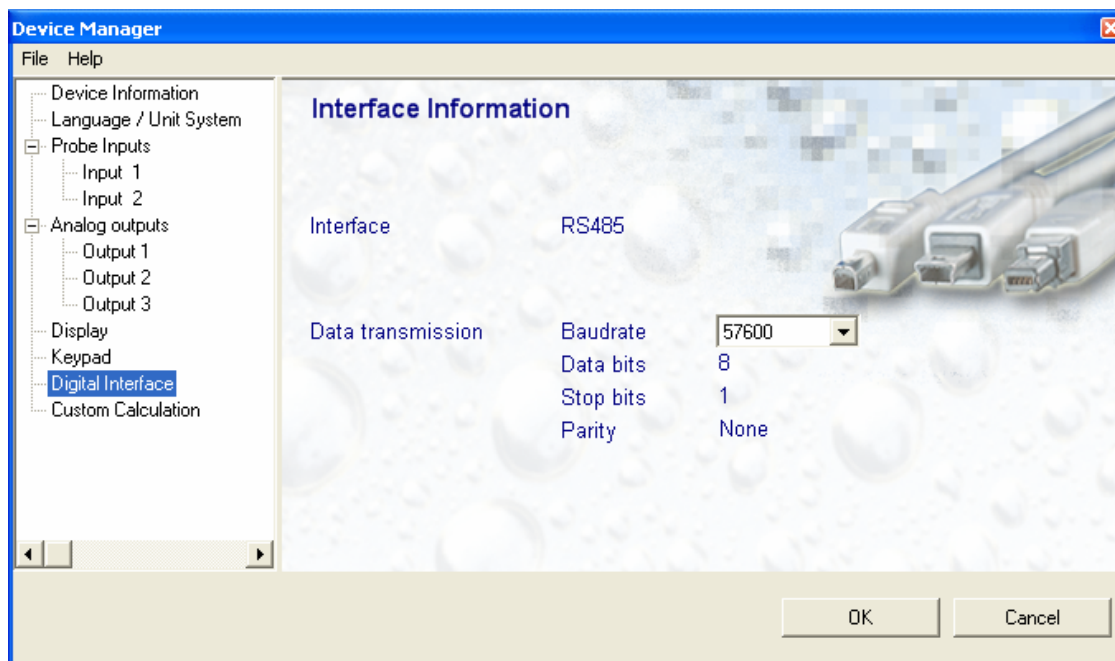


This form is used to select the RS-232 / RS-485 baud rate between 57600 and 19200. Other data cannot be changed from within the form. Generally, use a baud rate of 19200 when you plan on connecting more than 2 or 3 HygroFlex transmitters and similar devices to an RS-485 multi-drop.

IMPORTANT:

- To ensure backward compatibility with the ROTRONIC HW3 software, set the baud rate to 19200
- All instruments within the same RS-485 multi-drop should use the same baud rate. Communication within the RS-485 multi-drop will not work when instruments are configured with different baud rates.

18.9 Custom Calculation

The HygroFlex 3 allows defining a simple custom calculation such as the temperature difference between two probes connected to the HygroFlex.

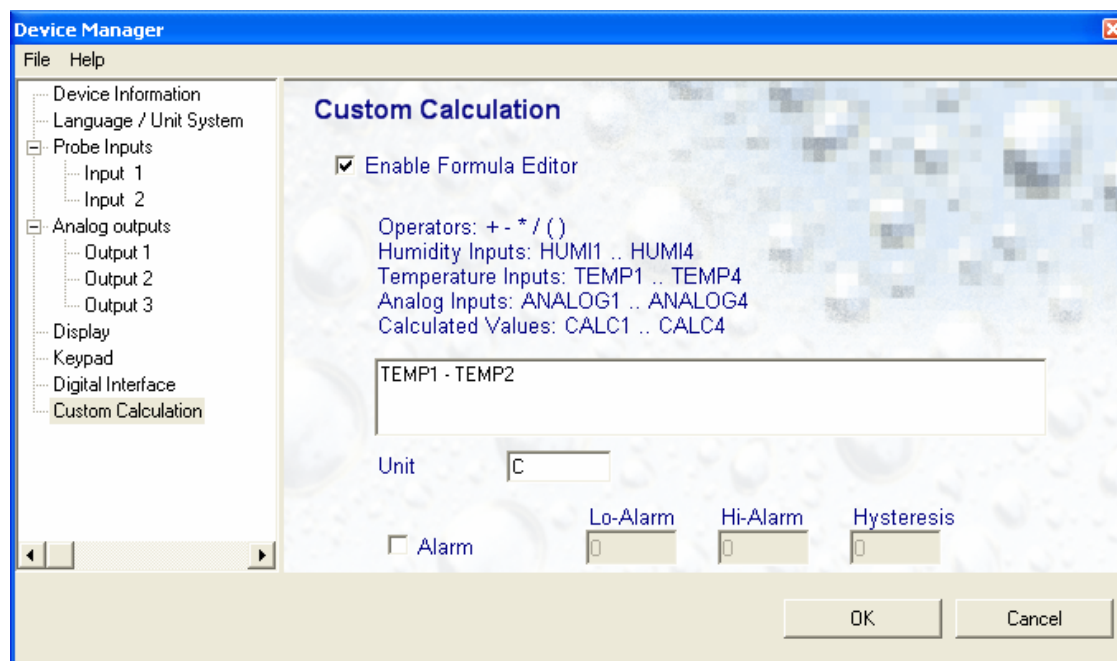

To define a custom calculation, click on the check box next to Enable Formula Editor.


Enter a custom formula in the text box. The parameters that can be used for a custom calculation and the operators are defined directly on the form. The above example illustrates a custom calculation consisting in the temperature difference between probe 1 and probe 2.

Specify the unit symbol to be shown on the HygroFlex display.

Alarm values can be entered in the same manner as for a HygroClip probe.

19 DEVICE MANAGER – HygroClip Digital Interface

The functions available for the HygroClip DI are illustrated below:

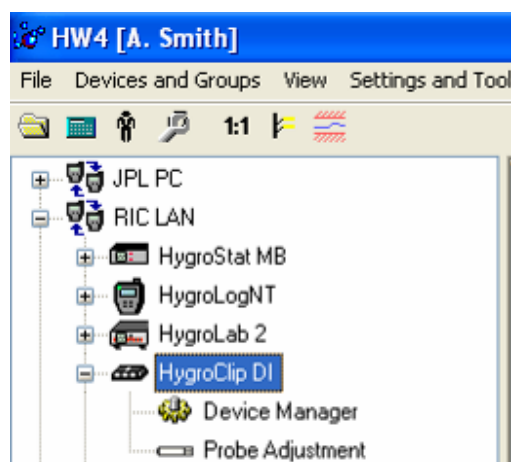

To open the Device manager form, click with the mouse on Device Manager. The Device Manager form is used to configure the HygroClip DI and to read device specific information. When Device Manager is started, it automatically interrogates the device and downloads its current configuration.



The forms that are available within the Device Manager form are listed in a tree located on the left pane of the form. To select a form, click on it with the left mouse button.

19.1 Device Manager Menu Bar

The Device manager menu bar is located at the top of the form.

File


- **Open:** opens the device configuration directory specified in the HW4 Settings Form - File Locations Tab - and displays all available device configuration files (extension .DAT). Any device configuration file that was previously saved can be opened to quickly configure an instrument. If so desired, any directory and any file type may be opened.
- **Save As:** saves the current configuration to a file (extension .DAT) in the device configuration directory specified in the HW4 Settings Form - File Locations Tab. If so desired, any directory and any file type may be specified.
- **Exit:** exits Device Manager

19.2 Device Information



- **Device name:** this text will be displayed next to the instrument model information. As far as possible use a unique device name.
- **RS485 address:** click on the underlined link to change the instrument address to be used in conjunction with an RS-485 network (multi-drop). Each network address should be unique and within the values of 0 to 63. Note: the default factory RS-485 address is 0. Unless necessary, do not manually modify this address. HW4 will automatically change the RS-485 address of the device, if so required.
- **Device Protection:** for a description of this function, see [Device Protection](#)

19.3 Language / Unit System

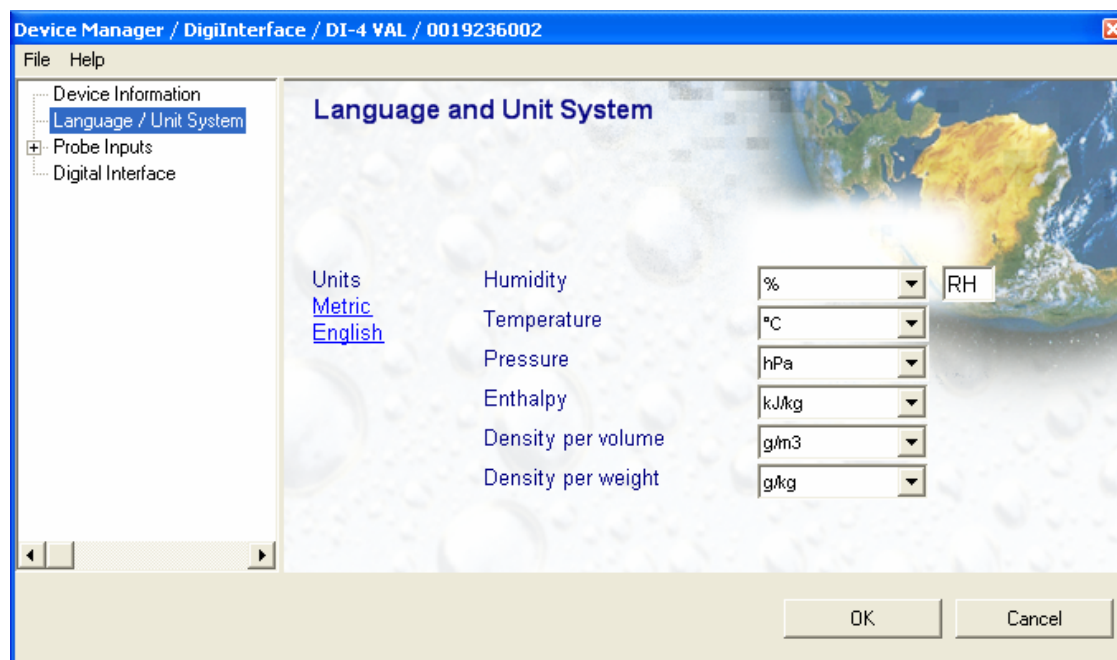

- **Units:** use the underlined links to the right of the form to globally change the unit system used by the instrument. If the unit system offers several choices such as In Hg or PSI for pressure, choose a unit by left clicking on the arrow to the right of each text box.


Humidity: the symbol to be used after the % symbol (RH) should be typed in the text box.

Note: the units selected for density per volume and density per weight also apply to vapor concentration and specific humidity.

IMPORTANT: Do not change the language / units settings while the PC is recording data.

19.4 Probe Inputs

Clicking on the plus sign to the left of “Probe Inputs” expands the tree and displays probe inputs 1 and 2.

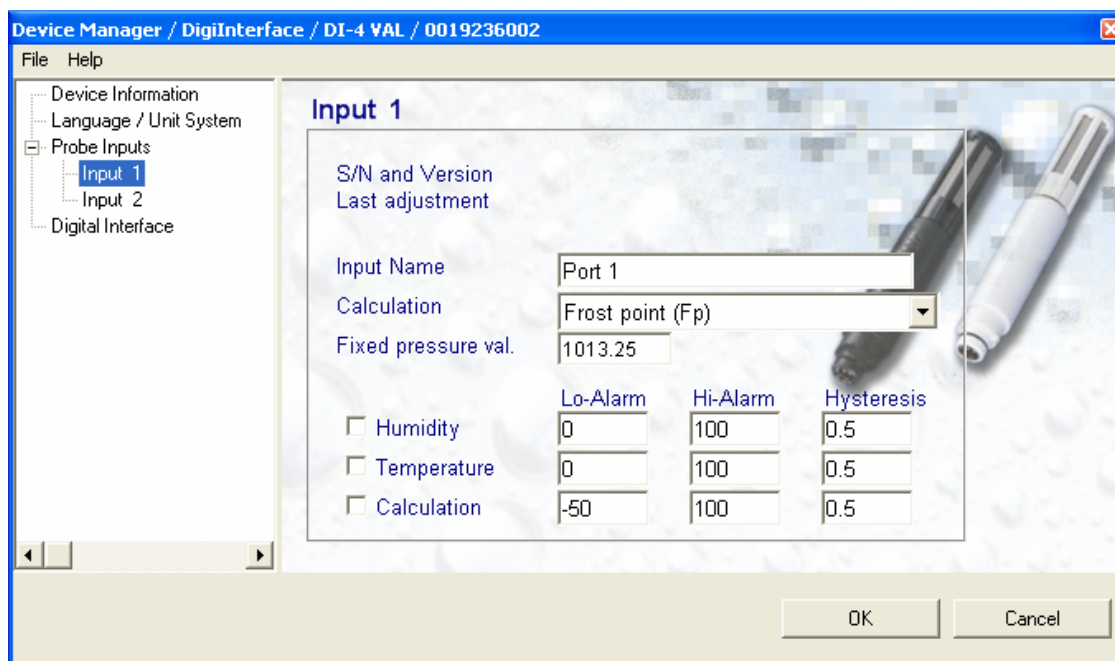


Note: the HygroClip DI accepts only HygroClip probes.

- **Input Name:** use a maximum of 12 characters
- **Calculation:** left click on the arrow to the right of the list box and select the parameter to be calculated by the HygroClip DI for this particular probe. This is also the calculated parameter that HW4 will display for this probe

Dew point (Dp) or Frost point (Fp): depending on the selection made here, HW4 will display either the symbol Dp or the symbol Fp. The symbol Fp indicates that any value below freezing is a frost point as opposed to a dew point. When selected, the symbol Fp is also displayed for values above freezing. This of course is to be understood as being the same as dew point.

- **Fixed pressure value:** enter here the fixed pressure value that will be used by the HygroClip DI. Note that this numerical value should be consistent with the pressure unit that was selected under Language / Unit System.
- **Alarm:** alarm conditions can be defined for humidity, temperature and the calculated parameter. Values that are below the low alarm value or above the high alarm value will trigger an alarm. A value can be specified for the alarm function hysteresis. This value is used for both the low and the high alarm. To trigger an alarm only in the event that a probe is missing, disconnected or not transmitting data, simply enable the alarm for one of the parameters and set the Hi and Lo values to the range of the probe (for example 0 and 100 for humidity). All versions of HW4 will show an out-of-limits value alarm in red on the monitor screen. HW4 Professional can also be configured (HW4 global settings - Alarm settings tab) to display an alarm table and generate a report whenever an out-of-limits condition occurs.

19.5 Digital Interface

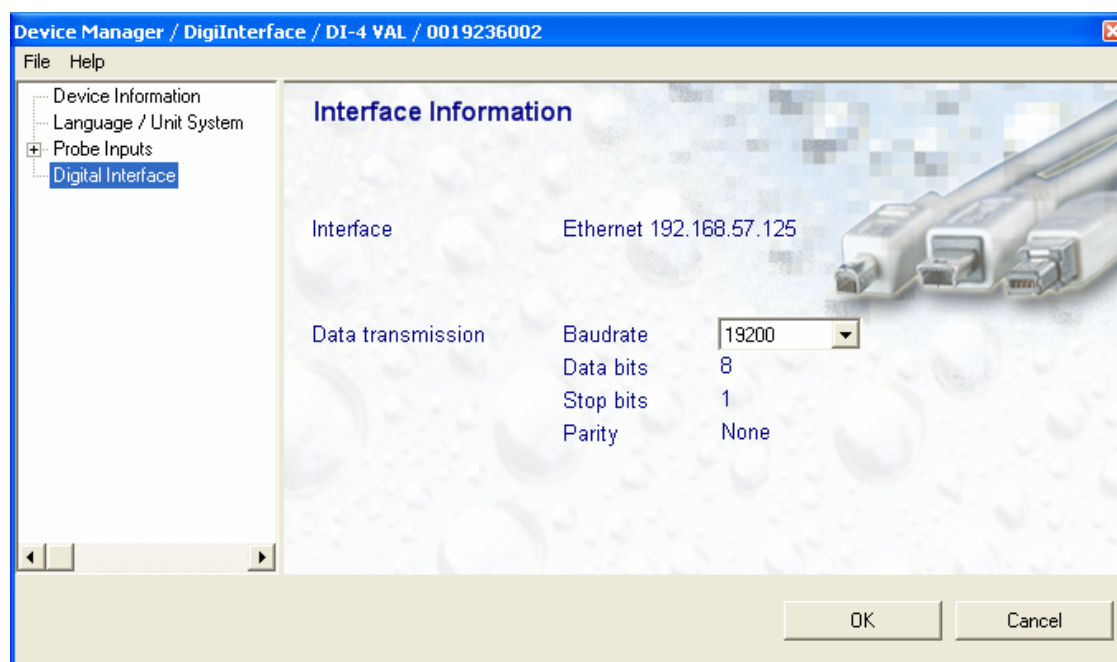

This form is used to select the RS-485 baud rate between 57600 (factory default) and 19200. Other data cannot be changed from within the form. Generally, leave the baud rate at 57600. Select 19200 when you plan on connecting the HygroClip DI and other devices such as the HygroFlex transmitter to an RS-485 multi-drop.


IMPORTANT:

All devices within an RS-485 multi-drop should use the same baud rate. Communication within the RS-485 multi-drop will not work when the devices are configured with different baud rates.

After changing the baud rate in Device Manager, you should also reflect the change in the configuration of the internal Digi International module used by the device to connect to the LAN. See: [Changing the baud rate of an Ethernet device](#)

20 DEVICE MANAGER – HygroClip Alarm

Device Manager is the only function available with the HygroClip Alarm (HCA). This is illustrated below:

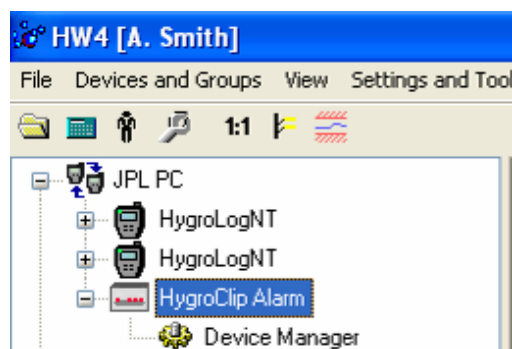

To open the Device manager form, click with the mouse on Device Manager. The Device Manager form is used to configure the HygroClip Alarm and to read device specific information. When Device Manager is started, it automatically interrogates the device and downloads its current configuration.



The forms that are available within the Device Manager form are listed in a tree located on the left pane of the form. To select a form, click on it with the left mouse button.

20.1 Device Manager Menu Bar

The Device manager menu bar is located at the top of the form.

File

- **Open:** opens the device configuration directory specified in the HW4 Settings Form - File Locations Tab - and displays all available device configuration files (extension .DAT). Any device configuration file that was previously saved can be opened to quickly configure an instrument. If so desired, any directory and any file type may be opened.
- **Save As:** saves the current configuration to a file (extension .DAT) in the device configuration directory specified in the HW4 Settings Form - File Locations Tab. If so desired, any directory and any file type may be specified.
- **Exit:** exits Device Manager

20.2 Device Information

Device Manager

File Help

Device Information

Device Type: HygroClip Alarm

Firmware Version: V2.1a

Serial Number: 0056258001

Device Name: HCA Test

RS-485 Address: 2 [Change RS address](#)

Device Time: 11/27/2006 12:57:04 [Synchronize with PC clock](#)

Device Protection: Disabled [Protect device](#)

☒ Alarm when Event Memory is Full [Open Event Viewer](#)
[Clear Event Memory](#)

OK Exit

- **Device name:** this text will be displayed next to the instrument model information. As far as possible use a unique device name.
- **RS485 address:** click on the underlined link to change the instrument address to be used in conjunction with an RS-485 network (multi-drop). Each network address should be unique and within the values of 0 to 63. Note: the default factory RS-485 address is 0. Unless necessary, do not manually modify this address. HW4 will automatically change the RS-485 address of the device, if so required.
- **Device Protection:** for a description of this function, see [Device Protection](#)


- **Open Event Viewer / Clear Event Memory:** click on one of the underlined links to view (and print) the events recorded by the device memory or to clear the memory of all recorded events. For more details, see: [HygroClip Alarm event memory and event viewer](#).
- **Alarm when Event Memory is Full:** check this box to have the HCA send an alarm to the PC whenever the event memory is full.

20.3 Logical Inputs

The screenshot shows the 'Device Manager' window with the 'Logical inputs' tab selected. The window has a menu bar with 'File' and 'Help'. On the left is a tree view with 'Device Information', 'Logical inputs' (selected), 'Relays', 'Programming', 'Digital Interface', and 'Buttons & Switches'. The main area is titled 'Logical inputs' and contains a 'Description' field with the text 'L-Inputs'. Below this are two checked checkboxes: 'L-Input 1' and 'L-Input 2'. To the right of these checkboxes is a section titled 'Logical Input Alarm' with two dropdown menus, both set to '1'. At the bottom right are 'OK' and 'Exit' buttons.


- **Description:** this text will be displayed by HW4 in the current values table.
- **Logical Input 1 / Logical Input 2:** put a check mark in the box next to Logical Input 1 and/or Logical Input 2 to enable monitoring of the logical input alarm status by HW4. The logical state that triggers an alarm can be defined under “Logical Input Alarm”.

20.4 Relays

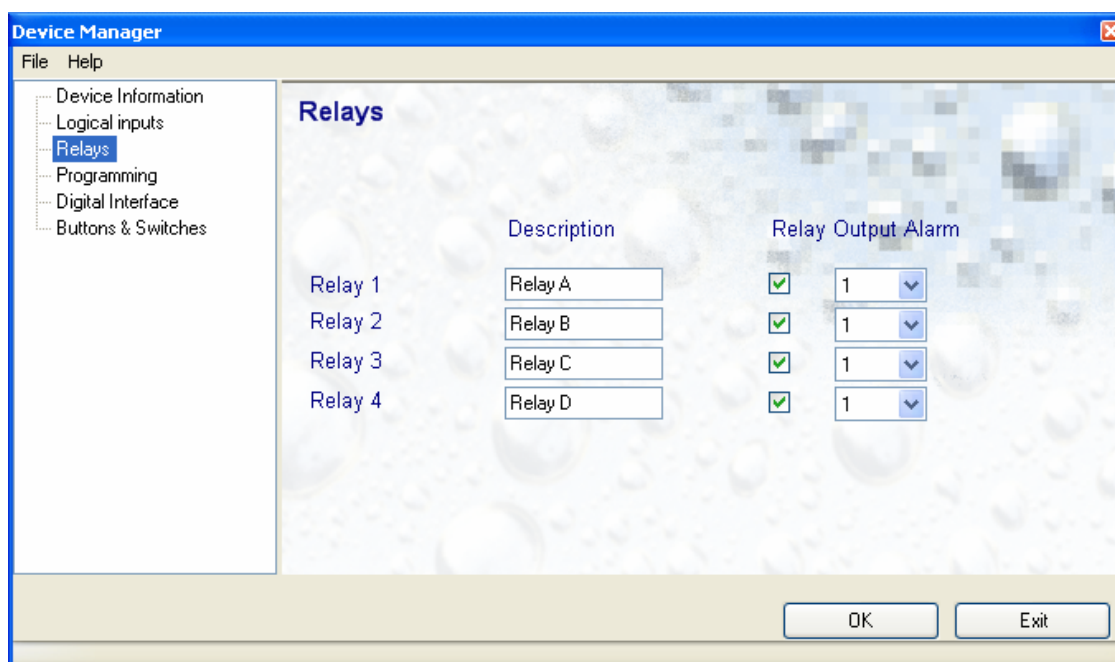

- **Description:** this text will be displayed by HW4 in the current values table.
- **Relay Output Alarm:** put a check mark in the box next to each relay to enable monitoring of the relay alarm status by HW4. The logical state that triggers an alarm can be defined under "Relay Output Alarm" (logical 1 means that the relay is energized).

20.5 Programming

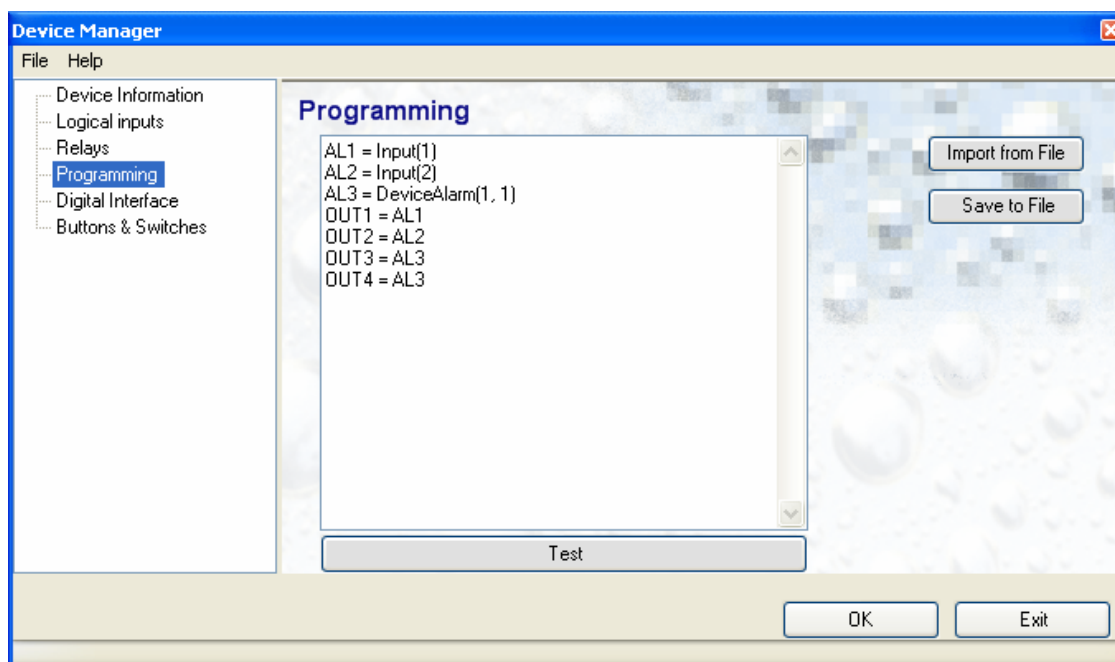


Program statements can be typed directly in the white area located under the heading "Programming". The functionality of this area is equivalent to a simple text editor.

For details regarding the HCA programming language, see separate HygroClip Alarm manual.


- **Import from File:** use this button to import a program from the PC (normally, programs are files with the TXT extension and are located in the ROTRONIC_HW4\SYS subdirectory).
- **Save to File:** use this button to save the current program to a file with the TXT extension (normally this file will be saved on the PC in the ROTRONIC_HW4\SYS subdirectory).

20.6 Digital Interface

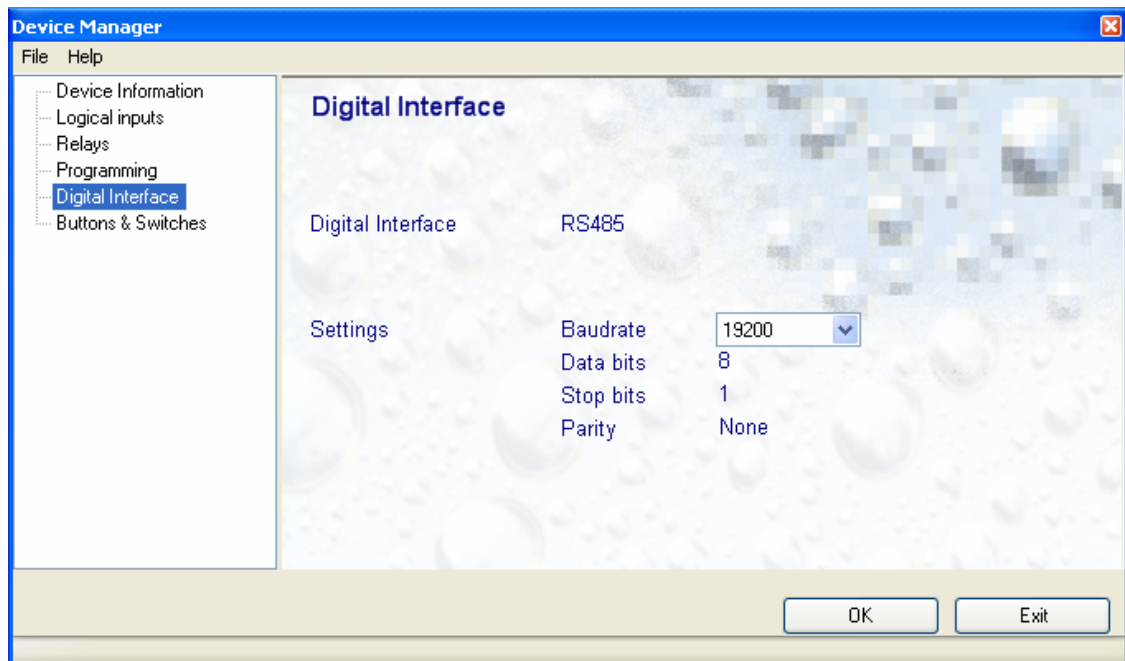

- **Baud rate:** use the box to select either 19200 (recommended) or 57600 (for use with HygroLog NT and HygroClip Interface only)

20.7 Buttons and Switches

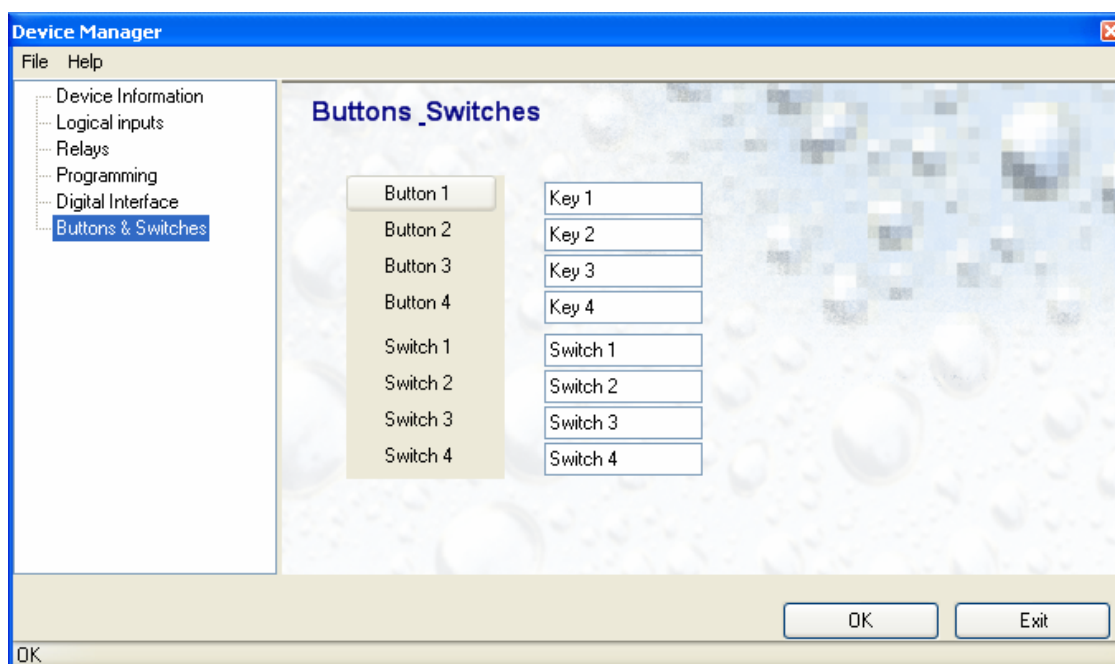


A virtual push button can be actuated by clicking on it with the mouse.

A virtual switch can be closed by clicking on it with the mouse as illustrated above for switch 1.

Click on the Exit button after clicking on a push button or changing the position of a switch.

21 DEVICE MANAGER – HygroStat MB

The functions available for the HygroStat MB are illustrated below:

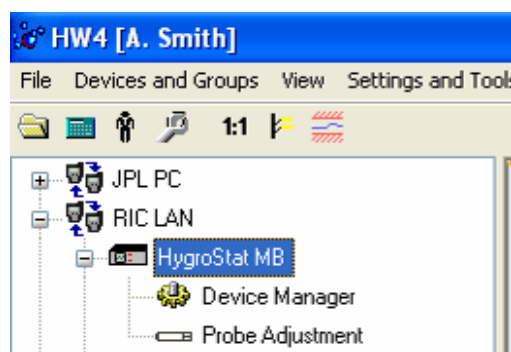

To open the Device manager form, click with the mouse on Device Manager. The Device Manager form is used to configure the HygroClip Alarm and to read device specific information. When Device Manager is started, it automatically interrogates the device and downloads its current configuration.



The forms that are available within the Device Manager form are listed in a tree located on the left pane of the form. To select a form, click on it with the left mouse button.

21.1 Device Manager Menu Bar

The Device manager menu bar is located at the top of the form.

File

- **Open:** opens the device configuration directory specified in the HW4 Settings Form - File Locations Tab - and displays all available device configuration files (extension .DAT). Any device configuration file that was previously saved can be opened to quickly configure an instrument. If so desired, any directory and any file type may be opened.
- **Save As:** saves the current configuration to a file (extension .DAT) in the device configuration directory specified in the HW4 Settings Form - File Locations Tab. If so desired, any directory and any file type may be specified.
- **Exit:** exits Device Manager

21.2 Device Information

Device Manager / HygroStat MB / 0087878001 / HygroStat MB

File Help

- Device Information
- Language / Unit System
- Probe Inputs
- Relays
- Display
- Digital Interface

Device Information

Device type: HygroStat MB

Firmware version: V1.0a

Serial number: 0087878001

Device Name: HygroStat MB
max. 30 characters

RS485 address: 5 [Change RS address](#)


Device Protection: Disabled [Protect device](#)

OK Cancel

- **Device name:** this text will be displayed next to the instrument model information. As far as possible use a unique device name.
- **RS485 address:** click on the underlined link to change the instrument address to be used in conjunction with an RS-485 network (multi-drop). Each network address should be unique and within the values of 0 to 63. Note: the default factory RS-485 address is 0. Unless necessary, do not manually modify this address. HW4 will automatically change the RS-485 address of the device, if so required.
- **Device Protection:** for a description of this function, see [Device Protection](#)

Note: the HygroStat MB internal clock is automatically synchronized with the HW4 PC clock when clicking on the Device Manager OK button.

21.3 Language / Unit System

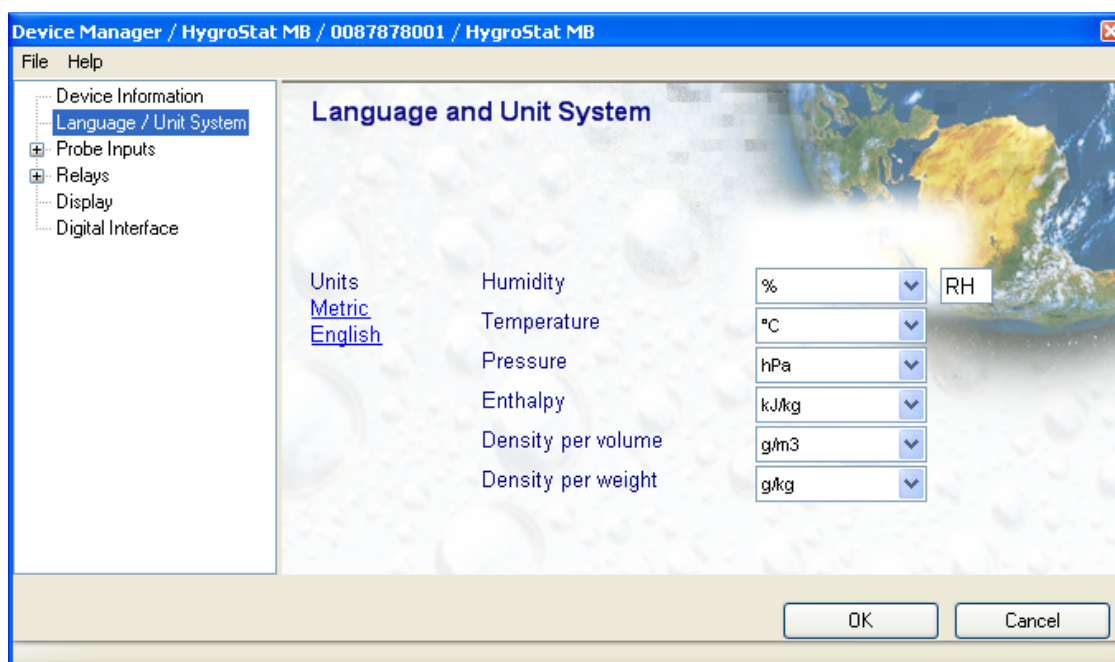

- **Units:** use the underlined links to the right of the form to globally change the unit system used by the instrument. If the unit system offers several choices such as In Hg or PSI for pressure, choose a unit by left clicking on the arrow to the right of each text box.

Humidity: the symbol to be used after the % symbol (RH) should be typed in the text box.

Note: the units selected for density per volume and density per weight also apply to vapor concentration and specific humidity.

IMPORTANT: Do not change the language / units settings while the PC is recording data.

21.4 Probe Inputs

Clicking on the plus sign to the left of “Probe Inputs” expands the tree and displays the probe input.

Device Manager / HygroStat MB / 0087878001 / HygroStat MB

File Help

Device Information
Language / Unit System
Probe Inputs
 Input 1
Relays
Display
Digital Interface

Input 1

S/N and Version 40516028 / V2.0
Last adjustment 12/5/2006 (87878001)

Input Name IW Pbe
Calculation Frost point (Fp)
Fixed pressure val. 1013.25

	Lo-Alarm	Hi-Alarm	Hysteresis
☐ Humidity	85	90	1
☐ Temperature	30	40	1
☐ Calculation	23	30	2

OK Cancel

Note: the HygroStat MB accepts only HygroClip probes.


- **Input Name:** use a maximum of 12 characters
- **Calculation:** left click on the arrow to the right of the list box and select the parameter to be calculated by the HygroStat MB. This is also the calculated parameter that HW4 will display for this probe

Dew point (Dp) or Frost point (Fp): depending on the selection made here, HW4 will display either the symbol Dp or the symbol Fp. The symbol Fp indicates that any value below freezing is a frost point as opposed to a dew point. When selected, the symbol Fp is also displayed for values above freezing. This of course is to be understood as being the same as dew point.

- **Fixed pressure value:** enter here the fixed pressure value that will be used by the HygroStat MB. Note that this numerical value should be consistent with the pressure unit that was selected under Language / Unit System.
- **Alarm:** These settings affect HW4 only and not the relays. Alarm conditions can be defined for humidity, temperature and the calculated parameter. Values that are below the low alarm value or above the high alarm value will trigger an alarm. A value can be specified for the alarm function hysteresis. This value is used for both the low and the high alarm. To trigger an alarm only in the event that a probe is missing, disconnected or not transmitting data, simply enable the alarm for one of the parameters and set the Hi and Lo values to the range of the probe (for example 0 and 100 for humidity). All versions of HW4 will show an out-of-limits value alarm in red on the monitor screen. HW4 Professional can also be configured (HW4 global settings - Alarm settings tab) to display an alarm table and generate a report whenever an out-of-limits condition occurs.


21.5 Relays

Clicking on the plus sign to the left of “Relays” expands the tree and displays Relay 1 and Relay 2.

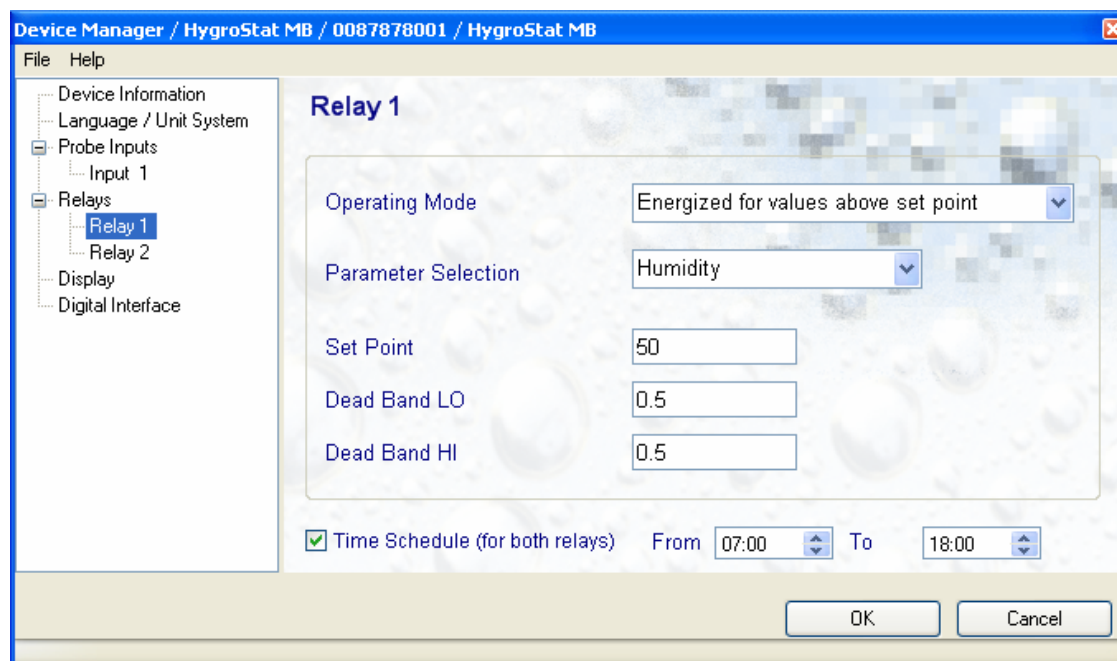


- **Operating Mode:** left click on the arrow to the right of the list box and select either “Energized for values above set point” or “Energized for values below set point”.
- **Parameter Selection:** left click on the arrow to the right of the list box and select the parameter to be monitored by the HygroStat MB relay: humidity, temperature or calculated parameter.
- **Set point:** enter here the value of the set point for the selected parameter.
- **Dead Band LO:** enter here the value of the low dead band. Depending on the operating mode, the relay is either energized or de-energized when the selected parameter is $\leq$ than the set point minus the value of the low dead band.
- **Dead Band HI:** enter here the value of the high dead band. Depending on the operating mode, the relay is either energized or de-energized when the selected parameter is $\geq$ than the set point plus the value of the high dead band.
- **Time Schedule:** if so desired, operation of both relays can be restricted to a time bracket.

21.6 Display

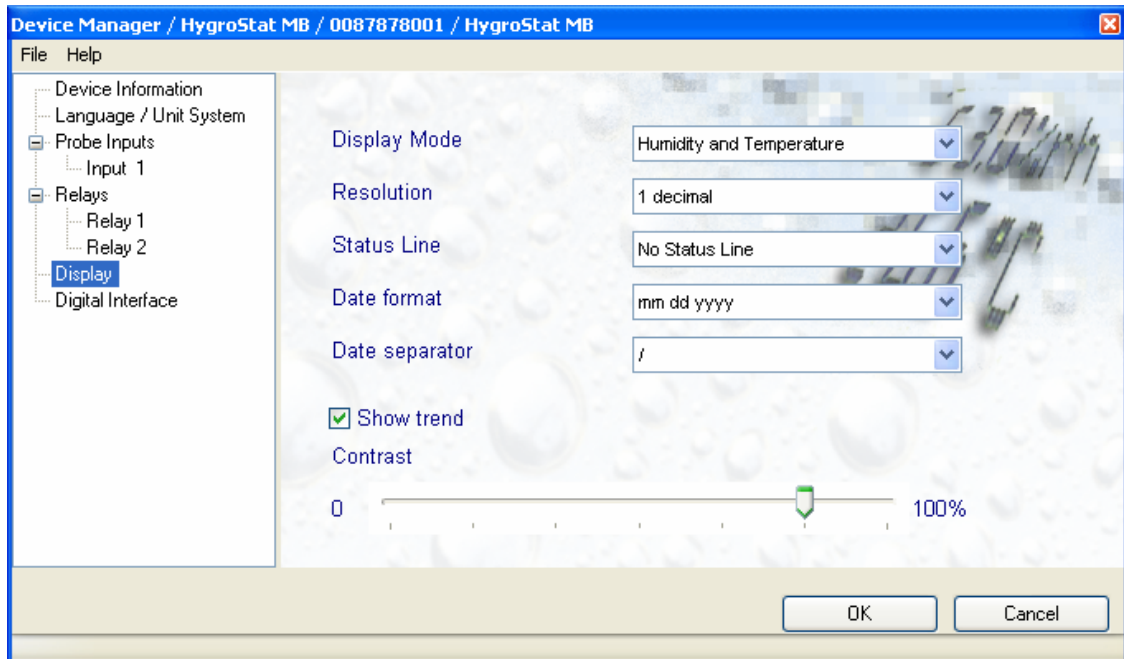

- **Display Mode:** left click on the arrow to the right of the list box and select the two parameters to be shown on the display of the HygroStat MB: humidity and temperature, humidity and calculated value, temperature and calculated value.
- **Resolution:** left click on the arrow to the right of the list box and select the number of decimals to be used on the display of the HygroStat MB: one, two or no decimal.
- **Status Line:** left click on the arrow to the right of the list box and select the information to be shown on the status line of the HygroStat MB display: no status line, time, date or input name.
- **Date Separator:** left click on the arrow to the right of the list box and select the separator symbol to be used when displaying the date.
- **Contrast:** drag the slider to adjust the contrast of the display

21.7 Digital Interface

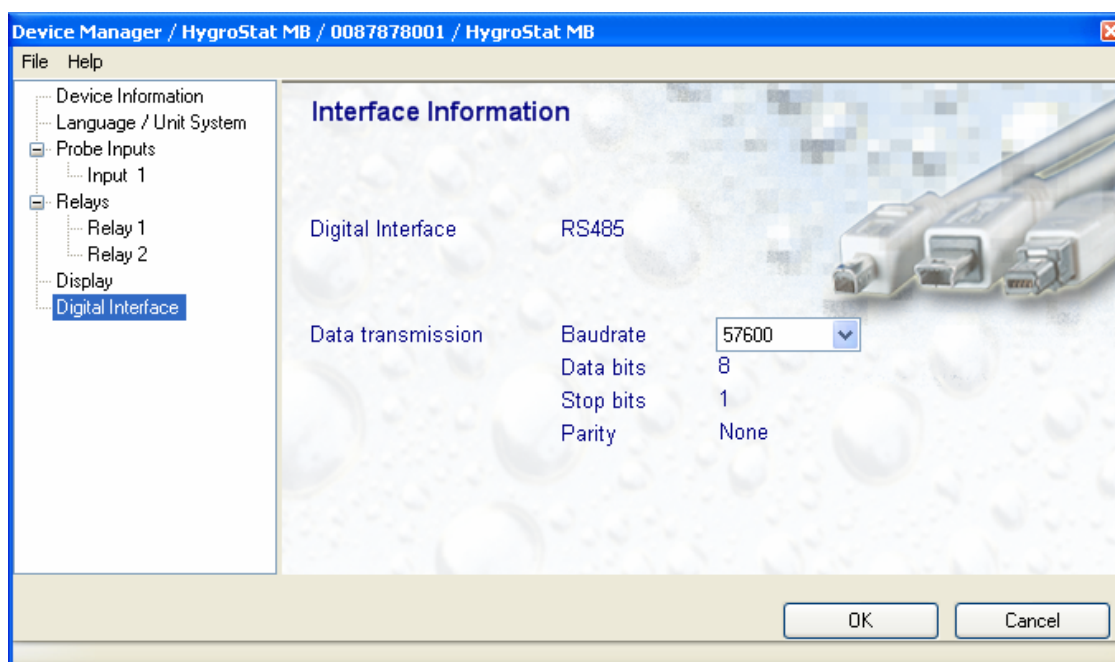

This form is used to select the RS-485 baud rate between 57600 (factory default) and 19200. Other data cannot be changed from within the form. Generally, leave the baud rate at 57600. Select 19200 when you plan on connecting the HygroStat MB and other devices such as the HygroFlex transmitter to an RS-485 multi-drop.

IMPORTANT: All devices within an RS-485 multi-drop should use the same baud rate. Communication within the RS-485 multi-drop will not work when the devices are configured with different baud rates.

22 DEVICE PROTECTION

This function is used to protect a device ¹ against malicious users. This is particularly useful when the device is exposed to the Internet. When a device is protected, the following functions are disabled:

- All Device Manager functions
- Log function programming (HygroLog NT)
- Probe adjustment
- Deletion of files from the device memory card (HygroLog NT)

¹ Availability of this function depends on the firmware version of the device.

Warning:

The following HW4 automatic functions are also disabled for the protected device:


- Automatic address change when adding new devices to a RS-485 multi-drop network
- Automatic synchronization with the HW4 PC clock

The procedure used to protect or unprotect a device is described, using the HygroClip DI as an example:


Select the device in the device tree and open Device Manager > Device Information:




To protect the device, click on the underlined link. HW4 opens the following form where a password can be entered (maximum 8 characters):

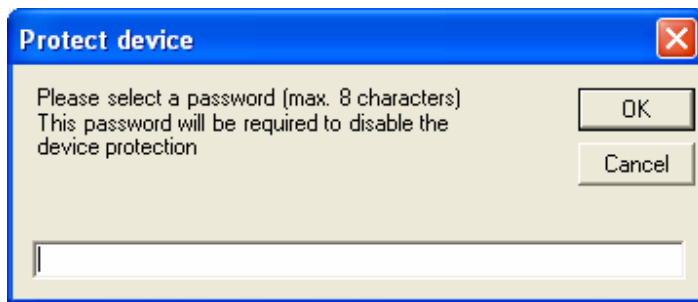

Enter the password and click on the OK button. HW4 confirms that the device is now protected:

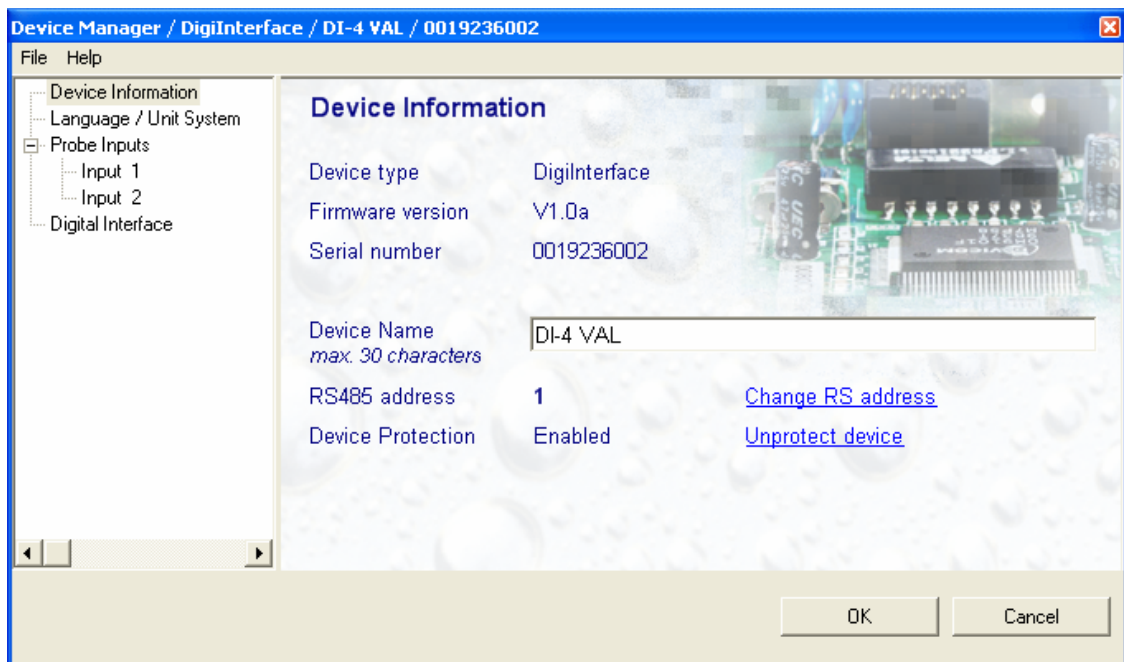


To unprotect the device and enable all functions, click on the underlined link. HW4 opens the password form:

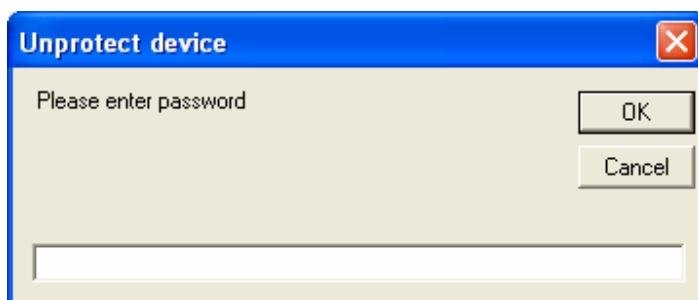

Enter the password and click on the OK button. HW4 confirms that the device is now unprotected.

FORGOT THE PASSWORD? - should you forget the password used to protect the device, remove power from the device (HygroLog NT: remove the battery). After restoring power to the device, you have about one minute to use the default password **!resume!** and unprotect the device (include the exclamation marks). After one minute the default password is not accepted.

23 DATA LOGGING – with the HygroLog NT

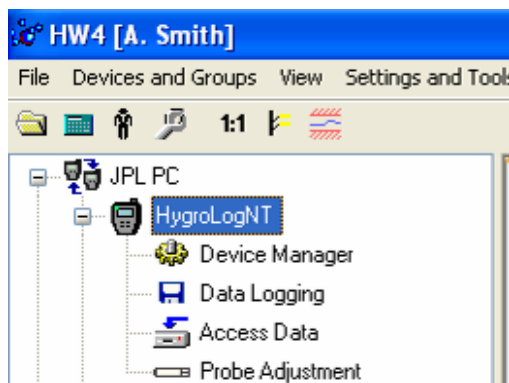


When HW4 has detected a HygroLog NT, the device appears in the left pane of the HW4 main screen. Expanding the device displays a list of the available device functions.

To select the Data Logging function, click on it with the left mouse button. HW4 opens the Data Logging form.

The Data Logging form is used to set the date and time to start and stop logging data for each individual input of the HygroLog NT and docking station. Other settings such as the log interval and file management are also defined with this form.

Clicking on the plus sign to the left of “Standard Inputs” or “Optional Inputs” expands the tree and displays a list of probe and logical inputs. Inputs are listed regardless of whether or not they do exist on the actual logger and docking station being used. Only those items selected in the View Tab are visible in Device Manager. When the list is either empty or incomplete, please exit Data Logging and check the selections made in the View Tab (see Right Pane: Device View Mode). To select an input, click on it with the left mouse button.

Clicking on Standard Inputs, displays information on the memory card and battery at the bottom of the form.

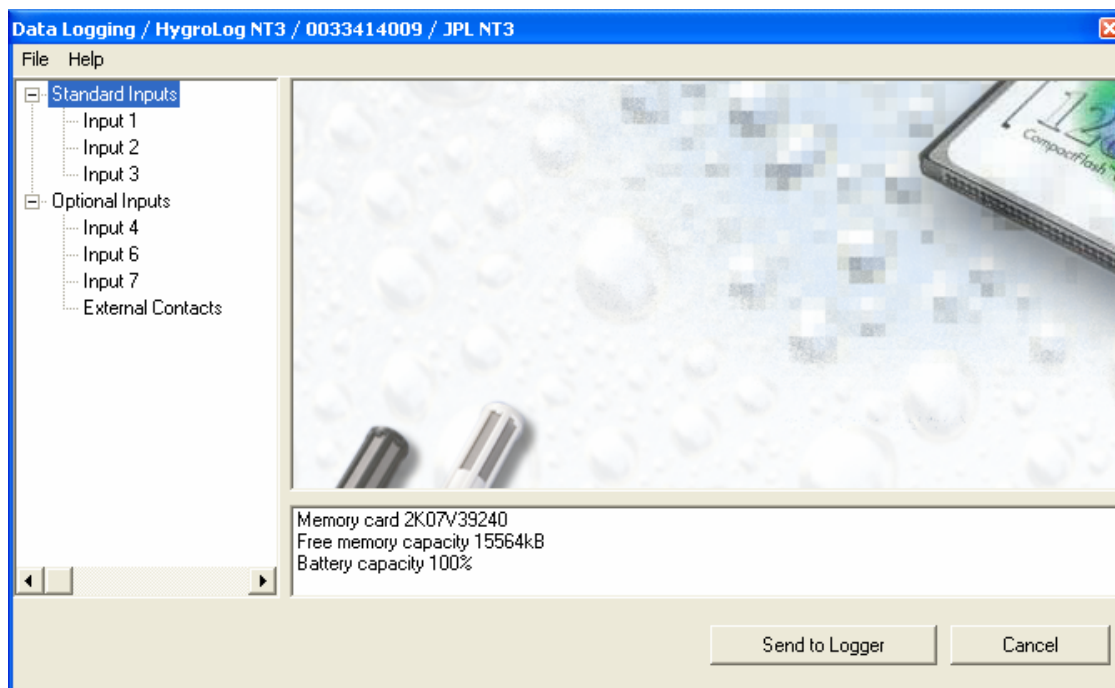

Note: the Data Logging form is not used to define the type of log file (text or binary). The type of log file is common to all probes and inputs and is defined from within Device Manager / Memory Card.

23.1 Menu Bar

The menu bar is located at the top of the form.

File

- **Exit:** closes the Data Logging form without making any changes.

Help

- **HW4 Help:** Opens HW4 Help
- **About HW4:** Displays the version number and ID number of HW4

23.2 Programming the HygroLog NT for data logging

To program the HygroLog NT for logging data, select first an input with the mouse.

- **Start:** enter the start date and the start time in each text box. To copy the current date and time from the PC, click on the “Right Now” blue link.

- **Stop:** enter the stop date and the stop time in each text box. To copy the current date and time from the PC, click on the “Right Now” blue link. **To stop logging immediately, click on this link, click on “Use these settings for all inputs” (only if you want to stop all other inputs) and click on the Send to Logger button.**

- **Log Interval:** use the mouse to highlight any of the following: hours, minutes or seconds. Use the up or the down arrow to set the log interval to the desired value.

- **New Log File:** it is good practice to limit the size of log files so as to avoid both the inconvenience and problems associated with having to download gigantic files from the logger. While the default setting is to close the file and start a new one after 200,000 data points, it often makes sense to start a new file sooner (for example every day). This decision depends

primarily on the log interval. The New Log File list box is used to tell the logger when to finish logging to the current file and when to start a new file. Click on the arrow located to the right of the list box to display a list of available options.

Note: the HygroLog NT creates a separate log file for each probe or input to be logged.

- **Description:** use up to 32 characters to enter a description / identification for the log file. This will be part of the file header.
- **Log:** check the box corresponding to each parameter to be logged (up to 3 parameters). In the case of an input with a single channel analog probe or with a pressure probe, you only need to check the top box. In the case of an RTD probe, you only need to check the Temperature box. In any case, no harm is done when the other boxes are also checked.


Notes:

- The selections made for one input can be transferred to all other inputs (except External Contacts) by clicking on **“Use these settings for all inputs”**. This includes the log file description. Alternatively, each probe or input can be individually programmed, using different settings.
- When a digital HygroClip probe is connected to the input, the serial number of the probe is displayed at the bottom of the form.
- When done, click with the mouse on the **Send to Logger** button. The HygroLog NT will start executing the log commands.
- The HygroLog NT automatically gives each log file a name comprised of the last 4 digits of the logger serial number, followed by the input number and a sequential run number.

23.2.1. Recording the status of external contacts (logical inputs)

Logical inputs are used to monitor the status of up to two external contacts, such as a door contact.

The HygroLog NT records the status of the logical inputs in a different manner as used for other input types. Recording is not based on a log interval. When a check mark is placed in one or two of the “Contact” boxes, data logging is triggered automatically each time that a logical input changes its value from 0 to 1 or from 1 to 0, and holds the new value for 5 or more seconds. In the situation where a logical input has either the value 0 or the value 1 during the entire recording time, the log file does not record any data.

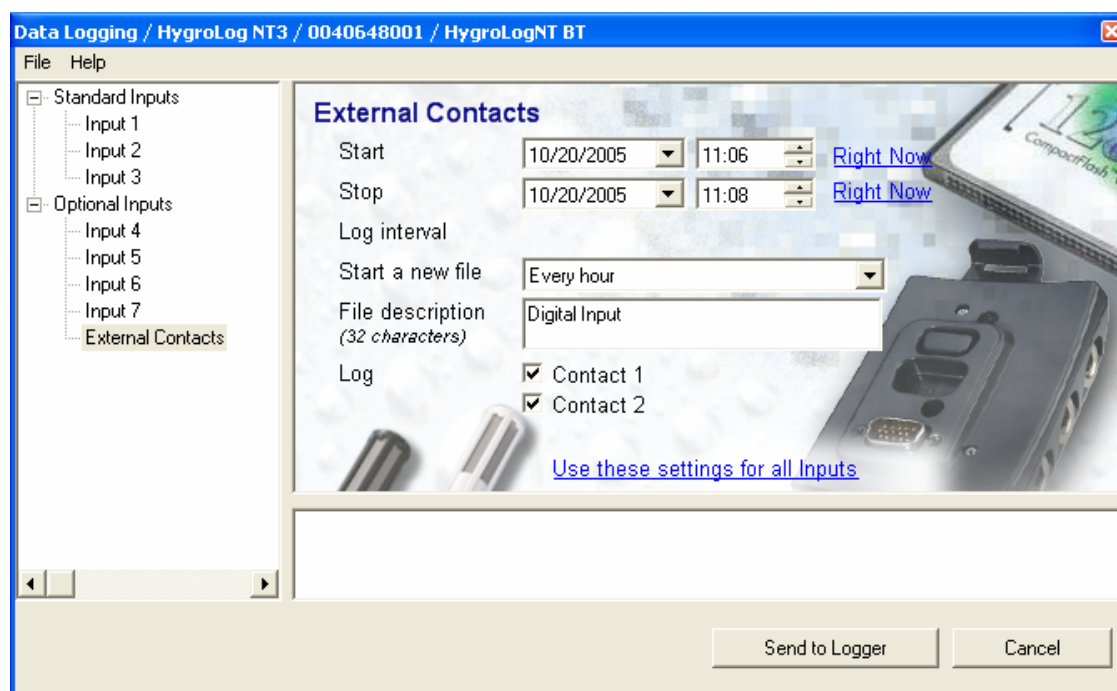

24 AUTOMATIC LOG FILE DOWNLOAD - HygroLog NT

Note: this function is available only with HW4 Professional

To access the automatic file download function, select the HygroLog NT in the device tree (left pane of the HW4 main screen). Select the Log to PC tab in the rights pane of the HW4 main screen.

The automatic log file function can be enabled and configured using the controls located at the bottom of the Log to PC tab:


Click with the mouse on the Start button. HW4 opens the following form:

Click with the mouse on the arrow located to the right of the box labeled Download Interval and select one of the available options from the drop down menu.

Click with the mouse on the arrow located to the right of the box labeled Downloaded Files and select one of the available options from the drop down menu.

To change the path where the files are to be downloaded, click with the mouse on the blue underlined link to the right of Path.

HW4 confirms that the function is enabled and displays the settings that were selected:

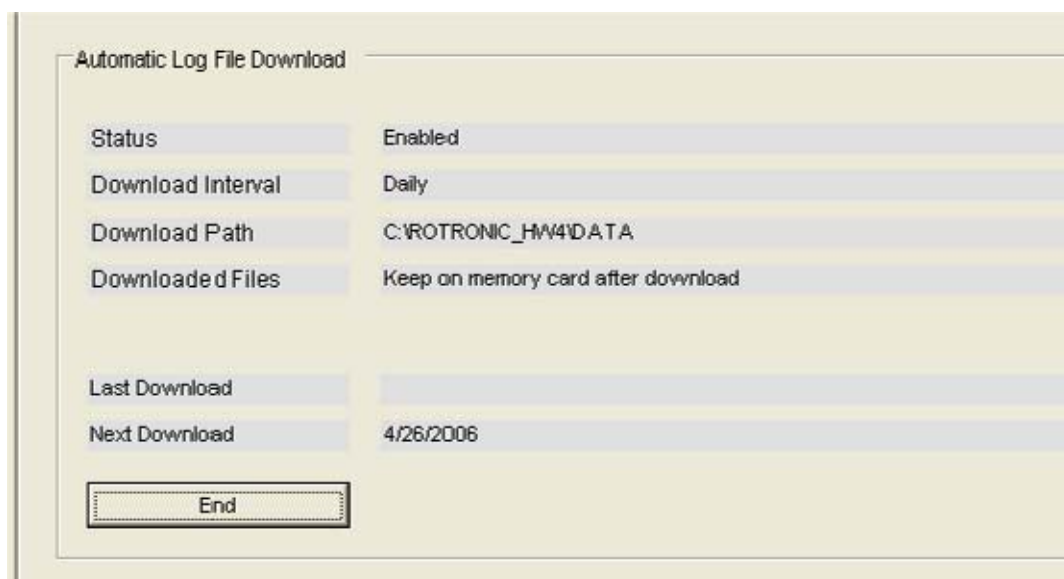


The screenshot shows a window titled "Automatic Log File Download". It contains a table of settings and a button at the bottom.

Automatic Log File Download	
Status	Enabled
Download Interval	Daily
Download Path	C:\ROTRONIC_HW4\DATA
Downloaded Files	Keep on memory card after download
Last Download	
Next Download	4/26/2006

End

After each file download, HW4 refreshes the Last Download and Next Download fields:

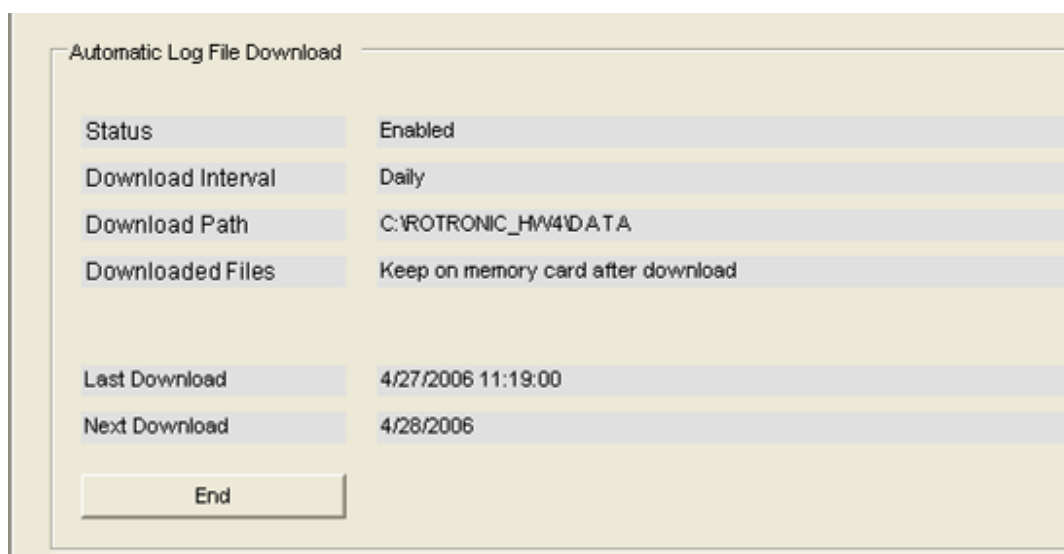

The screenshot shows the same window as before, but with updated values for the Last Download and Next Download fields.

Automatic Log File Download	
Status	Enabled
Download Interval	Daily
Download Path	C:\ROTRONIC_HW4\DATA
Downloaded Files	Keep on memory card after download
Last Download	4/27/2006 11:19:00
Next Download	4/28/2006

End

To disable the function, click with the mouse on the End button. The screen reverts to its default:

Automatic Log File Download

Status	Disabled
Download Interval	
Download Path	
Downloaded Files	
Last Download	
Next Download	

Start

25 DATA LOGGING - directly on the PC


The Log to PC tab is used to record data directly on the PC from any type of instrument that is selected in the Device Tree (left pane of the HW4 main screen). The Log to PC tab lists all input types recognized by HW4, regardless of whether they exist or not on the actual instrument. When the input list is either empty or incomplete, please go to the **View Tab** and check the selections. Only those items selected in the View Tab are visible in the Log to PC Tab.

Data from more than one input can be recorded, using a different file for each individual input.


Current Values	On-line Buffer	Log to PC	View	OPC Server		
Input	Description	Log Interval	Log File Management	Log File Name	File Location	
☐ 1	Probe 1	00:00:15	No file size limit			
☐ 2	Probe 2	00:00:15	No file size limit			
☐ 3	Probe 3	00:00:15	No file size limit			
☐ 4	Probe 4	00:00:15	No file size limit			
☐ 5	Probe 5	00:00:15	No file size limit			
☐ 6	Probe 6	00:00:15	No file size limit			
☐ 7	Probe 7	00:00:15	No file size limit			
☐ External Contacts	Contacts 1/2	00:00:15	No file size limit			
☐ Custom Calc.	Custom Calc.	00:00:15	No file size limit			

25.1.1.1. Start logging data to the PC

1. Select an instrument in the Device Tree. Click on the field "Log Interval" in the line corresponding to the input to be recorded. The following dialog box opens.

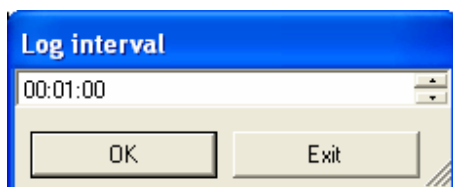


Click on the hours, minutes or seconds to highlight. In the following example, the seconds are highlighted. At this time, click on either the up or the down arrow to the right of the box to change the value.

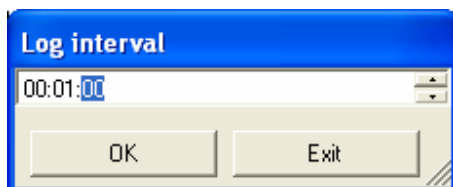

Click on OK when done.

2. Click on the field "Log File Management" in the line corresponding to the probe. The following dialog box opens.

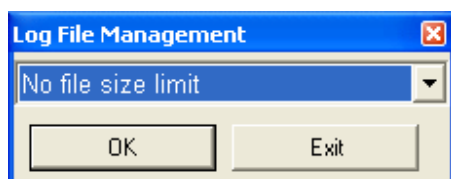


Click on the arrow to the right of the box to display a pull down menu of the available log modes:

- . No file size limit: log file of unlimited size (must be closed manually) ¹
- . Every hour: start a new file at the top of each hour.
- . Every day: start a new at the beginning of each day.
- . Every week: start a new file at the beginning of each week.
- . Every month: start a new file at the beginning of each month.

¹ From a practical point of view, avoid using the unlimited size mode. With any of the other selections, HW4 automatically starts a new file after closing the current file. The name of the new file is generated by HW4 by appending a 3-digit number to the original file name.

3. To start logging, click on the box located to the left of the line corresponding to the probe. The following box opens:

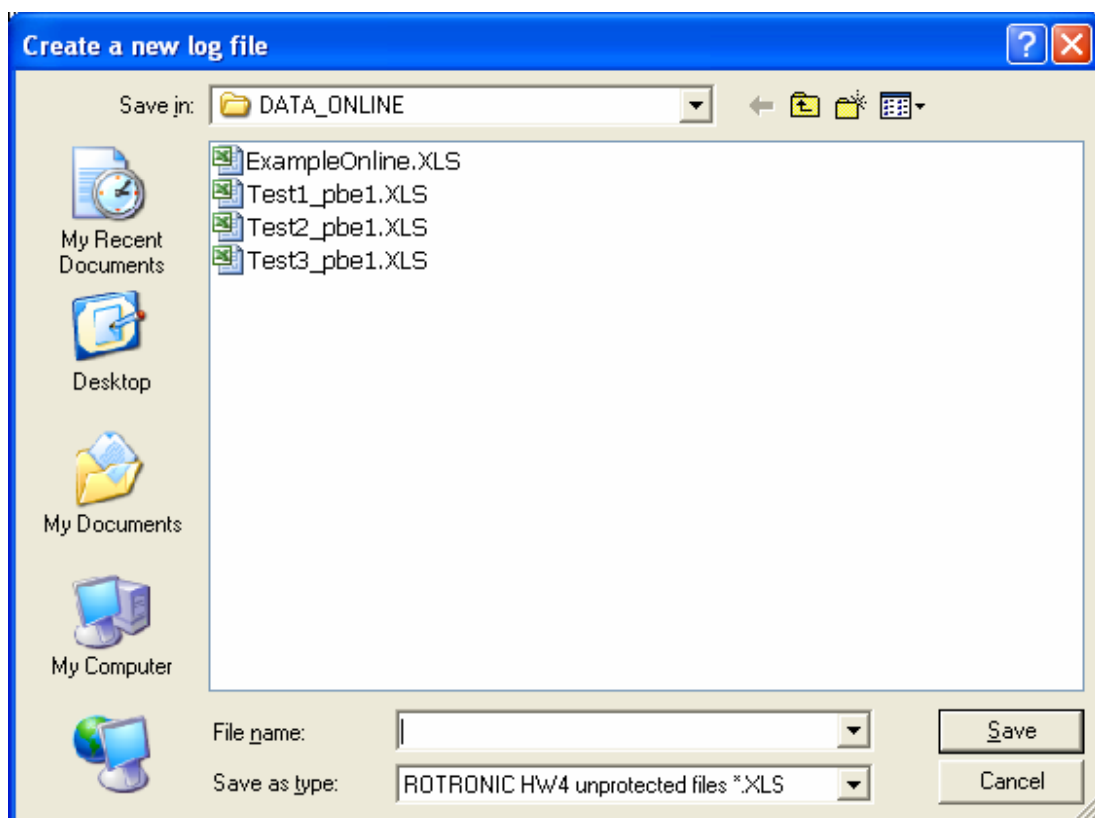

Enter the file name and select the file type: XLS for an unprotected data file or LOG for a protected (binary encoded) data file that meets FDA / GAMP requirements.

Note: if you attempt to use the name of an open or active file (logging in progress), you will get a warning from HW4.

Unprotected files (XLS) are not verified when they are opened with HW4. These files can also be opened and altered with an editor or with Microsoft Excel. Protected files (LOG) can be opened with HW4 but not with an editor or with Microsoft Excel. When opening a protected file, HW4 verifies that the file has not been tampered with. Note that protected binary files use less disk space than unprotected files.

When logging is active, a check mark appears in the box corresponding to the probe or probes.


Current Values On-line Buffer Log to PC View OPC Server					
Input	Description	Log Interval	Log File Management	Log File Name	File Location
☒ 1	Internal	00:00:15	No file size limit	Test_Pbe3.XLS	C:\Document
☐ 2	External 2	00:00:15	No file size limit		
☐ 3	External 3	00:00:15	No file size limit		

A diskette symbol also appears to the left of the probe in the Current Values tab.

Current Values On-line Buffer Log to PC View OPC Server					
!	#	Description	Temperature	Humidity	
	1	Internal	24.28°C -	56.02%RH -	
	2	External 2	--,--	--,--	

25.1.1.2. Stop logging data to the PC

In the Log to PC tab, click on the line corresponding to the probe. The following dialog box opens:

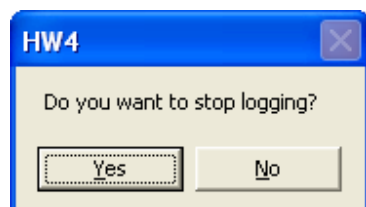

Click on Yes.


Note: if you attempt to exit HW4 while Log to PC is active, you will get a warning from HW4.

26 ACCESSING LOG FILES, PROTOCOLS AND EVENT FILES


26.1 HW4 main menu bar

All files generated by HW4 can be accessed from the HW4 main menu bar provided that they reside on the PC.

1. Select **File** and **Open**:

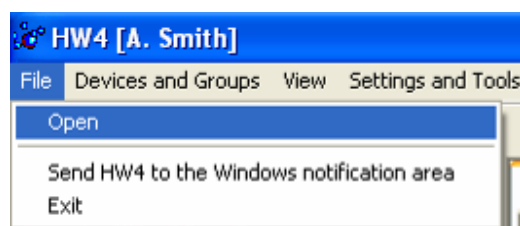

2. HW4 opens the following form

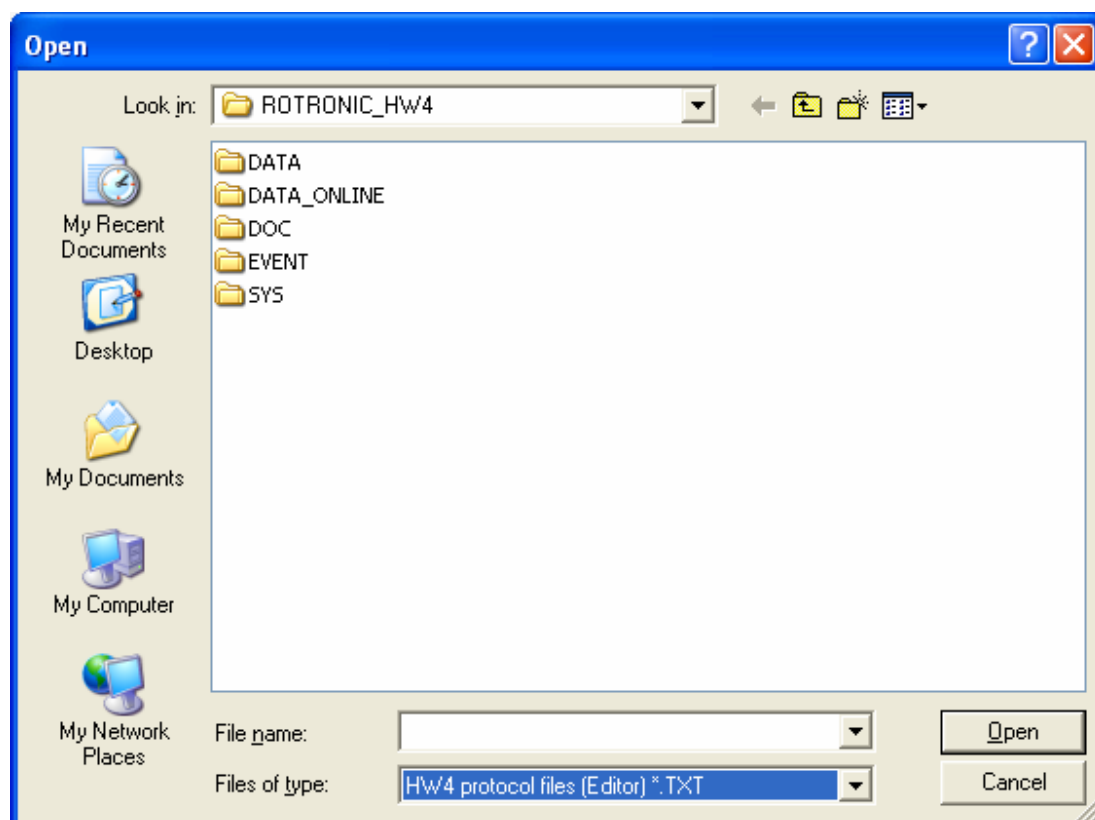



By default, the menu item File - Open (HW4 main menu bar) points to the directory C:\Documents and Settings\Windows User\Application Data\ROTRONIC_HW4. When the location of the files has been changed from the default location, you will have to manually look in the proper folder.

3. Select the folder where the type of file that you want to open is located:

- DATA: log files that have been copied or transferred (moved) from the data logger to the PC.
- DATA_ONLINE: log files directly created on the PC (Log to PC).
- DOC: protocol files (device configuration, logger programming and probe adjustment).
- EVENT: user event files as well as HW4 event files.
- SYS: frequently used device configuration files for future use with Device Manager

Example: Protocol files folder

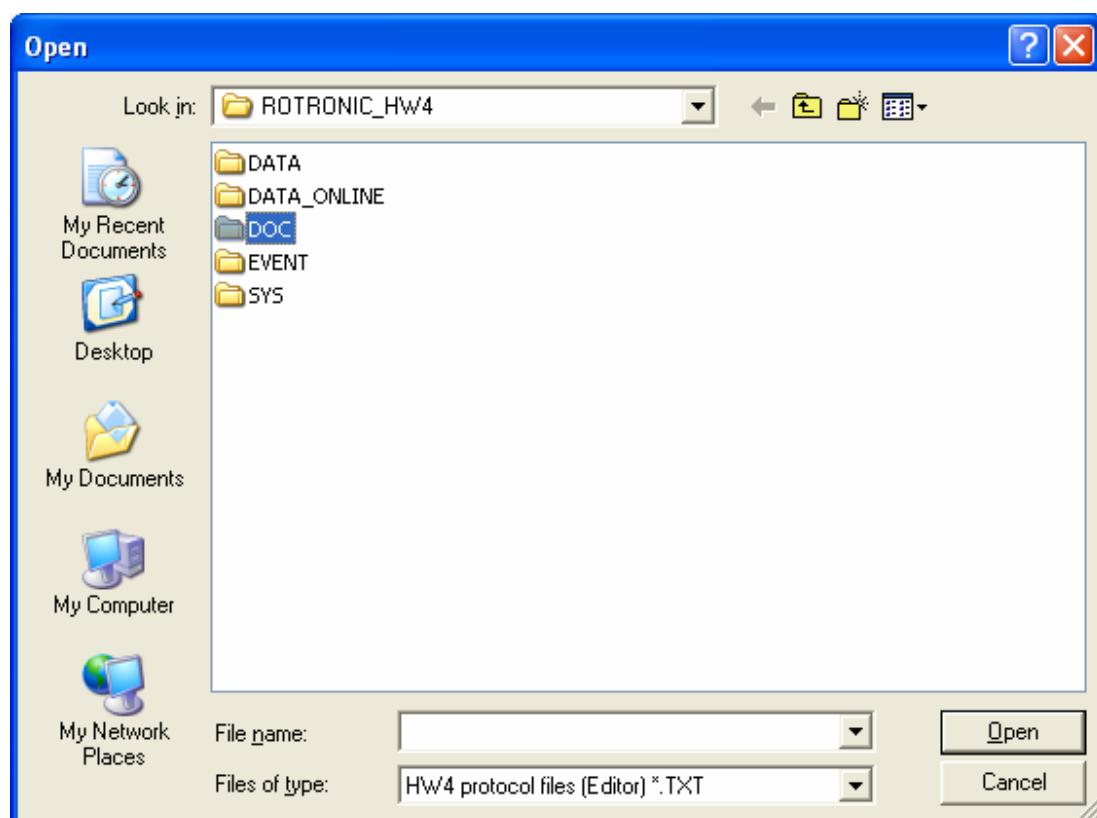

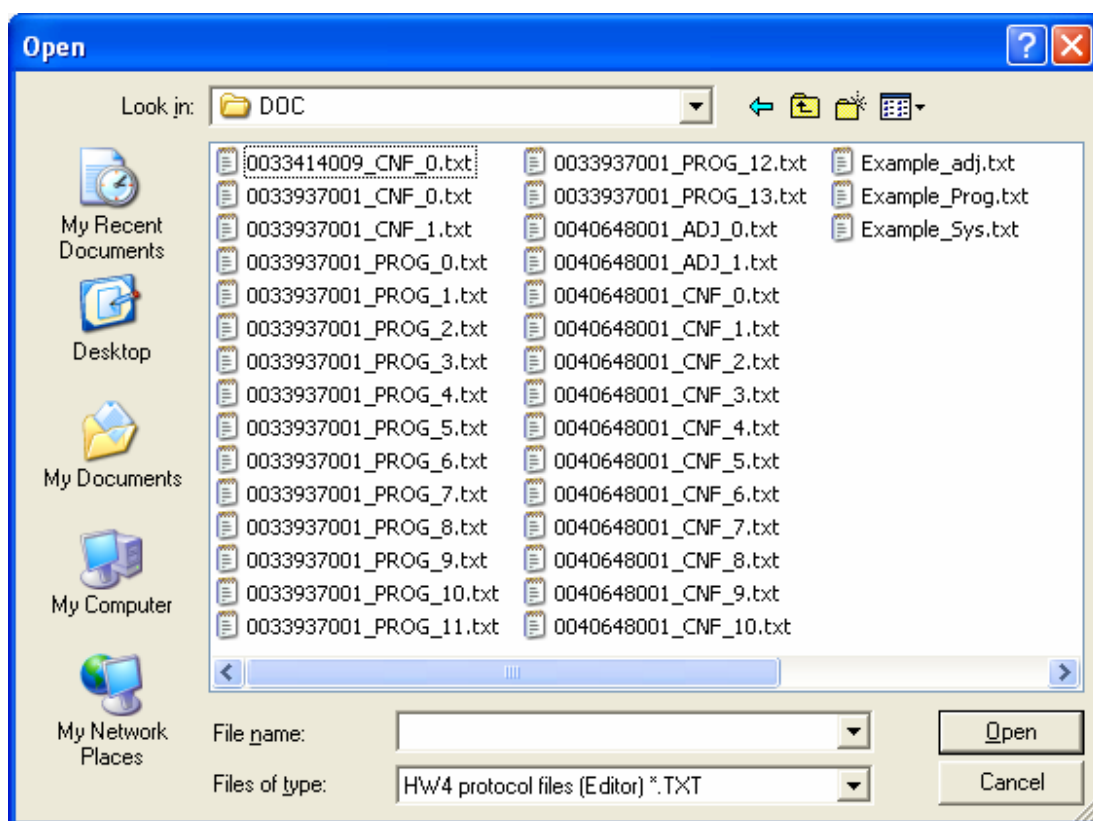

4. Select the proper file extension for the type of file you want to open, or select All files *.*.
If you do not see any file, check the folder and file extension selections.


5. Select a file with the mouse and click on the **“Open”** button.

Depending on the type of file, HW4 uses different file viewers.

Log files : see [View / Sign a Log File](#)

Protocols : see [View / Sign a Protocol](#)


26.2 Access Data (HygroLog NT)

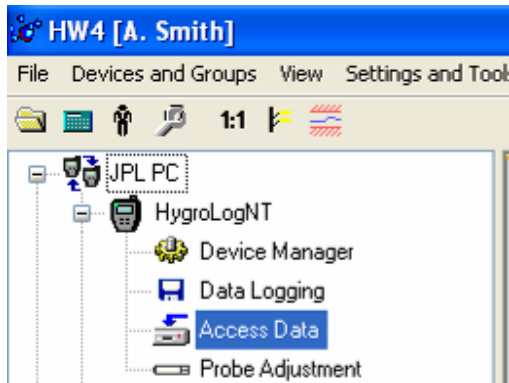


When HW4 has discovered a HygroLog NT, the device appears in the left pane of the HW4 main screen. Expanding the device displays a list of the available device functions.

To select the Access Data function, click on it with the left mouse button. HW4 opens the HW4 Explorer form.

26.2.1. HW4 Explorer

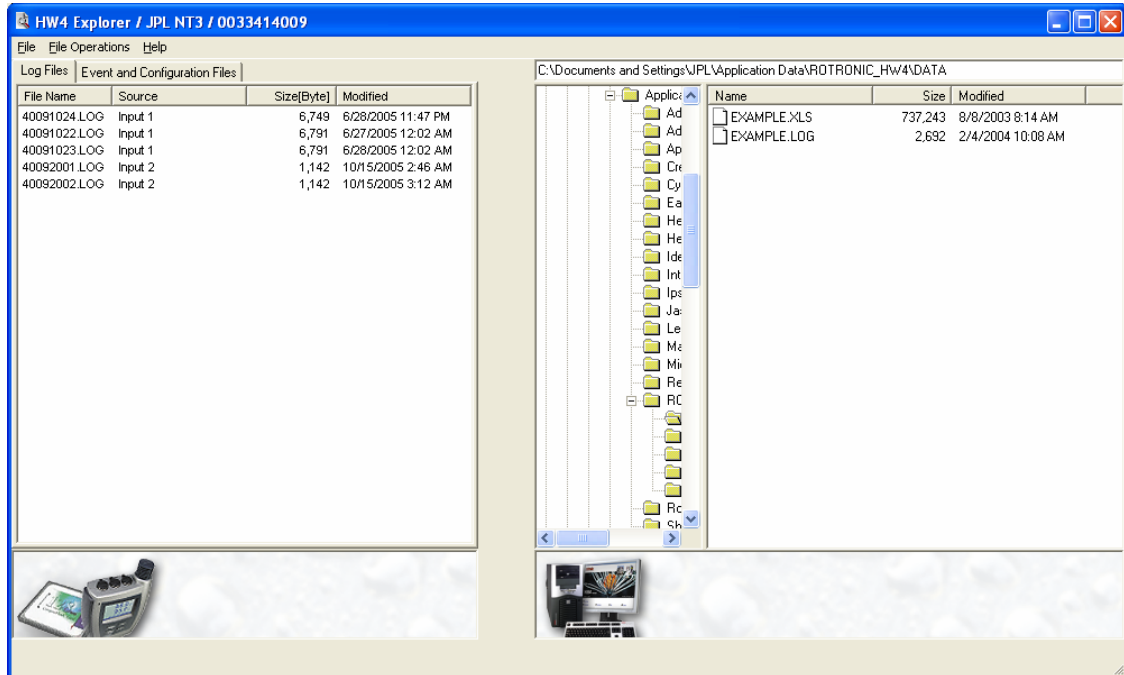


• **Left Pane:** the left pane of the HW4 Explorer form provides a list of the files currently present in the memory card of the HygroLog NT. Select the Log file tab to view log files. Select the Events and system data tab to view the logger event file (extension .EVT) and the logger configuration file (extension .DAT).

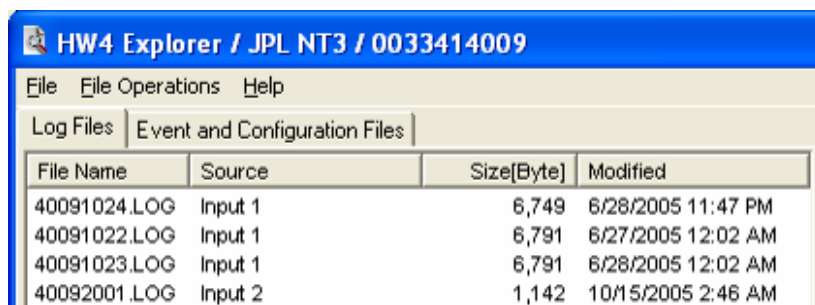Log files


The screenshot shows the 'Log Files' tab in the HW4 Explorer application. The window title is 'HW4 Explorer / JPL NT3 / 0033414009'. The menu bar includes 'File', 'File Operations', and 'Help'. The 'Log Files' tab is selected, and it displays a table with the following data:

File Name	Source	Size[Byte]	Modified
40091024.LOG	Input 1	6,749	6/28/2005 11:47 PM
40091022.LOG	Input 1	6,791	6/27/2005 12:02 AM
40091023.LOG	Input 1	6,791	6/28/2005 12:02 AM
40092001.LOG	Input 2	1,142	10/15/2005 2:46 AM

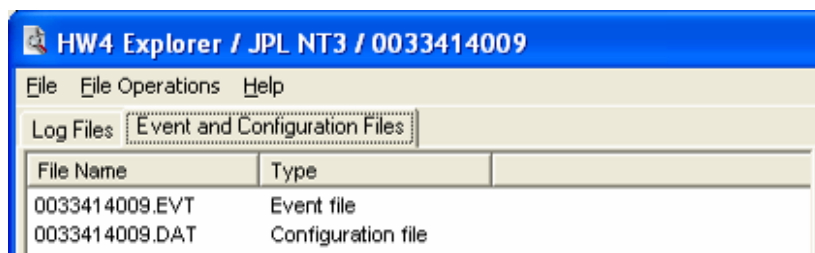
Logger event and logger configuration files


The screenshot shows the 'Event and Configuration Files' tab in the HW4 Explorer application. The window title is 'HW4 Explorer / JPL NT3 / 0033414009'. The menu bar includes 'File', 'File Operations', and 'Help'. The 'Event and Configuration Files' tab is selected, and it displays a table with the following data:

File Name	Type
0033414009.EVT	Event file
0033414009.DAT	Configuration file

A file with the extension .EVT and one with the extension .DAT are always present in the memory card. These files cannot be deleted with HW4.

• **Right pane:** the right pane of the HW4 Explorer form provides a list of the log files, protocols, event files and configuration files currently present on the PC. The location of these files depends on the directory that was specified for each type of file in HW4 Global Settings / File Locations.

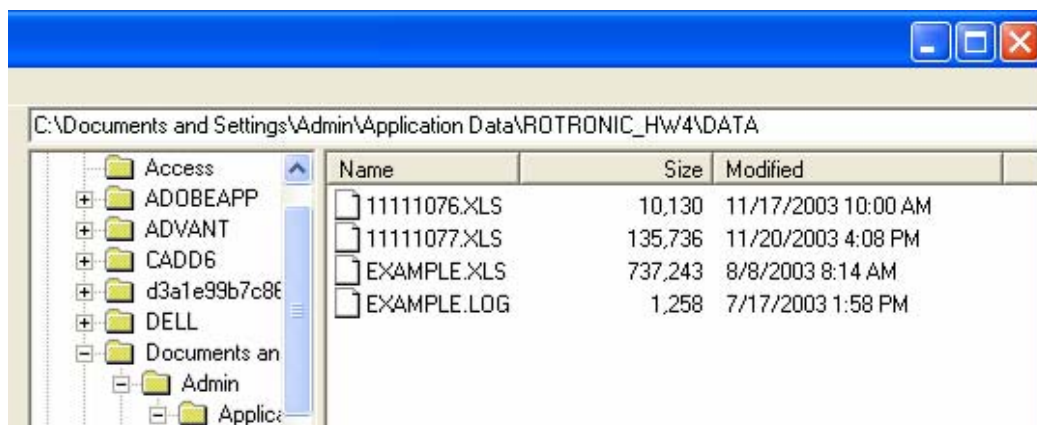

The screenshot shows a Windows Explorer window with the address bar set to 'C:\Documents and Settings\Admin\Application Data\ROTRONIC_HW4\DATA'. The left pane shows a tree view of the directory structure, including folders like 'Access', 'ADOBEAPP', 'ADVANT', 'CADD6', 'd3a1e99b7c86', 'DELL', 'Documents and Settings', 'Admin', and 'Applic'. The right pane displays a table of files with the following data:

Name	Size	Modified
11111076.XLS	10,130	11/17/2003 10:00 AM
11111077.XLS	135,736	11/20/2003 4:08 PM
EXAMPLE.XLS	737,243	8/8/2003 8:14 AM
EXAMPLE.LOG	1,258	7/17/2003 1:58 PM

Any file present in the memory card can be opened, copied or moved to the computer, or deleted from the memory card. The logger event and logger configuration files are also present in the internal EEPROM of the HygroLog NT. If not present, these two files will be automatically written to the memory card by the logger. Files that are on the computer can be either opened or deleted. These files cannot be copied or moved to the memory card.

26.2.2. Menu Bar

The menu bar is located at the top of the HW4 Explorer form.

File

- **Open:** opens any file that is highlighted either in the left pane of the Explorer (memory card) or in the right pane (computer).

Note: files located on the logger memory card are automatically copied to the PC before being opened.

- **Exit:** exits the HW4 Explorer

File Operations

Note: highlight the file first.

- **Move:** copies the file from the memory card to the computer and deletes the file from the memory card (except for event and configuration files).
- **Copy:** copies the file from the memory card to the computer and leaves the file in the memory card.
- **Delete:** deletes the file from the memory card or from the computer (except for event and configuration files located on the memory card). Several files may be selected at the same time for deletion.

Note: these functions are also directly available by right clicking with the mouse on the file.

Help


The Help menu consists of:

- **HW4 Help:** Opens HW4 Help
- **About HW4:** Displays the version number and ID number of HW4

26.2.3. Files located on the PC

In the right pane of HW4 Explorer, highlight the appropriate folder in C:\Documents and Settings\your Windows login name\ROTRONIC_HW4:

- DATA: this folder holds the log files that have been copied or transferred (moved) from the data logger to the PC.
- DATA_ONLINE: this folder holds the log files directly created on the PC (Log to PC).
- DOC: this folder holds protocol files (logger configuration, logger programming and probe adjustment).
- EVENT: this folder holds user event files as well as HW4 event files (used for troubleshooting software problems).
- SYS: this folder can be used to hold frequently used instrument configuration files for future use from within Device Manager

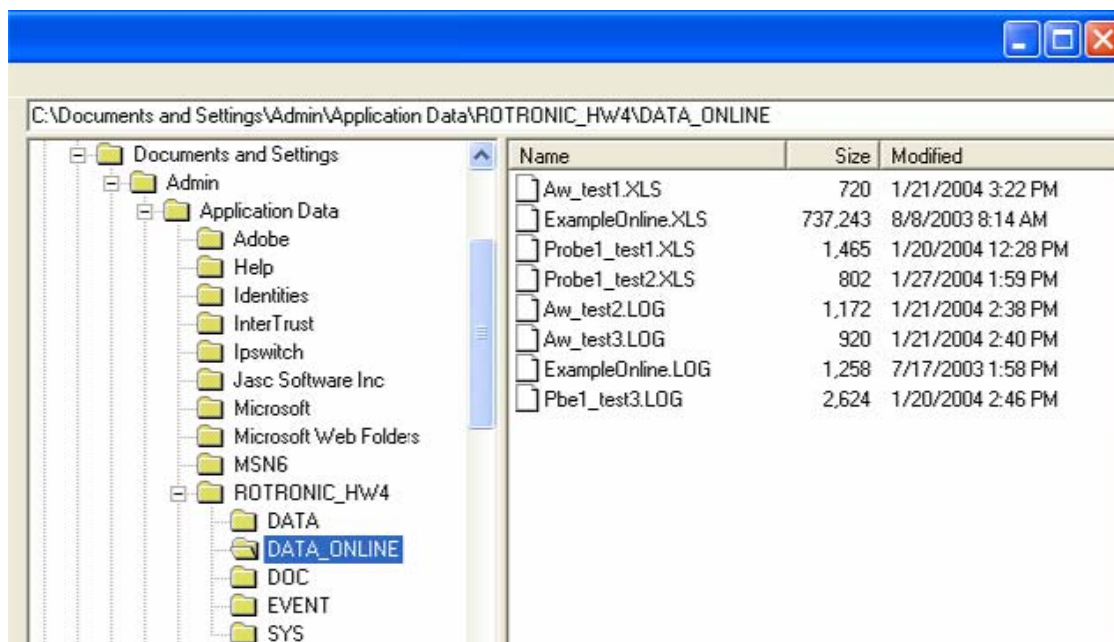


Highlight one of the files present in the folder. To open or delete the file, click on File Operations in the HW4 explorer menu bar or right click on the file to open a small menu next to the file name.

26.2.4. Files located on the logger

• Event and Configuration files

In the left pane of HW4 Explorer, select the Events and Configuration Files tab and right click on any file. This opens a small menu box next to the file name. Select with the mouse the desired file operation.

Example of a logger event file (file extension .EVT)

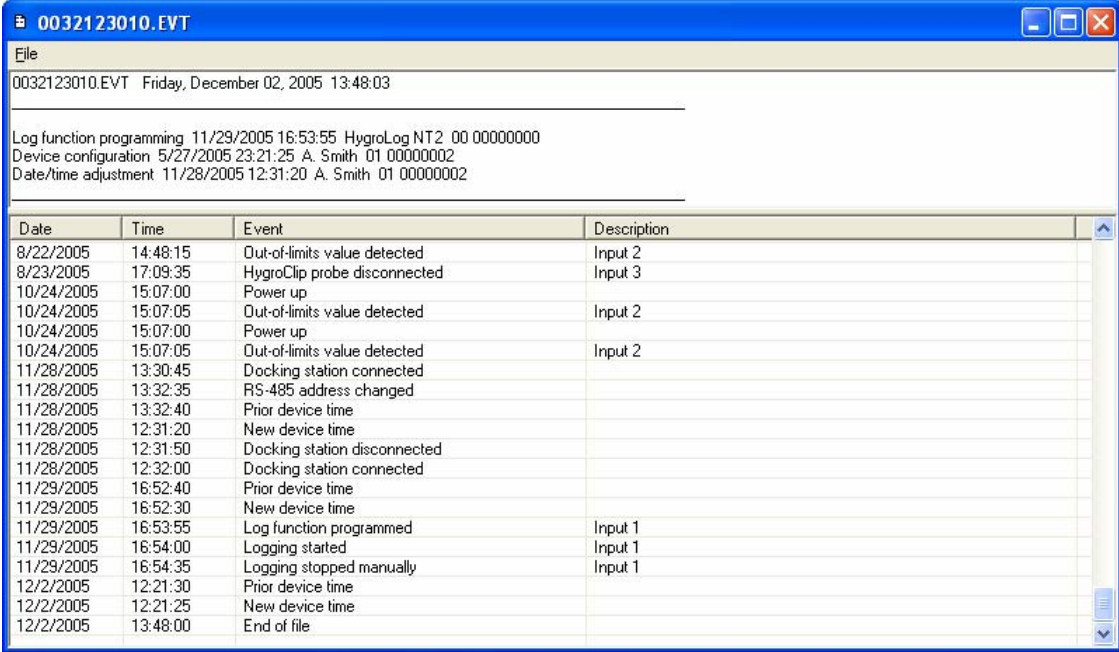

0032123010.EVT Friday, December 02, 2005 13:48:03

Log function programming 11/29/2005 16:53:55 HygroLog NT2 00 00000000
 Device configuration 5/27/2005 23:21:25 A. Smith 01 00000002
 Date/time adjustment 11/28/2005 12:31:20 A. Smith 01 00000002

Date	Time	Event	Description
8/22/2005	14:48:15	Out-of-limits value detected	Input 2
8/23/2005	17:09:35	HygroClip probe disconnected	Input 3
10/24/2005	15:07:00	Power up	
10/24/2005	15:07:05	Out-of-limits value detected	Input 2
10/24/2005	15:07:00	Power up	
10/24/2005	15:07:05	Out-of-limits value detected	Input 2
11/28/2005	13:30:45	Docking station connected	
11/28/2005	13:32:35	RS-485 address changed	
11/28/2005	13:32:40	Prior device time	
11/28/2005	12:31:20	New device time	
11/28/2005	12:31:50	Docking station disconnected	
11/28/2005	12:32:00	Docking station connected	
11/29/2005	16:52:40	Prior device time	
11/29/2005	16:52:30	New device time	
11/29/2005	16:53:55	Log function programmed	Input 1
11/29/2005	16:54:00	Logging started	Input 1
11/29/2005	16:54:35	Logging stopped manually	Input 1
12/2/2005	12:21:30	Prior device time	
12/2/2005	12:21:25	New device time	
12/2/2005	13:48:00	End of file	

The HygroLog NT keeps internally track of the last 170 operations, events and configuration changes in its internal EEPROM memory as well as in the .EVT file located on the memory card.

Example of a logger configuration file (file extension .DAT)

The contents of a logger configuration files are displayed using the Device Manager (for details, see Device Manager).

Device Manager / HygroLogNT / 111111111 / Rotronic HygroLogNT

File Settings Help

User Information

Name * John Doe [Company logo file](#)

Company * RIC [Display company logo](#)

Address 160 E. main Street

City/ZIP 11743 Huntington, NY

Country USA

Telephone (631) 427-3898

E-Mail doe@aol.com

☒ Include in log file

OK Cancel

Left pane tree view:

- User Information
- Device Information
- Language/Unit System
- Memory Card
- Probes
 - Probe1
 - Probe2
 - Probe3
- Display
- LEDs
- Keypad
- Sounds
- Power Source
- Interface
- Options
 - Probe4
 - Probe5
 - Probe6

• Log files

In the left pane of HW4 Explorer, select the Log Files tab and right click on the file. This opens a small menu box next to the file. Click with the mouse on Open. HW4 first downloads the file and makes a copy on the PC.


Depending on the type of file, HW4 uses different file viewers.

Log files : see [View / Sign a Log File](#)
 Protocols : see [View / Sign a Protocol](#)

27 VIEW / SIGN A PROTOCOL

HW4 uses the following viewer to display protocols on the monitor screen.

Example:

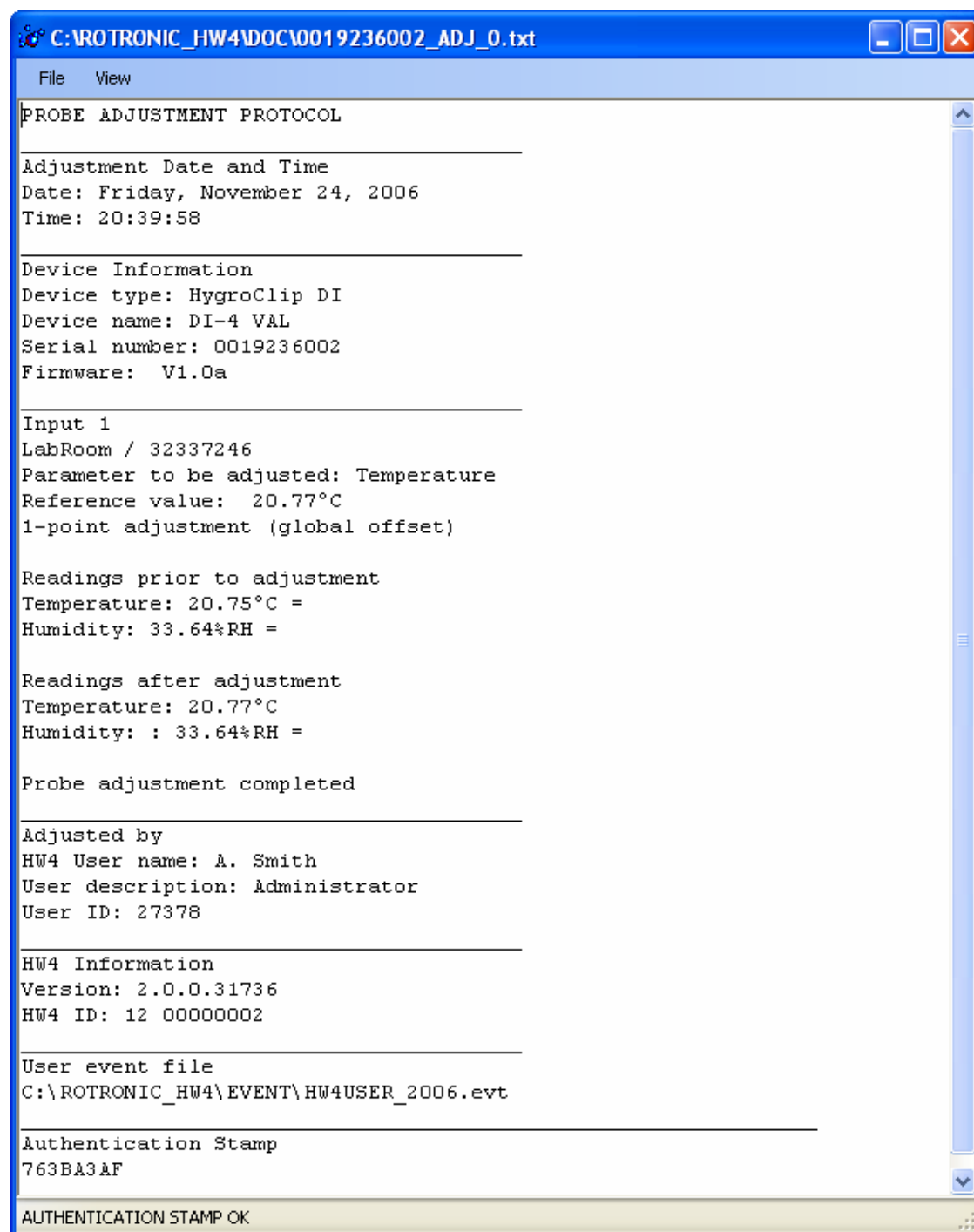

HW4 always checks the authentication stamp of a protocol when opening the file. The result of this verification appears on the bottom left corner of the form used to display the file contents:


AUTHENTICATION STAMP OK

WARNING: THIS IS NOT AN ORIGINAL HW4 FILE

27.1 Protocol viewer menu bar

File

- **Print:** opens the Windows printer form and prints the protocol on the specified printer.
- **Sign document:** use this menu item to add your HW4 user name and a comment text to the protocol. For authentication purposes you will need to enter your password.

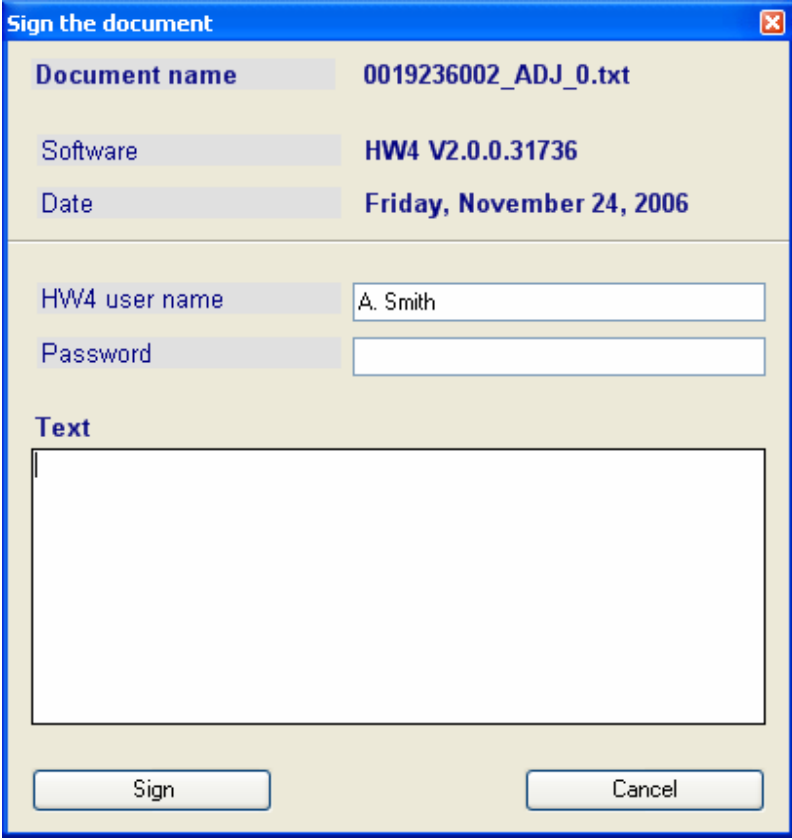


The image shows a Windows-style dialog box titled "Sign the document". It contains several fields and a text area. The fields are: "Document name" with value "0019236002_ADJ_0.txt", "Software" with value "HW4 V2.0.0.31736", and "Date" with value "Friday, November 24, 2006". Below these are "HW4 user name" with value "A. Smith" and "Password" (empty). At the bottom is a large "Text" area. There are "Sign" and "Cancel" buttons at the bottom right.

Document name	0019236002_ADJ_0.txt
Software	HW4 V2.0.0.31736
Date	Friday, November 24, 2006
HW4 user name	A. Smith
Password	
Text	
Sign	Cancel

- **Exit:** closes the viewer.

View

- **Document:** displays the entire document in the viewer.
- **Signature:** displays only the signature area of the document.

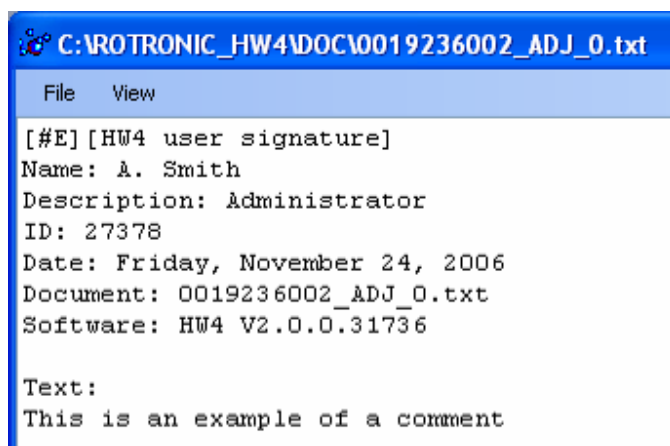

28 VIEW / SIGN A LOG FILE

View Data is the graphic module of HW4.

28.1 *Opening a log file with View Data*


A log file can be opened in View Data either from the HW4 Explorer or from the File menu in the HW4 main menu bar.

28.1.1. From HW4 Explorer


In HW4 Explorer, select either a log file located on the HygroLog NT (left pane - Log Files tab) or a log file located on the PC (right pane – DATA or DATA_ONLINE folder).

In the HW4 menu bar, click on the “File - Copy and Open” or “File - Copy” command. As an alternative right click with the mouse on the file and select “Open” from the small file menu.

When the file is located on the HygroLog NT, HW4 begins by making a copy of the file to the PC, prior to opening the file in View Data:

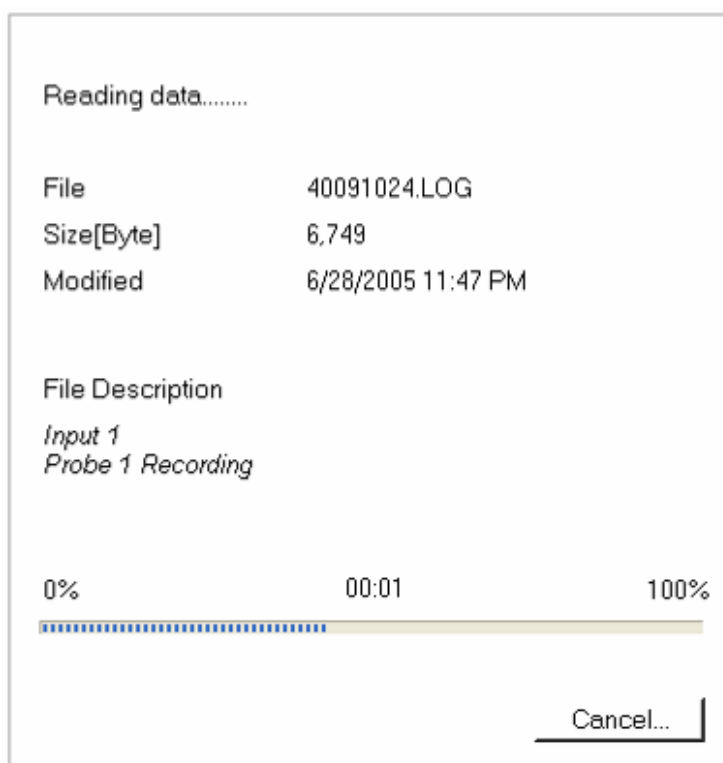

After copying the file to the PC, HW4 opens the file in View Data. When the file is already located on the PC, HW4 immediately opens the file in View Data:

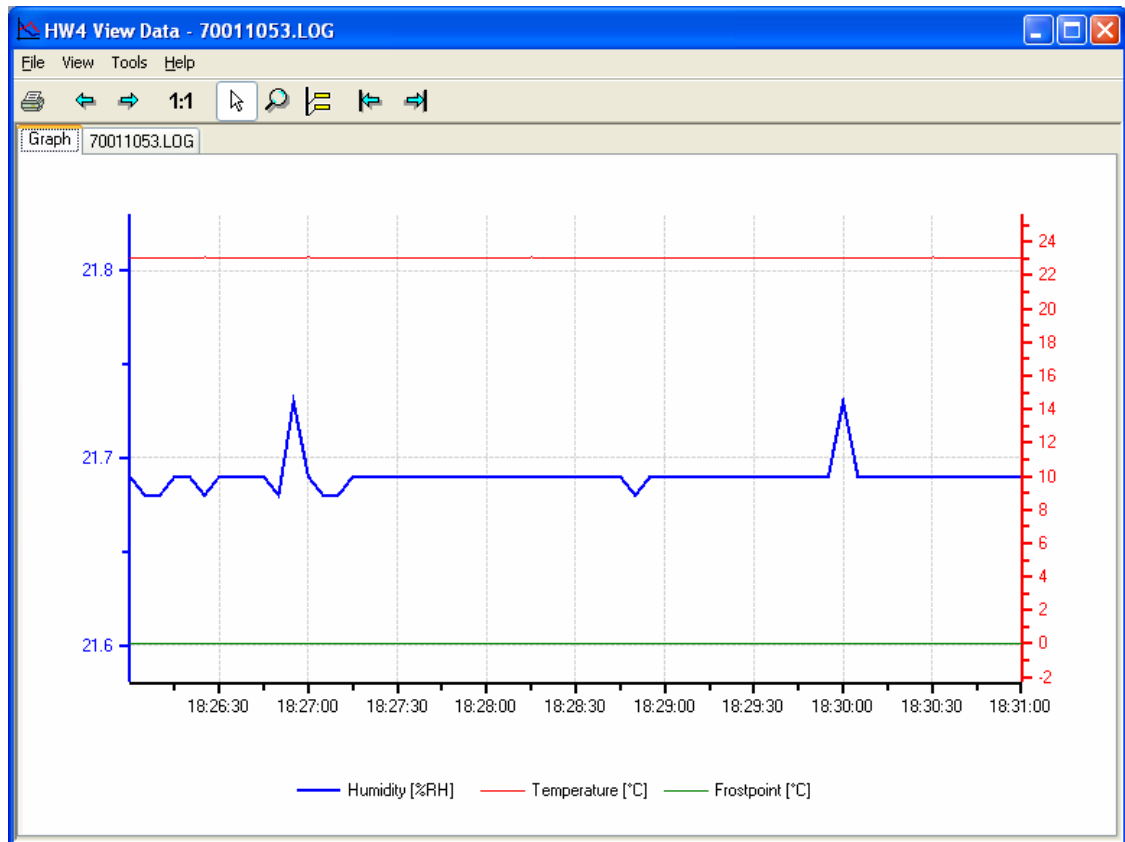


28.1.2. From the HW4 Main menu bar

Log files can be opened from View Data from the HW4 main menu bar provided that they reside on the PC.

Select **File** and **Open**:




HW4 opens the following form

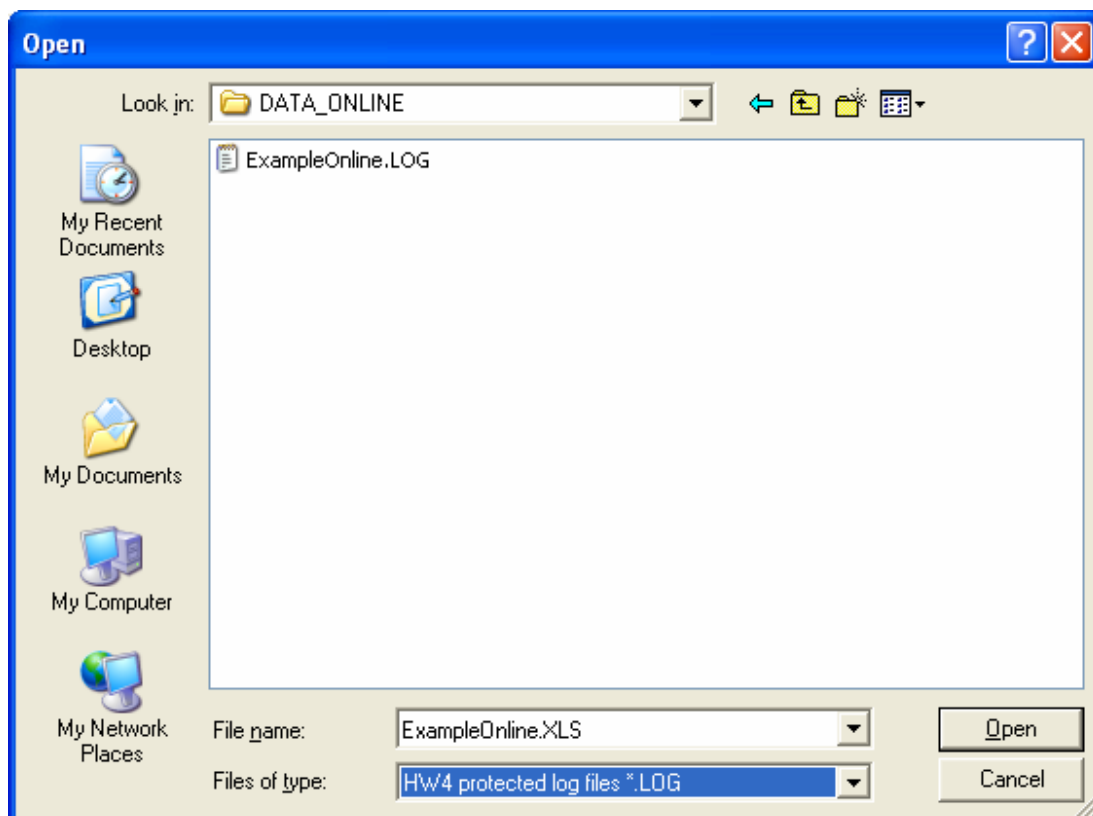

Select the folder where the file is located (DATA or DATA_ONLINE) as well as the type of log file (file with .LOG or .XLS extension). If you do not see the file, check the folder and file extension selections. Click on the **Open** button.

HW4 immediately opens the file in View Data:

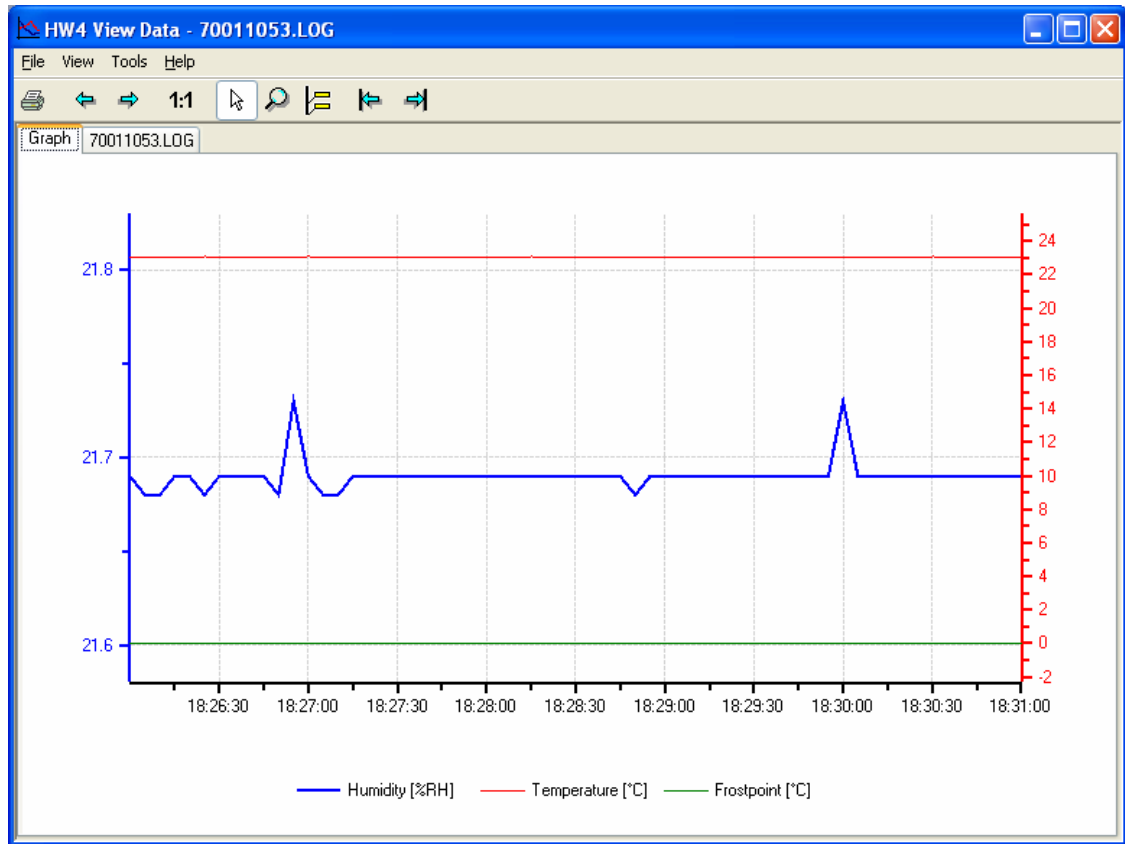

Click with the mouse on the tab with the file name to view the data in a table.

HW4 View Data - 70011053.LOG

File View Tools Help

Graph 70011053.LOG

[H][DEVICE NAME]
HyGen NT3

[P][PROCESS DESCRIPTION]
HygroGen / HygroClip 35/80/10/35

[S][DEVICE INFORMATION]
Release = V1.2b
Serial Number = 0034277001
Filename = 70011053.LOG
Interval = 00:00:05
Probe = V2.0/33220033
Name = Port 1

Num...	Date	Time	Humidity %RH	Temperature °C	Frostpoint °C	
4	11/24/2006	18:26:15	21.69	23.05	0.01	
5	11/24/2006	18:26:20	21.69	23.02	-0.02	
6	11/24/2006	18:26:25	21.68	23.07	0.02	
7	11/24/2006	18:26:30	21.69	23.03	-0.01	
8	11/24/2006	18:26:35	21.69	23.05	0.01	
9	11/24/2006	18:26:40	21.69	23.02	-0.02	
10	11/24/2006	18:26:45	21.69	23.05	0.01	
11	11/24/2006	18:26:50	21.68	23.03	-0.01	
12	11/24/2006	18:26:55	21.73	22.99	-0.02	
13	11/24/2006	18:27:00	21.69	23.07	0.03	
14	11/24/2006	18:27:05	21.68	23.03	-0.01	
15	11/24/2006	18:27:10	21.68	23.03	-0.01	
16	11/24/2006	18:27:15	21.69	23.06	0.02	
17	11/24/2006	18:27:20	21.69	23.05	0.01	
18	11/24/2006	18:27:25	21.69	23.02	-0.02	
19	11/24/2006	18:27:30	21.69	23.05	0.01	
20	11/24/2006	18:27:35	21.69	23.02	-0.02	

Click with the mouse on the Graph tab to view the data again in a graph.


28.2 View Data menu bar

28.2.1. File

- **Overlay another file:** this menu item is available only with HW4 Professional and allows bringing several log files into a single graph (the maximum number of files depends on the size of the files). Note that the files to be overlaid must be present on the PC (if necessary, begin by copying the desired files from the logger to the PC). The files being overlaid can be a mixture of the LOG and XLS format.

Note: Files recorded with the HygroLog NT can be started all at the exact same time. By contrast, this is not possible when recording several inputs on the PC as illustrated in the picture below.

The **Save As** command offers the possibility of saving the overlaid data under a new file name. The new file is always in the XLS format and it separately retains the file header and the data from each original file. The new file is not Write protected. After overlaying two or more files, it sometimes makes sense to view the data with a common time origin. For instructions, see **Options – Overlay synchronization**.

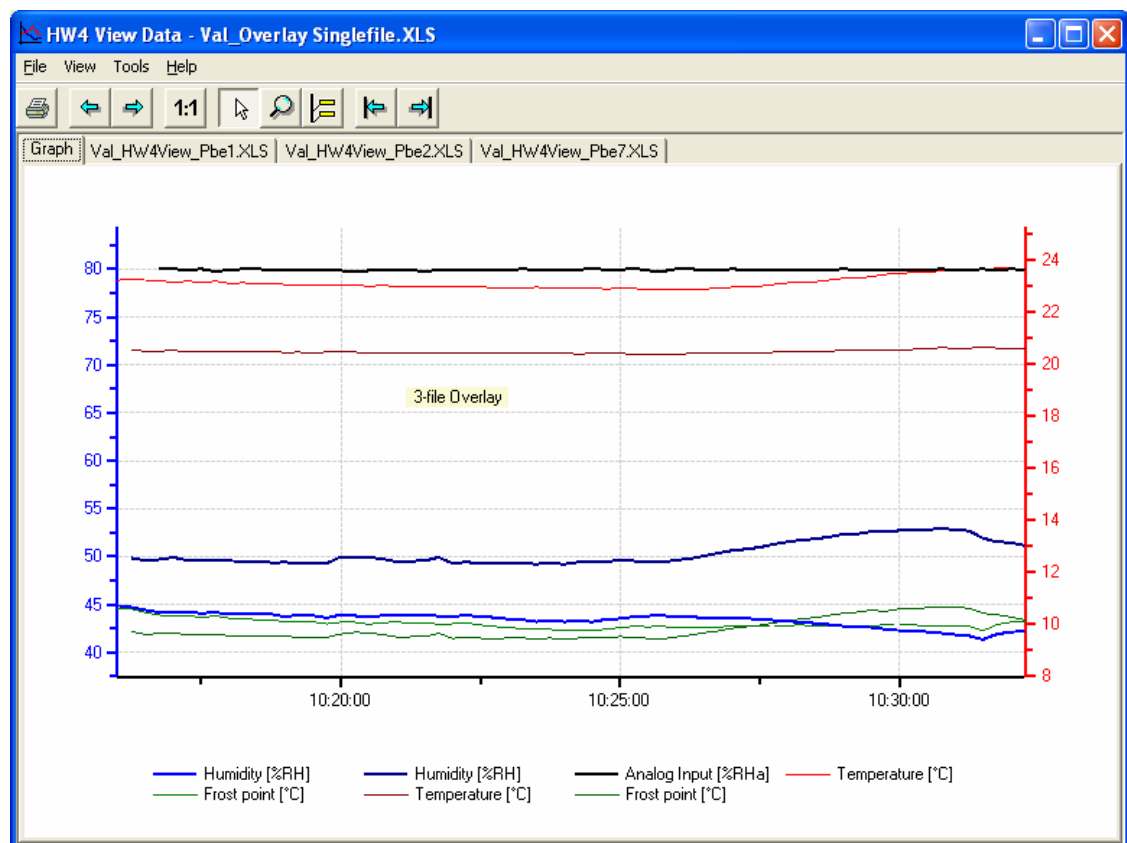


- **Save:** Saves a file located on the PC to its original location under the original file name. This command will not work when the file is Write protected and therefore cannot be used with the original files created by HW4.

- **Save As:** saves a copy of the file to any location on the computer. The file is always saved in the .XLS format. This command can be used to do the following:
 - Convert a binary log file (LOG) to the text format (XLS) and save the converted file to disk, under a new name (the original LOG file is still available)
 - Save a write protected file under a new name after adding notes, text boxes, etc.
 - Save the results of overlaying several files under a new file name. The new file separately retains the file header and data of each original file.
- **Print:** the effect of this menu item depends which tab (Graph or File Name/Data Table) is selected in HW4 View Data. The command opens the Windows printer form and prints the graph or table on the specified printer.

Graph: the log file name is added to the top of the graph when the graph is being printed.


Data table: the file header is always printed before the table. It is possible to select with the mouse a number of lines from the data table, prior to printing. Depending on the selection made within the Windows printer form, either all of the data or just the selected lines will be printed.
- **Sign document:** use this menu item to add your HW4 user name and a comment text to the log file. For authentication purposes you will need to enter your password.

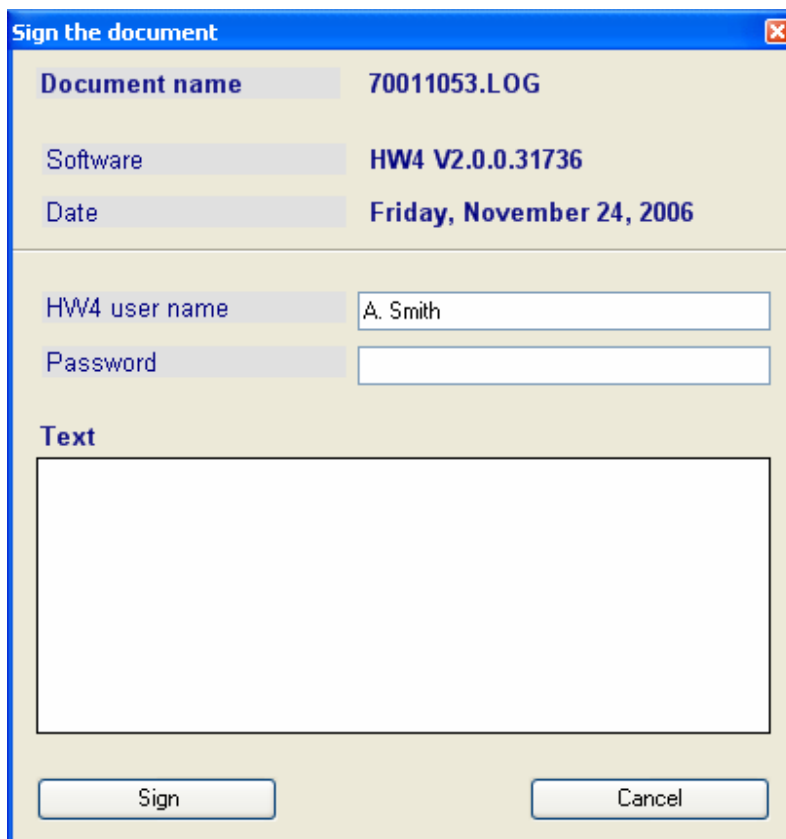

The image shows a Windows-style dialog box titled "Sign the document" with a close button (X) in the top right corner. The dialog has a light beige background and a blue border. It contains several fields and buttons:

- Document name:** 70011053.LOG
- Software:** HW4 V2.0.0.31736
- Date:** Friday, November 24, 2006
- HW4 user name:** A. Smith (text input field)
- Password:** (empty text input field)
- Text:** A large empty text area for a comment.
- Buttons:** "Sign" and "Cancel" at the bottom.

Signing a log files creates an additional tab in View Data:

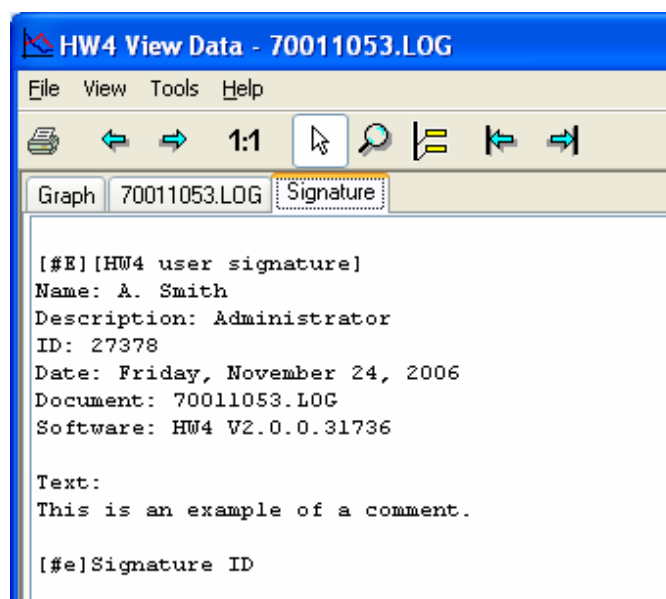

- **Exit:** exits HW4 View Data

28.2.2. View

- **Initial View**¹: returns the graph to the view used by HW4 View Data when opening the file
- **Previous View**¹: after using the magnifying glass or zooming on a graph area several times, returns the graph to the previous view²
- **Next View**¹: changes the graph to the next view when the magnifying glass or zoom has been used several times²
- **Graph time window**¹: used to select the time window corresponding to the graph X axis. The available selections appear in a submenu: Hour, Day, Week, Month or Year. For example, selecting Day, makes the X axis correspond to a 24-hour time window.
- **Next time window**¹: used to increment the origin of a time window (X axis) by one time window unit. For example, if the graph is displaying a week worth of data, this command changes the graph to the following week.
- **Previous time window**¹: used to decrement the origin of a time window (X axis) by one time window unit. For example, if the graph is displaying a week worth of data, this command changes the graph to the previous week.
- **Toolbar**: shows or hides the toolbar (the toolbar shortcut buttons are used only in graph view and duplicate a number of items of the View menu).
- **Notes**: shows or hides the Notes tab. This tab allows extensive notes to be associated with the log file. To retain these notes, a write protected file should be saved under a new name

¹ this item is grayed out when the Graph tab is not selected

² this item is context sensitive and is grayed out when not applicable

28.2.3. Tools


● **Merge the overlays**¹: when several log files have been opened using the HW4 View Data File menu and Overlay another file, this menu item merges the files into a single file (single data table and a single graph). Merging the overlays should be used only when there is no overlapping time between the log files (otherwise the merge results will not make sense). The new file created with this command can be saved in .XLS format, under a new name, using **File – Save As**.

¹ This menu item is available only with HW4 Professional


Example:

Two log files are overlaid into a single graph, prior to being merged into a single file.

Note: the recording time of the two files does not overlap.

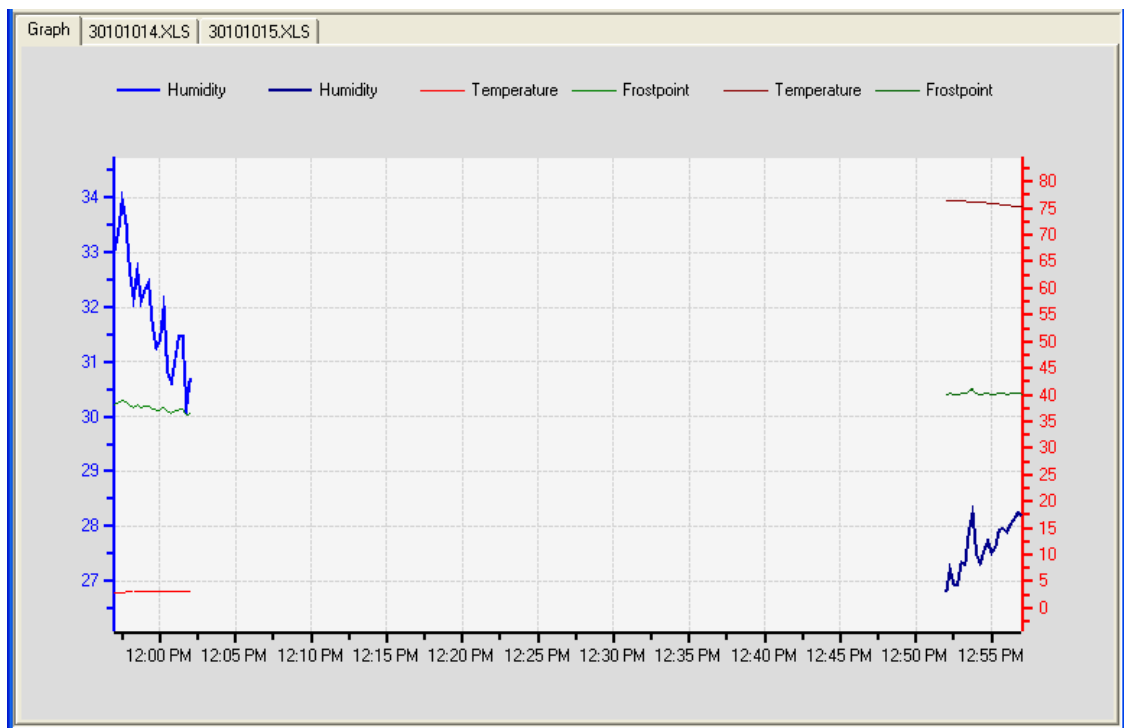

The results of merging data from the same two log files into a single file:

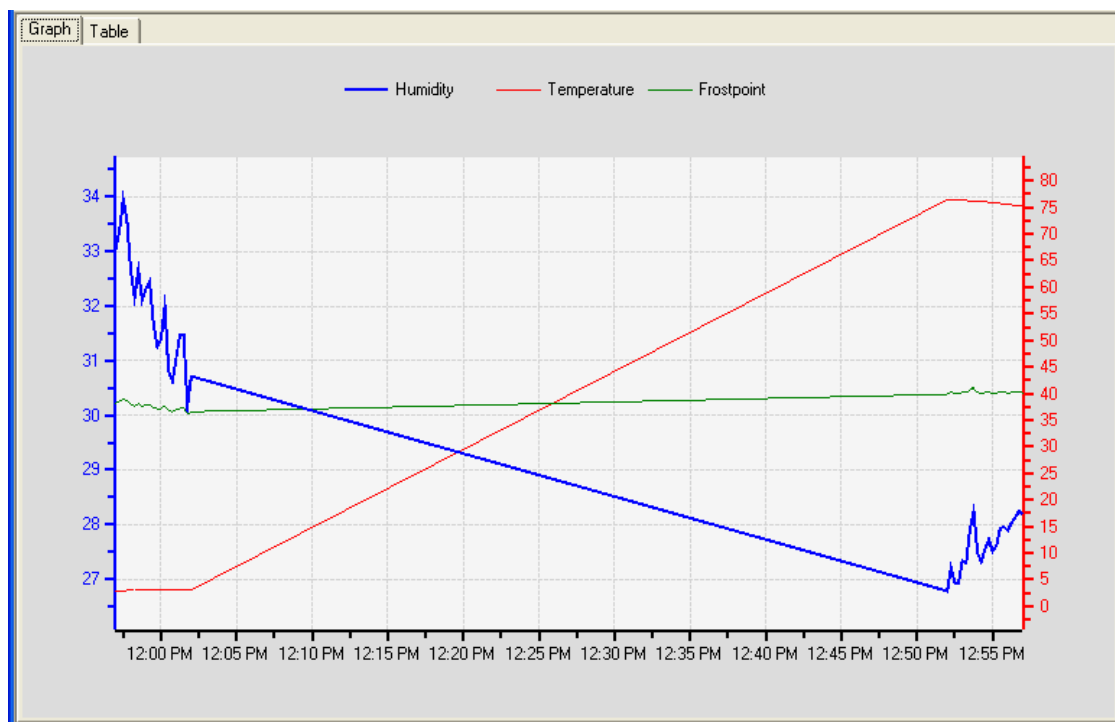

The two tabs with the original file names have been replaced with a single tab labeled Table. Data points that were initially separated by a time interval are now joined.

Note: The header of the merged file is the header of the first file that was opened. The headers of the other files are lost. For example, if the other files have alarm settings that are not the same as the first file, this information will be lost.

• **Overlay synchronization**¹: when the data from different files was not all recorded starting at the exact same time, it is sometimes useful to use a common time origin for viewing view all traces after overlaying the files in a single graph. The effect of this menu item is to shift some of the traces data along the graph time (X) axis. The actual time corresponding to a number of data points is lost as far as the graph is concerned, but this can be reversed.

¹ This menu item is available only with HW4 Professional

This command opens the following sub-menu:

- Earliest time
- Hour
- Day
- Week
- Month
- Year
- Undo


Earliest time: causes HW4 to use the earliest time found in the log files as a common time origin for all data in the graph.

Hour, Day, Week, Month or Year: causes HW4 to use the top of an hour (such as 11:00) or the beginning of a day, week, month or year as the common time origin for all data. Using any of these selections makes sense only when the log data spans a sufficient length of time. For example, do not select Day if the log data does not cover more than a day.


Undo: reverses the current synchronization and must be used before applying a different synchronization method.

Example:

1) Initial result after overlaying two files started at a different time

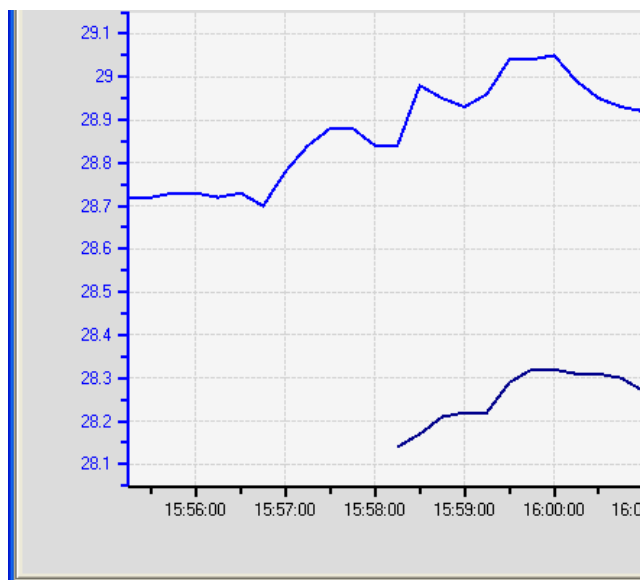


2) After synchronizing the files based on “earliest time”, both traces start at the exact same time.

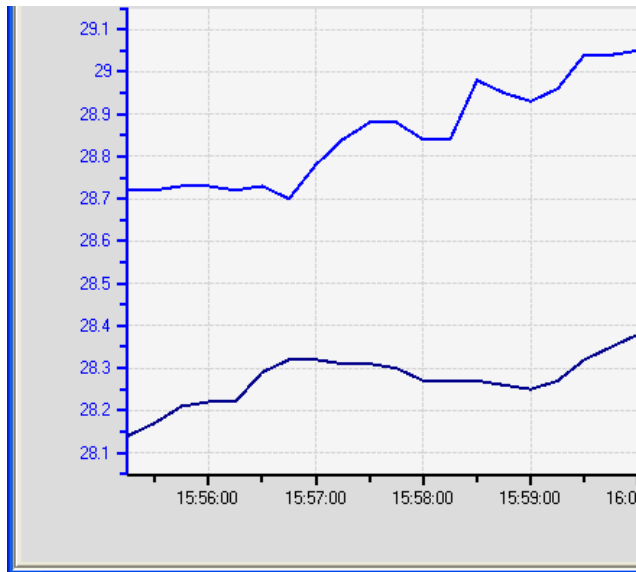

• **Statistical data:** this menu item opens a new tab in HW4 View Data that displays the number of data points, minimum, maximum, average and standard deviation for each parameter that was recorded. A sub-menu offers the choice of displaying statistical information only for the data currently shown on the graph (for example after zooming) or for an entire log file. When the graph consists of more than one log file (overlays), all of the data is being used regardless of its origin.

Mean kinetic temperature: this additional data element is given only for temperature and uses the value of activation energy entered in HW4 Global Settings – Language / Unit System tab.

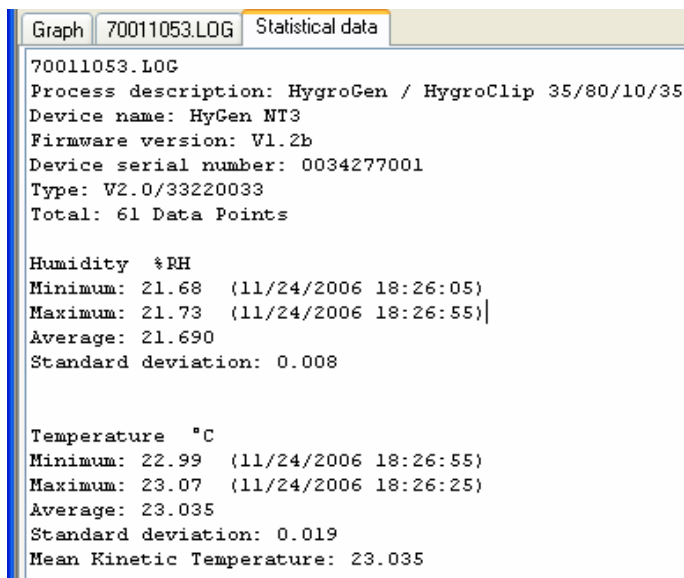

When high and low alarm values have been programmed for an input, each individual occurrence of an alarm condition is reported in a separate block (Alarm 1, Alarm2, etc.). A separate block is also used for each recorded parameter, even when alarm conditions occur simultaneously. Contiguous data points which are out-of-limits values are reported as one block. Statistical information is provided for each block.

```
Alarm 1
Humidity
Number of data points: 9
From: 10/18/2005 14:54:00
To: 10/18/2005 14:56:00
Minimum: 47.75 (10/18/2005 14:56:00)
Maximum: 88.71 (10/18/2005 14:54:30)
Average: 71.410
Standard deviation: 14.922

Alarm 2
Temperature
Number of data points: 4
From: 10/18/2005 14:54:00
To: 10/18/2005 14:54:45
Minimum: 25.98 (10/18/2005 14:54:00)
Maximum: 28.59 (10/18/2005 14:54:30)
Average: 27.483
Standard deviation: 1.093
```


- **Graph settings:** this menu item opens the HW4 Global Settings form with Graph Settings Tab selected. For details, see Settings and Tools, HW4 global settings.

Note: any change made here to the graph settings will affect all graphs, including the on-line graph

28.2.4. Help

- **HW4 Help:** opens HW4 Help.
- **About HW4:** displays the version number and ID number of HW4.


28.3 View Data Toolbar

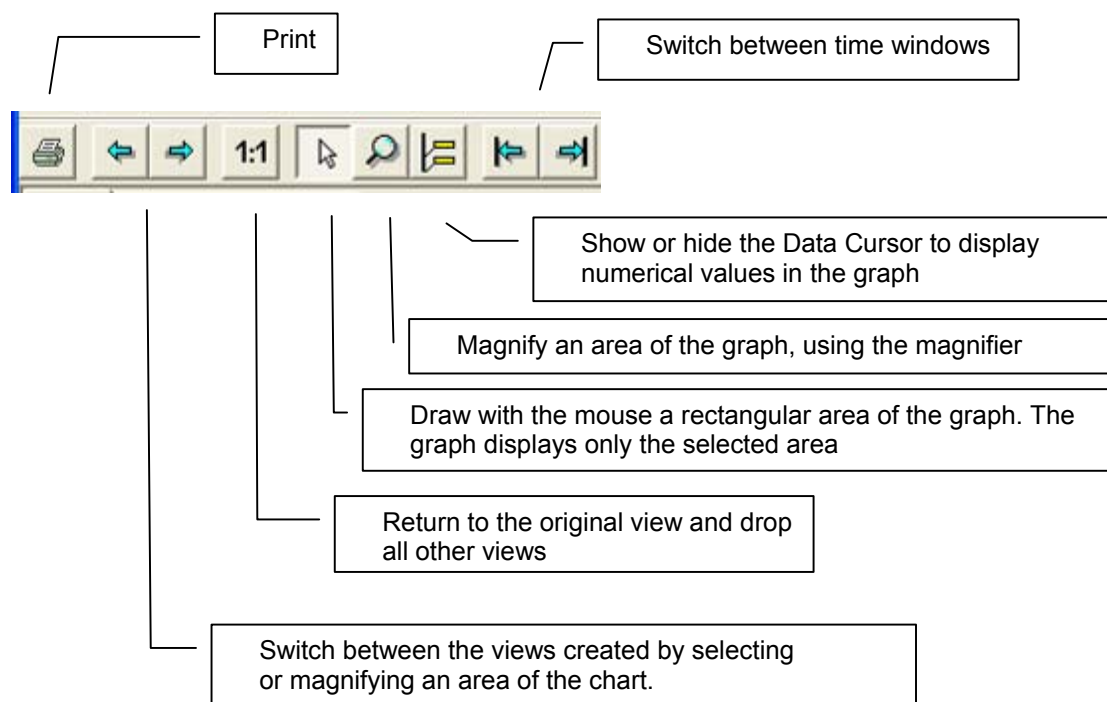

28.4 Graph Tab

28.4.1. Using the Data Cursor


In the toolbar, left click on the Data Cursor button. The Data Cursor appears as a vertical bar on the graph. Labels appear next to the data cursor to display numerical values.

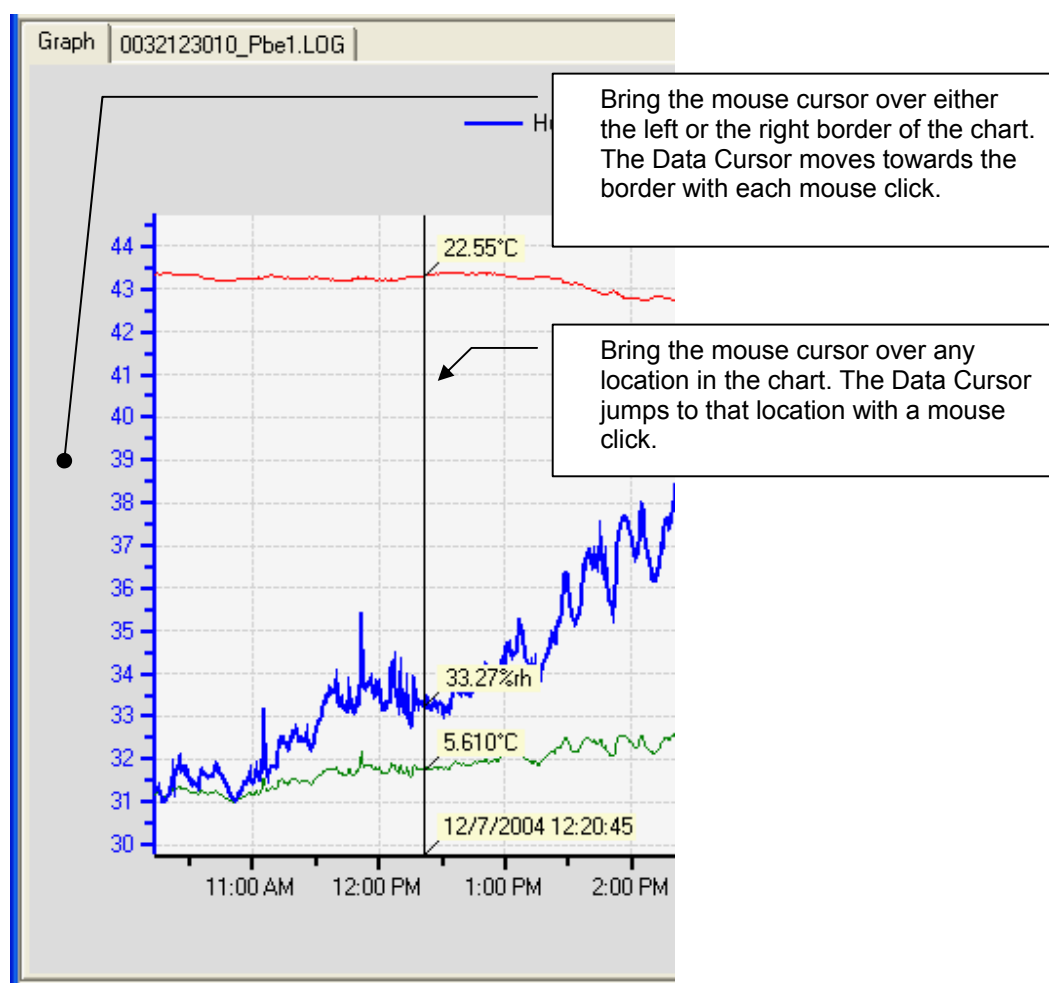


28.4.2. Displaying hidden data labels

When two traces on the chart are close one to the other, it may happen that the data labels attached to the Data Cursor are on top one of the other. Left clicking with the mouse on any visible label makes the label move 90° clockwise and can be used to uncover a hidden label.

a) The temperature data label is hidden behind the humidity data label

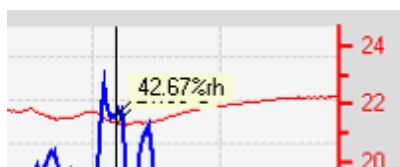

b) Left click with the mouse on the humidity data label. The temperature data label becomes visible.

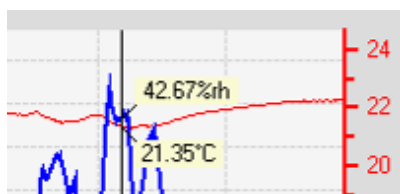

28.4.3. Selecting which traces are displayed by the graph

Right clicking with the mouse over any area of the graph opens the following menu:

- Insert text box (see Adding text notes directly in the graph, further down)
- Show all traces: displays all available traces on the graph

To select a trace, bring the mouse cursor over one of the chart traces. The cursor changes to a hand. Right click with the mouse to select the trace. The color of the trace changes and the following menu opens:

- Hide this trace
- Show this trace only
- Show alarm band: shows on the graph the alarm settings (if any) for the selected trace. The alarm band can be displayed only for one trace at a time. The alarm band appears as two shaded areas of the same color as the trace. To hide, deselect the menu item.
- Attach text/data box (see further down)
- Show all traces


Note: showing the alarm band for one of the traces usually results in a different scale.

28.4.4. Attaching a text / data box to a trace

Bring the mouse cursor close to a trace. Right click with the mouse when the cursor changes to a hand. Doing this selects the trace and also opens a small menu box.


The menu item Attach Text/Data Box creates a text / data box that is attached to a specific location of the selected trace. Upon selecting this menu item, the mouse cursor changes to a cross. By default, HW4 fills the text box with data corresponding to the location of the mouse cursor. Move the mouse cursor (cross) to any location on the trace. The data in the box changes as the box moves along the trace.

To attach the text box, left click with the mouse when the cursor is at the desired location. The mouse cursor jumps inside the text box and the text inside the box is highlighted (edit mode). Left click with the mouse to edit the contents of the box. Use the keyboard arrow keys to move the cursor to a location where text is to be inserted or added. When done, left click with the mouse anywhere in the graph.

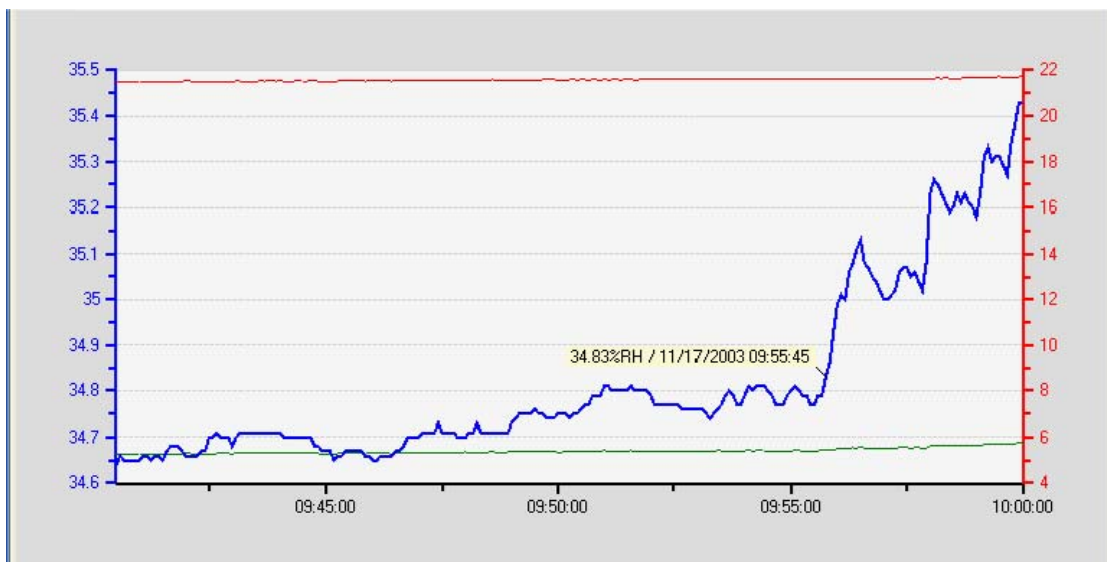

Bringing the mouse cursor over the text box makes the cursor change into a hand. At that time, right clicking with the mouse opens another menu:

- Edit text box: use to add text after the numerical data or to replace the numerical data with text
- Delete text box
- Orientation: use one of the available options to place the text box at a convenient location

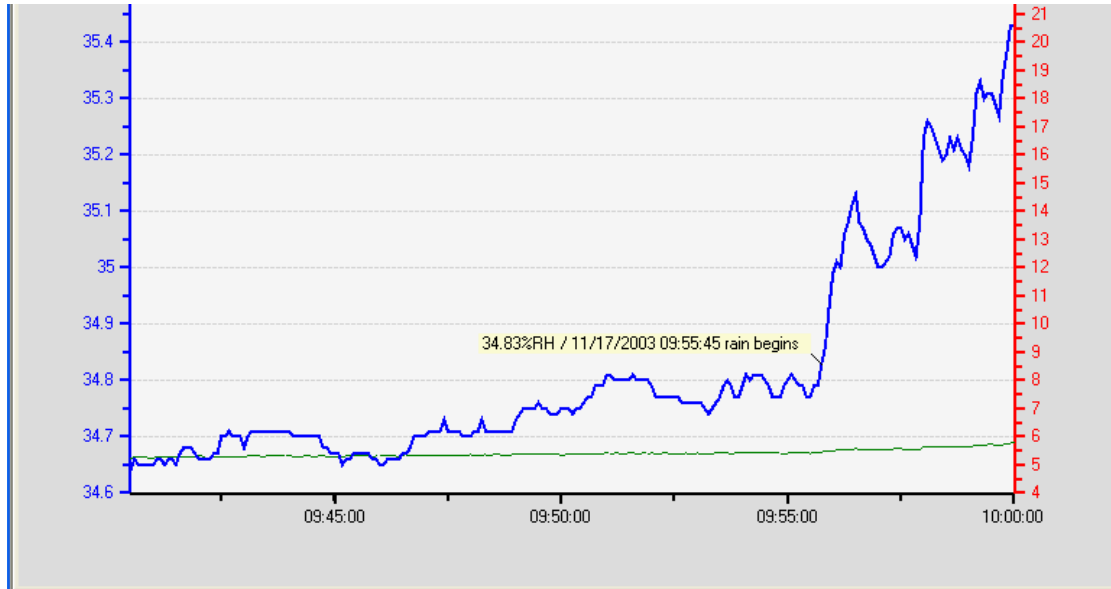

28.4.5. Adding a text note directly in the graph

To create a text note that is not attached to a specific trace, use the following two steps:

a) Create and position the text box: bring the mouse cursor to the location of the graph where you want to add a text note. Right click with the mouse. This opens a small menu. In this menu, select Insert Text Box. The mouse cursor changes to a cross and a textbox is created with the text: "Enter text here". Move the cross to the desired location of the text box and left click with the mouse.

b) Enter the text: right click with the mouse on the text box. This opens a small menu. Select Edit Text Box and type the text in the blue area. When done, left click with the mouse

28.4.6. Deleting a text note or text/data box

Make sure that the Data Cursor is not visible. Bring the mouse cursor over the text/data box. The cursor changes to a hand (the hand appears only when the data cursor is not visible). Right click with the mouse. This opens a small menu. Select Delete Text Box and left click with the mouse

28.5 Data Table / File Name Tab

Clicking on the file name tab, displays both the contents of the file header (top of the form) and the data in the form of a table (bottom of the form):

Graph: 0032123010_Pbe2.LOG

[#U]User
RIC
A. Smith
160 E. Main Street
11743
Huntington
(631) 427-3898
USA

[#I]Device name
RIC Demo2

[#P]Process type

Number	Date	Time	Humidity	Temperature	Frost point
1322	12/7/2004	15:45:45	39.44	20.11	5.90
1323	12/7/2004	15:46:00	39.54	20.07	5.90
1324	12/7/2004	15:46:15	39.67	20.07	5.95
1325	12/7/2004	15:46:30	39.64	20.04	5.91
1326	12/7/2004	15:46:45	39.71	20.01	5.91
1327	12/7/2004	15:47:00	39.78	20.00	5.93
1328	12/7/2004	15:47:15	39.84	19.97	5.92
1329	12/7/2004	15:47:30	39.86	19.97	5.93
1330	12/7/2004	15:47:45	39.85	19.97	5.92
1331	12/7/2004	15:48:00	39.92	19.97	5.95
1332	12/7/2004	15:48:15	39.97	19.91	5.92
1333	12/7/2004	15:48:30	40.04	19.90	5.93
1334	12/7/2004	15:48:45	40.07	19.90	5.94
1335	12/7/2004	15:49:00	40.13	19.91	5.97
1336	12/7/2004	15:49:15	40.07	19.89	5.93
1337	12/7/2004	15:49:30	40.07	19.89	5.93
1338	12/7/2004	15:49:45	40.10	19.89	5.94

Values that correspond to an alarm condition (if any was specified for the probe in Device Manager) are shown over a red background.

28.6 Working with log files and graphs (How To)

[Superimpose several log files within a single graph](#)

[Merge several log files into a single file](#)

[Use a common time origin when displaying log files within a single graph](#)

[Display statistical data](#)

[View numerical data directly within a graph \(Data Cursor and Data Labels\)](#)

[Display a data label that is hidden behind another label](#)

[Attach a text or data box to a trace within a graph](#)

[Add a text note to a graph](#)

[Delete a text note or a text/data box](#)

[Select which trace to display in a graph](#)

[Display in a graph the alarm band corresponding to a trace](#)

[Make a copy of a file in LOG format, converted to the text format \(XLS\)](#)

[Sign a log file](#)

29 PROBE ADJUSTMENT



When HW4 has discovered a HygroLog NT or any other type of instrument, the device appears in the left pane of the HW4 main screen. Expanding the device displays a list of the available device functions.


To select the Probe Adjustment function, click on it with the left mouse button. HW4 opens the Probe Adjustment form.

The Probe Adjustment function is used to adjust against reference humidity and temperature conditions any HygroClip digital probe connected to any instrument or docking station. When several probes are connected to the instrument or docking station, each probe can be adjusted individually or all probes can be simultaneously adjusted.

HW4 cannot be used to adjust ROTRONIC analog probes or probes from third parties.

IMPORTANT:

- Do not interrupt the adjustment process while HW4 and the instrument are communicating as this may give unexpected results.
- Adjustment of a probe is not recommended while data is being logged.

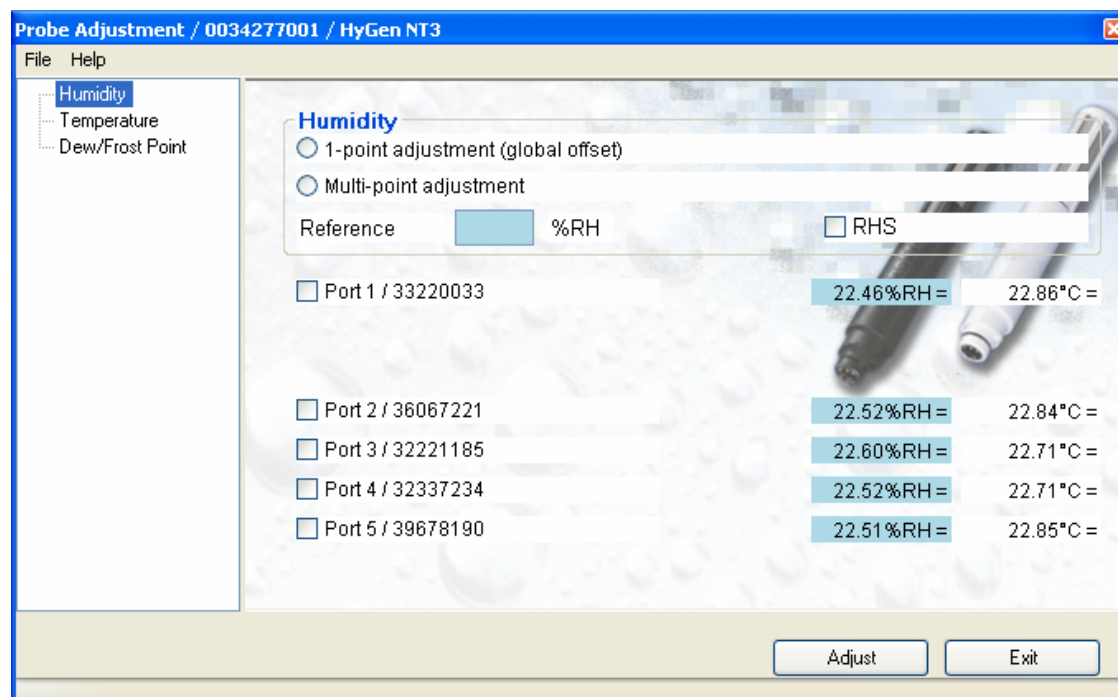

29.1 Overview

The HW4 probe adjustment function matches closely the probe adjustment function built in the firmware of our instruments. In fact, HW4 does very little other than to collect instructions and reference values from the user and to pass these to the instrument. The instrument is where all the necessary calculations are made to adjust a probe. The exception is the adjustment of relative humidity against a reference dew or frost point. When this type of adjustment is selected, HW4 calculates the reference relative humidity based on the value of the dew or frost point and on the temperature measured by the probe.

IMPORTANT:

- **Be sure that the reference environment is stable and that the probe to be adjusted has equilibrated with the environment.**
- **Always verify that the probe temperature indication is accurate. If necessary, adjust temperature prior to any humidity adjustment.**
- HW4 can be used to do either a 1-point probe adjustment (global offset) or a multi-point adjustment (offset, slope and linearization) against a known reference environment.
During a multi-point adjustment, it is essential to observe the prescribed sequence since each adjustment point influences the next.

The Probe Adjustment form displays a list of the instrument probe inputs. Unused inputs can be hidden to simplify the appearance of the form (to hide or show probe inputs, select the instrument in the device tree and use the View tab in the right pane of the HW4 main screen).

29.2 Menu Bar

The menu bar is located at the top of the Probe Adjustment form.

File

- **Exit:** closes the Probe Adjustment form

Help

- **HW4 Help:** Opens HW4 Help
- **About HW4:** Displays the version number and ID number of HW4

29.3 Temperature adjustment

- **Select Temperature** (click with the mouse in the left pane of the Probe Adjustment form)

Probe Adjustment / 0034277001 / HyGen NT3

File Help

Humidity
Temperature
 Dew/Frost Point

Temperature

☐ 1-point adjustment (global offset)
☐ Multi-point adjustment

Reference

☐ Port 1 / 33220033	22.47%RH =	22.89°C =
☐ Port 2 / 36067221	22.53%RH =	22.84°C =
☐ Port 3 / 32221185	22.61%RH =	22.74°C =
☐ Port 4 / 32337234	22.53%RH =	22.72°C =
☐ Port 5 / 39678190	22.53%RH =	22.85°C =

Adjust Exit

The form displays the humidity and temperature data from each probe connected to the instrument. The temperature data is shown on a blue background indicating that this parameter will be adjusted.

The data is updated with each polling interval. If enabled in HW4 Global Settings – View tab, a trend indicator is shown next to each parameter.

- **Select the type of adjustment**

Click with the mouse on one of the radio buttons located to the left of 1-point adjustment and Multi-point adjustment.

- **Enter the reference temperature value**

Click with the mouse on the text box labeled Reference and type in the reference temperature value. Regarding the reference value, please observe the following:

1-point temperature adjustment

For a 1-point adjustment, use a reference value that is within the limits allowed by the firmware version of your instrument (consult your instrument manual). Depending on the firmware version, the limits for the temperature reference are either $\geq -20^{\circ}\text{C}$ / -4°F and $< 40^{\circ}\text{C}$ / 104°F or for older versions $\geq -20^{\circ}\text{C}$ / -4°F and $< 70^{\circ}\text{C}$ / 158°F .

Note: this is the only adjustment available in the case of RTD probes used with any of the following HygroLog NT docking stations: DS-PT-2, DS-PT-4 and DS-PT-4-WL.

Multi-point (2-point) temperature adjustment

Be sure to observe the following sequence: Adjust the probe first at the low value (T-low). As far as possible, use a reference temperature of 20°C / 68°F¹. The second adjustment value (T-high) should be above the minimum required by the firmware version of your instrument (consult your instrument manual). Depending on the firmware version, this minimum is either ≥ 40°C / 104°F or for older versions ≥ 70 °C / 158 °F.

¹ Based on the high temperature value, the instrument firmware computes the gain using 20°C as a fixed pivot point. Because of this, T-low should be as close as possible to 20°C / 68°F so as not to introduce an error.

Note: when the reference value used for the second adjustment point is lower than the minimum programmed in the instrument firmware, the result of the adjustment is the same as with a 1-point adjustment.

- **Select the probe or probes to be adjusted**

In the list of probes, use the mouse to check the box to the left of the probe or probes to be adjusted.

☒	Port 2 / 36067221	22.54%RH =	22.87°C =
☒	Port 3 / 32221185	22.63%RH =	22.77°C =
☒	Port 4 / 32337234	22.54%RH =	22.70°C =
☒	Port 5 / 39678190	22.54%RH =	22.85°C =

Adjust Exit

- **Click on the Adjust button**

29.4 Adjustment against a reference relative humidity value

- **Select Humidity** (click with the mouse in the left pane of the Probe Adjustment form)

Probe Adjustment / 0034277001 / HyGen NT3

File Help

Humidity
Temperature
Dew/Frost Point

Humidity

☐ 1-point adjustment (global offset)
☐ Multi-point adjustment

Reference %RH ☐ RHS

☐ Port 1 / 33220033 22.48%RH = 22.89°C =

☐ Port 2 / 36067221 22.54%RH = 22.84°C =

☐ Port 3 / 32221185 22.64%RH = 22.71°C =

☐ Port 4 / 32337234 22.54%RH = 22.71°C =

☐ Port 5 / 39678190 22.54%RH = 22.81°C =

Adjust Exit

The form displays the humidity and temperature data from each probe connected to the instrument. The relative humidity data is shown on a blue background indicating that this parameter will be adjusted.

The data is updated with each polling interval. If enabled in HW4 Global Settings – View tab, a trend indicator is shown next to each parameter.

- **Select the type of adjustment**

Click with the mouse on one of the radio buttons located to the left of 1-point adjustment and Multi-point adjustment.

- **Enter the reference relative humidity value**

Click with the mouse on the text box labeled Reference and type in the reference relative humidity value. Regarding the reference value, please observe the following:

1-point humidity adjustment

For a 1-point adjustment, you may use any a reference value.

Multi-point humidity adjustment (2, 3 or 4 points)

In the case of a multi-point adjustment, be sure to observe the following sequence as each adjustment point influences the next adjustment point:

When adjusting against the ROTRONIC humidity standards (2, 3 or 4 points), always follow the sequence 35 %RH, 80 %RH, 10 %RH or 5 %RH, 0.5 %RH.

When using a reference other than the ROTRONIC humidity standards, use reference humidity values that are within the following brackets:

35 %RH	: used to compute the adjustment offset ¹
>55 %RH	: used to compute the adjustment gain
>1 %RH...≤25 %RH	: sensor linearity adjustment
≤ 1%RH	: sensor linearity adjustment

¹ Based on the second humidity adjustment, the instrument firmware computes the gain using 35 %RH as a fixed pivot point. Because of this, the first reference value should be as close as possible to 35 %RH so as not to introduce an error.

Note: The action taken by the instrument firmware is based on comparing the reference humidity value with the brackets programmed in the firmware. For example, using a reference value between 25 %RH and 55 %RH for the second adjustment point has the same effect as doing a 1-point adjustment.

- **When using a ROTRONIC humidity standard as the reference, check the box labeled RHS**

When the RHS box is checked, the effect of temperature on the standard is automatically compensated and no further correction is required. Information on the effect of temperature on each standard is provided on the cover of each box of standards. The value of the standards is not affected by barometric pressure (altitude or site elevation).

- **Select the probe or probes to be adjusted**

In the list of probes, use the mouse to check the box to the left of the probe or probes to be adjusted.

☒ Port 2 / 36067221	22.54%RH =	22.87°C =
☒ Port 3 / 32221185	22.64%RH =	22.74°C =
☒ Port 4 / 32337234	22.54%RH =	22.72°C =
☒ Port 5 / 39678190	22.54%RH =	22.81°C =

- **Click on the Adjust button**

29.5 Adjustment against a reference chilled mirror

- **Select Dew/Frost Point** (click with the mouse in the left pane of the Probe Adjustment form)

Port	Adjustment (°C)	%RH	Calculated Temp (°C)
Port 1 / 33220033	0.35°C	22.49%RH =	22.86°C =
Port 2 / 36067221	0.37°C	22.55%RH =	22.84°C =
Port 3 / 32221185	0.32°C	22.65%RH =	22.71°C =
Port 4 / 32337234	0.26°C	22.55%RH =	22.71°C =
Port 5 / 39678190	0.37°C	22.55%RH =	22.84°C =

The form displays the dew/frost point, humidity and temperature data from each probe connected to the instrument. The dew or frost point data is shown on a blue background indicating that this parameter will be adjusted.

The data is updated with each polling interval. If enabled in HW4 Global Settings – View tab, a trend indicator is shown next to each parameter.

- **Select the type of adjustment**

Click with the mouse on one of the radio buttons located to the left of 1-point adjustment and Multi-point adjustment.

- **Enter the reference dew or frost point**

Click with the mouse on the text box labeled Reference and type in the reference dew or frost point. Regarding the reference value, please observe the following:

1-point adjustment against a chilled mirror

For a 1-point adjustment, you may use any a reference value.

Multi-point adjustment against a chilled mirror (2, 3 or 4 points)

In the case of a multi-point adjustment, be sure to observe the following sequence as each adjustment point influences the next adjustment point:

When using a chilled mirror instrument as the reference, use reference values that correspond to relative humidity values within the following brackets:

35 %RH	: used to compute the adjustment offset ¹
>55 %RH	: used to compute the adjustment gain
>1 %RH...≤25 %RH	: sensor linearity adjustment
≤ 1%RH	: sensor linearity adjustment

¹ Based on the second humidity adjustment, the instrument firmware computes the gain using 35 %RH as a fixed pivot point. Because of this, the first reference value should be as close as possible to 35 %RH so as not to introduce an error.

Note: The action taken by the instrument firmware is based on comparing the reference humidity value with the brackets programmed in the firmware. For example, using a reference value between 25 %RH and 55 %RH for the second adjustment point has the same effect as doing a 1-point adjustment.

• Specify the calculation basis for the reference %RH

When the radio button labeled Dp is checked, HW4 calculates the reference humidity value, based on the assumption that the chilled mirror provides a dew point, even at values below freezing.

When the radio button labeled Dp/Fp is checked, HW4 calculates the reference humidity value, based on the assumption that the chilled mirror provides a frost point at values below freezing.

In both cases, HW4 uses the temperature read by each individual probe to calculate the resulting reference relative humidity.

• Select the probe or probes to be adjusted

In the list of probes, use the mouse to check the box to the left of the probe or probes to be adjusted.

☒ Port 2 / 36067221	0.38°C	22.57%RH =	22.84°C =
☒ Port 3 / 32221185	0.33°C	22.65%RH =	22.72°C =
☒ Port 4 / 32337234	0.26°C	22.57%RH =	22.70°C =
☒ Port 5 / 39678190	0.38°C	22.56%RH =	22.84°C =

• Click on the Adjust button

30 ALARM INDICATION & REPORTING - Overview

30.1 Standard alarm notification (all HW4 versions)

When a device detects an alarm condition, all versions of HW4 give a visual notification on the PC monitor. Device alarm conditions are defined when configuring the device with Device Manager.

30.1.1. Device tree

Any individual device reporting an alarm condition appears in red in the device tree. Device groups also appear in red when a member of the group reports an alarm condition.




30.1.2. Right Pane - Current Values Tab

Both in Device View and in Group View, values that are out of limits appear in red.

Current Values On-line Buffer Log to PC View OPC Tags						
!	#	Source	Description	Humidity	Temperature	Calculation
	2	HygroClip Pbe	External 2	69.05%RH +	23.95°C =	17.93°C +
		Logical Input	Contacts 1/2			
		Relay Status	Relay Status			

Note: limit values for each individual device and input are defined in Device Manager

Other types of alarms are also shown in red on the screen as illustrated below:

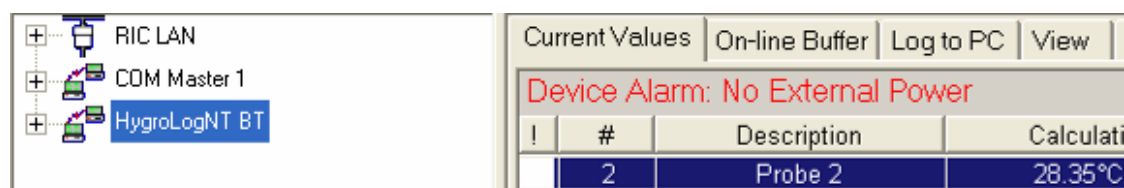

30.2 Optional alarm notification (HW4 Professional)

In addition to visualizing alarm conditions directly on the PC monitor, HW4 Global Settings can be configured so that the PC will take any combination of the following actions whenever an alarm condition occurs:

a) Give an acoustic alarm (a single beep)

HW4 Global Settings Form

View General Language/Unit System File Locations Graph Settings **Alarm Settings** Events

Trigger an alarm for the following:

Condition	Priority	Time delay [hh:mm:ss]
☒ No communication with device	Low	00:00:30
☒ Measured values out of limits or N/A	High	00:00:00
☐ HygroLog NT errors	High	00:00:00
☐ Unavailable log-to-PC or HW4 path	High	
☐ HW4 software errors	High	
☐ Alarm Test	Low	Interval: Daily

[Run an alarm test now](#)

Time limits for the alarm function

☐ Enable workdays only [Mo-Fr]

☐ Enable during limited hours Enable time [HH:mm] 00:00 Disable time [HH:mm] 00:00

When an alarm is triggered:

Action	For priority	Low	Medium	High
☐ Send e-mail		☒	☒	☒
☐ Print report		☒	☒	☒
☒ Give acoustic signal on PC (once)		☒	☒	☒
☐ Display the alarm table				

E-mail settings

Recipient: plafarie@optonline.net

☐ MAPI

☒ SMTP: mail.optonline.net

OK Cancel Help

b) Automatically display the Alarm Table

Alarm Table

File View Options

	Start Time	End Time	Device	Device name	Input name	Alarm Description
	11/24/2006 8:38 PM		HygroLab 2	HyGen Lab 2	HG Ctrl	No communication with probe

The red exclamation mark denotes an alarm condition that is still ongoing.

c) Print a report and/or send an e-mail.

When an alarm is triggered, it is either printed and / or an email is sent. The printout or e-mail reproduces the information provided in the alarm table for this particular alarm.

30.3 Alarm Reporting (HW4 Professional)

HW4 offers several ways of globally reviewing current and past alarms

30.3.1. Viewing the Alarm Table

The alarm table can be viewed at any time by clicking on Settings and Tools in the HW4 main menu bar and selecting View Alarm Table.

Note: alarm conditions that occur simultaneously are reported separately in the alarm table. For example, when both the value of temperature and the value of humidity measured by the same probe correspond to an alarm condition, the alarm table displays two separate lines, one for temperature and the other for humidity.

By default, the alarm table contains only alarms that occurred during the current HW4 session. The alarm table offers the choice of viewing past events and of filtering the alarm events in several ways.

See also [Alarm table](#).

30.3.2. Alarm summary report

An alarm summary report can be printed by selecting File and Print Summary Report from the Alarm Table menu bar. This report reproduces the current contents of the alarm table. Consequently, the contents of the report depend on the selections made as to which events and which dates to display in the Alarm Table.

30.4 Additional sources of alarm data


30.4.1. Log files

Measurement data, and other information such as the status of logical inputs, can be recorded on the PC and / or locally on a data logger. HW4 offers the possibility of automatically starting a new log file every hour, every day, etc. Within Device Manager, alarm conditions can be defined for any measured or calculated parameter or logical input.

A log file (PC or data logger) can be displayed by HW4 in the form of either a chart or a table. In both cases, values that correspond to an alarm condition can be made clearly visible on the PC screen:

- When data is viewed as a graph, any alarm band that may be associated with an individual parameter can be displayed directly on the graph (limited to one parameter at a time).
- When data is viewed as table, values that are out of limits appear over a red background (however, the red background cannot be printed).

Summarized (statistical) data can also be displayed on the PC screen and printed. As illustrated below, each individual occurrence of an alarm condition is clearly reported in the statistical data as separate block (one parameter and one occurrence per block).

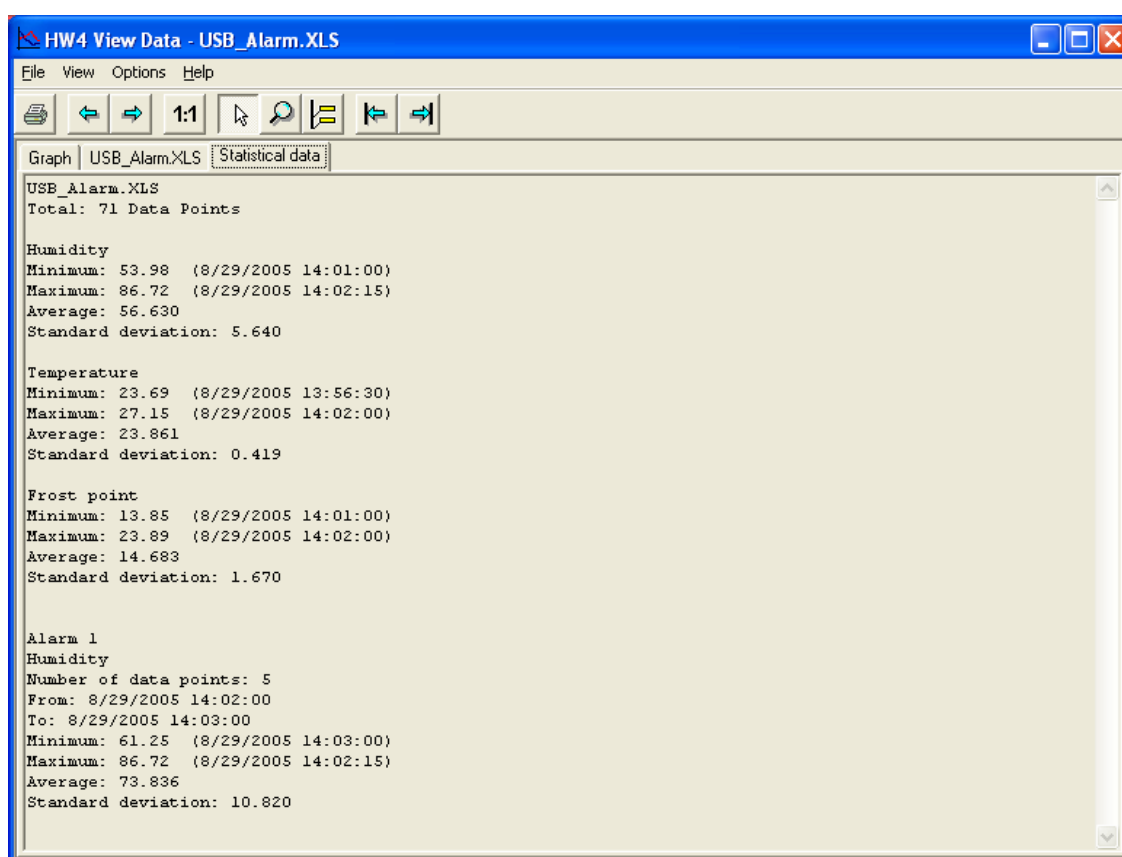

HW4 offers the possibility of printing log file data either as a graph or as a table. Statistical data can be printed separately.

30.4.2. HygroLog NT Logger event file

The HygroLog NT maintains a file that records the date, time and description of all logger events, including any alarm condition that may have been defined when configuring the logger with Device Manager. This file is split between the internal memory of the logger (up to 170 events) and the flash memory card (unlimited number of events). All past logger events are available to HW4 as long as the flash memory card is not removed from the logger.

HW4 offers the possibility of displaying and printing this file.


30.4.3. HygroClip Alarm event memory and event viewer

The HygroClip Alarm (HCA) features an internal memory that can retain up to 680 events. When the memory is full, the oldest event is erased whenever a new event is recorded. The following events are automatically recorded by the HCA:

- Logical input alarm
- End of logical input alarm
- Relay output alarm
- End of relay output alarm
- Device power up
- Watchdog overflow
- Configuration written to device
- Date/Time prior to being changed
- Date/Time after being changed
- RS-485 address changed
- Event memory erased


Accessing the events memory requires the HCA to be connected via RS-485 and a master device to a PC with the HW4 software. Using HW4, the Event Viewer is accessed or the event memory is cleared from the HCA device manager.

Event Viewer: the Event Viewer displays the most current events.

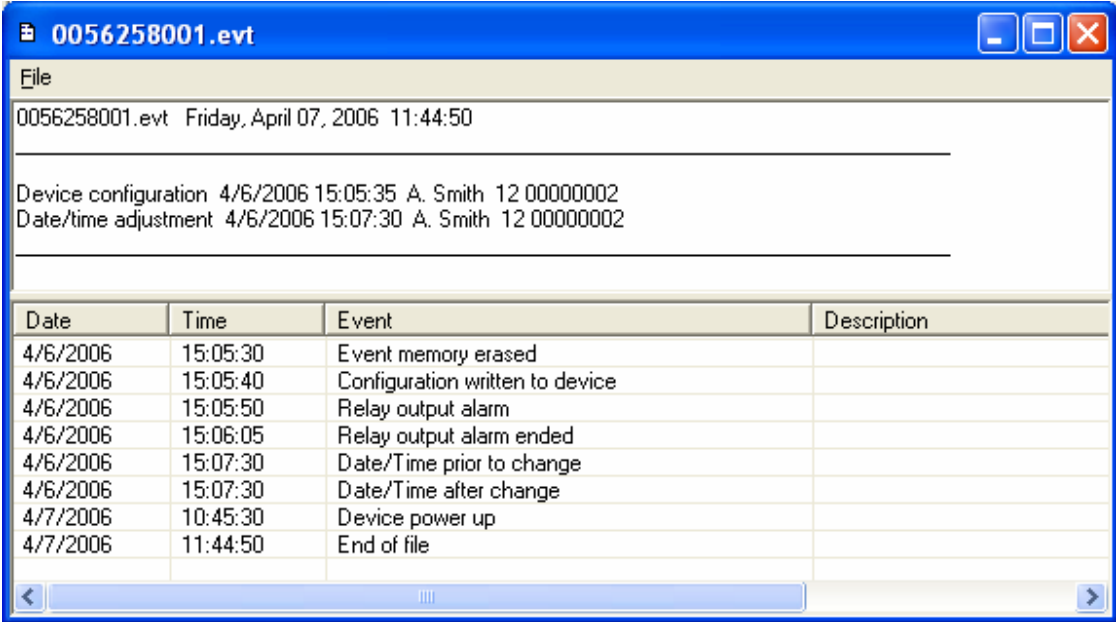


Date	Time	Event	Description
4/6/2006	15:05:30	Event memory erased	
4/6/2006	15:05:40	Configuration written to device	
4/6/2006	15:05:50	Relay output alarm	
4/6/2006	15:06:05	Relay output alarm ended	
4/6/2006	15:07:30	Date/Time prior to change	
4/6/2006	15:07:30	Date/Time after change	
4/7/2006	10:45:30	Device power up	
4/7/2006	11:44:50	End of file	

File: the file item is located on the Event Viewer menu bar:

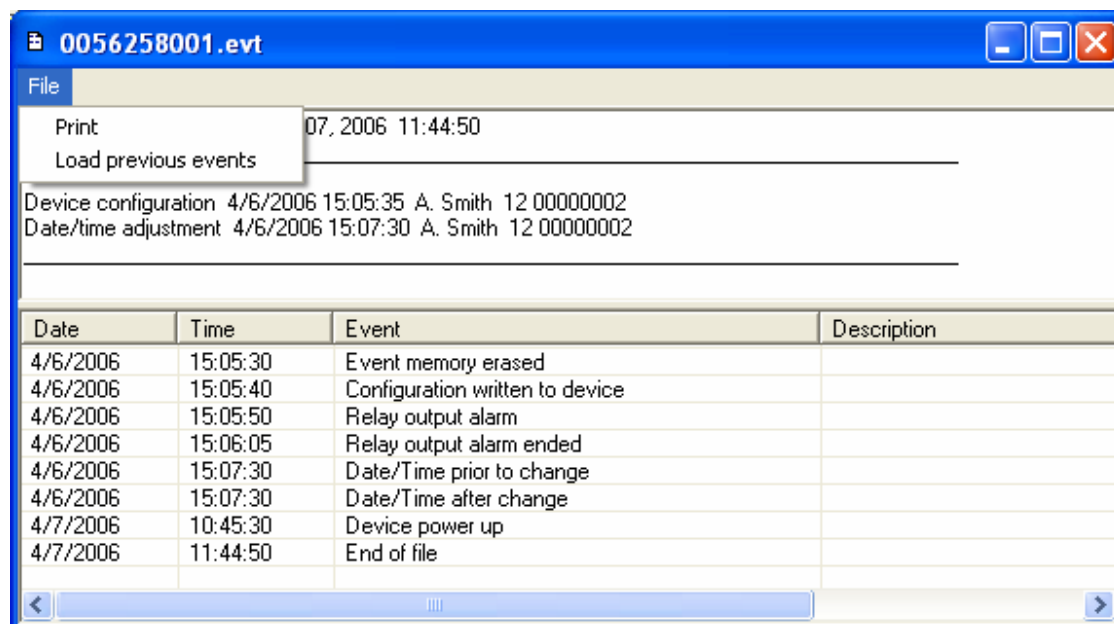

Print: prints the current contents of the event viewer to a printer or to a file.

Load previous events: opens a date selection box and allows loading in the viewer events from a previous day (this operation can be repeated several times).

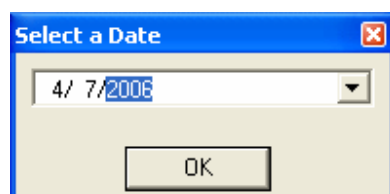

31 FDA / GAMP COMPLIANCE (HW4 Professional)

31.1 *Required settings and selections*

The following settings and selections are required in order to comply with FDA / GAMP requirements regarding electronic records, electronic signatures (data integrity) and the tracking of software problems.

- **Main Screen Menu Bar – Users and Passwords:** create at least one user with administrative rights.
- **HW4 Global Settings – General Tab:** Enable system monitoring (tracking of software problems)
- **HW4 Global Settings – Events Tab:** Enable authentication stamps and enable the monitoring of user events.
- **HW4 Global Settings – Events Tab:** enable protocols (see Record keeping by HW4)
- **Device Manager – User Information Form:** put a check mark in the box labeled “include in log file”
- **Device Manager – Keypad Form:** disable the MENU key of the HygroLog NT to prevent unauthorized or undocumented operations.
- **Log file type:** select the file extension LOG for the log files (measurement data) recorded either with HW4 on the PC or with the HygroLog NT. For the HygroLog NT, the file type is selected from Device Manager – Memory Card. For the files recorded by HW4 to the PC, the selection is done in the Log to PC tab on the Main Screen at the time the file name is entered.

31.2 *Electronic records*

In compliance with FDA / GAMP requirements regarding electronic records, electronic signatures (data integrity) and the tracking of software problems, HW4 maintains a number of event files and protocols. To effectively comply with FDA/GAMP, both types of record must be enabled in **HW4 Global Settings**. Details on the event files and protocols are provided in the “**Record keeping**” section of this manual.

31.3 Log File Format

Both HW4 (direct data logging to the PC) and the HygroLog NT offer the choice of two different types of file format to record the measured data. Both file types have two main sections: the file header and the measurement data.

- **Binary Files (LOG):** the header section of this type of file can be read with a regular text editor. As opposed to this, the data section is in binary format and cannot be read with a text editor or imported into a program such as Microsoft Excel. Both the header and data sections are protected against alterations. If the file contents are somehow modified, HW4 will display an error message when trying to open the file.

For maximum protection of the recorded data and to comply with FDA/GAMP, use the LOG file type.

- **Text Files (XLS):** these files are entirely in text format and can be read with a regular text editor. This type of file is easily imported into Microsoft Excel. Like all the other files created by HW4, files with the XLS extension are saved with the "Read Only" attribute. This attribute provides a protection against inadvertent file operations such as file delete, file move and saving the file under the same name and location (eventually after altering the file contents). Since it is possible to remove the "Read Only" attribute, this attribute does not provide protection against intentional alterations. **HW4 does not detect alterations to a file with the XLS extension.**

As an additional protection, HW4 keeps track of the date and time when a log file is created or copied to the PC. This information is kept in the protected user event file and can be compared with the file creation / file modification date and time recorded by Windows.

Note: the HygroLog NT automatically gives each log file a name comprised of the last 4 digits of the logger serial number, followed by the input number and a sequential run number.

32 RECORD KEEPING - Overview (HW4 Professional)

In compliance with FDA / GAMP requirements regarding electronic records, electronic signatures (data integrity) and the tracking of software problems, HW4 maintains the following records when enabled to do so in HW4 Global Settings:

32.1 Event Files

Both HW4 and the HygroLog NT maintain a number of files to keep track of major events. With the exception of the HW4 event file (extension ERR), these files are protected. Altering the contents of these files will cause any line that has been modified to appear on a red background when the file is opened with HW4. By default event files are located in the directory **C:\Documents and Settings\your Windows login name\Application Data\ROTRONIC_HW4\EVENT**.

32.1.1. HW4 events

A text file with the extension ERR is generated whenever HW4 encounters a software execution problem. The HW4 Event files are not directly useful to the HW4 user. These files are meant to be sent to the manufacturer and hold information useful to ROTRONIC for troubleshooting software problems. The file name consists of the text ErrorHW4Event followed by the date and time.

Example: ErrorHW4Event2005_9_18_11_58_31.ERR

To enable HW4 to keep track of software problems, place a check mark in the box labeled **Enable System Monitoring**, located in **HW4 Global Settings – General tab**.

32.1.2. User events

A user event file with the extension EVT is created by HW4 as soon as the very first session begins.

The file name consists of the text HW4USER followed by the year.

Example: HW4USER_2005.EVT

To the purpose of limiting the size of individual files, HW4 starts a new user event file at the beginning of each calendar year. Data from several years can still be reviewed on the screen and printed using the user events form.

To enable HW4 to keep track of user events, place a check mark in the box labeled **Enable Monitoring of User Events** located in **HW4 Global Settings – Events tab**.

A separate record with date and time is entered in the user event file for each event listed in the table below:

Event File Text	Event
Start HW4 ID#	HW4 ID # started
User logon	A legitimate user has logged on correctly
Log function programmed	The log function of input # has been programmed
Protocol generated	Protocol # has been generated
Authentication stamp	Checksum of protocol #
Configuration written to device	New or existing configuration written to device #
Error during probe adjustment	The probe connected to input # could not be adjusted due to an error
Probe adjusted	The probe connected to input # was adjusted
Device adjusted to PC date and time	The date and time of device # was adjusted to the PC date and time
New user created	A new user was created
Old RS address	RS-485 address of device # prior to change
New RS address	RS-485 address of device # after the change

Event File Text	Event
Language file uploaded to device	A new or existing internal language file was uploaded to device #
File deleted from PC	File # was deleted from the PC
File deleted from logger	File # was deleted from the logger memory card
User data / rights changed	User data or rights changed
Exit HW4	Exit HW4
Device deleted	Device # was deleted from the device tree
User deleted	User # was deleted
Log-to-PC started	Data recording to the PC started for input #
Log-to-PC stopped	Data recording to the PC stopped for input #
Log file path and name	Log file path and name
HW4 Global Settings saved	The HW4 Global Settings were either changed or written over
Logger data file saved to PC	Logger data file # was saved on the PC
Logger data file downloaded	Data file # was read from the logger
Logger configuration saved to PC	The logger configuration file # was saved to the PC
Logger configuration downloaded	The configuration of logger # has been read
Logger event file saved to PC	Logger event file # has been saved on the PC
Logger event file downloaded	The event file of logger # has been read
Authentication stamp of language file not valid	The HW4 language file is not a factory original
Unauthorized user tried to log on	An unauthorized user tried to log on
User could not log on (password error)	An invalid password was entered
User log off	The current user logged off

32.1.3. Logger events

The HygroLog NT maintains an internal event file with the extension EVT. No particular configuration is required to enable this feature.

Part of the procedure to ensure conformity to the FDA requirements is to disable the MENU key on the HygroLog NT keypad. The MENU key is the only one that can be used to make changes to the logger. Therefore, the events recorded by the HygroLog NT are normally the result of an interaction with the HW4 software. When the MENU key is not disabled, a limited number of events are recorded in the logger event file and no entries are made in the logger event file header. The two tables below provide a list of the events tracked by the HygroLog NT.

The logger event file is split between the internal memory of the logger (up to 170 events) and the flash memory card (practically unlimited number of events). HW4 offers the possibility of downloading, opening and printing the entire file contents. All past logger events are available to HW4 as long as the flash memory card is not removed from the logger. The serial number of the logger is used as the file name.

Example: 1111111111.EVT

The logger event file consists of a file header and a file body. The file header provides the following information:

- Most recent programming of the log function: programming date and time, user and HW4 product ID
- Most recent device configuration: date and time, user and HW4 product ID
- Most recent adjustment to the PC date and time: date and time, user and HW4 product ID

An individual record with date and time is entered in the event file body for each event listed in the table below:

Event File Text	Event (interaction with HW4)
Power up	The battery was inserted or the instrument powered down and up after an internal problem
Watchdog overflow	Internal instrument error (normally should not occur)
Writing device configuration	New or existing configuration written to logger
Memory card removed	While recording data, the logger could not find the memory card. The data was written to the EEPROM and are not yet lost
Memory card full	No free memory space on the memory card
New memory card inserted	The memory card was replaced while data was being logged. The data will be split between two different memory cards. Possibly, some data has been lost.
Humidity adjusted	Humidity adjustment of the probe connected to input #
Temperature adjusted	Temperature adjustment of the probe connected to input #
Logging started	Start recording data from input #
Logging stopped manually	Data recording of input # ended before the programmed stop time

Event File Text	Event (interaction with HW4)
Logging ended automatically	Data recording of input # ended at the programmed stop time
Out-of-limits value detected	An out-of-limits value was newly detected on input #
Battery almost empty	Battery voltage dropped below 6.5V
Battery empty	The battery is empty and the logger has powered itself off (keeping power up could result in erroneous data or loss of data)
Beginning accumulator charge	Starting to charge the rechargeable battery
Accumulator charge ended	Rechargeable battery full
MFG command	<i>Reserved for the factory</i>
Lost data, memory card not ready	While recording data, the logger could not find the memory card. The data could not be written to the EEPROM, and was lost
HygroClip probe connected ¹	A HygroClip probe was connected to input #
HygroClip probe disconnected ¹	A HygroClip probe was removed from input # or the input can no longer communicate with the probe
External power connected	The A/C adapter was connected and is being powered
External power removed or faulty	The A/C adapter was disconnected, or failed, or is not being powered
Device time changed / adjusted	Device date and time prior and after adjustment (up to firmware v1.1d only)
Prior device time	Device date and time prior to change (firmware v1.2a and up)
New device time	Device date and time after change (firmware v1.2a and up)
RS-485 address changed	The RS-485 address was changed
EEPROM erased	<i>Reserved for the factory</i>
Docking station disconnected ²	The docking station was disconnected or there is no longer any communication with it
Docking station connected ²	A docking station was connected and communication was established
Logger language file downloaded	A different internal language file was loaded or the same file was loaded again
Log function programmed	The log function has been programmed for input #
Event file deleted	<i>Reserved for the factory</i>

¹ connection / removal of analog probes is not recorded

² only when the docking station has internal electronics

Event File Text	Event (triggered from the Keypad)
Humidity adjusted	Humidity adjustment of the probe connected to input #
Temperature adjusted	Temperature adjustment of the probe connected to input #
Logging started	Start recording data from input #
Logging stopped manually	Data recording of input # ended before the programmed stop time
Prior device time	Device date and time prior to change (firmware v1.2a and up)
New device time	Device date and time after change (firmware v1.2a and up)
Log function programmed	The log function has been programmed for input #

32.2 Protocols

Protocols are text files (extension .txt or .doc) that can be generated, opened and printed from within HW4. By default, protocols are located in the directory **C:\Documents and Settings\your Windows login name\Application Data\ROTRONIC_HW4\DOC**

Protocols are enabled by placing a check mark either in the box labeled **“Enable Protocols”** or in the box **“Display Protocols”** both located in **HW4 Global Settings – Events tab**. When protocols are enabled, HW4 automatically creates the following types of file:

- **HW4 configuration protocol:** a file is created whenever the HW4 Global Settings are changed.
- **Instrument configuration protocol:** a file is created whenever the Device Manager function of HW4 is used to configure a device / instrument
- **Log function programming protocol:** a file is created whenever the Data Logging function of HW4 is used to program a data logger

Note: HW4 does not create a protocol when data is being logged directly to the PC. However, a record is entered in the user event file whenever log-to-PC is started or stopped. This record includes the date, time, location and name of the log file.

- **Probe adjustment protocol:** a file is created whenever the Probe Adjustment function of HW4 is used to adjust a probe connected to an instrument

In addition to a detailed description of each event, protocol files include the following information:

- Date and time
- HW4 current user
- HW4 version and ID number (product key)

To prevent tampering, HW4 Professional adds an authentication stamp at the end of each protocol. To have HW4 automatically validate the authentication stamp whenever a protocol file is opened from within HW4, place a check mark in the box labeled **“Authentication stamp”** box located in **HW4 Global Settings – Events tab**.

When the tracking of user events is enabled:

- Whether protocols are enabled or not, HW4 enters a record of the event in the user event file.
- When protocols are enabled, HW4 enters an additional record in the user event file for cross-reference purposes. This record includes the location and name of the corresponding protocol file. Similarly, the protocol text file includes a reference to the user event file (location and file name).

32.2.1. Protocol file names

32.2.1.1. HW4 configuration protocols

The file name for this type of protocol consists of the text HW4_Extras, followed by a sequential number generated by HW4.

Example of an HW4 configuration protocol file name: HW4_Extras_1.txt

32.2.1.2. Instrument configuration protocols

The file name for this type of protocol consists of the serial number of the instrument followed by the letters CNF and a sequential number generated by HW4.

Example of a configuration protocol file name: 111111111_CNF_3.txt

32.2.1.3. Log function programming protocols

The file name for this type of protocol consists of the serial number of the instrument followed by the letters PROG and a sequential number generated by HW4

Example of a log function programming protocol file name: 111111111_PROG_0.txt

32.2.1.4. Probe adjustment protocols

The file name for this type of protocol consists of the serial number of the instrument followed by the letters ADJ and a sequential number generated by HW4


Example of a probe adjustment protocol file name: 111111111_ADJ_2.txt

32.3 Cross referencing protocols and event files

By cross referencing the date and time of log files - or the date, time and authentication stamp of the event files and protocols, it is possible to keep track of events and to identify which user was logged in HW4 at the time of the event.

Example: cross referencing a probe adjustment protocol with the user event file

a) Open the probe adjustment protocol:

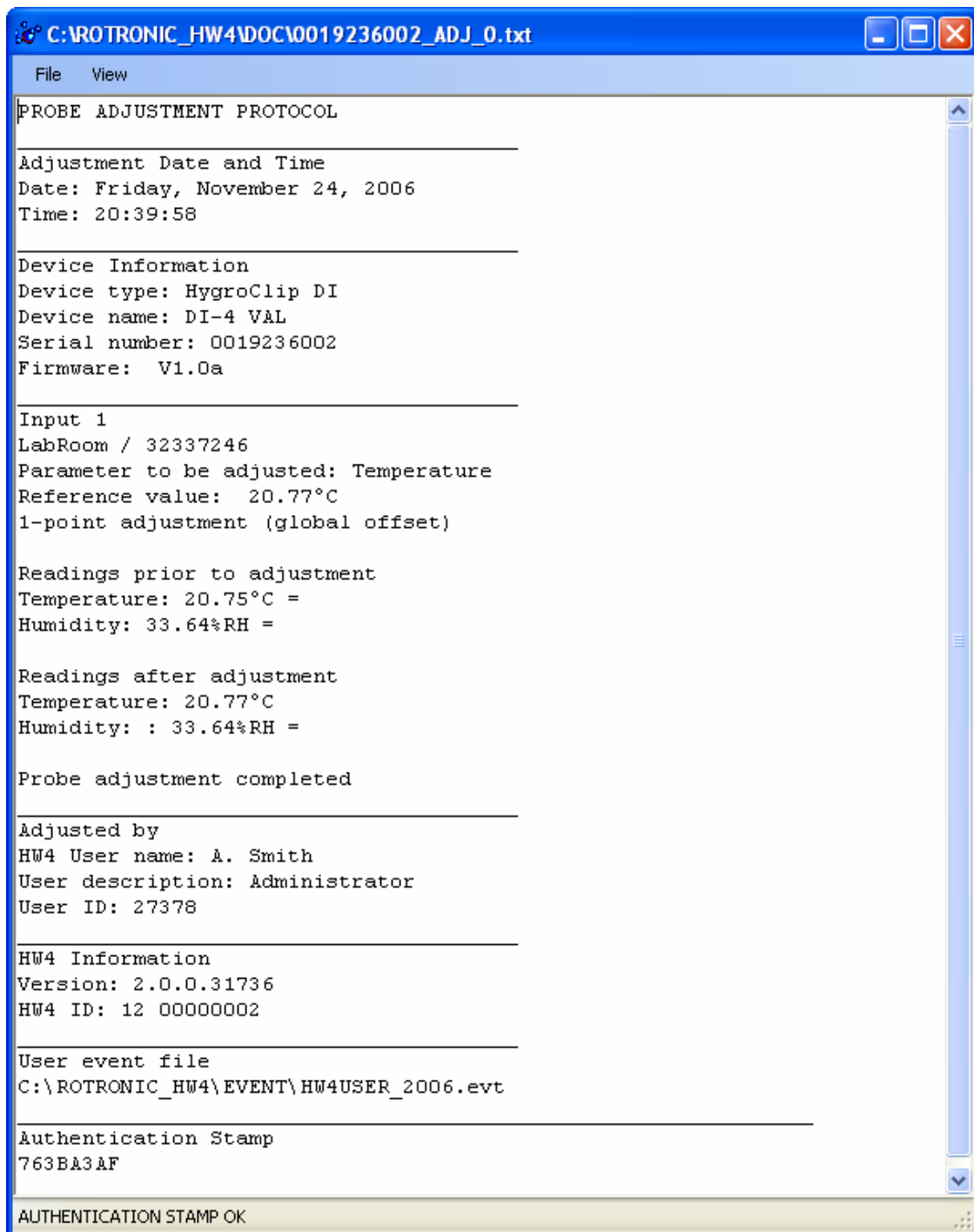

The following information is located at the end of the probe adjustment protocol:

- User information
- HW4 information
- User event file: path and name
- Authentication stamp

b) Open the user event file referenced in the probe adjustment protocol:

HW4USER_2006.evt			
File			
Date	Time	Event	Description
11/24/2006	20:37:55	Start HW4 ID#	12 00000002 V2.0.0.31736
11/24/2006	20:38:00	User Log on	A. Smith
11/24/2006	20:38:00	Authentication Stamp of language file not valid	C:\Program Files\HW4\Language\English.txt
11/24/2006	20:38:50	Error during probe adjustment	HygroLogNT 0033414009
11/24/2006	20:38:50	Protocol generated	C:\ROTRONIC_HW4\DOC\0033414009_ADJ_0.txt
11/24/2006	20:38:50	Authentication Stamp of protocol	528EB2D1
11/24/2006	20:39:15	Error during probe adjustment	HygroLogNT 0033414009
11/24/2006	20:39:15	Protocol generated	C:\ROTRONIC_HW4\DOC\0033414009_ADJ_1.txt
11/24/2006	20:39:15	Authentication Stamp of protocol	C1BB0E27
11/24/2006	20:39:55	Probe adjusted	HygroClip DI 0019236002
11/24/2006	20:39:55	Protocol generated	C:\ROTRONIC_HW4\DOC\0019236002_ADJ_0.txt
11/24/2006	20:39:55	Authentication Stamp of protocol	763BA3AF

Look for the event “Probe adjusted” (in this example, line 10 of the user event file). Since there could be several such events, make use of the date and time of the probe adjustment protocol to find the correct event. Compare the following information with the information located at the end of the probe adjustment protocol: user name, description and group

Lines 10 and 11 of the user event file report that a probe connected to the HygroClip DI S/N 0019236002 was adjusted and that a protocol was generated. The date and time should match the date and time of the probe adjustment protocol. The authentication stamp on line 12 of the user event file should match the authentication stamp of the probe adjustment protocol.

33 FILE PROTECTION

The data and record keeping files generated by HW4 are of two types: binary files and text files. Binary files use the file extension LOG while text files use a variety of file extensions: TXT, XLS, EVT, ALR and ERR.

All data and record keeping files created by HW4 are saved with the “Read Only” attribute. This provides a protection against inadvertent file operations such as file delete, file move and saving the file under the same name and location (eventually after altering the file contents). Since it is possible to remove the “Read Only” attribute, this attribute does not by itself provide protection against intentional alterations.

33.1 Authentication stamp

HW4 makes use of one or more authentication stamps to verify that a file has not been altered. With the exception of files with the XLS and ERR file extensions all data and record keeping files generated by HW4 are stamped. HW4 will display an error message when opening a stamped file that has been altered.

- **Binary Files (LOG):** Both the header and data sections are protected by an authentication stamp against alterations.
- **Text Files (TXT):** The file extension TXT is used for all protocols generated by HW4. These files are protected by a single authentication stamp.
- **Text Files (EVT):** The file extension EVT is used for all event files (user and device events). Each individual file entry is protected by an authentication stamp.
- **Text Files (ALR):** The file extension EVT is used for the file holding the HW4 alarm records captured by HW4. Each individual file entry is protected by an authentication stamp.

33.2 Overview

File extension	Usage	Protected
LOG	Log files (recommended format)	YES
XLS	Log files (alternate format)	NO
ERR	Service files used for trouble shooting	NO
TXT	Protocols	YES
EVT	User and device events	YES
ALR	Alarm records	YES

34 RELOCATING THE HW4 USER FOLDER

WARNING:

Relocating the HW4 User Folder may require re-entering your HW4 ID Number, all HW4 users, the HW4 Global Settings and may also require a new search for the devices connected to the PC.

Prior to relocating the HW4 User Folder, be sure to have a record of your HW4 ID Number. When HW4 is running, this information can be found by clicking on Help – About ROTRONIC HW4 in the HW4 main menu bar.

During the initial startup, HW4 automatically creates a user folder named **ROTRONIC_HW4** folder and several subfolders.

The ROTRONIC_HW4 folder holds the HW4 configuration data. The subfolders DATA, DATA_ONLINE, SYS, DOC and EVENT hold the log files (measurement data), protocols, event files and alarm records.

By default HW4 creates the User Folder on drive C of the PC starting HW4 for the first time and uses the following path:

C:\Documents and Settings\Windows User\Application Data

Windows User is the Windows user account that was logged in at the time HW4 was installed.

Windows automatically gives *Windows User* **exclusive access** to the directory C:\Documents and Settings\Windows User and subfolders.

If for any reason the default path for the HW4 user folder is not either not convenient or suitable, HW4 can be configured to use a different path ¹.

¹ Example: in order to share the contents of the **ROTRONIC_HW4** folder between several workstations, the folder can be located on a Windows 2003 server. The server drive where ROTRONIC_HW4 is located should be mapped to all workstations that will use HW4 and sufficient permissions to the ROTRONIC_HW4 folder and subfolders must be given by the server administrator to the different workstation users.

IMPORTANT:

- Do not change the name of the HW4 user folder. This name must always be ROTRONIC_HW4.
- You can change freely the name of any of the ROTRONIC_HW4 subfolders

Two different methods are available to change the path of the HW4 user folder:

34.1 Method 1 (using Windows Explorer)


Note: Method 1 can be used just after installing HW4 and prior to starting HW4 for the first time or at any time thereafter when HW4 is not running.

Basic assumptions:

- HW4 is already installed on the PC or workstation, in the folder C:\Program Files\HW4


- **Step 1:**

Open the directory C:\Program Files\HW4\tools (created during the installation of HW4).

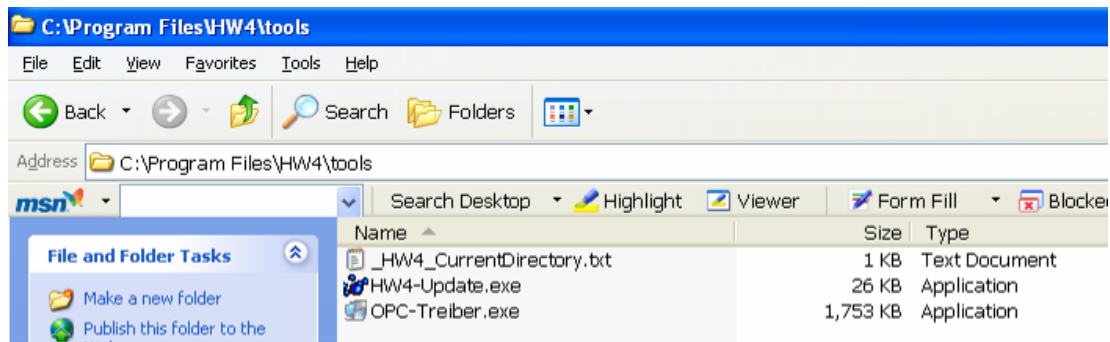


- **Step 2:**

Select the text file **_HW4_CurrentDirectory.txt** and copy it in the same directory.

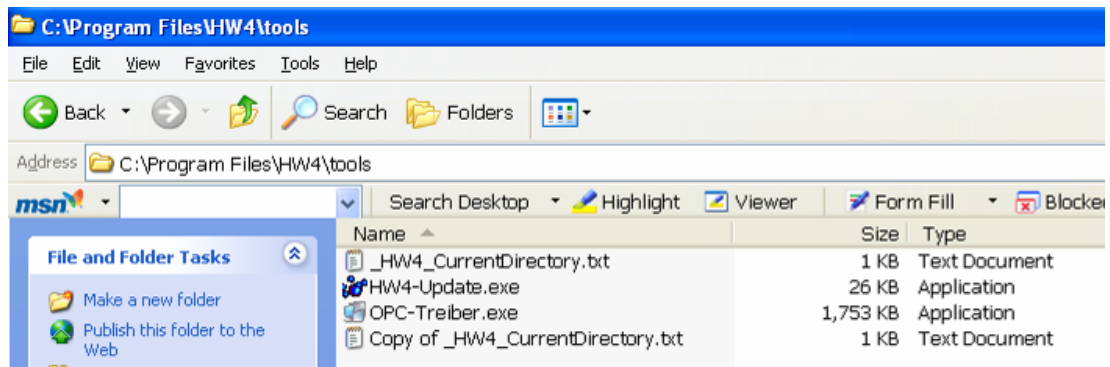


- **Step 3:**

Right click on **Copy of _HW4_CurrentDirectory.txt** and select **Properties**. Remove any Read-only attribute that may be present.

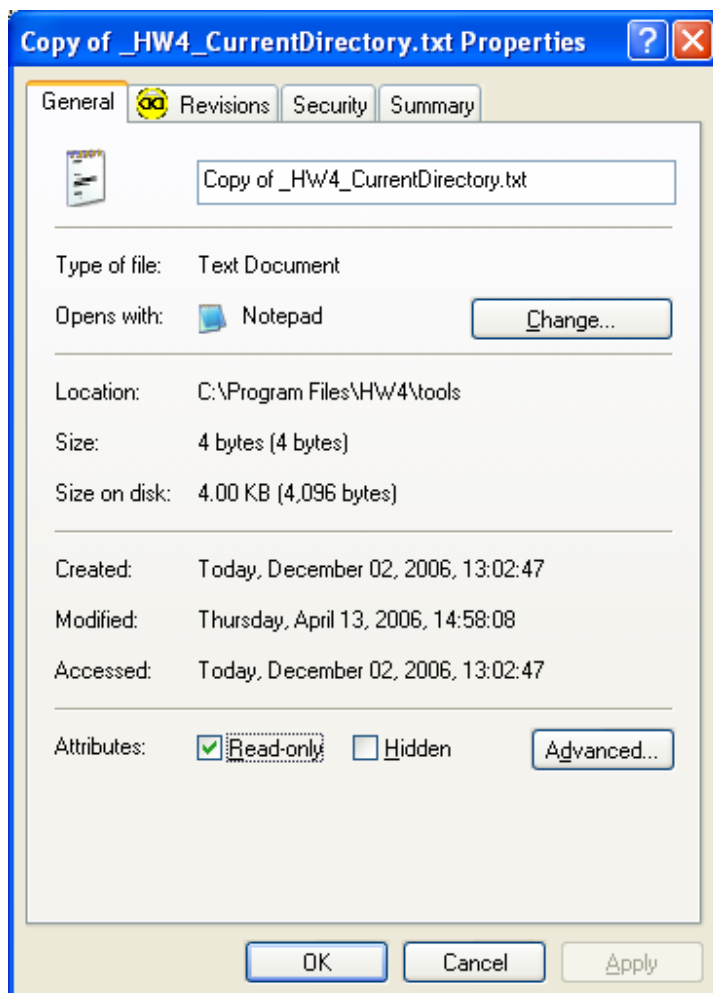

- **Step 4:**


Rename the file copy to **HW4_CurrentDirectory.txt** (no leading underline).

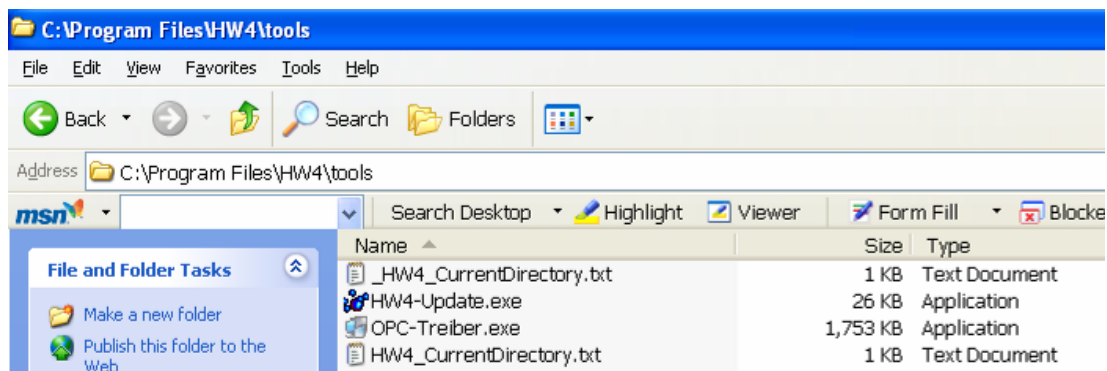


• Step 5:

Open **HW4_CurrentDirectory.txt** with Notepad. Delete the existing text line and replace it with the new path for the HW4 User Folder (**do not put a backslash at the end of the line**).


Examples:

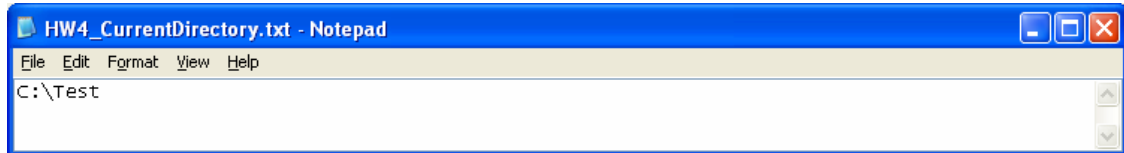

End result: C:\Test\ROTRONIC_HW4

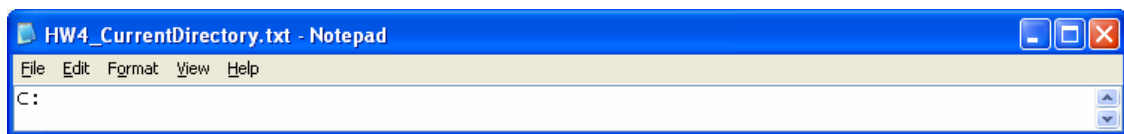

End result: C:\ROTRONIC_HW4

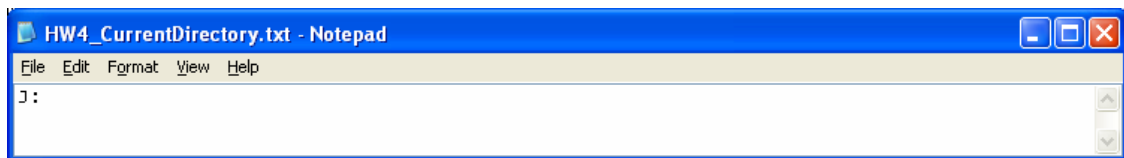

End result: J:\ROTRONIC_HW4 (J: is a mapped server drive)

Save and close the file.

• Step 6:

HW4 will automatically create a new ROTRONIC_HW4 and subfolders at the desired location.

IMPORTANT:

- To avoid the inconvenience of a new initial start-up proceed as explained in [Retrieving your previous settings and other data](#).
- If you do not want to retrieve your settings and data, start HW4 and proceed as described in [Initial Start-up](#).
- To make HW4 use again the default location of the ROTRONIC_HW4 folder (C:\Documents and Settings\Windows User\Application Data), delete the file **HW4_CurrentDirectory.txt**.

34.2 Method 2 (after the initial HW4 start-up)

- **Basic assumptions:**

- HW4 is already installed on the PC and has been started at least one time

- **Step 1:**

In the HW4 main menu bar select Settings and Tools > HW4 Global Settings > File locations.

Note: after the initial start-up, the default location of the HW4 user folder and subfolders is as follows:

HW4 Global Settings Form

View | General | Language/Unit System | **File Locations** | Graph Settings | Alarm Settings | Events

Measurement Data	
Logger Data Files	C:\Documents and Settings\admin.ROTRONIC\Application Data\ROTRONIC HW4\DATA
Log to PC Files	C:\Documents and Settings\admin.ROTRONIC\Application Data\ROTRONIC HW4\DATA ONLINE

Configuration Data	
User Folder Path	C:\Documents and Settings\admin.ROTRONIC\Application Data
Device Configuration Files	C:\Documents and Settings\admin.ROTRONIC\Application Data\ROTRONIC HW4\SYS

Protocols and Events	
Protocols	C:\Documents and Settings\admin.ROTRONIC\Application Data\ROTRONIC HW4\DOC
Device Event Files	C:\Documents and Settings\admin.ROTRONIC\Application Data\ROTRONIC HW4\EVENT
Alarm Records	C:\Program Files\HW4\EVENT
HW4 Event Files	C:\Documents and Settings\admin.ROTRONIC\Application Data\ROTRONIC HW4\EVENT
User Event Files	C:\Documents and Settings\admin.ROTRONIC\Application Data\ROTRONIC HW4\EVENT

OK Cancel Help

- **Step 2:** Click on the blue link next to User Folder Path. HW4 issues the following warning:


HW4

WARNING: changing the HW4.ini path will require re-entering the HW4 license data, HW4 users and HW4 Global Settings and will also require a new search for devices connected to the PC. Proceed ?


OK Cancel

Click on the OK button.

HW4 opens the following box:

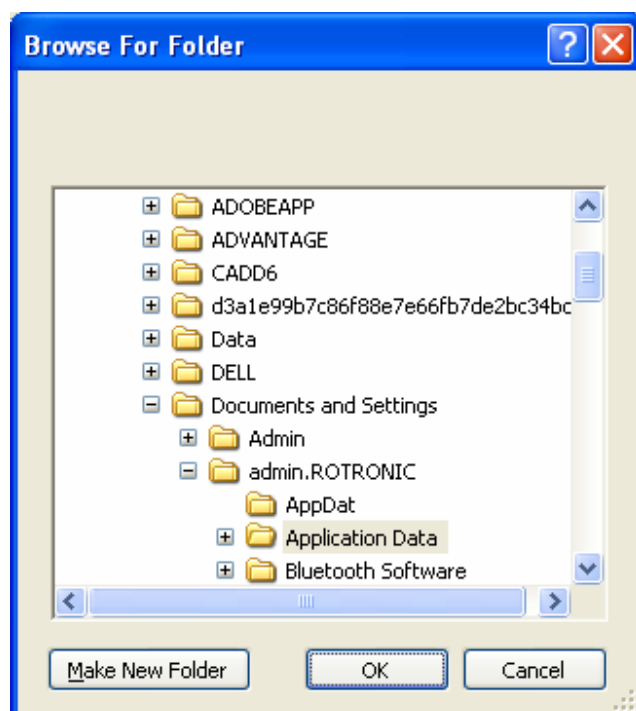


- **Step 3:** Make the selection pointing HW4 to the desired location of the HW4 User Folder (in this example the root directory of drive C)

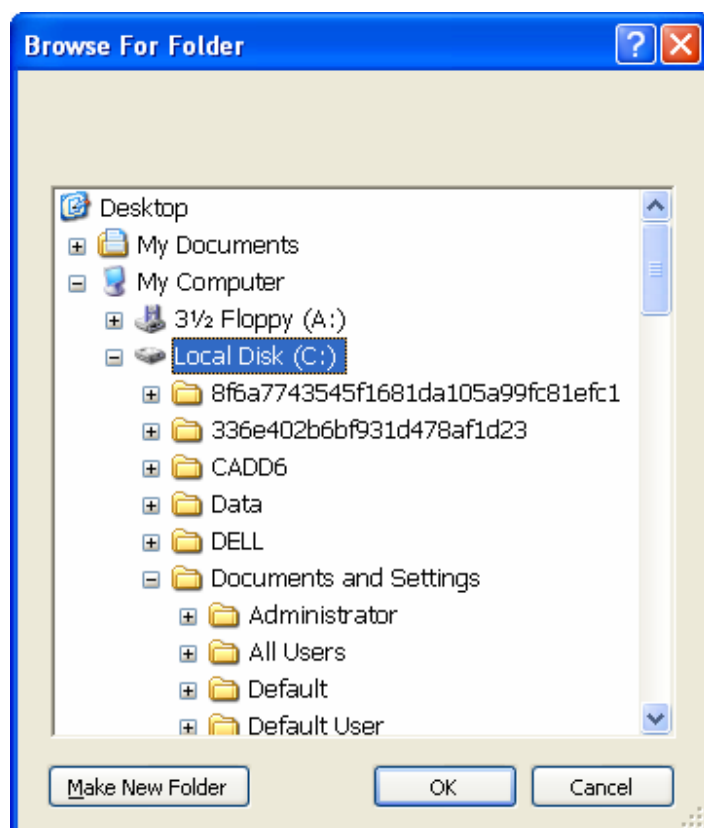


Click on the OK button.

The following message box appears. Click on the OK button:

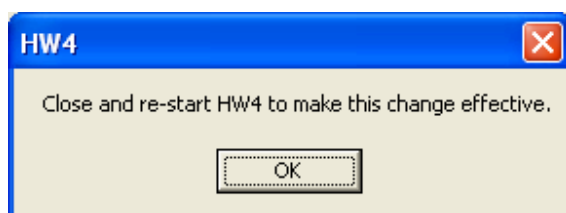

The HW4 Global settings form changes as follows:

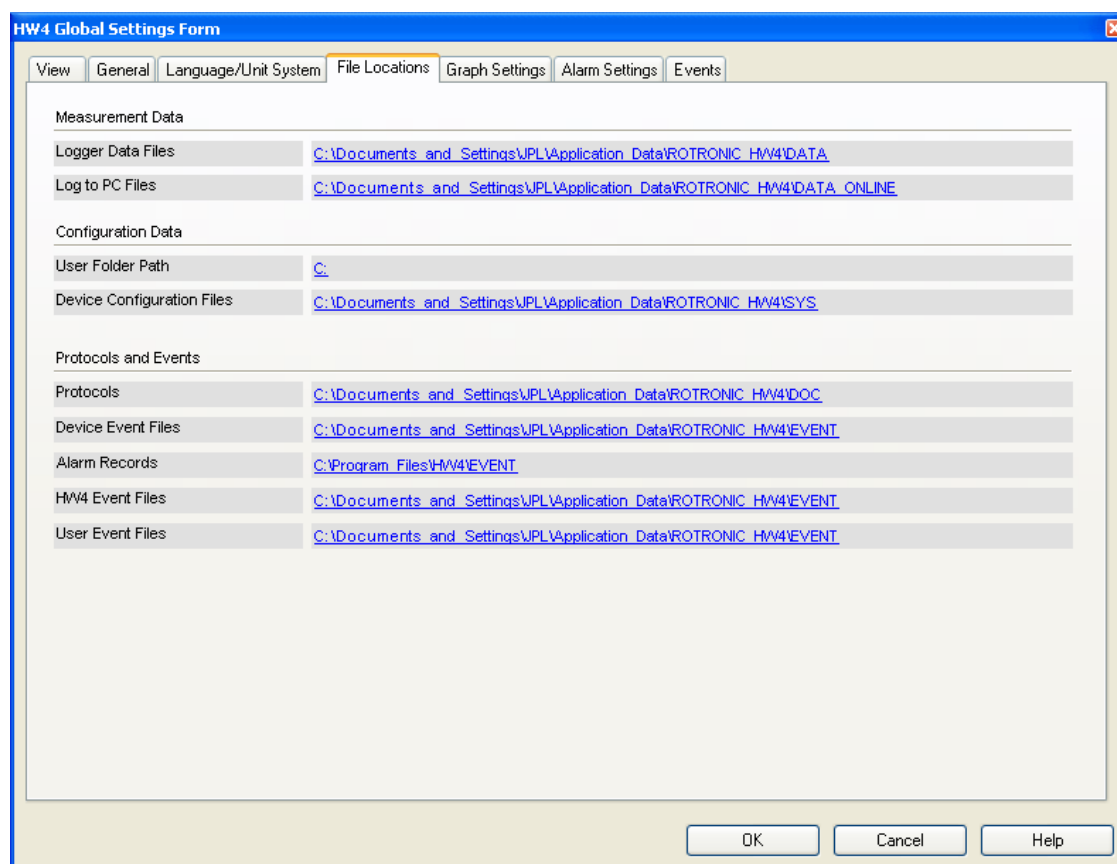

Click on the OK button at the bottom of the HW4 Global settings form

• **Step 4:**

IMPORTANT: to avoid the inconvenience of a new initial start-up proceed as explained in [Retrieving your previous settings and other data](#).

If you do not want to retrieve your settings and data , start HW4 and proceed as described in [Initial Start-up](#).

HW4 will automatically create a new ROTRONIC_HW4 and subfolders at the desired location.

34.3 Retrieving your previous settings and other data

After relocating the HW4 User Folder, close but do not restart HW4.

- After closing HW4, open My Computer. Select C: > Documents and Settings > *Windows User* > Application Data. Open the Application Data folder.

Windows User is the Windows user account that was logged in at the time HW4 was installed.

Select the folder named ROTRONIC_HW4. Assuming that the new location of the HW4 User Folder is the root directory of drive C, either copy or move ROTRONIC_HW4 to C:\ - **do not rename the folder**.

Note: you may have to change the Windows permissions to the folder and subfolder (Security Tab in folder properties).

- Start HW4 again.

In the HW4 main menu bar select Settings and Tools > HW4 Global Settings > File locations:

The screenshot shows the 'HW4 Global Settings Form' with the 'File Locations' tab selected. The form is organized into three sections: Measurement Data, Configuration Data, and Protocols and Events. Each section contains a list of file paths.

Section	Item	Path
Measurement Data	Logger Data Files	C:\Documents and Settings\JPL\Application Data\ROTRONIC_HW4\DATA
	Log to PC Files	C:\Documents and Settings\JPL\Application Data\ROTRONIC_HW4\DATA_ONLINE
Configuration Data	User Folder Path	C:
	Device Configuration Files	C:\Documents and Settings\JPL\Application Data\ROTRONIC_HW4\SYS
Protocols and Events	Protocols	C:\Documents and Settings\JPL\Application Data\ROTRONIC_HW4\DOC
	Device Event Files	C:\Documents and Settings\JPL\Application Data\ROTRONIC_HW4\EVENT
	Alarm Records	C:\Program Files\HW4\EVENT
	HW4 Event Files	C:\Documents and Settings\JPL\Application Data\ROTRONIC_HW4\EVENT
	User Event Files	C:\Documents and Settings\JPL\Application Data\ROTRONIC_HW4\EVENT

At the bottom of the form are three buttons: OK, Cancel, and Help.

- Change the path of the different subfolders used to hold data, configuration, protocol and event files. To change a path, click on the corresponding blue link.

The screenshot shows the 'HW4 Global Settings Form' with the 'File Locations' tab selected. The form is organized into three sections: Measurement Data, Configuration Data, and Protocols and Events. Each section contains a list of file paths, most of which are blue links indicating they can be edited.

Section	Setting	Value
Measurement Data	Logger Data Files	C:\ROTRONIC_HW4\DATA
	Log to PC Files	C:\ROTRONIC_HW4\DATA_ONLINE
Configuration Data	User Folder Path	C:
	Device Configuration Files	C:\ROTRONIC_HW4\SYS
Protocols and Events	Protocols	C:\ROTRONIC_HW4\DOC
	Device Event Files	C:\ROTRONIC_HW4\EVENT
	Alarm Records	C:\ROTRONIC_HW4\ALARM
	HW4 Event Files	C:\ROTRONIC_HW4\EVENT
	User Event Files	C:\ROTRONIC_HW4\EVENT

At the bottom of the form are three buttons: OK, Cancel, and Help.

- Click on the OK button at the bottom of the form.
- Close and restart HW4

35 CHANGING THE BAUD RATE OF AN ETHERNET DEVICE

When an Ethernet device ¹ is used as the master device in a RS-485 multi-drop network, it may be necessary to change the factory default baud rate because all devices in the RS-485 multi-drop must use the same baud rate. See: [Connecting and discovering additional devices](#)

¹ Examples: HygroLog NT docking station, HygroClip DI

Changing the baud rate of an Ethernet device is a two step process:

- (1) Change the baud rate in HW4 > Device Manager > Interface
- (2) Change the baud rate using the Web Interface of the internal Digi International module: Configuration > Serial Port > Basic Serial Settings

The Web Server of the Digi International module can be accessed either directly from Internet Explorer (type in the IP address of the module and click on Go) or from the HW4 main menu bar > Settings and Tools > Ethernet Configuration Tool (Digi Device Discovery). When Digi Device Discovery is open, right click on the device with the mouse and select “**Open web interface**”.

In the case of a device with a wired connection, no user name or password are required to open the web server. By contrast, a device with a wireless connection (WI-ME) requires a user name and password. The ROTRONIC factory defaults are as follows:

- User name: **rotronic**
- Password : **wlan**

Note: if for any reason, you had to reset the WI-ME module to the original Digi International manufacturer defaults, the user name reverts to **root** and the password **dbps**.

The home page of the Digi module web server is as follows:

Home

Configuration

- Network
- Serial Port
- GPIO
- Alarms
- Security
- System

Management

- Serial Ports
- Connections

Administration

- File Management
- Backup/Restore
- Update Firmware
- Factory Default Settings
- System Information
- Reboot

[Logout](#)

Home

Getting Started

Tutorial Not sure what to do next? This Tutorial can help.

System Summary

Model: Digi Connect ME

IP Address: 192.168.57.152

MAC Address: 00:40:9D:24:60:3D

Description: JPL Test1

Contact: JPL Test 1

Location: JPL Test 1

User Interfaces

Web Interface (Default)

Custom Interface:

The left hand side of the web server home page displays a menu. Left click with the mouse on Serial Port to gain access to the Serial Port Configuration page:

Serial Port Configuration - HygroLog NT

▼ Port Profile Settings

Current Port Profile: **TCP Sockets** [Change Profile...](#)
The TCP Sockets Profile allows a serial device to communicate over a TCP r

TCP Server Settings

Connect directly to the serial device using the following TCP ports on the n

☒ Enable Telnet access using TCP Port:

☒ Enable Raw TCP access using TCP Port:

☒ Enable Secure Socket access using TCP Port:

TCP Client Settings

Automatically establish bi-directional TCP connections between the serial d

☐ Automatically establish TCP connections

Establish connection under one of the following conditions:

☒ Always connect and maintain connection

☐ Connect when data is present on the serial line

Match string:

☐ Strip string before sending

☐ Connect when DCD (Data Carrier Detect) line goes high

☐ Connect when DSR (Data Set Ready) line goes high

Establish connection to the following network service:

IP Address:

Service:

TCP Port:

► Basic Serial Settings

► Advanced Serial Settings

Click with the mouse on Basic Serial Settings.



Digi Connect ME Configuration and Managemen

[Home](#)

Configuration

[Network](#)

[Serial Port](#)

[GPIO](#)

[Alarms](#)

[Security](#)

[System](#)

Management

[Serial Ports](#)

[Connections](#)

Administration

[File Management](#)

[Backup/Restore](#)

[Update Firmware](#)

[Factory Default Settings](#)

[System Information](#)

[Reboot](#)

Serial Port Configuration - HygroLog NT

► Port Profile Settings

▼ Basic Serial Settings

Description:

Baud Rate: ▼

Data Bits: ▼

Parity: ▼

Stop Bits: ▼

Flow Control: ▼

► Advanced Serial Settings

Set the baud rate to match the value entered in Device Manager > Interface.

Warning: exercise caution when using the device web server:

- Be sure to use settings that are compatible with your LAN.
- The menu item "Serial Port" on the home page allows the configuration of the TCP ports used by the device. No configuration should be required. Whatever you do, "Enable raw TCP access using TCP port" should always be checked and the port number should always be 2101.


36 CONCURRENT HW4 SESSIONS ON DIFFERENT WORKSTATIONS

HW4 can be installed on a number of workstations connected to a file server and the HW4 User Folder (ROTRONIC_HW4) can be relocated to a file server drive that is mapped to all workstations. The purpose of doing this is to be able to use the same data, protocols, events and alarms on all workstations.

IMPORTANT: Installing the HW4 software on the file server is not recommended and may not work.

LIMITATIONS: Running concurrent sessions of HW4, each on a different workstation is subject to the following limitations.

- The HW4 user who opens the first session (or Master session) can use HW4 with all the rights given to that user by the HW4 administrator.
- When an additional session starts on one of the other workstations, HW4 issues the following warning:

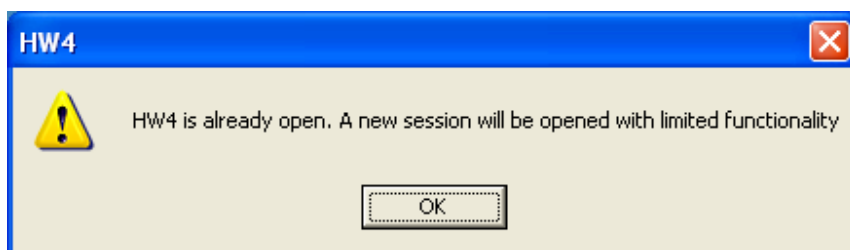

After clicking on the OK button, HW4 starts the new session. The user of this session has only minimum rights (essentially limited to viewing the measurement data), regardless of the rights given to this user by the HW4 administrator. This limitation remains in effect until the end of the session, even when the original Master session has been ended.

- The new session is not tracked by the HW4 user event file (user login, etc.).
- The Master session is the only one to write to the HW4.ini file (HW4 main screen configuration, etc.). In the main, this happens only at the time when this session is closed
- When running concurrent sessions of HW4, the live memory of each workstation holds the information that was in the HW4.ini file at the time the workstation opened the HW4 session. No complete update of this file takes place during a session. This means that workstations that are running the additional HW4 sessions are essentially limited to viewing measurement data and that they do not necessarily have current information regarding any of the changes made by the Master first session or any activity of that session.
- The HW4 edition is automatically the same as the edition used to open the Master session.
- Both the online buffer and online graph are local to each workstation and can be locally cleared without affecting the other workstations.
- The only devices that can be seen simultaneously by several workstations are devices connected either directly to the LAN (Ethernet) or indirectly connected to the LAN by means of a RS-485 multi-drop. Connections via Internet are similar to a LAN connection.

- Conflicts may result when the same device is simultaneously polled by several workstations. Using a longer polling interval may make such conflicts less frequent. When there is a conflict, a device icon may temporarily appear in the device tree with a red cross over it. This does not necessarily mean that communication with the device was lost and it may simply show that the device is busy. For example, a HygroLog NT that is in the process of uploading a log file to the server is likely to appear busy to a number of workstations. False alarms may be avoided by using an appropriate time delay for the alarm type “no communication with device”.


CONFIGURATION: After installing HW4 on the different workstations, each workstation can be configured to make HW4 use a common HW4 User Folder (ROTRONIC_HW4) that is located on a drive of the file server that is mapped to all workstations.

See [Relocating the HW4 User Folder](#)

It is important to distinguish between a Windows user and a HW4 user. Do not confuse HW4 rights with the permissions to a Windows folder. The administrator of the Windows 2003 server should grant sufficient permissions to the HW4 user folder and subfolders (Folder Properties – Security tab) so that all workstations have access to the same files.

NOTE: The file HW4RunningStatus.dat is located in the HW4 User Folder (ROTRONIC_HW4). This file holds the IP address of the workstation that is (or was) running the Master session. This information may be used to identify the location of this workstation.

HW4RunningStatus.dat can be opened with Notepad:

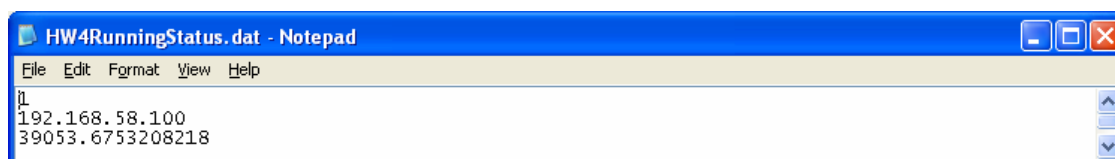


The first line is either a 1 or a 0 (1 means the Master session is still open / 0 means that a new Master session can be started).

37 MULTIPLE HW4 SESSIONS ON THE SAME PC

We do not recommend running simultaneously multiple sessions of HW4 on the same PC or workstation. Doing so may give unpredictable results.

Nevertheless, if you so desire, HW4 will let you open more than one session on a PC or workstation and will issue the following warning:

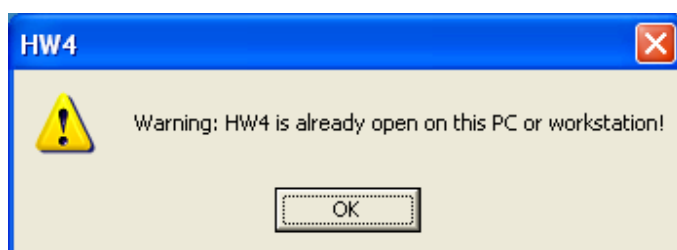

38 STARTING HW4 AUTOMATICALLY with Windows

HW4 can be started automatically by Windows, for example after re-starting the PC. To do this a few steps are required.

- Configure Windows so that it can be started without requiring any keystroke or mouse click.
- During the initial start-up HW4 creates the HW4 User Folder following folder in C:\Documents and Settings\Windows User\Application Data\ROTRONIC_HW4.

Relocate this folder so that it can be accessed by all Windows users.

See: [Relocating the HW4 User Folder](#)


- In the HW4 main menu bar, select Settings and Tools > HW4 Global Settings > General. Use the mouse to place a check mark in the box labeled "Auto-Start Mode".

See: [General Tab](#)

- In Windows, right click on Start and select Explore all Users.

Select Documents and Settings > All Users > Start Menu > Programs. In Programs, open the folder Startup. Right click with the mouse on a blank spot within the folder view and select New > Shortcut.

This opens the create shortcut wizard:

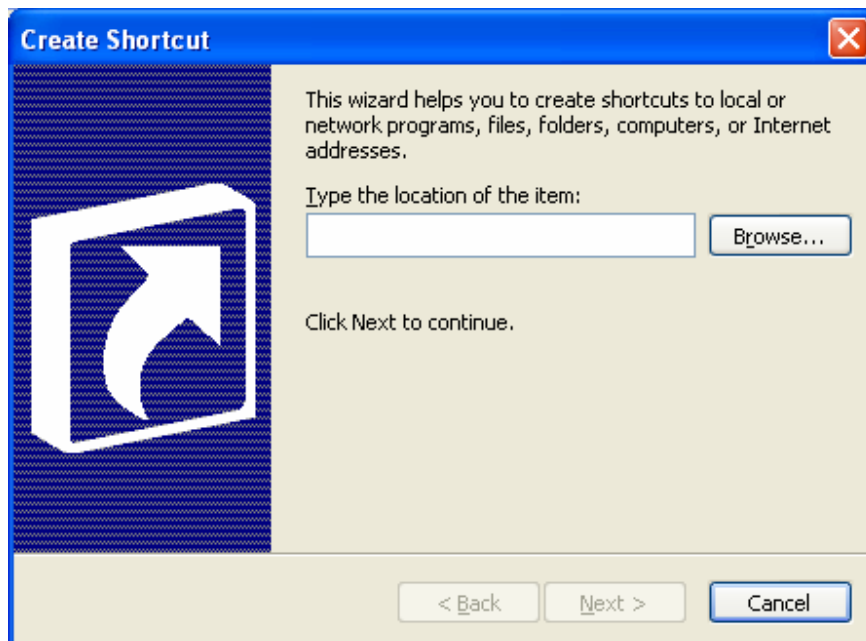

Use the Browse button to locate HW4.exe in C:\Program Files\HW4.

Click on the OK button to save the changes. Because the shortcut is in the Startup folder in All Users, HW4 will start automatically each time that Windows starts up.

39 PROBE ADJUSTMENT BASICS

39.1 *Temperature*

Note: the stability of the Pt100 RTD sensor used to measure temperature is such that temperature adjustment in the field is seldom required. In order to be able to correctly evaluate the accuracy of the temperature measurements provided by the probe, you should be able to meet the following requirements:

- The probe and the reference thermometer should be both ventilated with the same stream of air. Any dust filter cartridge used to protect the sensors should be removed from the probe. If the probe has a protective slotted cap, this should be left on the probe.
- Air velocity at the sensor should be within the limits of 200 to 500 ft/minute (1 to 2.5 m/sec). Any comparison between two instruments at a velocity under 200 ft/minute may not be valid. Air velocity above 500 ft/minute may damage the unprotected humidity sensor.
- The temperature of the air stream should be practically constant.

39.2 *Humidity*

- Temperature stability is the single most important requirement. If possible, adjust the probe at room temperature (18 to 25°C). Room temperature should be stable to $\pm 0.25^{\circ}\text{C}$ or better during the period of time required for each adjustment. Do not adjust close to an air vent or a heater, in direct exposure to sun rays, etc.

If using a humidity generator to calibrate the probe, make sure that the probe is as fully immersed in the generator as possible to minimize temperature effects.

ROTRONIC provides easy-to-use, certified humidity standards (RHS) for users who do not have access to a humidity generator. To use these standards, you will need a calibration device that is suitable for the probe.

39.2.1. Calibration device

A calibration device is a small airtight container that fits on the probe and seals around the humidity sensor. During calibration, known reference humidity is produced inside the calibration device by means of a humidity standard (usually an aqueous salt solution). To find out which model of calibration device is suitable for a specific probe, please consult the probe documentation.



39.2.2. ROTRONIC certified humidity standards

The ROTRONIC certified standards are available in boxes of 5 glass ampoules of the same value, which can be stored indefinitely. Standards in the range of 5 to 95 %RH are non-saturated aqueous salt solutions that are precisely titrated at our factory for the right concentration. The 0 %RH humidity standard is made of small granules of a highly porous ceramic that have been dried at a high temperature. A Material Safety Data Sheet is available for each standard. Since most standards are a salt solution, parts which have come in contact with the liquid should be cleaned after each use.

Each box of standards comes with a certificate that provides statistical information on the manufacturing batch of the standard. Information on the effect of temperature on each standard is provided on the cover of each box of standard. When calibrating with the HW4 software (RHS box checked), the effect of temperature on the standards is compensated by the software and no further correction is required. The value of the standards is not affected by barometric pressure (site elevation).



Using the ROTRONIC humidity standards

- Install the calibration device on the probe so that the receptacle (or solution holder) is under the probe. Check for a tight fit and remove the receptacle from the calibration device. Place one fiber disc (each box of standards includes 5 discs) in the receptacle of the calibration device. The purpose of this disc is to prevent accidental spilling of the solution inside the calibration device or on the humidity sensor.
- Tap the top of the ampoule so that all liquid drops to the bottom of the ampoule. Snap off top and empty contents on fiber disc. Since the ampoule is made of glass, exercise proper caution (gloves, safety glasses) when snapping off the top.
- Put the receptacle back on the calibration device and make sure that the solution does not come in contact with the sensor: The solution inside the calibration device should never be on top of the probe sensors.
- Allow at least 60 minutes insuring that the calibration device, the solution and the sensor are in a state of equilibrium. This is verified by monitoring the values being displayed.
- After adjusting the probe, remove the receptacle from the calibration device. Throw away the wet disc (non reusable). Thoroughly wash and dry the receptacle, removing all traces of the humidity standard.

40 BASIC TCP/IP CONCEPTS

40.1 *Compatibility requirements*

In order to use a device / docking station on an LAN, the device (also called a host) must to be configured so as to be compatible with the LAN.

- **IP address** (Network ID and Host ID)

An IP address consists of 4 numbers (1 to 3 digits each) separated with dots. On any network, each IP address should be unique.

IP addresses are classified as follows:

Class A: **1**.x.x.x to **127**.x.x.x

Class B: **128.0**.x.x to **191.255**.x.x

Class C: 192.0.0.x to 223.255.255.x

In the above, x is a number between 0 and 255.

Within an IP address, the network ID is defined depending on the class of the IP address and is indicated above in bold. The remainder of the address is the host ID.

Typically a local area network (LAN) uses class C IP addresses. In a class C address, the first 3 numbers define the network ID and therefore, the maximum number of devices on the network is limited to 254 (host IDs of all zeros or all ones are not allowed).

In order for two devices in a network to communicate together, their IP address must contain the same network ID.

- **Gateway Address**

The gateway address used by the host should be the LAN IP address of the router used to interface with other networks (e.g. Internet) or sub-networks. When there is no router, you may try any unused IP address on the network.

- **Subnet Mask**

The subnet mask is used to differentiate the network ID portion of an IP address from the host ID.

Class A: **255.0.0.0**

Class B: **255.255.0.0**

Class C: **255.255.255.0**

Note: if more network ID's are required, the last number of the subnet mask can be used but this reduces the maximum number of hosts available for each network ID.

In order for two devices in the same network to communicate together, their subnet mask must be the same.

- **TCP port**

TCP ports are used to specify how requests are sent to a device and how a device listens to requests. HW4 makes exclusive use of port 2101 to communicate with devices on the LAN. If using a ROTRONIC Ethernet docking station, this port number is already configured and should not be changed.

40.2 DHCP

A TCP device / docking station may be configured to enable DHCP (Dynamic Host Configuration Protocol). This makes sense only if a DHCP server such as a router is present on the network.

With DHCP enabled, the DHCP server senses when a new device (or host) is connected to the network and automatically assigns / configures the device with the next available IP address. The device is also automatically configured with the network subnet mask and gateway address. Using DHCP is simple, guarantees compatibility with the network and eliminates the possibility of conflicting IP addresses.

An IP address assigned using DHCP is also known as a dynamic address. As long as DHCP remains enabled on the host, the host IP address is subject to change. This is not desirable when using HW4.

Once the dynamic address of the host is known, and communication with the host can be established, we recommend to change the host IP address to an unused network address and to disable DHCP on the host. In this manner, the new host IP address is a static address, not subject to change.

40.3 MAC Address

Each individual Ethernet device is given by its manufacturer a unique MAC address that cannot be changed by the user. In a way the MAC address is akin to a social security number.

41 WATER ACTIVITY MEASUREMENT WITH HW4

41.1 *Water Activity: definition and applications*

41.1.1. Definitions

The moisture content of a product can be defined as the percentage weight of water in relation to the dry weight of the product.

Products in which moisture can be present can be classified in two categories: hygroscopic and non hygroscopic. Examples of hygroscopic materials are salts, vegetal fibers, most metal oxides, many polymers, etc. Examples of non hygroscopic products are metal powders, glass granules, etc.

Regarding the moisture content of a product, we define static equilibrium as a set of conditions under which the product does not exchange any moisture with its environment. Under conditions of static equilibrium, the moisture content of a hygroscopic product depends on the nature of the product and also on the two following factors:

- (a) the partial pressure of water vapor in the immediate environment of the product
- (b) the temperature of the product

If the moisture content of a product is not dependent on both these factors, the product is not hygroscopic.

Hygroscopic products may absorb water in different ways: sorption with formation of a hydrate, binding by surface energy, diffusion of water molecules in the material structure, capillary condensation, formation of a solution, etc. Depending on the absorption process, water is bound to the product with more or less strength. Moisture content can include both an immobilized part (e.g. water of hydration) and an active part.

Water activity A_w (or equilibrium relative humidity %ERH) measures the vapor pressure generated by the moisture present in a hygroscopic product.

$A_w = p / p_s$ and %ERH = 100 x A_w , where:

p : partial pressure of water vapor at the surface of the product

p_s : saturation pressure, or the partial pressure of water vapor above pure water at the product temperature

Water activity reflects the active part of moisture content or the part which, under normal circumstances, can be exchanged between the product and its environment.

Water activity is usually defined under static conditions of equilibrium. Under such conditions, the partial pressure of water vapor (p) at the surface of the product is equal to the partial pressure of water vapor in the immediate environment of the product. Any exchange of moisture between the product and its environment is driven by a difference between these two partial pressures.

Finally, water vapor can also be present in a gas or gas mixture. The relative humidity of a gas is defined as %RH = 100 x p/p_s , where (p) is the partial pressure of the water vapor present in the gas mixture and (p_s) is the saturation pressure, or the partial pressure of water vapor above pure water at the temperature of the gas.

41.1.2. A_w and Temperature

Both water activity (materials) and relative humidity (gases) are referred to the saturation pressure (ps) or partial pressure of water vapor above pure water:

$$A_w = p / p_s$$

$$\%RH = 100 \times p/p_s$$

The saturation pressure (p_s) is strongly dependent on temperature. At normal room temperature, (p_s) increases by about 6.2% for a 1°C increase in temperature. In an open environment that is not saturated with water vapor, the partial pressure of water vapor (p) does not change with temperature. In a closed environment, (p) changes proportionally to the °K temperature (°K temperature = °C temperature + 273.16). At normal room temperature, the change in (p) caused by a small change in the °C temperature is practically negligible. Because (p) does not change with temperature while (p_s) does, the relative humidity of a gas (%RH = 100 x p/p_s) is strongly temperature dependent. At 95 %RH and room temperature, an increase of 1°C results in a relative humidity decrease of about 6 %RH. At 50%RH, the same temperature increase causes relative humidity to decrease by about 3 %RH.

The water activity of most hygroscopic products is not as strongly dependent on temperature. At room conditions, research data typically shows that water activity varies only by roughly 0.0005 to 0.005 A_w (0.05 to 0.5 %RH) when temperature changes by 1°C.

This is explained by the fact that the partial pressure (p) at the surface of a hygroscopic product does vary with temperature. Above most hygroscopic products, the magnitude of the change in the partial pressure of water vapor (p) with temperature is similar (but not exactly equal) to the magnitude of the change of the saturation pressure (p_s) above pure water.

In summary, a change in temperature causes the partial pressure of water vapor above a hygroscopic product to change. At the same time, the partial pressure in the air above the product is practically unchanged. It follows that any change in the temperature of a hygroscopic product automatically causes the product to exchange moisture with the air (or gas) that surrounds it. Moisture is exchanged until the partial water vapor pressure at the surface of the product and in the air is equal. When measuring water activity, it is essential to keep temperature as constant as possible.

41.1.3. Applications

The active part of moisture content and, therefore, water activity, provide better information than the total moisture content regarding the micro-biological, chemical and enzymatic stability of perishable products such as foods and seeds. For similar reasons, water activity is equally relevant in the pharmaceutical industry where it provides useful information regarding the cohesion of tablets and pills, or the adherence of coatings. Water activity can be directly compared with the relative humidity of the ambient air to prevent dimensional changes in a product (paper, photographic film), to prevent hygroscopic powders (powdered sugar, salt) from caking or turning into a solid block, etc.

Water activity can be used with some products (mostly synthetic products) as a means of indirectly measuring the total moisture content. This requires developing sorption isotherms to this purpose. Sorption isotherms are graphs that provide the relationship between water activity and moisture content at constant temperature. For most natural products, repeatable sorption isotherms cannot be reliably developed and water activity should be regarded as separate from moisture content.

41.2 Instruments and probes for measuring water activity

In theory, any ROTRONIC instrument and probe from can be used together with HW4 to measure water activity. Usually product samples are being measured, typically with a probe such as the ROTONIC model AW-DIO. This probe can be connected without adapter cable to the HygroPalm and HygroLab indicators, the HygroLog NT data logger and to the HygroClip DI Interface.

41.3 Water activity measurement modes in HW4

The AwE / AwQuick tab is available only with HW4 Professional with AwQuick. HW4 features two modes for measuring water activity and the name of the tab depends on which mode has been selected.

• AwE mode

In this mode HW4 monitors the stability of both temperature and humidity. The measurement is automatically ended as soon as both humidity and temperature reach equilibrium. The natural (or static) equilibration of most products typically requires from 45 to 60 minutes and can take as long as a couple of hours.

• AwQuick mode:

In the AwQuick mode, HW4 uses an algorithm to project the full equilibrium value (water activity) of the measured product. The measurement is automatically ended and typically requires about 5 minutes.

HW4 performs the following tasks:

- 1) The value of the humidity signal is constantly monitored
- 2) The stability of the temperature signal is constantly monitored
- 3) After an initial period of time (dwell time), HW4 uses the humidity data to project the end value of the equilibration process (water activity). The measurement ends automatically as soon as the projected Aw value is stable.

With the dwell time set to 4 minutes, measurements typically require about 5 minutes. When temperature conditions are stable (both at the product and probe), the measurement obtained with the AwQuick mode is generally within ± 0.005 aw of the measurement that would be obtained by waiting for full equilibration (AwE mode).

41.4 Mode selection and settings

To select the mode and to enter or modify the settings to be used in each mode , click on the Settings button located on the AwE / AwQuick tab.

The screenshot shows the 'HW4 Global Settings Form' with the 'AwQuick Mode' tab selected. The form contains two sections: 'AwQuick Mode' and 'AwE Mode'. The 'AwQuick Mode' section has a radio button and two input fields: 'Dwell Time [min]' with a value of 4, and 'Temperature stability [°C/min]' with a value of 0.01. The 'AwE Mode' section has a radio button and two input fields: 'Humidity stability [Aw/min]' with a value of 0.0001, and 'Temperature stability [°C/min]' with a value of 0.01. At the bottom right, there are three buttons: 'OK', 'Cancel', and 'Help'.

Mode	Parameter	Value
AwQuick Mode	Dwell Time [min]	4
	Temperature stability [°C/min]	0.01
AwE Mode	Humidity stability [Aw/min]	0.0001
	Temperature stability [°C/min]	0.01

Note: the values shown above are the most frequently used.

IMPORTANT: Both the selected mode and its settings apply globally to all instruments and probes.

41.5 Using the AwE mode

- 1) Select the instrument to be used in the device tree (left pane of the HW4 main screen)
- 2) In the right pane of the HW4 main screen, select the AwE / AwQuick tab (the label of this tab changes depending on which mode is currently selected)
- 3) If necessary, click on the settings button and select the AwE mode in the form. Click on the form OK button and verify that the tab label reads AwE Mode
- 4) To start the measurement, click on the input corresponding to the desired probe. A check mark appears on the corresponding box and the Status column changes to "Running". If so desired, you can click on each of the inputs available with the instrument selected in the device tree.

Note: only those inputs that are selected in the View tab appear in the AwE / AwQuick tab.

Current Values Online Buffer Log to PC View AwE Mode							
Input	Input Name	Probe S/N	Humidity [Aw]	Temperature [°C]	Elapsed Time	Water Activity [Aw]	Status
☒ 1	Port 1	33220033	0.176 <=>	24.13 <=>	00:00:05		Running
☐ 4	Port 2	36067221	0.177	24.13			
☐ 5	Port 3	32221185	0.178	23.99			
☐ 6	Port 4	32337234	0.178	23.97			
☐ 7	Port 5	39678190	0.176	24.13			
Settings Generate report							

When HW4 detects equilibrium condition, the measurement ends automatically and the background color of the line corresponding to the probe turns to green. The status column changes to "End Result".

Current Values Online Buffer Log to PC View AwE Mode							
Input	Input Name	Probe S/N	Humidity [Aw]	Temperature [°C]	Elapsed Time	Water Activity [Aw]	Status
☒ 1	Port 1	33220033	0.177 <=>	24.08 <=>	00:06:55	0.177	End Result
☐ 4	Port 2	36067221	0.177	24.06			
☐ 5	Port 3	32221185	0.178	23.95			
☐ 6	Port 4	32337234	0.178	23.96			
☐ 7	Port 5	39678190	0.176	24.05			
Settings Generate report							

All data on the line with the green background remains frozen until the process is ended by clicking on the input box with the mouse to remove the check mark.

41.6 Using the AwQuick mode

- 1) Select the instrument to be used in the device tree (left pane of the HW4 main screen)
- 2) In the right pane of the HW4 main screen, select the AwE / AwQuick tab (the label of this tab changes depending on which mode is currently selected)
- 3) If necessary, click on the settings button and select the AwE mode in the form. Click on the form OK button and verify that the tab label reads AwE Mode
- 4) To start the measurement, click on the input corresponding to the desired probe. A check mark appears on the corresponding box and the Status column changes to "Dwell Phase" and to "Running" at the end of the dwell time. If so desired, you can click on each of the inputs available with the instrument selected in the device tree.

Note: only those inputs that are selected in the View tab appear in the AwE / AwQuick tab.

Current Values Online Buffer Log to PC View AwQuick Mode							
Input	Input Name	Probe S/N	Humidity [Aw]	Temperature [°C]	Elapsed Time	Water Activity [Aw]	Status
☒ 1	Port 1	33220033	0.298	23.12	00:00:10		Dwell Phase
☒ 4	Port 2	36067221	0.299	23.08	00:00:10		Dwell Phase
☐ 5	Port 3	32221185	0.301	22.99			
☐ 6	Port 4	32337234	0.299	22.95			
☐ 7	Port 5	39678190	0.299	23.10			
Settings Generate report							


When HW4 generates a stable projection, the measurement ends automatically and the background color of the line corresponding to the probe input turns to green. The status column changes to "End Result".

Current Values Online Buffer Log to PC View AwQuick Mode							
Input	Input Name	Probe S/N	Humidity [Aw]	Temperature [°C]	Elapsed Time	Water Activity [Aw]	Status
☒ 1	Port 1	33220033	0.298 <=>	23.09 <=>	00:05:10	0.297	End Result
☒ 4	Port 2	36067221	0.299 <=>	23.10 <=>	00:04:20	0.299	End Result
☐ 5	Port 3	32221185	0.301	22.98			
☐ 6	Port 4	32337234	0.299	22.92			
☐ 7	Port 5	39678190	0.299	23.10			
Settings Generate report							


The data on the line with the green background remains frozen until the process is ended by clicking on the input box with the mouse to remove the check mark.

41.7 Water activity measurement report

Click on the “Generate Report” button to generate a measurement report. HW4 automatically saves the report in the subfolder ROTRONIC_HW4\DOC.as a text file with the name AW-Quick followed by an automatically generated number.

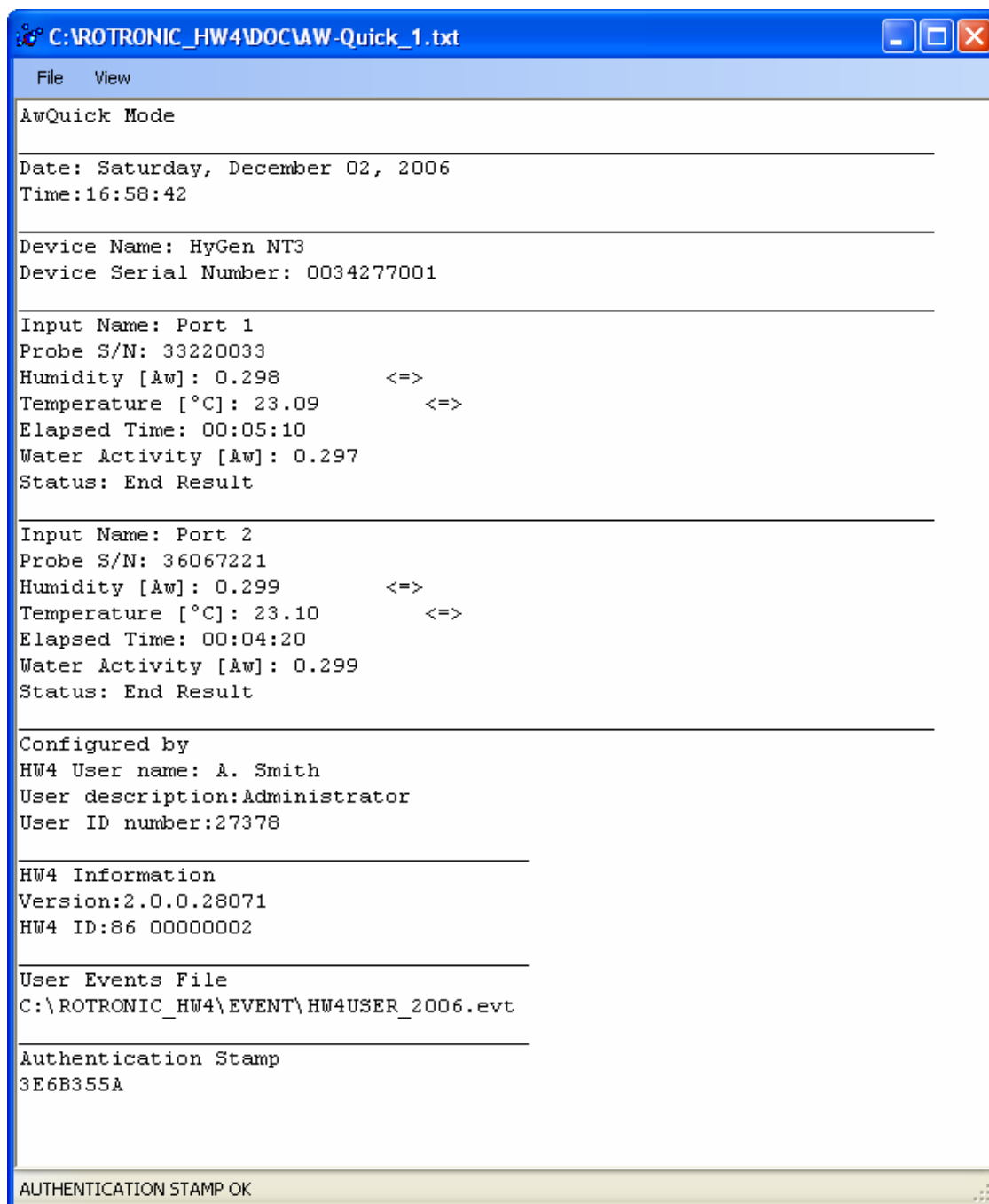

Use the report menu bar to sign and / or print the report.

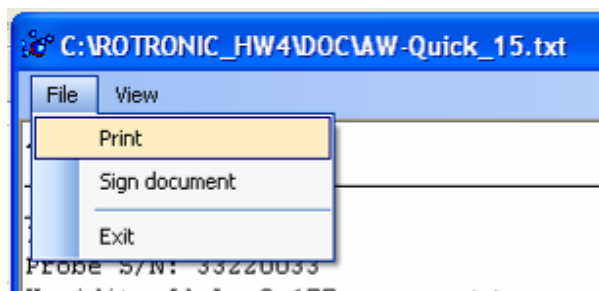

41.7.1. Signing the report

Select “Sign document” from the file menu to add your HW4 user name and a comment text to the report. For authentication purposes you will need to enter your password.

A screenshot of a dialog box titled 'Sign the document'. The dialog has a blue title bar with a close button. Inside, there are several labeled fields: 'Document name' with the value 'AW-Quick_1.txt', 'Software' with 'HW4 V2.0.0.28071', and 'Date' with 'Saturday, December 02, 2006'. Below these are two input fields: 'HW4 user name' containing 'A. Smith' and 'Password' containing a series of asterisks. At the bottom, there is a 'Text' label above a large text area containing the text 'This is a test of the AwQuick mode'. At the very bottom of the dialog are two buttons: 'Sign' and 'Cancel'.

A new section is inserted at the bottom of the report:

