

Keysight Technologies CX3300 Series Device Current Waveform Analyzer

CX3322A 2-Channel Model
CX3324A 4-Channel Model

Quick Operation
Guide

Notices

Copyright Notice

© Keysight Technologies 2018

No part of this manual may be reproduced in any form or by any means (including electronic storage and retrieval or translation into a foreign language) without prior agreement and written consent from Keysight Technologies as governed by United States and international copyright laws.

Manual Part Number

CX3300-90300

Edition

Edition 1, June 2018

Published by:

Keysight Technologies International
Japan G.K.
9-1, Takakura-cho, Hachioji-shi, Tokyo
192-0033 Japan

Technology Licenses

The hardware and/or software described in this document are furnished under a license and may be used or copied only in accordance with the terms of such license.

Declaration of Conformity

Declarations of Conformity for this product and for other Keysight products may be downloaded from the Web. Go to <http://www.keysight.com/go/conformity> and click on "Declarations of Conformity." You can then search by product number to find the latest Declaration of Conformity.

U.S. Government Rights

The Software is "commercial computer software," as defined by Federal Acquisition Regulation ("FAR") 2.101. Pursuant to FAR 12.212 and 27.405-3 and Department of Defense FAR Supplement ("DFARS") 227.7202, the U.S. government acquires commercial computer software under the same terms by which the software is customarily provided to the public. Accordingly, Keysight provides the Software to U.S. government customers under its standard commercial license, which is embodied in its End User License Agreement (EULA), a copy of which can be found at <http://www.keysight.com/find/sweula>. The license set forth in the EULA represents the exclusive authority by which the U.S. government may use, modify, distribute, or disclose the Software. The EULA and the license set forth therein, does not require or permit, among other things, that Keysight: (1) Furnish technical information related to commercial computer software or commercial computer software documentation that is not customarily provided to the public; or (2) Relinquish to, or otherwise provide, the government rights in excess of these rights customarily provided to the public to use, modify, reproduce, release, perform, display, or disclose commercial computer software or commercial computer software documentation. No additional government requirements beyond those set forth in the EULA shall apply, except to the extent that those terms, rights, or licenses are explicitly required from all providers of commercial computer software pursuant to the FAR and the DFARS and are set forth specifically in writing elsewhere in the EULA. Keysight shall be under no obligation to update, revise or otherwise modify the Software. With respect to any technical data as defined by FAR 2.101, pursuant to FAR 12.211 and 27.404.2 and DFARS 227.7102, the U.S. government acquires no greater than Limited Rights as defined in FAR 27.401 or DFARS 227.7103-5 (c), as applicable in any technical data.

Warranty

THE MATERIAL CONTAINED IN THIS DOCUMENT IS PROVIDED "AS IS," AND IS SUBJECT TO BEING CHANGED, WITHOUT NOTICE, IN FUTURE EDITIONS. FURTHER, TO THE MAXIMUM EXTENT PERMITTED BY APPLICABLE LAW, KEYSIGHT DISCLAIMS ALL WARRANTIES, EITHER EXPRESS OR IMPLIED, WITH REGARD TO THIS MANUAL AND ANY INFORMATION CONTAINED HEREIN, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. KEYSIGHT SHALL NOT BE LIABLE FOR ERRORS OR FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES IN CONNECTION WITH THE FURNISHING, USE, OR PERFORMANCE OF THIS DOCUMENT OR OF ANY INFORMATION CONTAINED HEREIN. SHOULD KEYSIGHT AND THE USER HAVE A SEPARATE WRITTEN AGREEMENT WITH WARRANTY TERMS COVERING THE MATERIAL IN THIS DOCUMENT THAT CONFLICT WITH THESE TERMS, THE WARRANTY TERMS IN THE SEPARATE AGREEMENT SHALL CONTROL.

Safety Information

CAUTION

A CAUTION notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in damage to the product or loss of important data. Do not proceed beyond a CAUTION notice until the indicated conditions are fully understood and met.

WARNING

A WARNING notice denotes a hazard. It calls attention to an operating procedure, practice, or the like that, if not correctly performed or adhered to, could result in personal injury or death. Do not proceed beyond a WARNING notice until the indicated conditions are fully understood and met.

Contents

1. File

- How to save the data
- How to load the data
- What types of data formats are supported on the CX3300 ?
- How to save rarely occurring waveforms automatically

2. Setup

- How to change the horizontal scale and position
- How to change the sampling rate
- How to change the memory depth
- How to change the 14-bit or 16-bit resolution
- How to change the vertical scale and position
- How to change the current range
- How to set the bandwidth limit
- How to minimize the offset current
- How to perform a Skew adjustment
- What determines the effective bandwidth of the instrument ?
- How to zoom in on the waveform
- How to use the waveform memory
- How to overlay multiple acquisitions
- How to display the secondary Y axis

3. Trigger

- How to set the trigger

4. Measure

- How to apply a Measurement Function
- How to apply the A-B marker
- How to place pointers on waveforms
- How to add an Annotation

5. Analyze

- How to apply the post-filter
- How to perform the FFT (Fast Fourier Transform)
- How to analyze the waveform in the specified area with an Area marker
- How to apply the Power and Current Profiler
- How to apply the XY analysis

6. Utilities


- How to change the measurement units (Clone to Ampere Hour) in Area analysis
- How to perform User calibration
- How to check the firmware revisions
- How to simulate the 50 MHz or 100 MHz bandwidth option on the mainframe with the 200 MHz option
- Where to find the Online help
- How to display the menu bar

Keysight Technologies CX3300 Current Waveform Analyzer


Quick Operation Guide

How to save the waveform data


- How to save the Waveform data in CSV format.

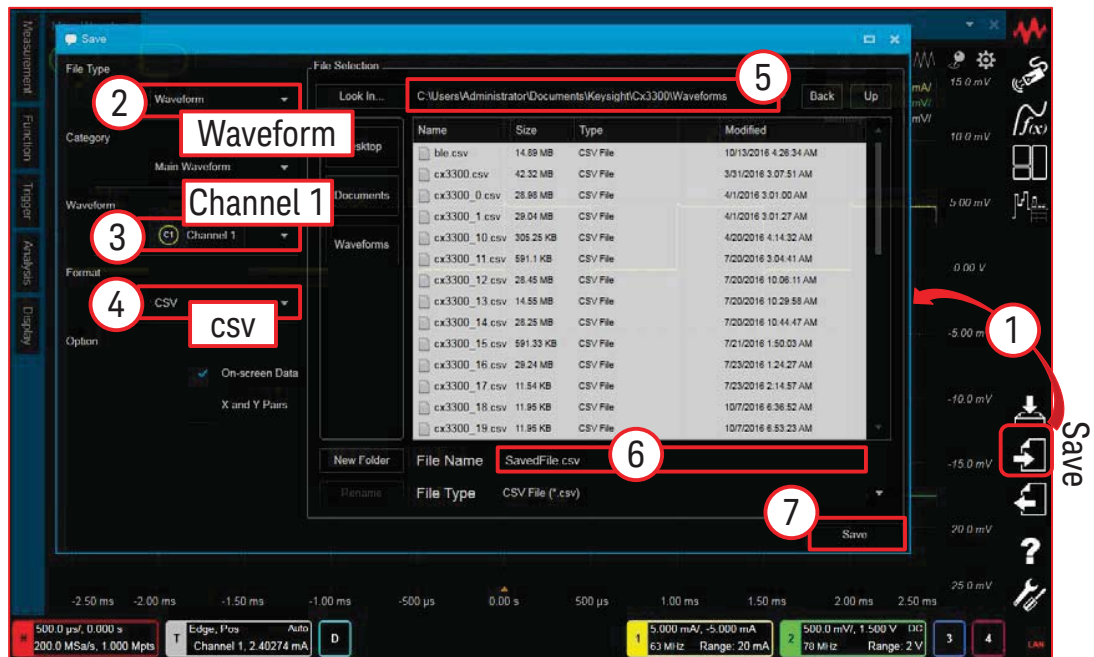


- How to save the Waveform FFT data in CSV format.

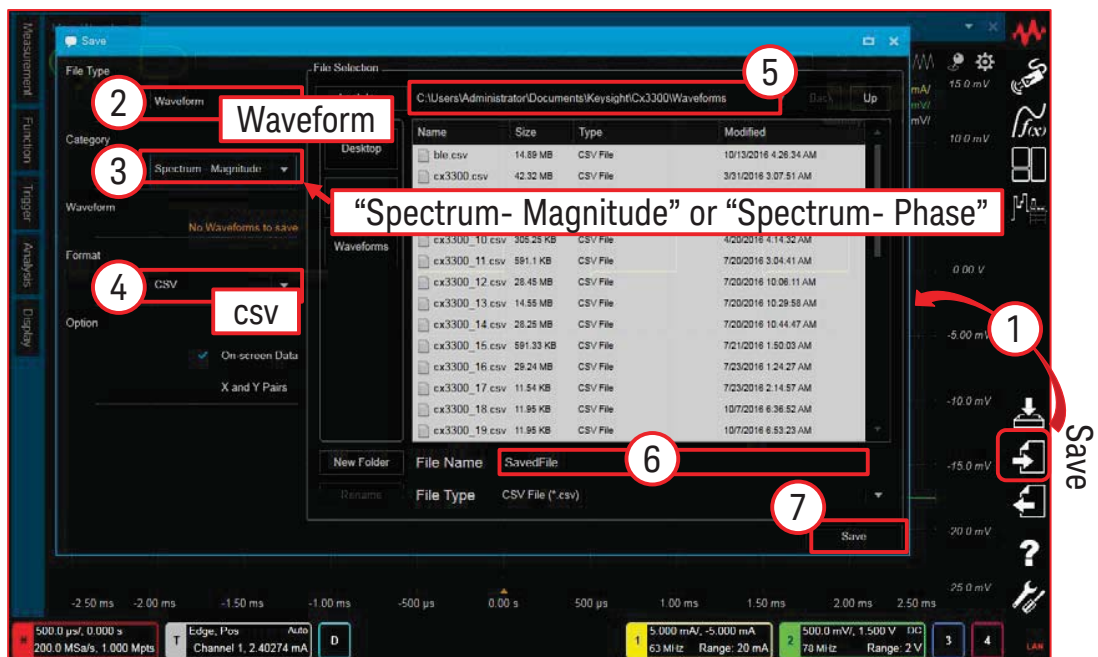

- How to save the Window Image data (Format : PNG / JPEG / Bitmap).

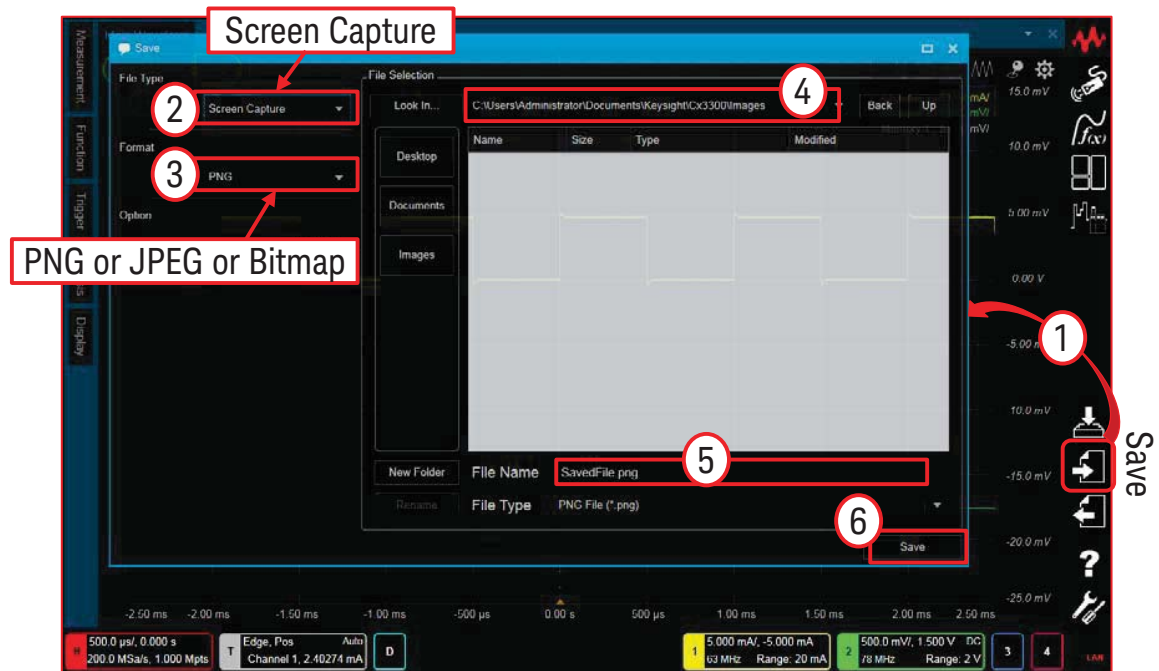

- How to save all setting data and waveform data (File Type : Composite).

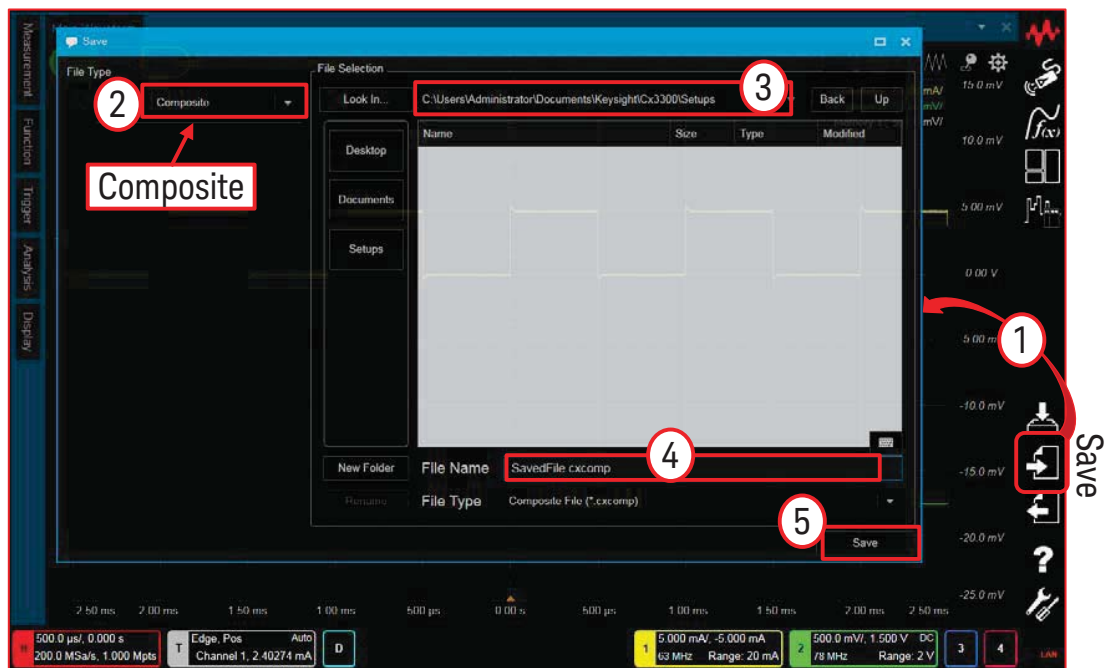


Note:

- For further details about the File Type and Format, please refer to “What types of data formats are supported on the CX3300 ?”.

Quick Operation Guide

How to load the data

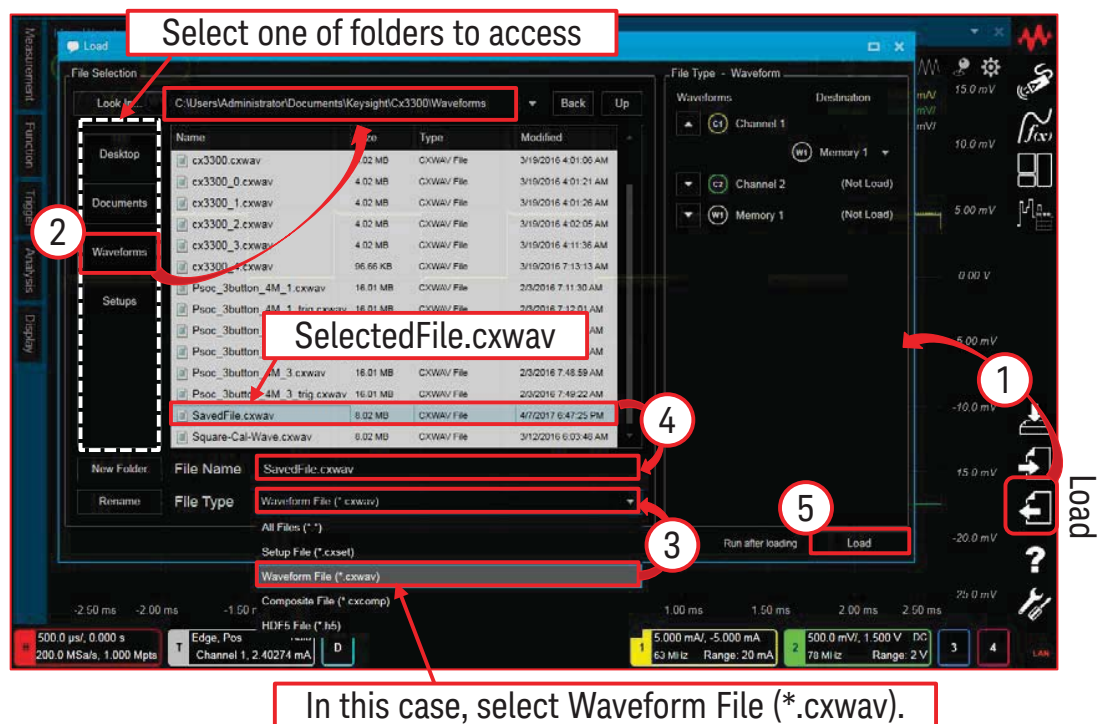

Note:

- For further details about the File Type and Format, please refer to “What types of data formats are supported on the CX3300?”.

Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

What types of data formats are supported on the CX3300 ?

The CX3300 supports several File Types and Formats when the setup data or the waveform data is saved or loaded.

[Detail of the File Type and Format]

File Type	Format	Save		Load	
		Setup	Data	Setup	Data
Waveform (Waveform data)	CSV	No	Yes	No	No
	TSV				No
	Waveform				Yes
	HDF5				Yes
Report (Outputs the measurement and analysis result report in the dedicated style)	XPS	Acquisition and channel setups on the summary bar	Measurement and analysis result	No	
	CSV				
	TSV				
	Text				
Composite (All setup data and waveform data)	-	Yes	Yes	Yes	Yes
Screen Capture (Image of the CX3300 full screen)	PNG	The images on the screen		No	
	JPEG				
	Bitmap				
Setup (All setup data or trigger setup data only)	-	Yes	No	Yes	No

Keysight Technologies


CX3300 Current Waveform Analyzer

Quick Operation Guide



How to save rarely occurring waveforms automatically

"Auto Save" can be used to automatically save a current waveform if it satisfies specific conditions.

1. Select "Utilities > Configuration > Auto Save".

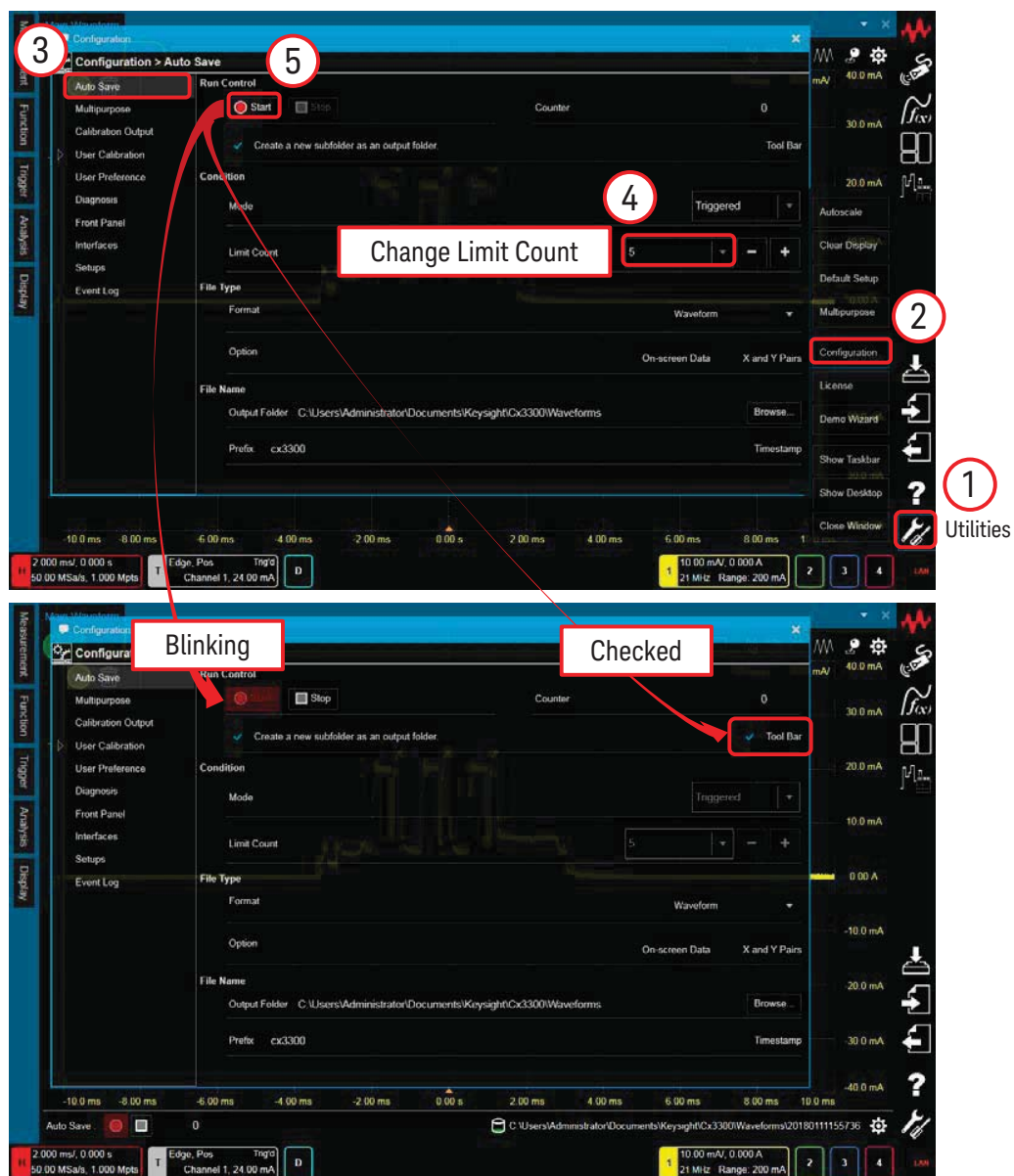

2. Push the “Auto Save” (red circle) button on the lower left of the display. The waveforms consistent with the trigger condition are stored up until the Limit Count set on the “Configuration > Auto Save” window.

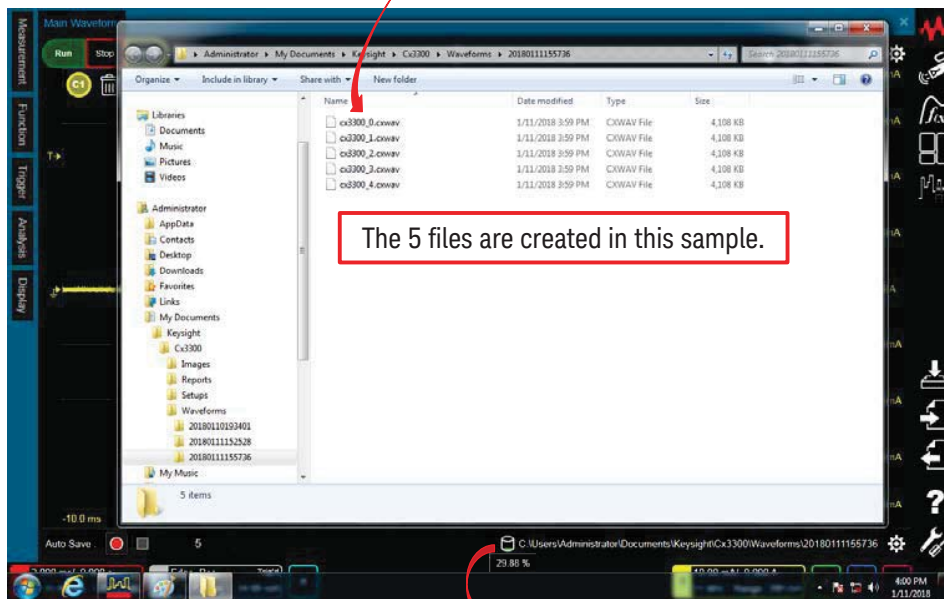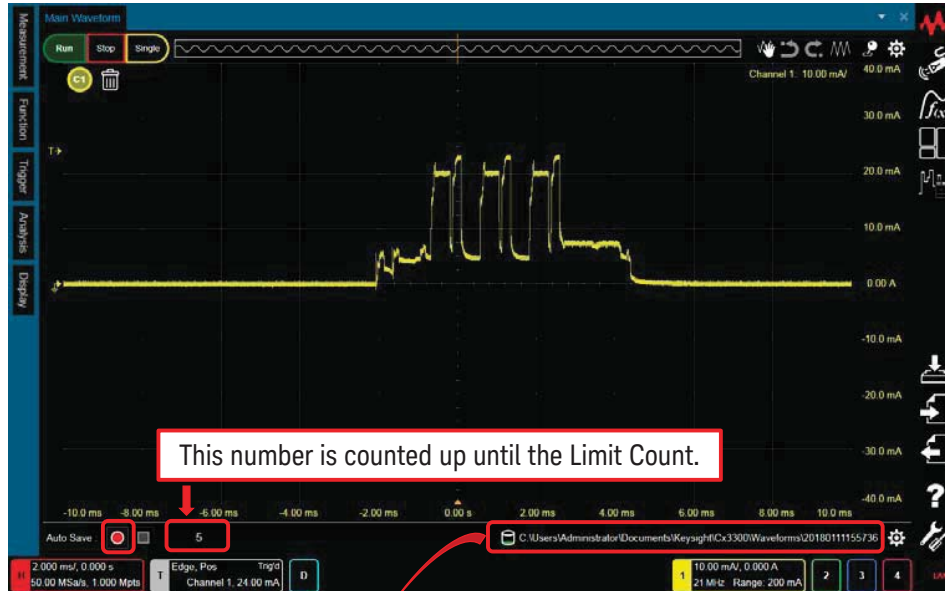

The icon shows the disk usage. When the free space becomes 10% or less, the color changes from green to red.

Note:

- During the processing for previously acquired waveforms, the data is not saved even if the specified condition is satisfied.
- Depending on the File Type and Limit Count, the size or number of the files to be saved increases.


Keysight Technologies CX3300 Current Waveform Analyzer

Quick Operation Guide

How to change the horizontal scale and position

You can change the horizontal scale and position using any of the following methods:

- Using the knobs on the Horizontal Controls
- Using the Horizontal mini dialog box
- Using the X Axis dialog box for fine control.

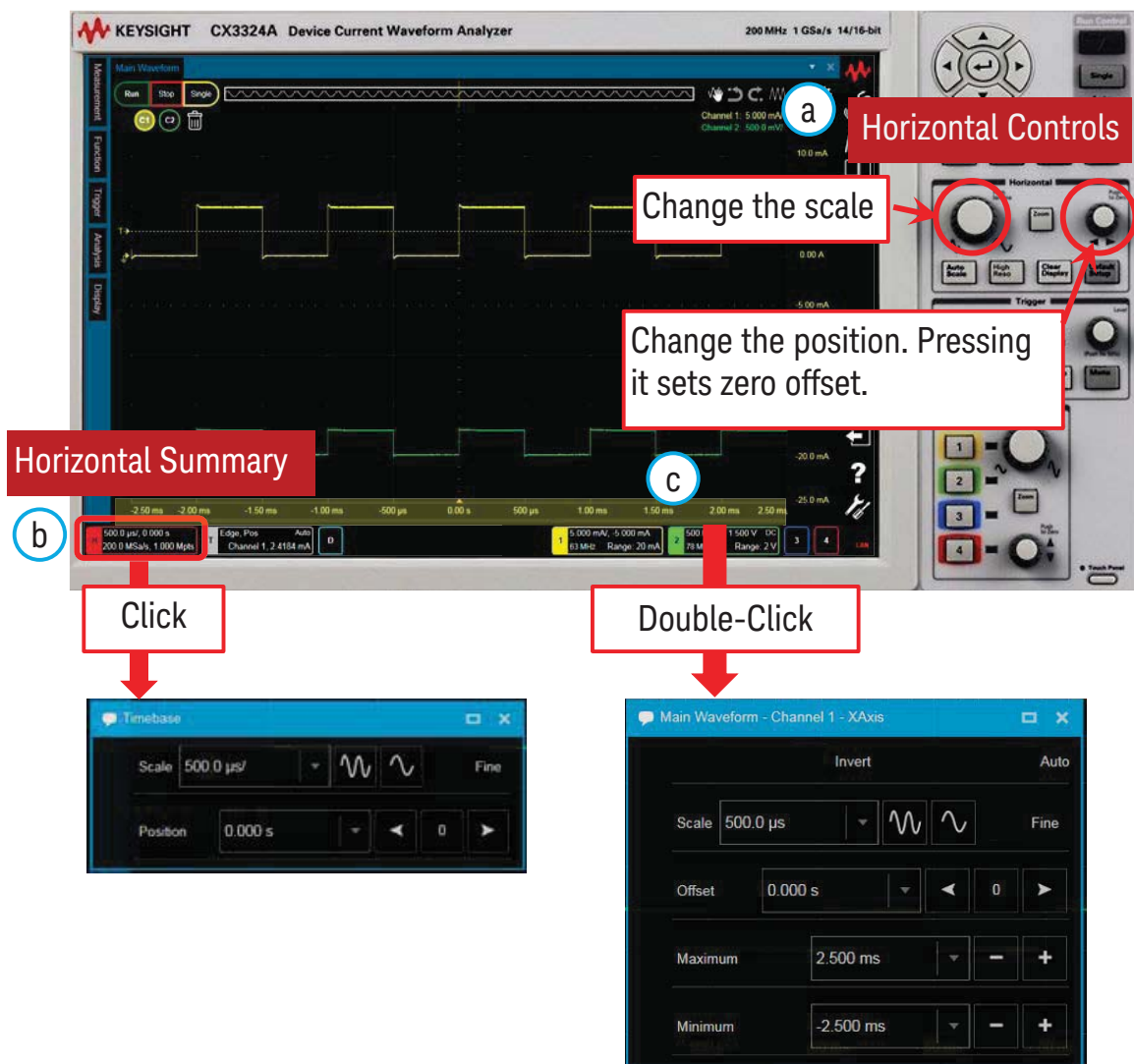

Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to change the sampling rate

1. Click on the Horizontal summary to open the Horizontal mini dialog box.
2. Un-check "Automatic".
3. Set the "Sampling Rate".



Note:

- If "Auto" is enabled, a suitable sampling rate is selected automatically by the horizontal scale setting controlled by the Horizontal scale control knob.


Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to change the memory depth

1. Click on “Horizontal summary” to open the “Timebase” mini dialog box.
2. Click on the “□” at the upper right of the mini dialog box to open the Setting dialog box.
3. Select “Acquisition”.
4. Un-check “Automatic”.
5. Set the “Memory Depth”.

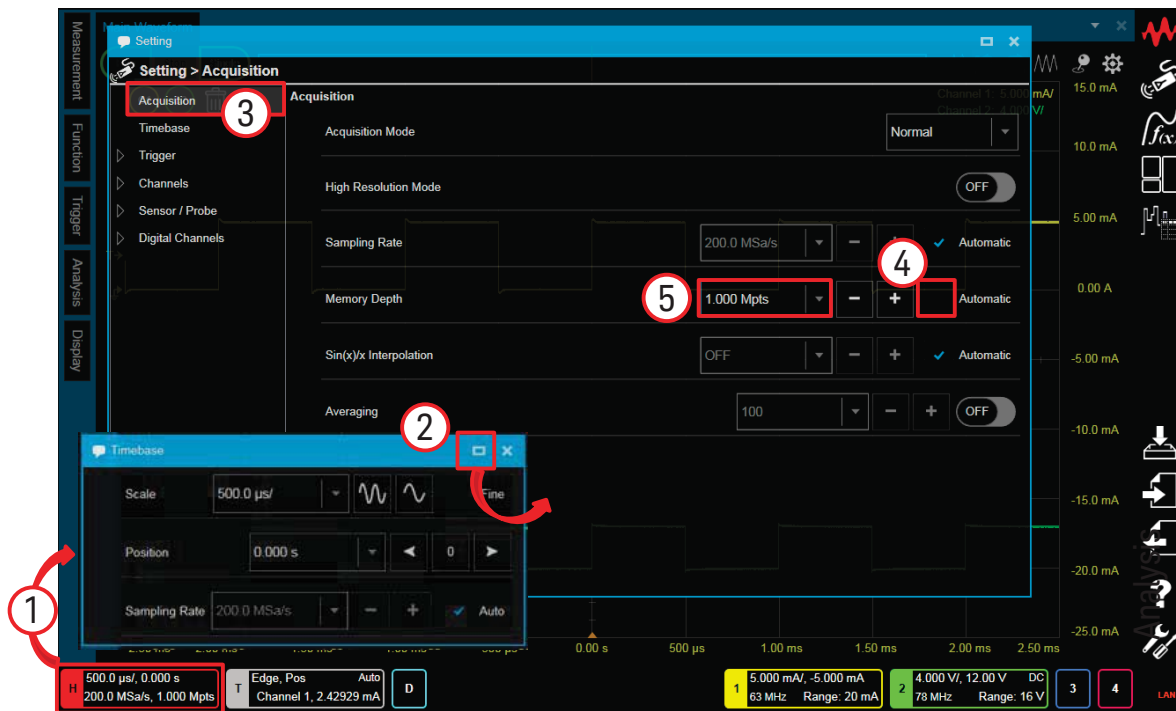


Horizontal Summary

Keysight Technologies CX3300 Current Waveform Analyzer

Quick Operation Guide

How to change the 14-bit or 16-bit resolution.

Change a 14-bit (high speed mode) or a 16-bit (high resolution mode) using the [High Reso] key.

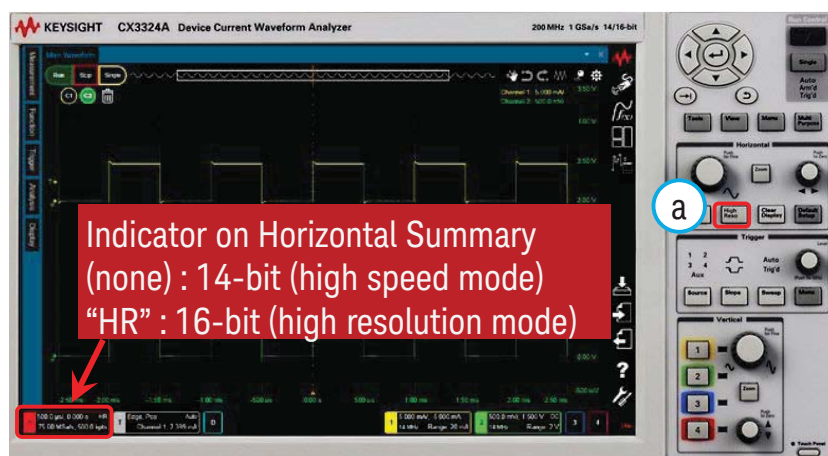

Note:

The maximum sampling rate and bandwidth is automatically changed by the mode.

Resolution	Max. Sampling	Max. Bandwidth
14-bit (high speed mode)	1 G Sa/s	200 MHz
16-bit (high resolution mode)	75 Msa/s	15 MHz

Keysight Technologies CX3300 Current Waveform Analyzer


Quick Operation Guide

How to change the vertical scale and position

1. How to select the channel.

Select the target waveform by either:

- Pressing [1], [2], [3] or [4] for the analog input channels
- OR
- Clicking on the icons on the Waveform tray for the analog input channels or function waveform.

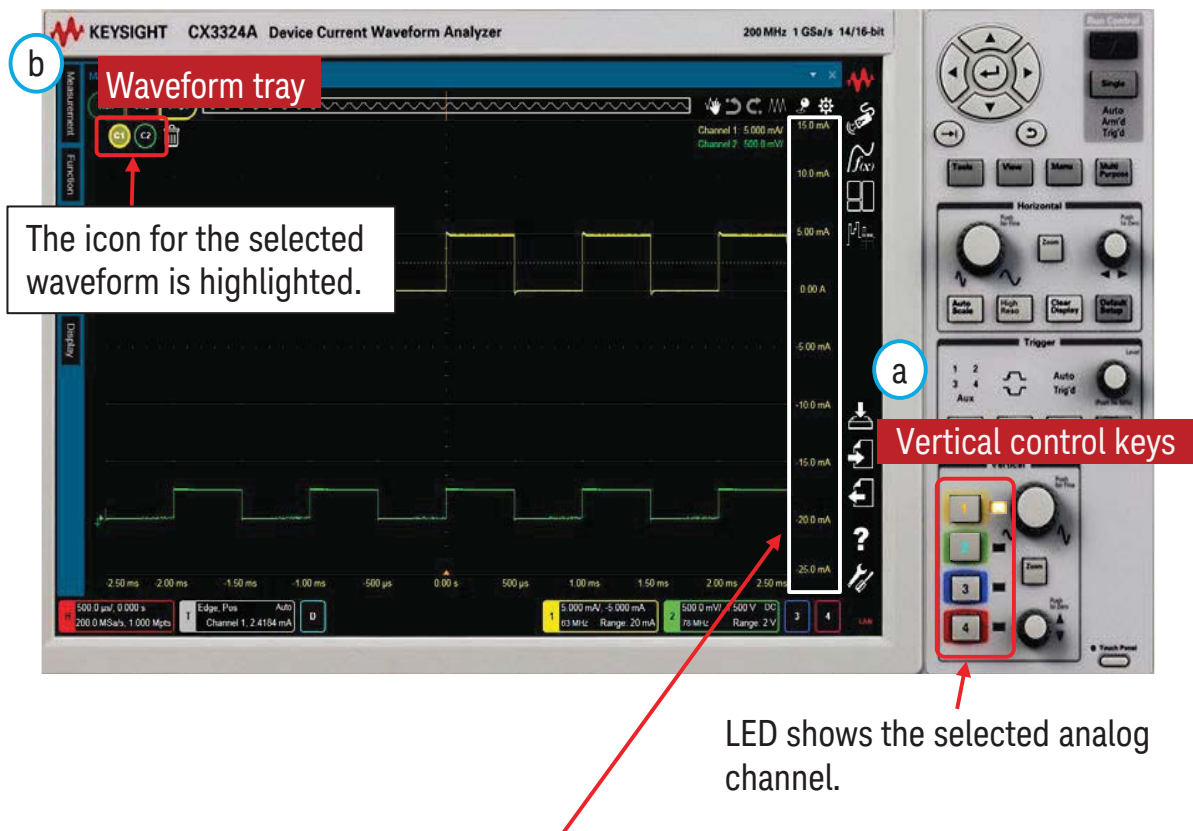


The color of the Y axis is the same color as the selected waveform.

2. How to Change the horizontal scale and position

You can change the horizontal scale and position using any of the following

methods:

- Using the knobs on Vertical Controls
- Using the Vertical mini dialog box
- Using the Y Axis dialog box for fine control.

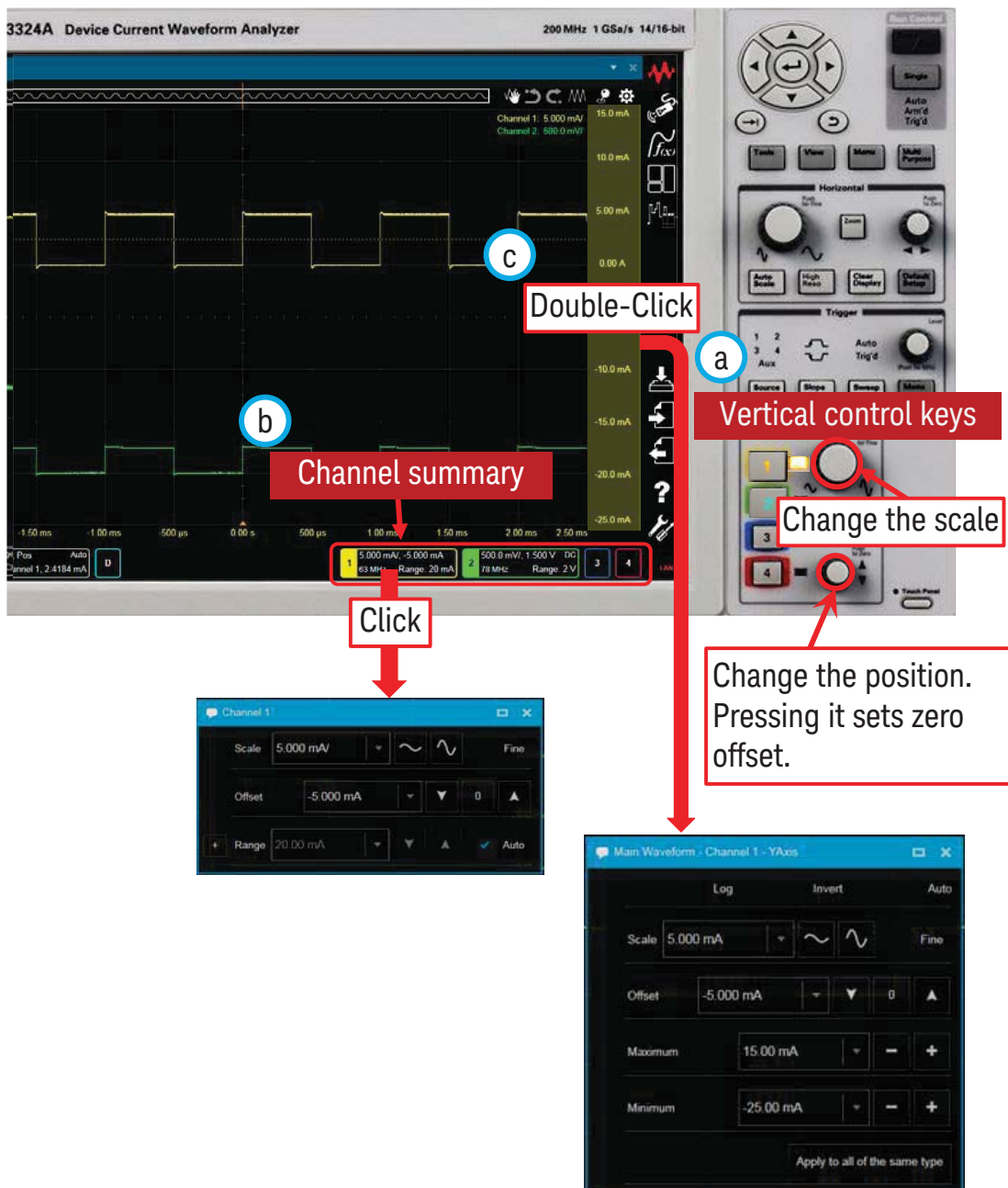

Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to change the current range

1. Click on the Channel summary to open the Channel mini dialog box.
2. Un-check "Auto".
3. Select the "Range".



Click



Note:

- If "Auto" is enabled, a suitable range is selected automatically by the vertical scale setting controlled by the Vertical scale control knob.
- The measurement bandwidth and the input resistor are automatically changed by changing the range.


Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to set the bandwidth limit

1. Click on the Channel summary to open the Channel mini dialog box.
2. Click on the "□" at the upper right of the window to open the Setting dialog box.
3. Select "ON".
4. Set the bandwidth limit.

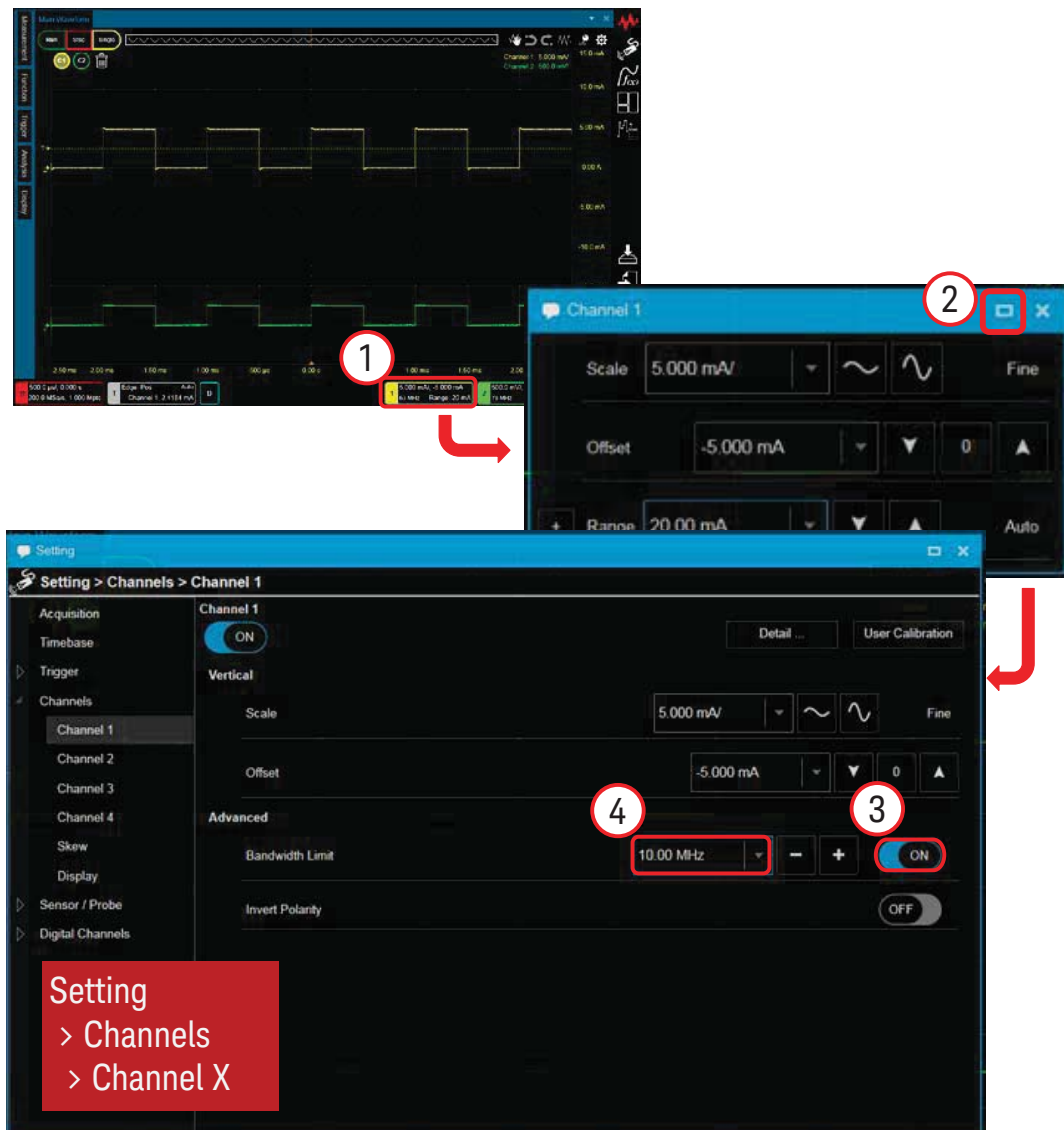


Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to minimize the offset current

1. Perform "User Calibration". The offset current of the current sensor will be minimized.
2. Perform a actual measurement and fix the connection and the CX3300 settings such as the current range.
3. Perform a "Null" measurement.

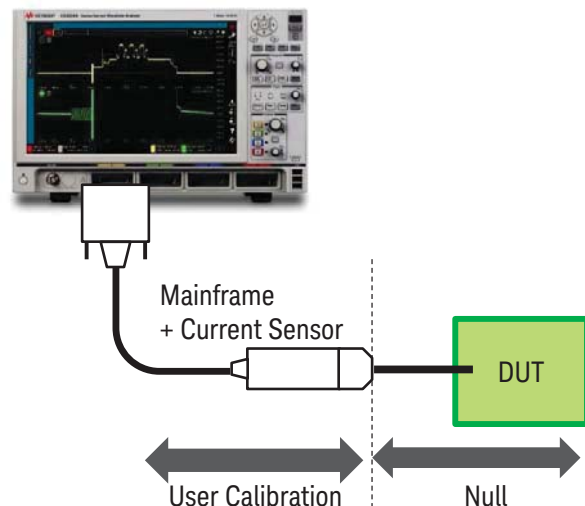


3-1. Make a connection to prevent the current flow into the current sensor.

Note for connection:


- Connect either the (+) or (-) terminal of the current sensor to the stable voltage on the DUT.
- Connect the ground lead to the circuit common of the DUT.

Connection Example

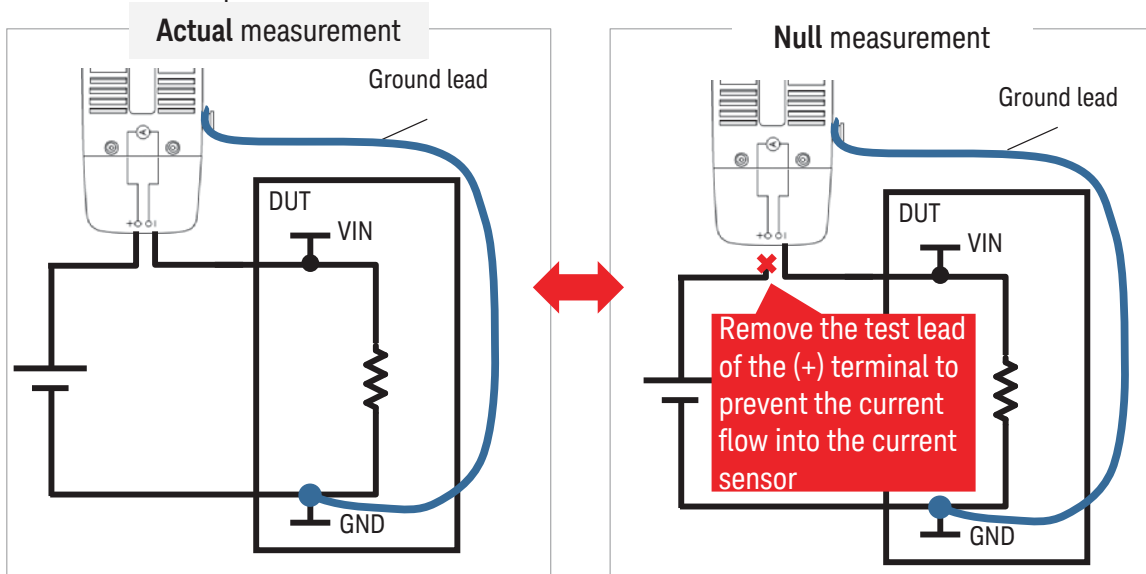

3-2. Null measurement

1. Click on the Channel summary to open the Channel mini dialog box.
2. Turn "ON" the Null measurement.
3. Press "Measure".

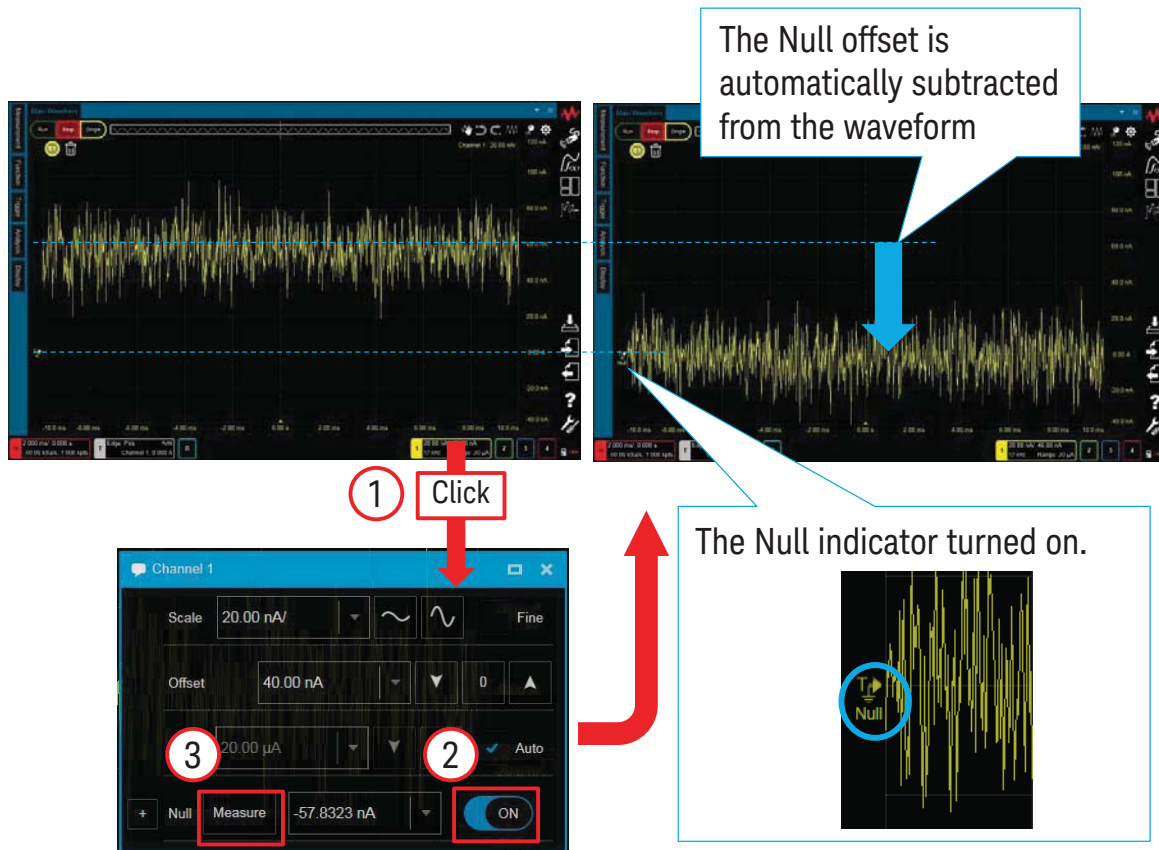

4. Restore the connection then make an actual measurement again.

Note:

The Null offset is set to zero when:

- The current sensor is disconnected
- The User Calibration is performed
- The CX3300 is closed or turned off.

Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to perform a Skew adjustment


The skew between the channels on the measurement device can be adjusted.

1. Preparation

- 1-1. Connect the Current Sensor (CX1101A/CX1102A/CX1103A) or the Passive Probe Interface Adapter (CX1151A) to the device.
- 1-2. Execute the measurement.

2. Settings

- 2-1. Click on the "Setting" icon.
- 2-2. The "Setting" dialog box opens.
- 2-3. Click on "Setting>Channels>Skew".
- 2-4. The Skew setup panel opens and the waveforms for the effective channels are shown behind the "Setting" dialog box.
- 2-5. Adjust the skew of the other channels to the triggered channel while observing the waveforms behind the "Setting" dialog box.

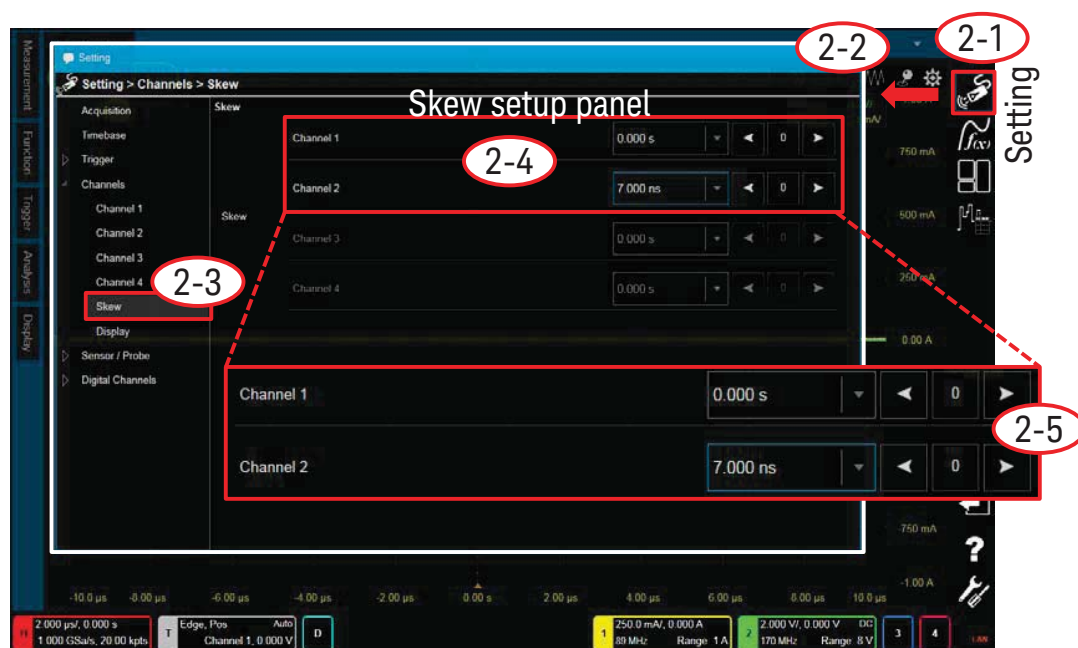


Keysight Technologies CX3300 Current Waveform Analyzer

Quick Operation Guide

What determines the effective bandwidth of the instrument?

The effective bandwidth of the instrument is determined by the following settings and displayed in each “Channel summary” on the front panel.

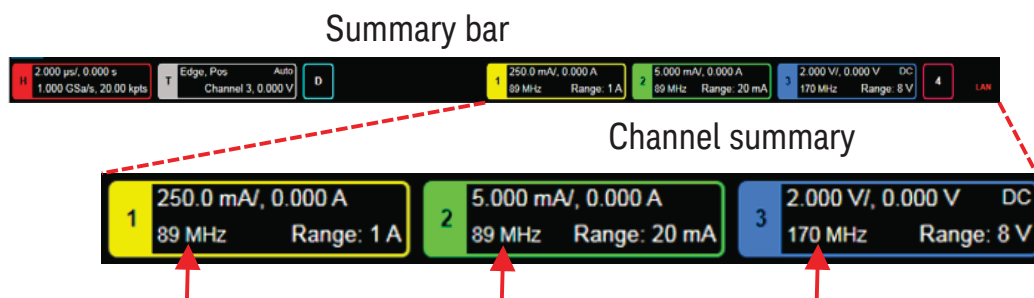

- Bandwidth options of the mainframe
- ADC (14-bit: High speed mode or 16-bit: High resolution mode)
For details of settings, please refer to “*How to change the 14-bit or 16-bit resolution?*”.
- Measurement range
For details of settings, please refer to “*How to change the current range?*”.
- Sampling rate
For details of settings, please refer to “*How to change the sampling rate?*”.
- Bandwidth Limit
For details of settings, please refer to “*How to set the bandwidth limit?*”.

Tips: In addition, the connection type and the DUT can limit the actual measurement bandwidth at DUT. Please refer to the application note, “7 Hints for Precise Current Measurements with the CX3300 Series Device Current Waveform Analyzer” (Literature number: 5992-2118) for more information.

Keysight Technologies CX3300 Current Waveform Analyzer

Quick Operation Guide


How to zoom in on the waveform

The CX3300 has two kinds of zoom in functions:

- a. Zoom at single display
- b. "Anywhere" Zoom.

a. Zoom at single display


- a-1. Click on the "Move" icon to change the 'Zoom' icon.
- a-2. Drag the cursor until the displayed box covers the area of interest.

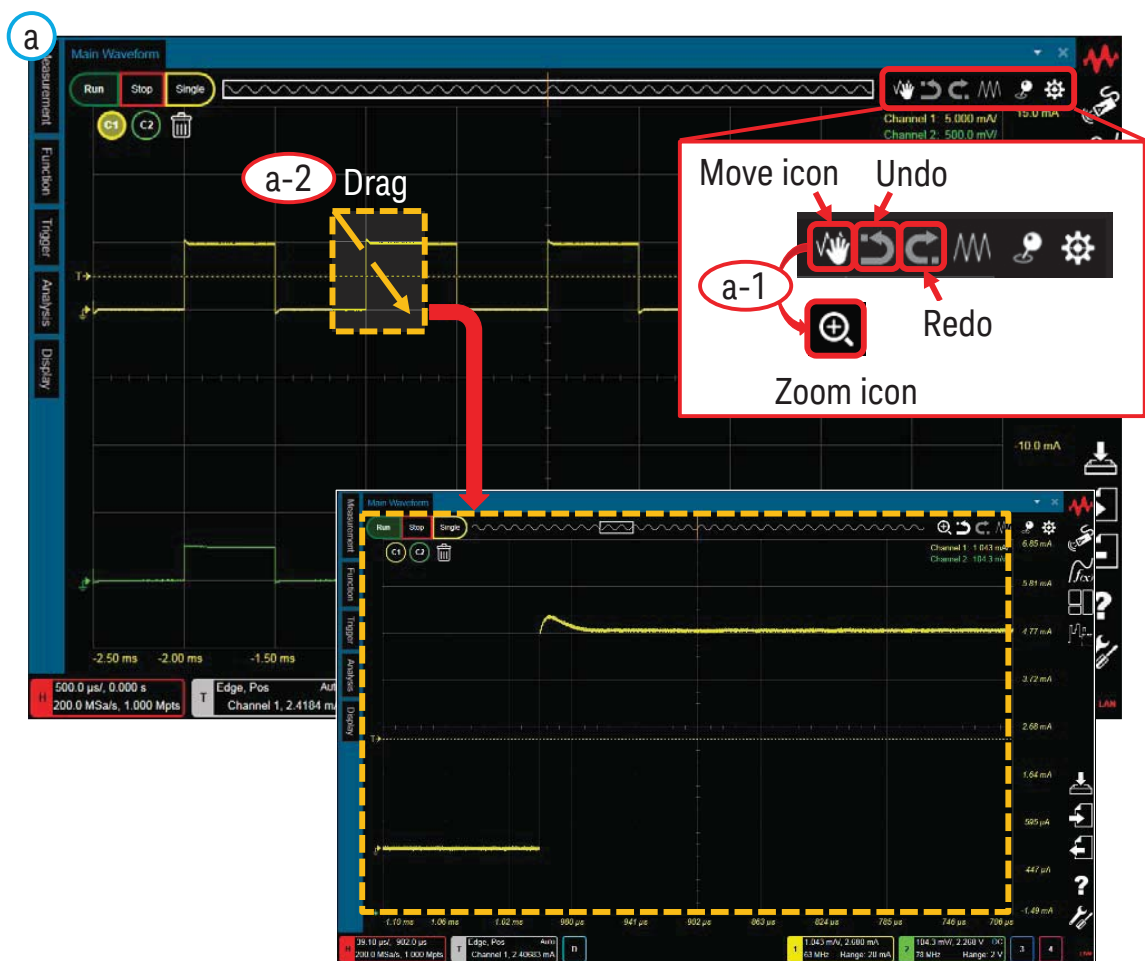


Tips: If you want to cancel the zoom, click the Undo icon.

Tips: By using the Scale Reference function, you can change the base point of the Zoom.

- Utilities > Configuration > User Preference > Scale Reference

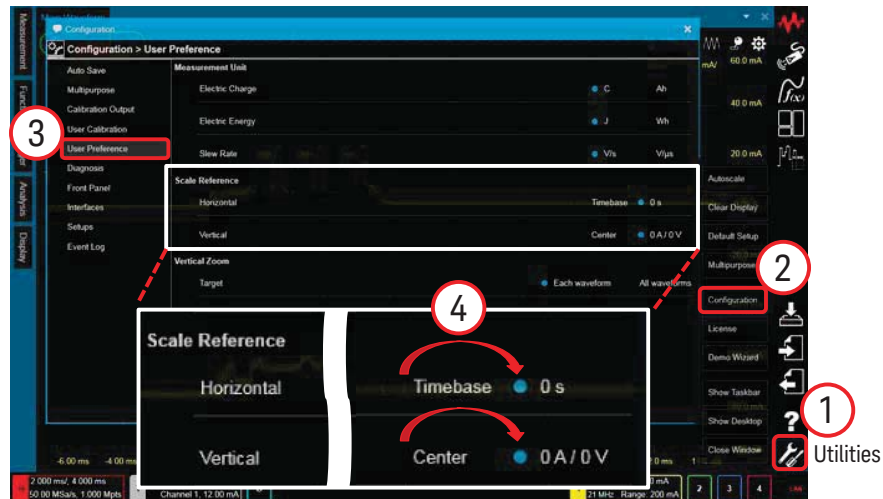

[Scale Reference > Horizontal Sample]

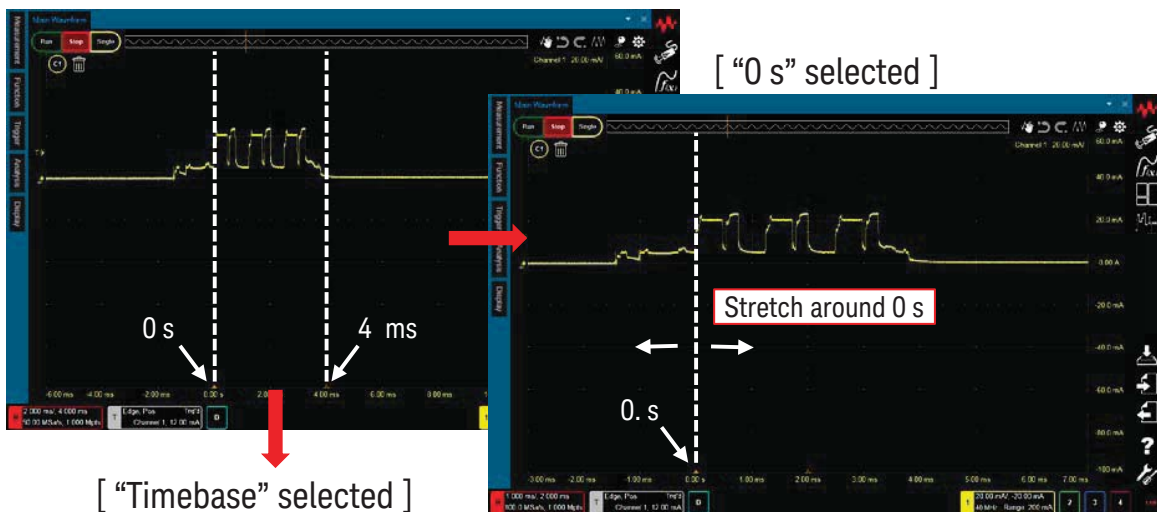

Note:

- The setting is also retained in CX3300 power cycle and [Default Setup].
- It returns to the factory condition with [Factory Default]¹.


1. [Factory Default] can be selected on the "Utilities > Configuration > Setups" window.

b. "Anywhere" Zoom


b-1. Click on the "Waveforms" icon.

b-2. Click on the "1 & zoom" icon.

b-3. Drag the boundary line of the area.

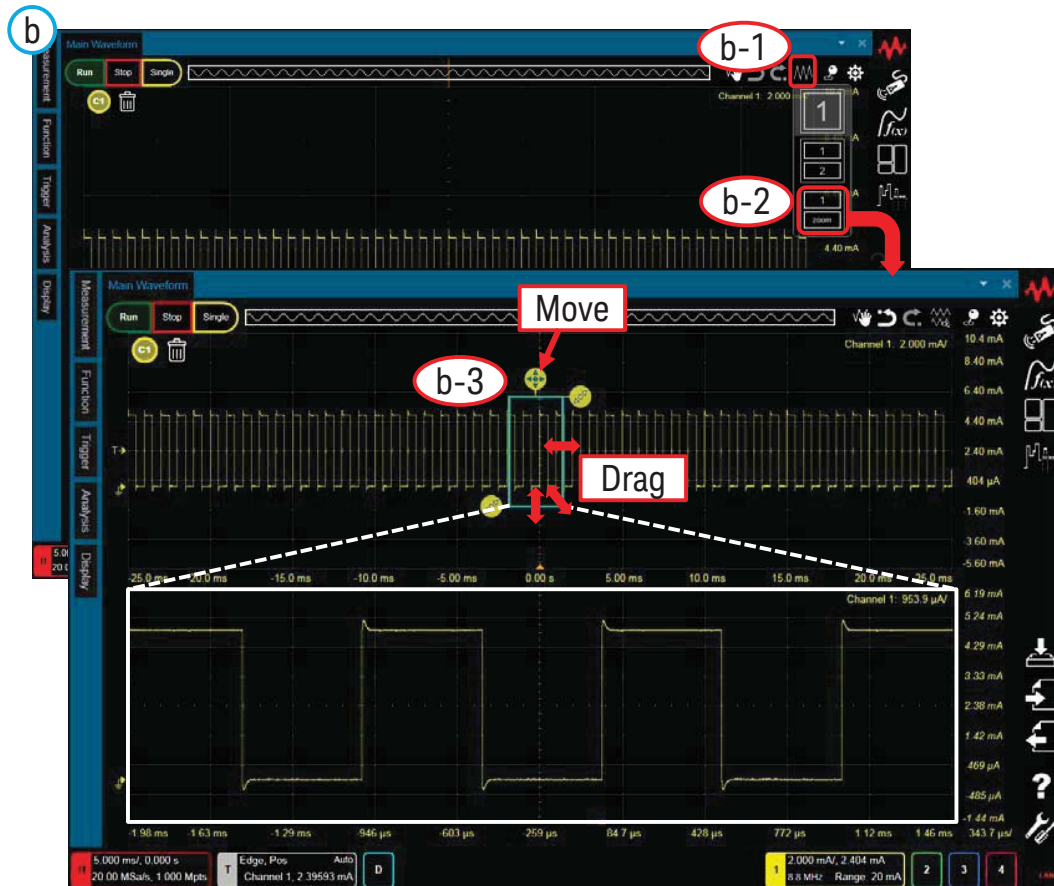

Tips: You can also specify the scale and position on the Axis dialog.

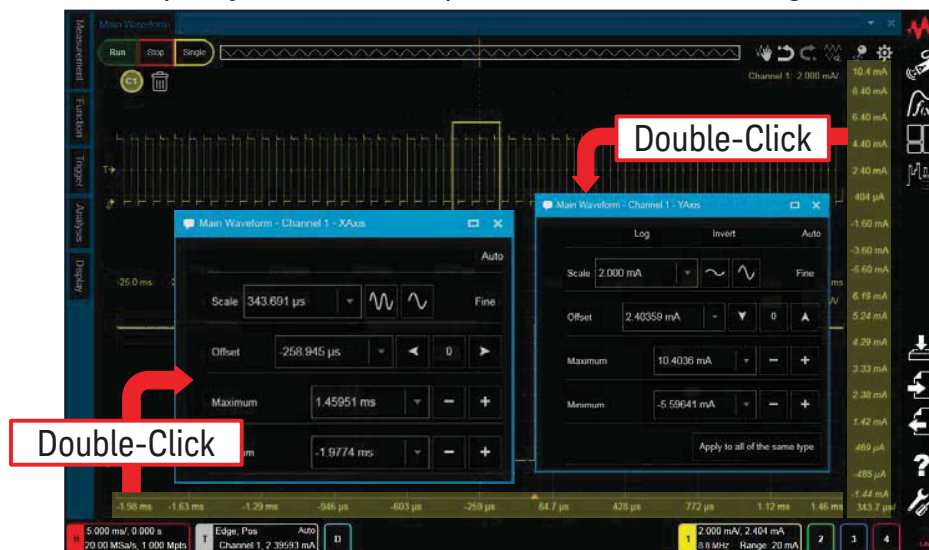


Keysight Technologies CX3300 Current Waveform Analyzer

Quick Operation Guide

How to use the waveform memory

1. Click on the “Analysis” icon to open the “Analysis” dialog box.
2. Select “Waveform Memories”.
3. Select “ON”.
4. Select “Memory X” (X is a number. X=1 in this case.)
5. Press [Load]. The waveform is loaded from Memory X.
6. Change the color, if necessary.

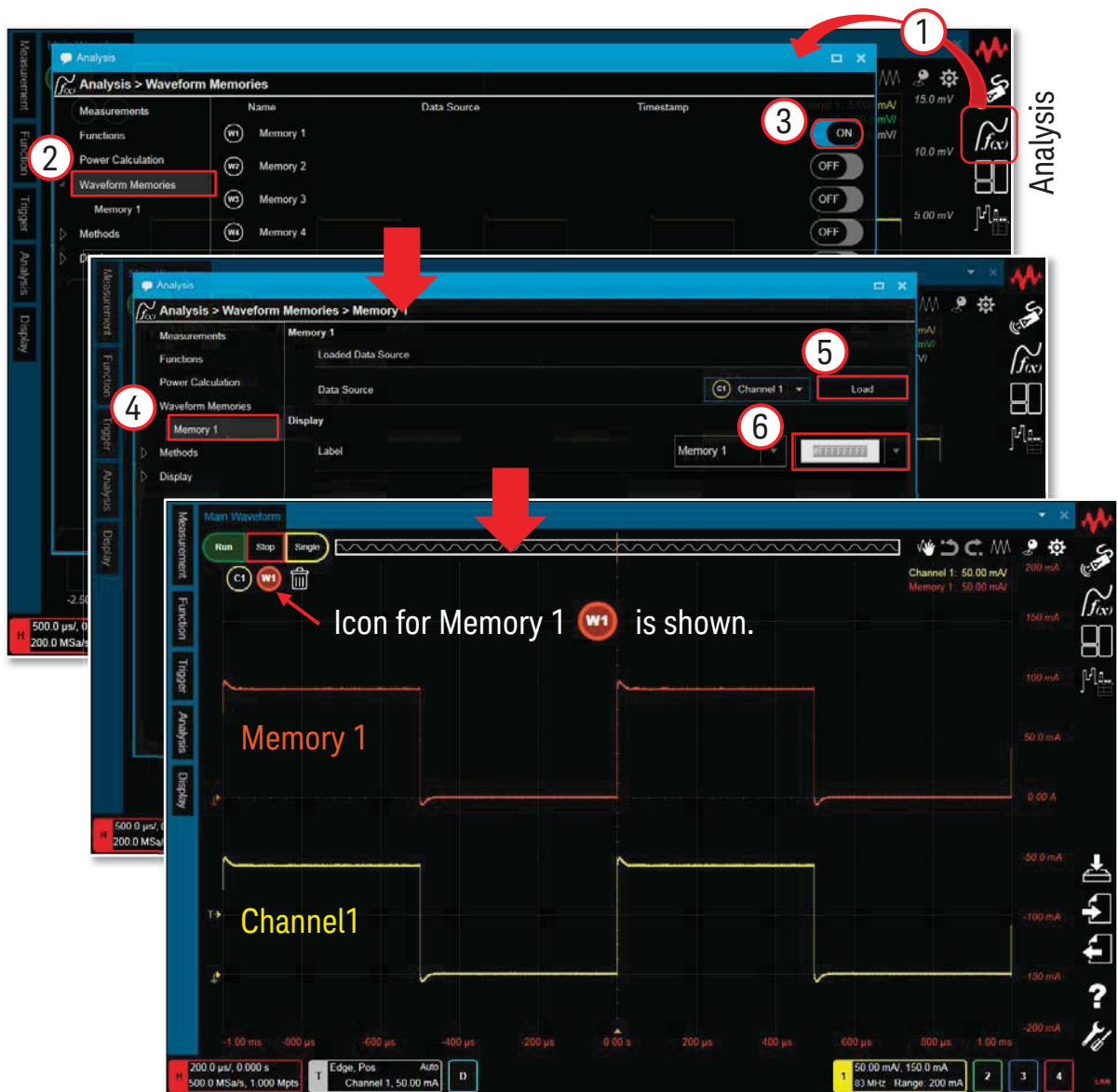

Keysight Technologies CX3300 Current Waveform Analyzer


Quick Operation Guide

How to overlay multiple acquisitions

1. ON : Display > Persistence > Infinity
 - 1-1. Click the “Display” icon on the Tool Palette.
 - 1-2. Select “Persistence” on the Display tool palette.
 - 1-3. Drag and drop the “Infinity” icon on the Display tool palette onto the area of interest on the waveform.



2. OFF : Display > Persistence > Normal
 - 2-1. Click the “Display” icon on the Tool Palette.
 - 2-2. Select “Persistence” on the Display tool palette.
 - 2-3. Drag and drop the “Normal” icon on the Display tool palette onto the area of interest on the waveform.

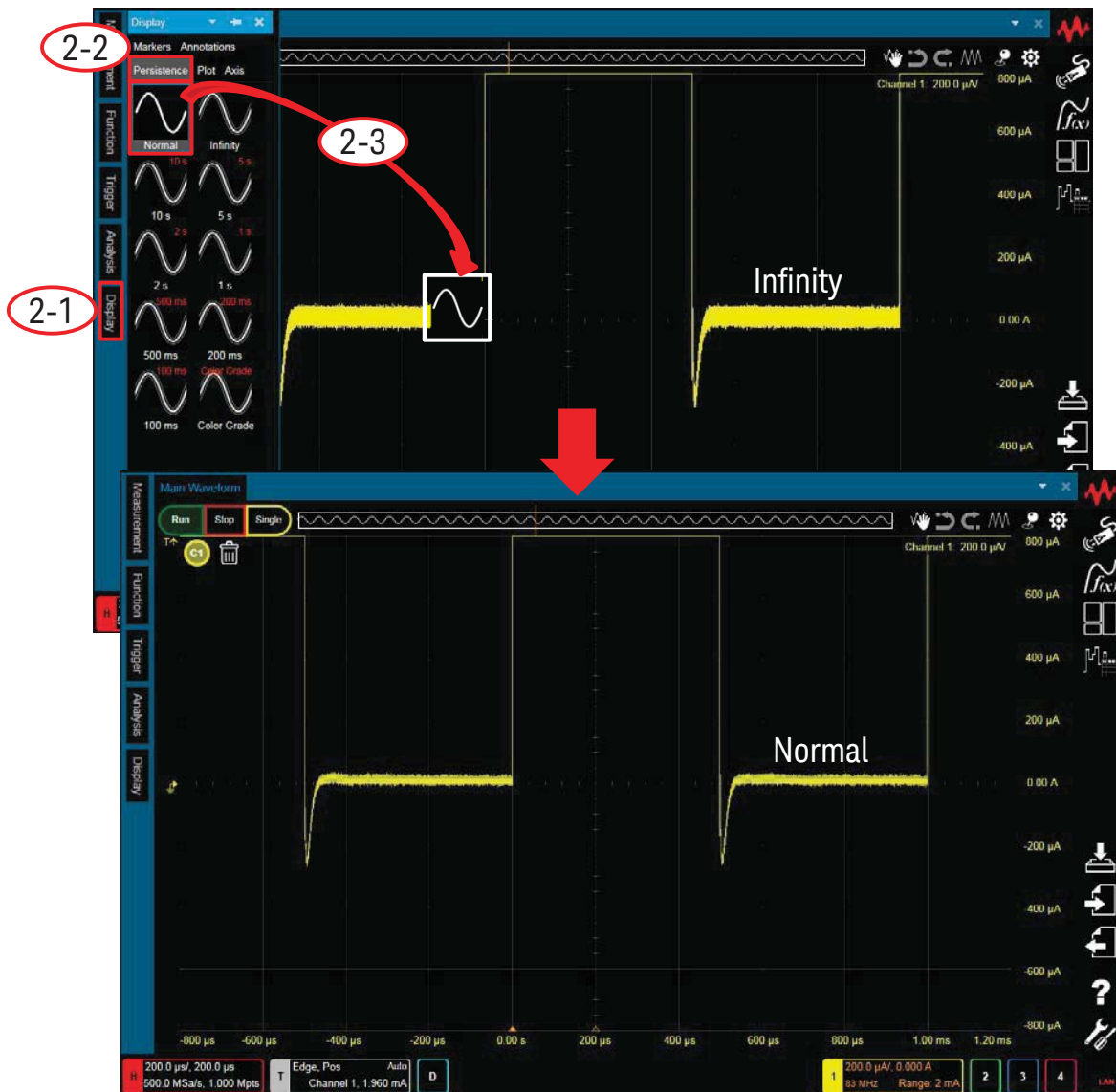


Keysight Technologies CX3300 Current Waveform Analyzer

Quick Operation Guide


How to display the secondary Y axis

1. Select “Axis” in the Display Tool Pallet. Then, drag and drop the “X1, Y1 and Y2” icon on the waveform.

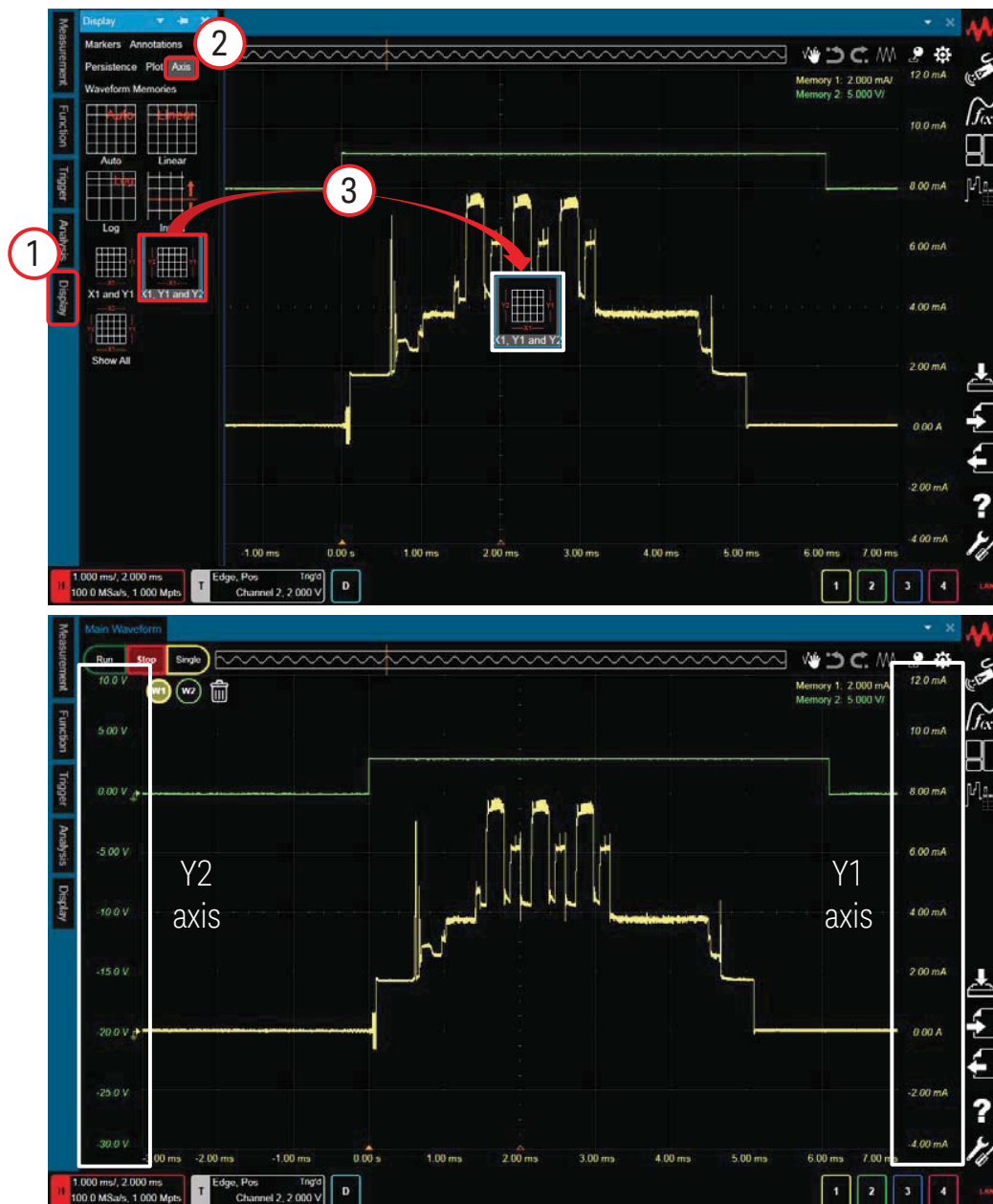


2. How to delete the secondary Y Axis.
There are 2 ways to delete the secondary Y Axis.

2-1. Drag and drop the “Display > Axis > X1 and Y 1” icon.

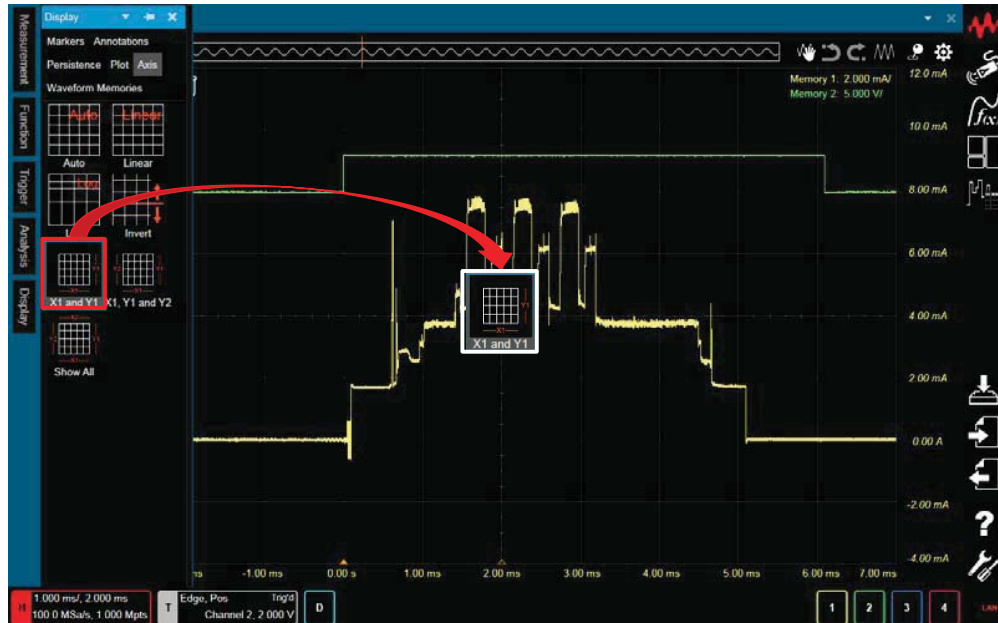


2-2. Double-click somewhere on the screen and switch OFF the “Secondary Y Axis” on the “Plot Properties” window.

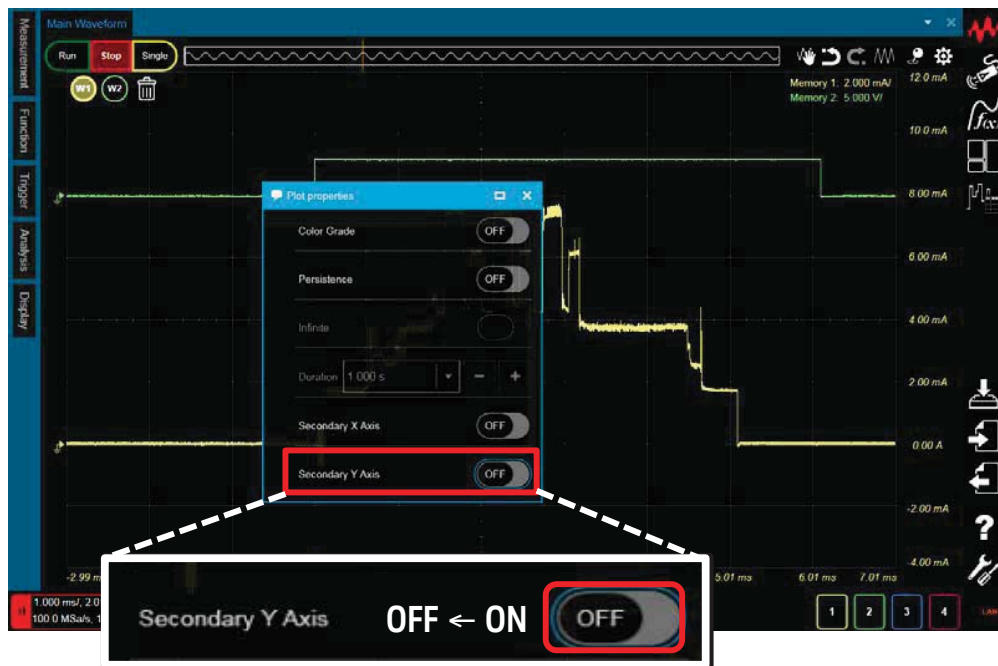


Quick Operation Guide

How to set the Trigger

1. Press [Sweep] to set the Trigger to ON.
2. Select the Trigger Type. (Default: “Edge”)
 - 2-1. Select “Trigger” on the Tool Palette
 - 2-2. Select the Trigger Type (“Edge” in this case) from the Trigger Tool Palette
 - 2-3. Drag-and-drop the icon onto the target waveform. The color of the icon shows the target waveform.

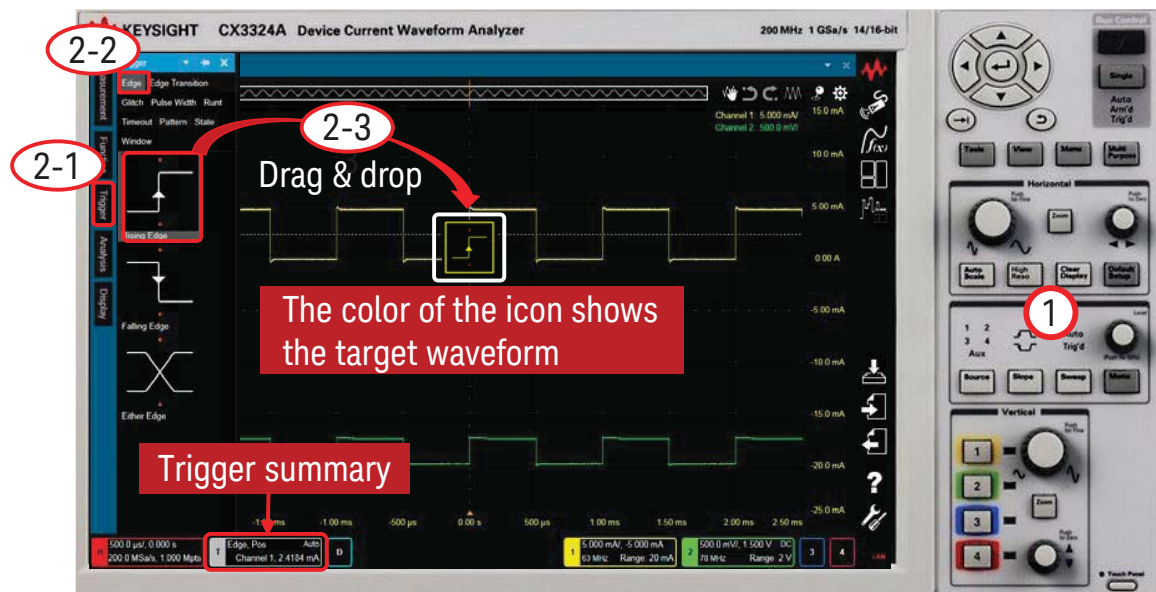

3. Change the trigger source, slope and level using either:
 - a) The Trigger Controls
 - OR
 - b) The Trigger Mini dialog box.
4. For more detailed conditions, open the “Setting > Trigger dialog box” by either:
 - a) Pressing the [Menu] button
 - OR
 - b) Clicking the “□” at the upper right of the Trigger mini dialog box.

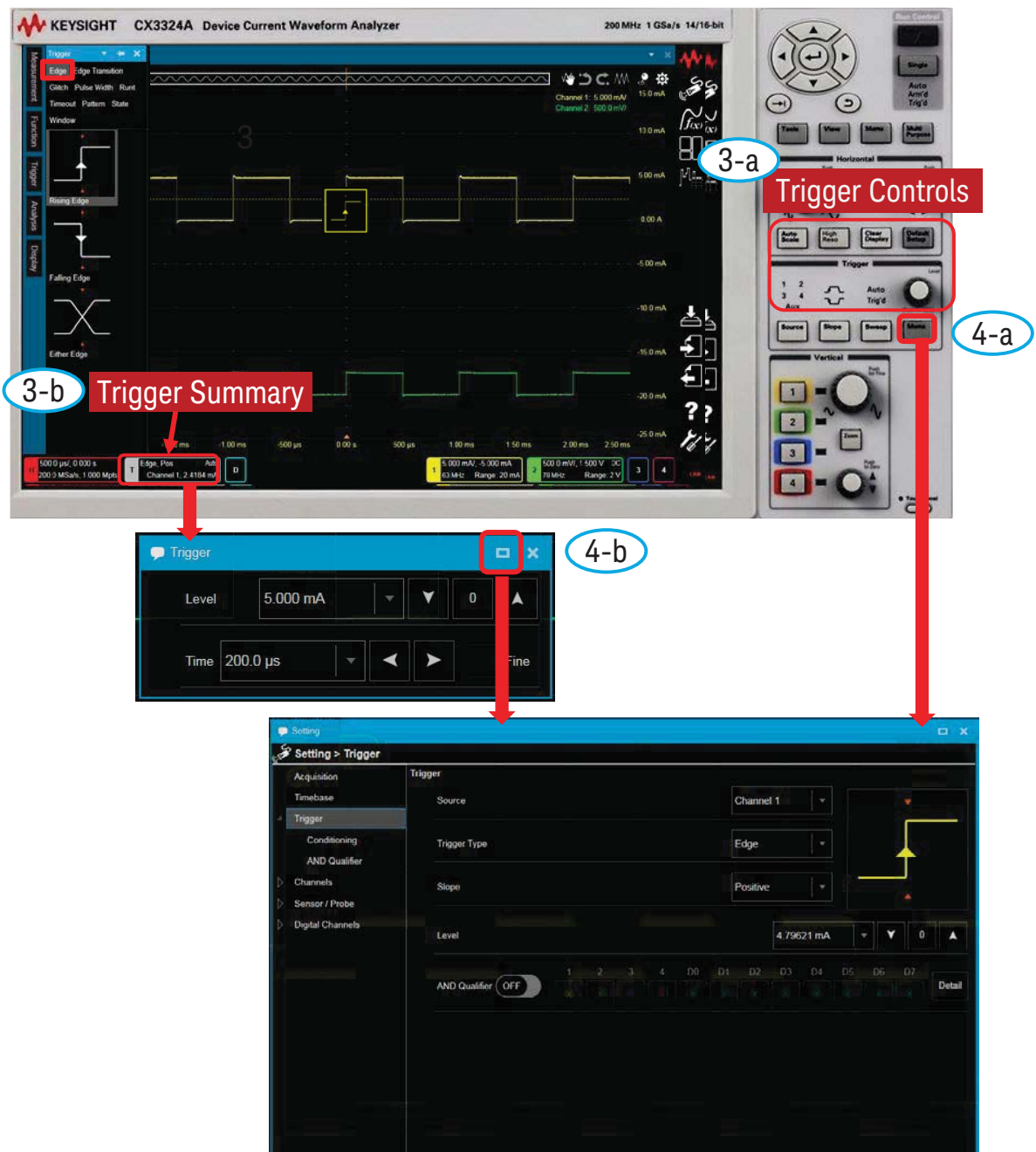


Keysight Technologies

CX3300 Current Waveform Analyzer


Quick Operation Guide

How to apply a Measurement Function

1. Apply a Measurement Function (“+width” in this case)
 - 1-1. Click the “Measurement” icon on the Tool Palette.
 - 1-2. Drag and drop the “+width” icon onto the target. The color of the icon shows the target waveform.

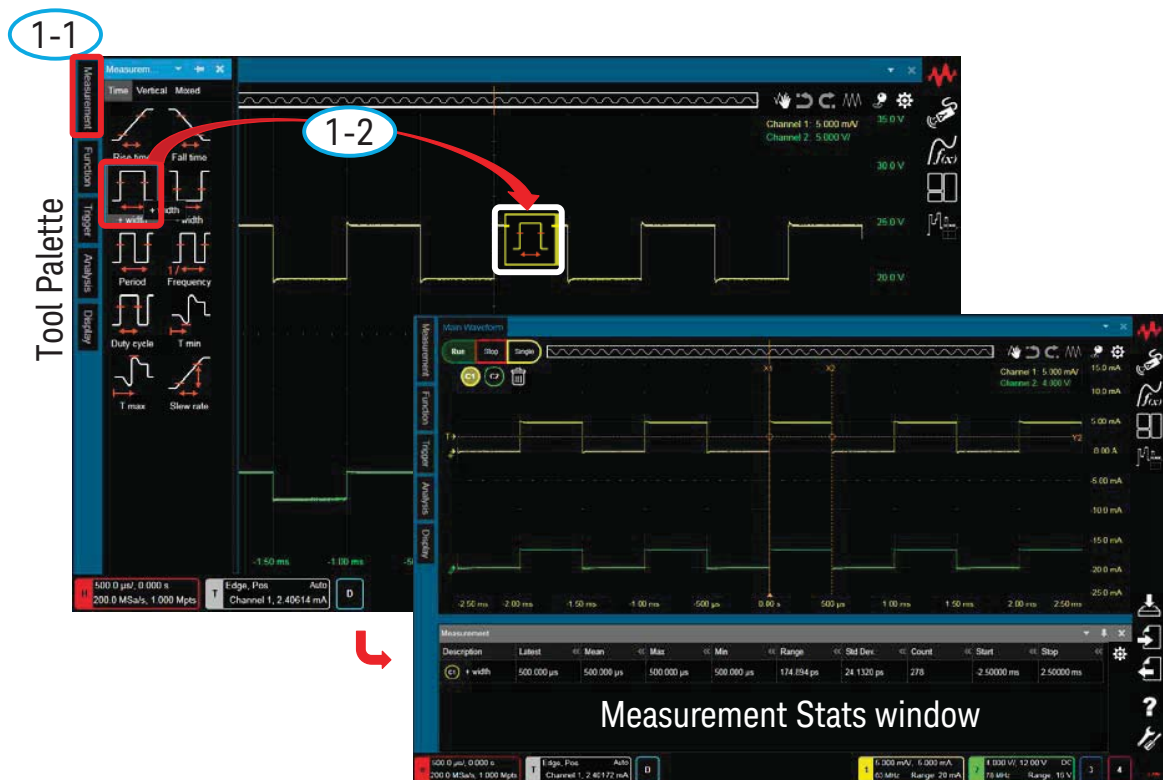


- For more detailed conditions, open the Measurements pane by clicking the Setting icon on the Measurement Stats window and change the settings.

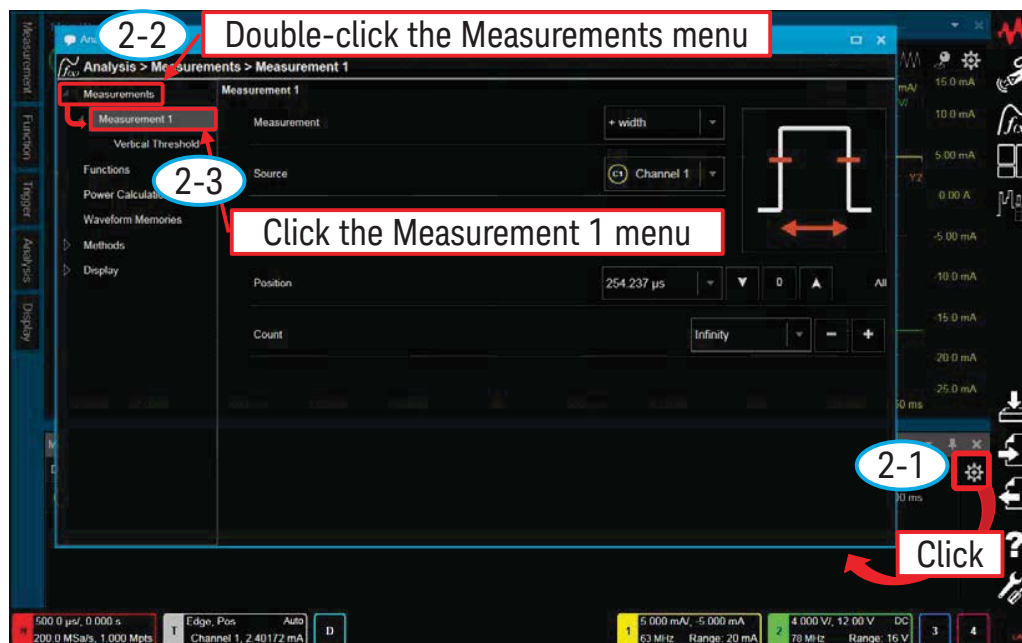

You can also access the Measurement pane via Side Bar > Analysis > Measurement.

- Turn the measurement function off.

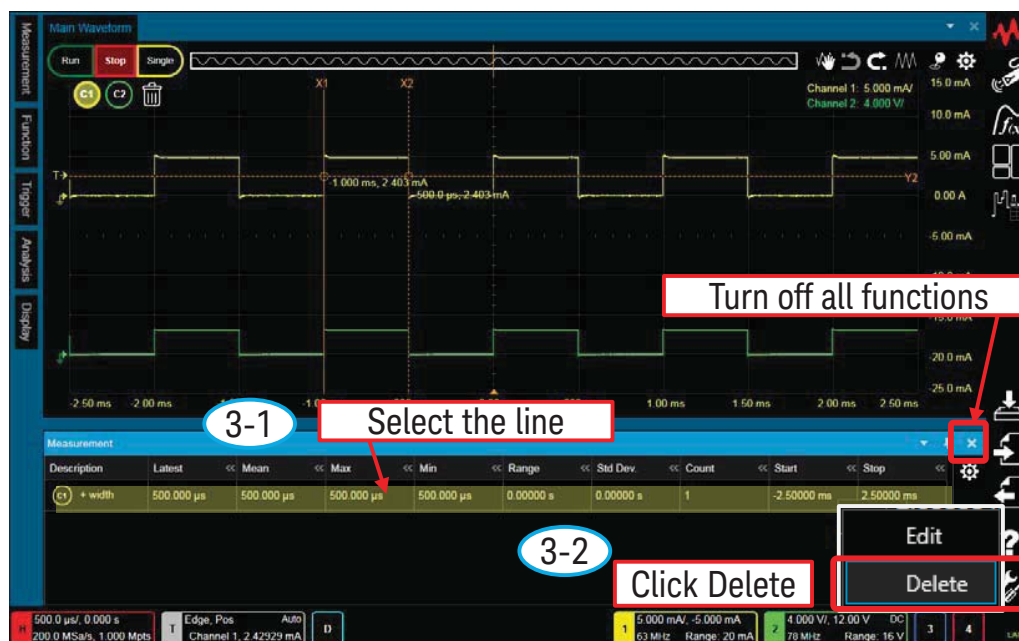

You can switch the measurement function ON and OFF on the Measurement pane as described in procedure 2.


Keysight Technologies

CX3300 Current Waveform Analyzer


Quick Operation Guide

How to apply the A-B marker


1. Turn the Marker ON.
 - 1-1. Click “Display” on the Tool Palette.
 - 1-2. Click “Markers”.
 - 1-3. Drag and drop the A-B marker icon onto the target. The color of the icon shows the target waveform.
 - 1-4. Move X1 and X2 to the target point.

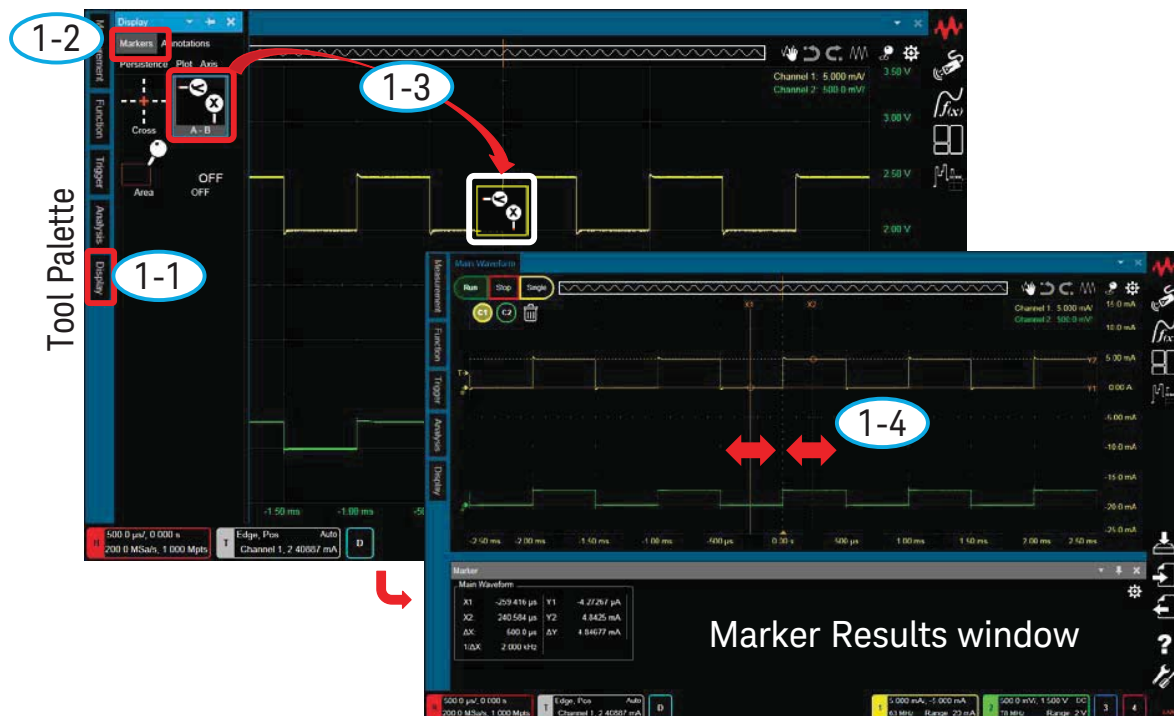


Tips: You can switch the Marker Results window ON and OFF via Side Bar > Layout > Marker Results.

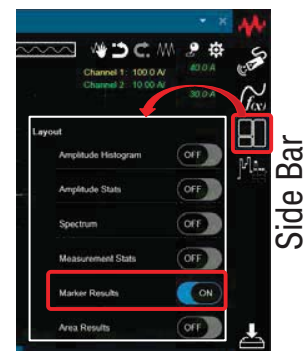

2. For more detailed conditions, open the Markers dialog box by clicking the Setting icon on the Marker Results window and change the settings.

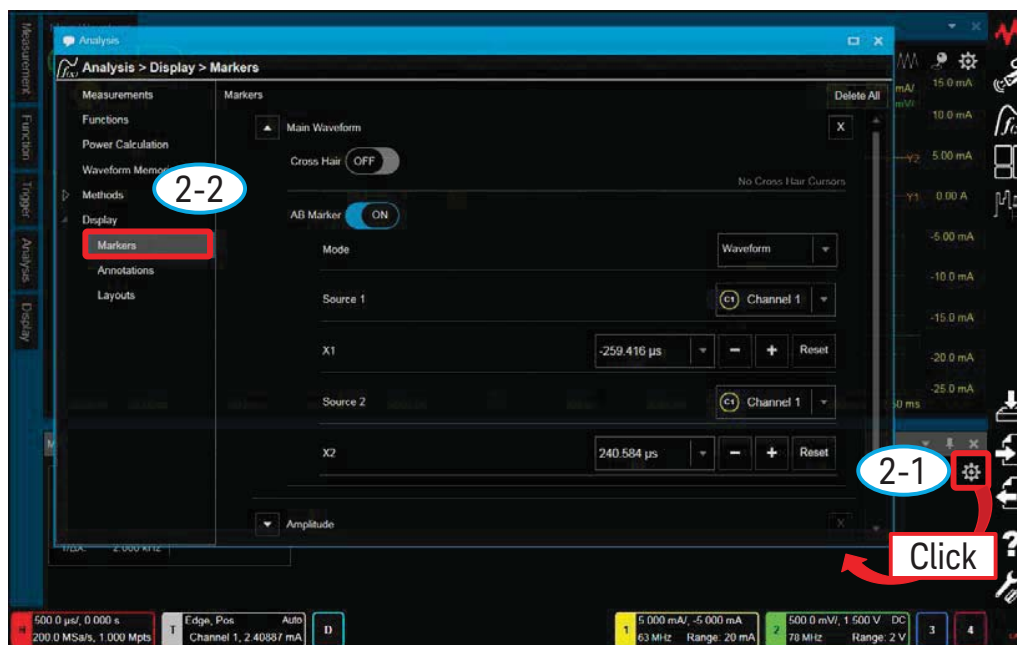

3. Turn the Marker OFF.
 - 3-1. Click “Display” on the Tool Palette.
 - 3-2. Click “Markers”.
 - 3-3. Drag and drop the OFF icon onto the target.

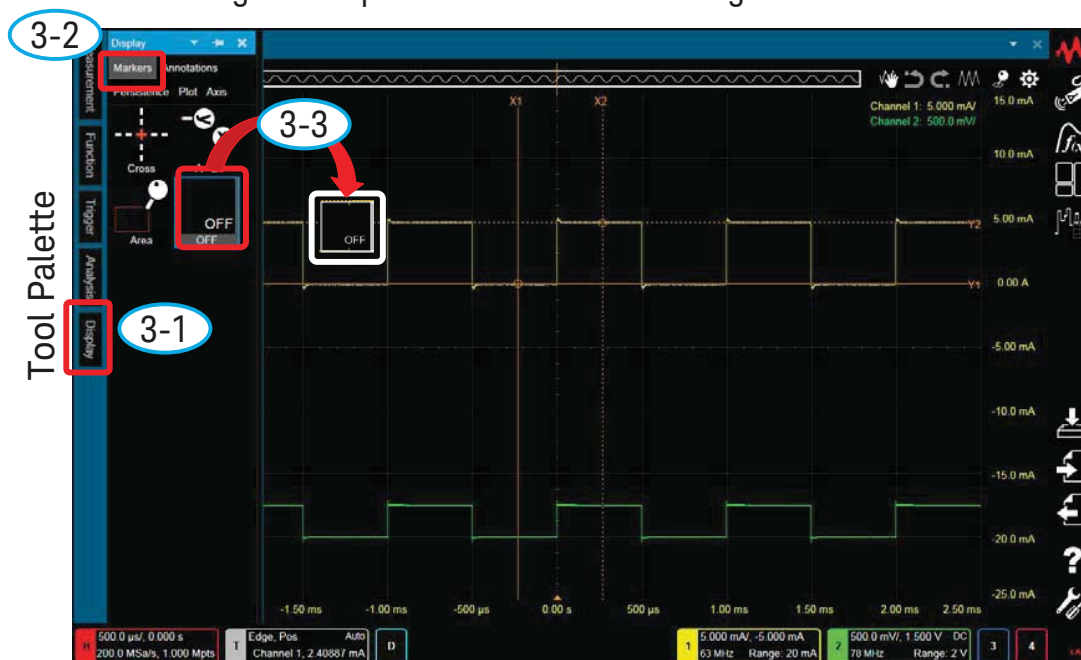


Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide


How to place pointers on waveforms

1. Select "Markers" in the Display Tool Pallet. Then, drag and drop the "Cross" icon on the point that you want to monitor.

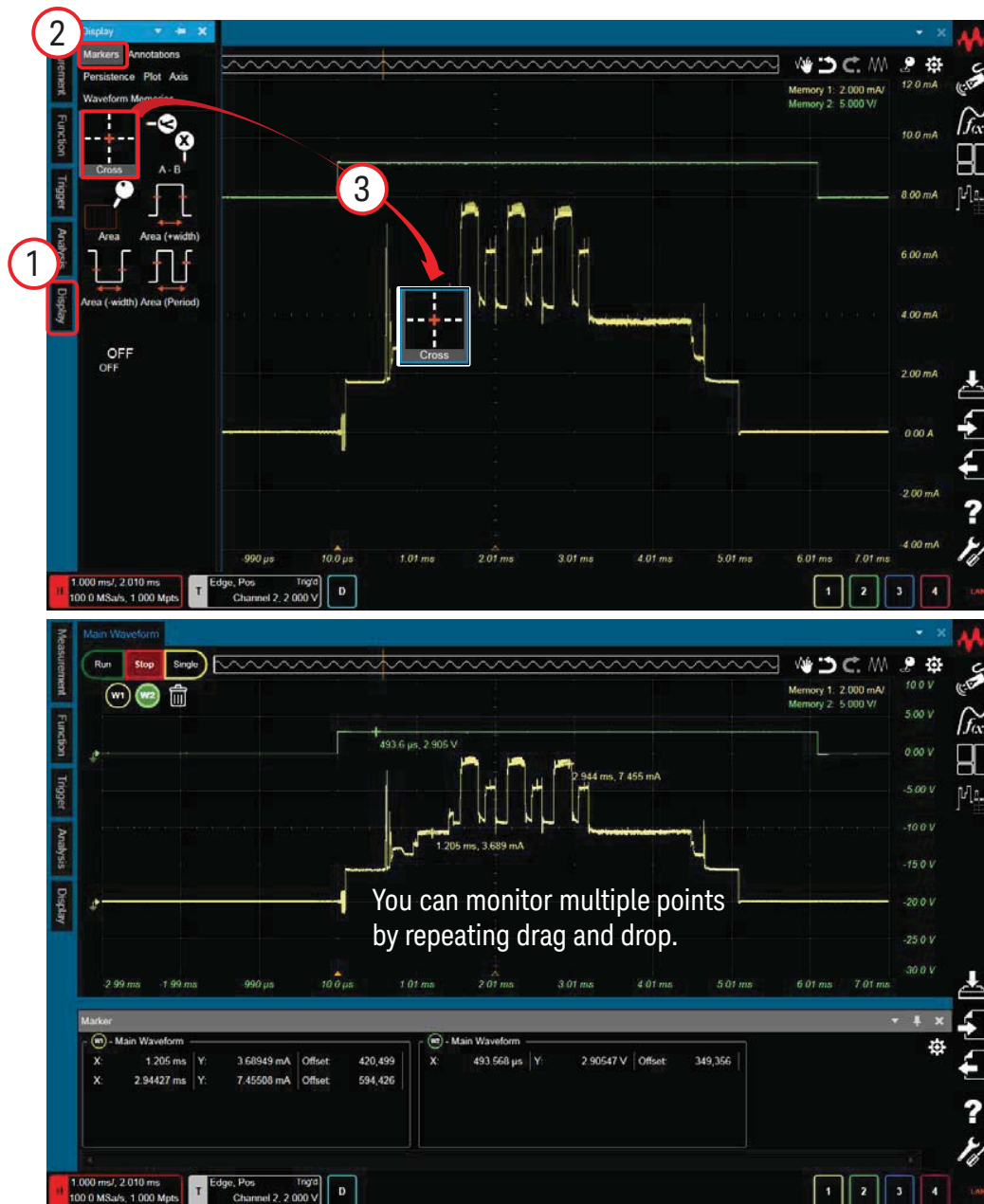


2. How to delete the Cross Hair Marker.
There are 3 ways to delete the cross hair marker.




2-1. Right-click the mouse on the cross hair marker and choose "Delete".

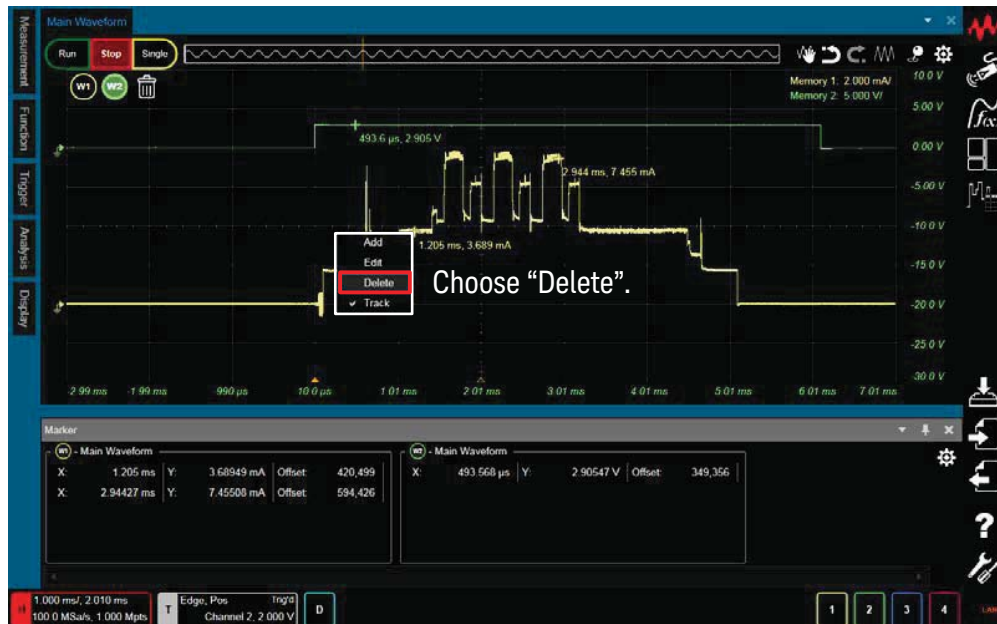

2-2. Click the gear icon on the upper right corner of the main screen.

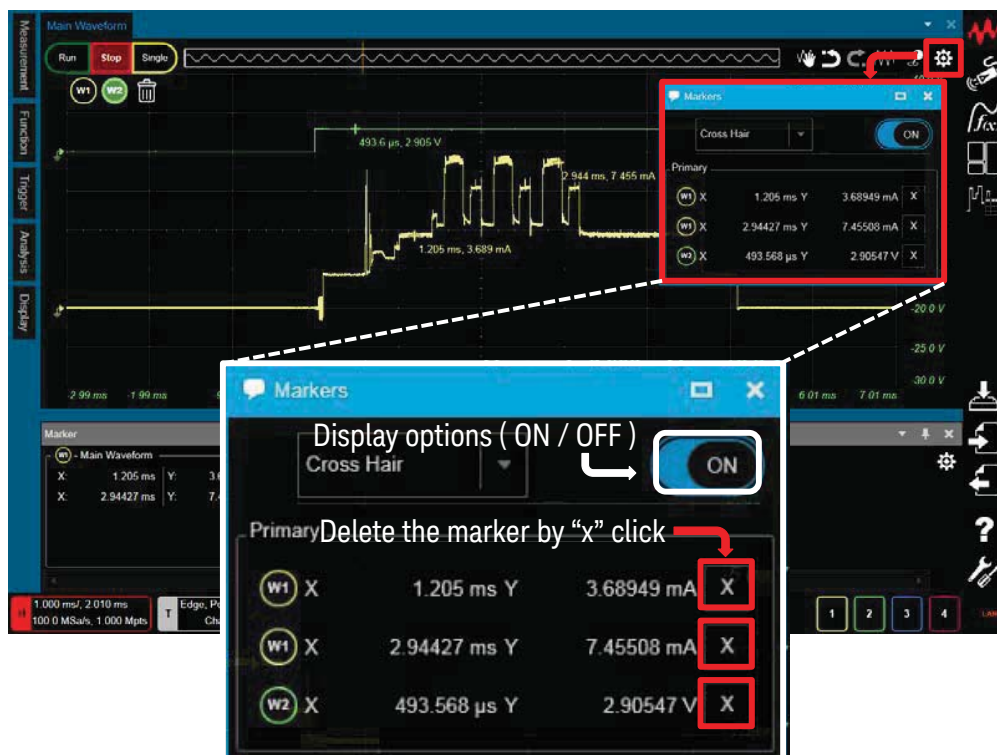



2-3. Click the gear icon on the upper right corner of the “Marker” screen.

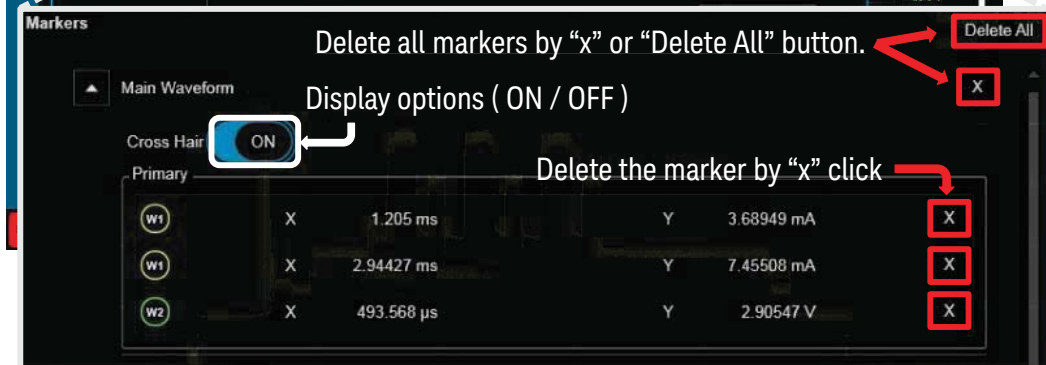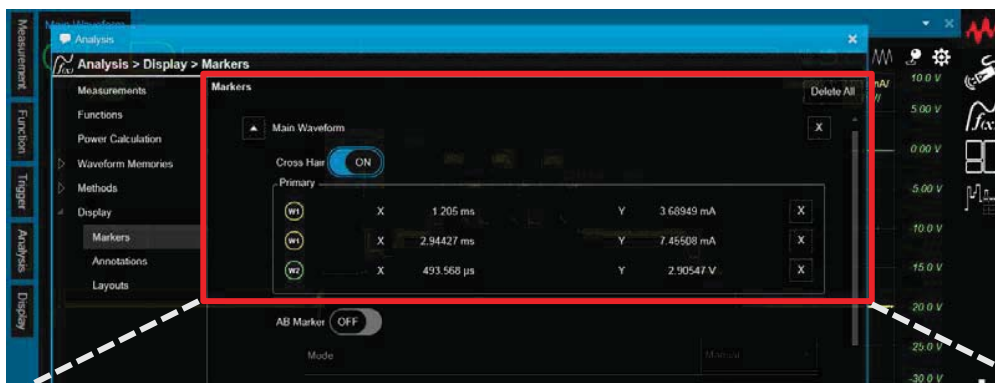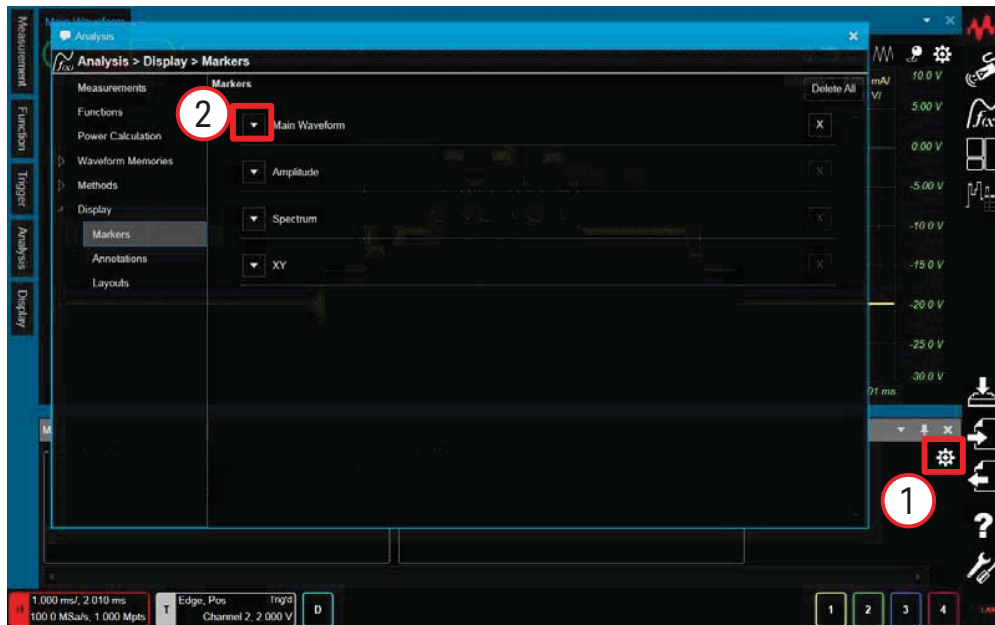

Tips :


- You can select whether or not the cross hair marker is tracked as follows:

Right-click the mouse on the cross hair marker.

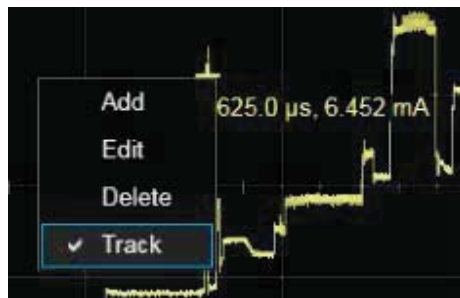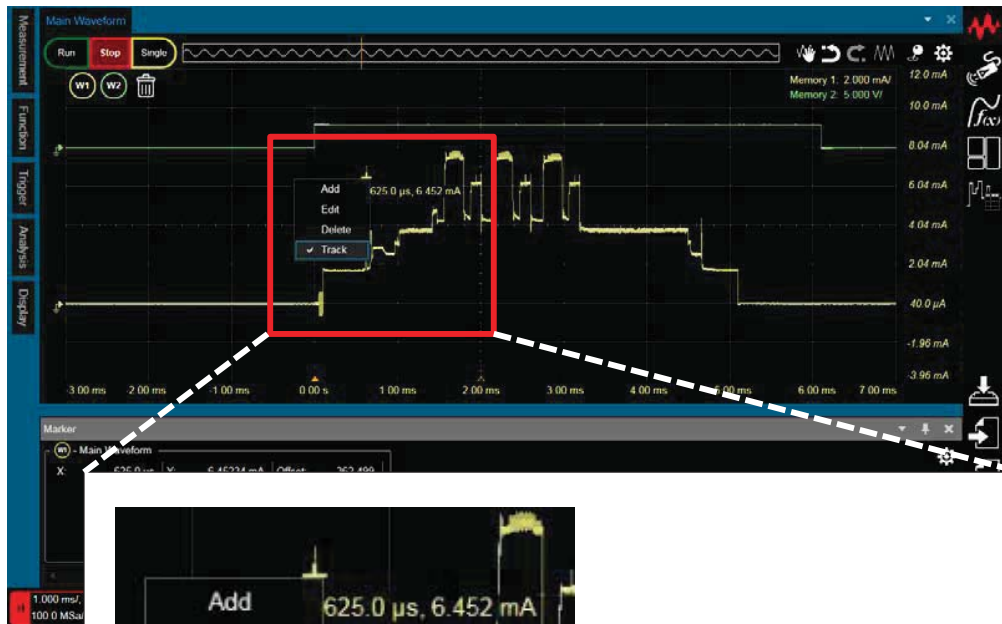

This cross hair marker can be tracked on the waveform in this mode.

Choose "Track".

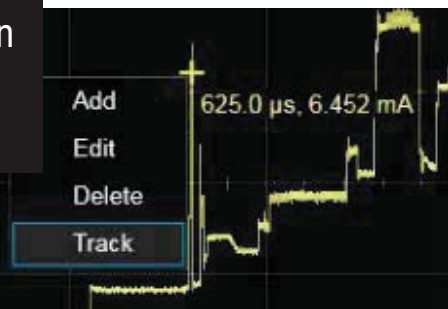

This cross hair marker can NOT be tracked on the waveform in this mode.


Keysight Technologies

CX3300 Current Waveform Analyzer


Quick Operation Guide

How to add an Annotation

1. Drag and drop the Display icon onto the target waveform.
 - 1-1. Click the “Display” icon on the Tool Palette.
 - 1-2. Click the “Annotations” tab.
 - 1-3. Drag and drop the “Display” icon on the Display tool palette onto the target waveform.

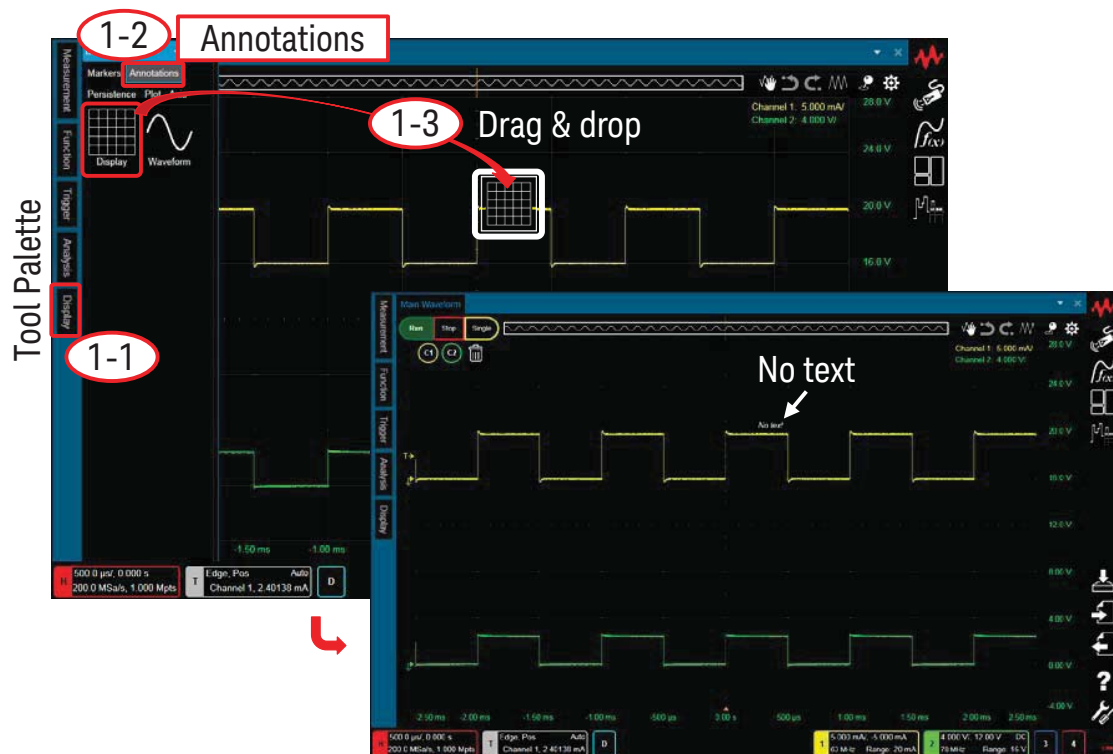


2. How to edit the Annotation
 - 2-1. Select Edit by right clicking on “No text”.
 - 2-2. You can edit the Annotation.

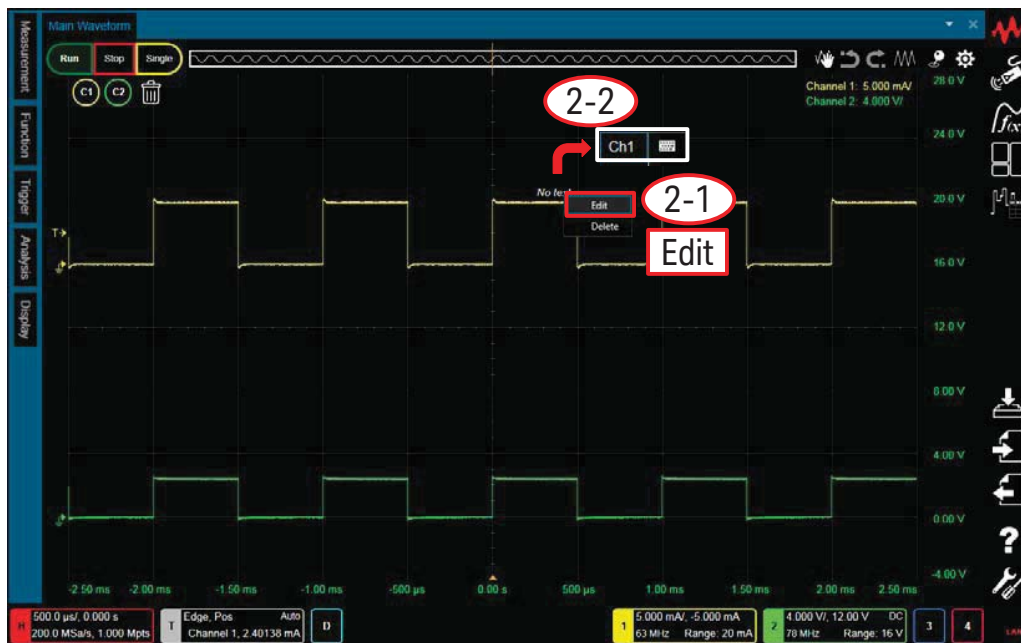

3. How to delete the Annotation

You can delete the Annotation by right clicking on it, as shown with “Ch1” in the example below, and selecting Delete.

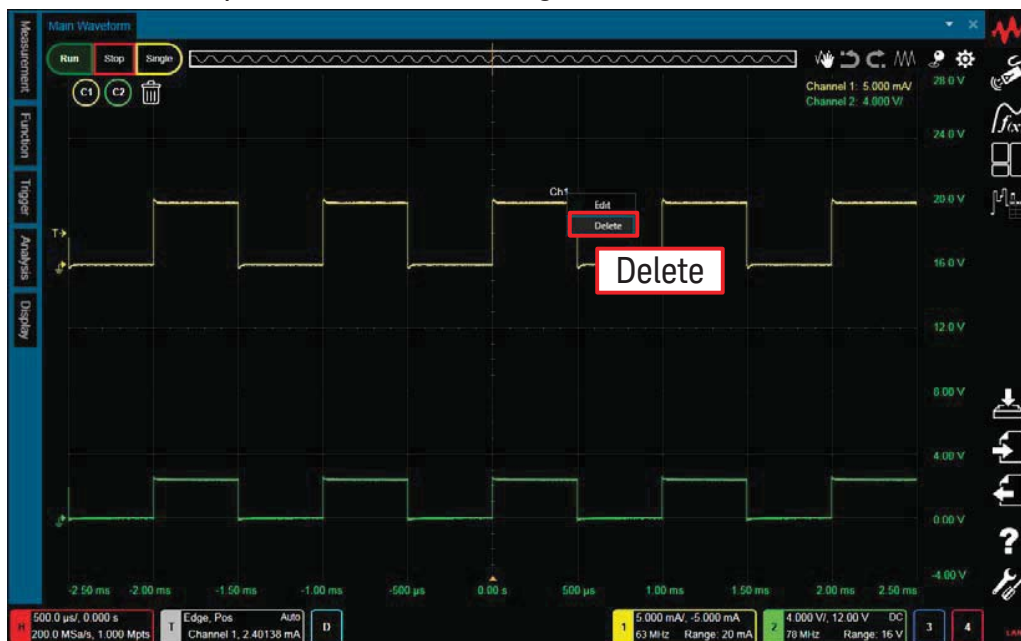


Keysight Technologies CX3300 Current Waveform Analyzer


Quick Operation Guide

How to apply the post-filter

1. How to apply the post-filter
 - 1-1. Click the “Function” icon on the Tool Palette.
 - 1-2. Select “Filter” on the Function tool palette.
 - 1-3. Drag and drop the “Smooth 100pts” icon on the Function tool palette onto the area of interest on the waveform.

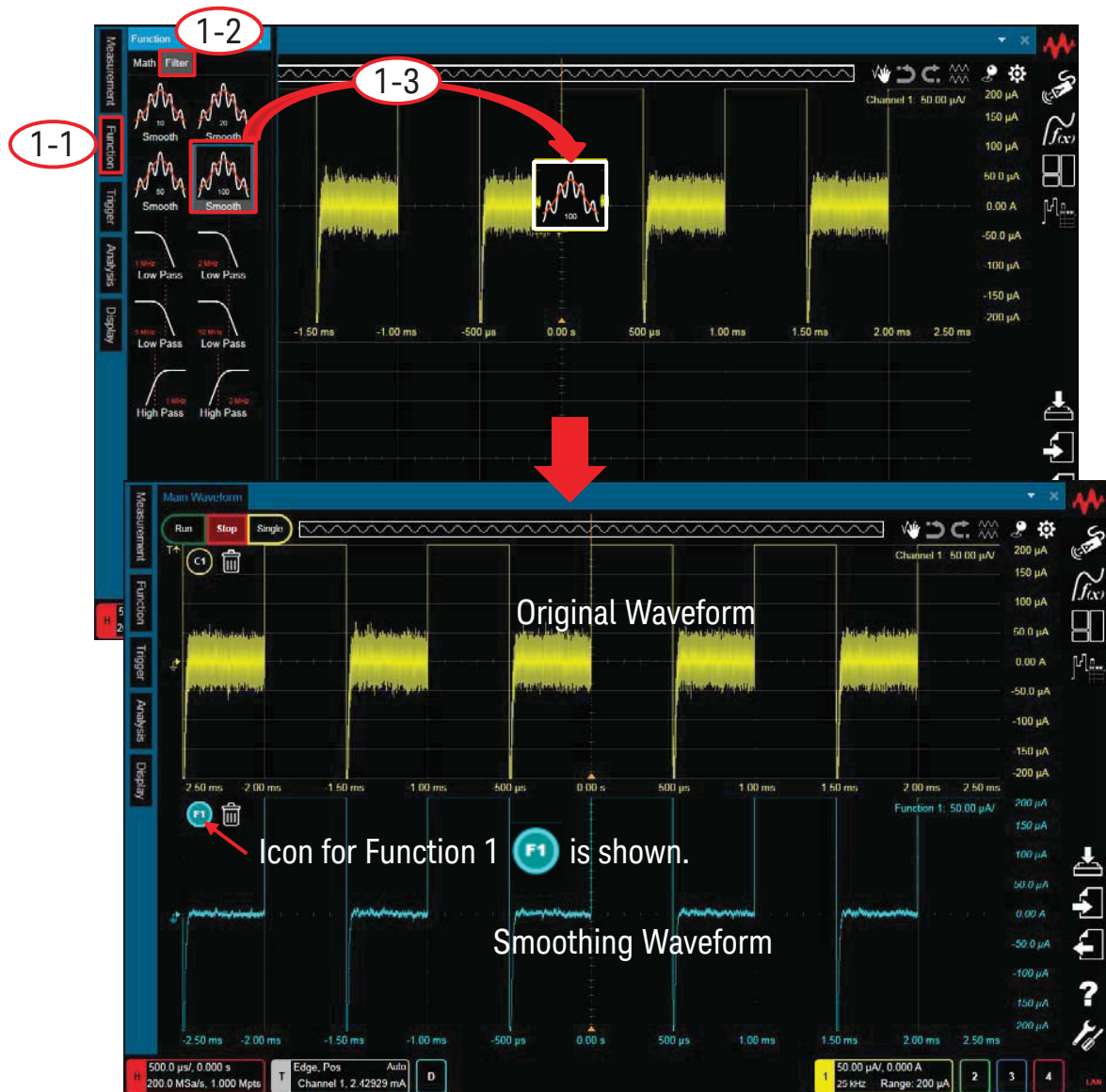

2. How to change the settings
 - 2-1. Click on the “Analysis” icon to open the Analysis dialog box.
 - 2-2. Select “Function 1”.
 - 2-3. You can change the Smoothing Points.

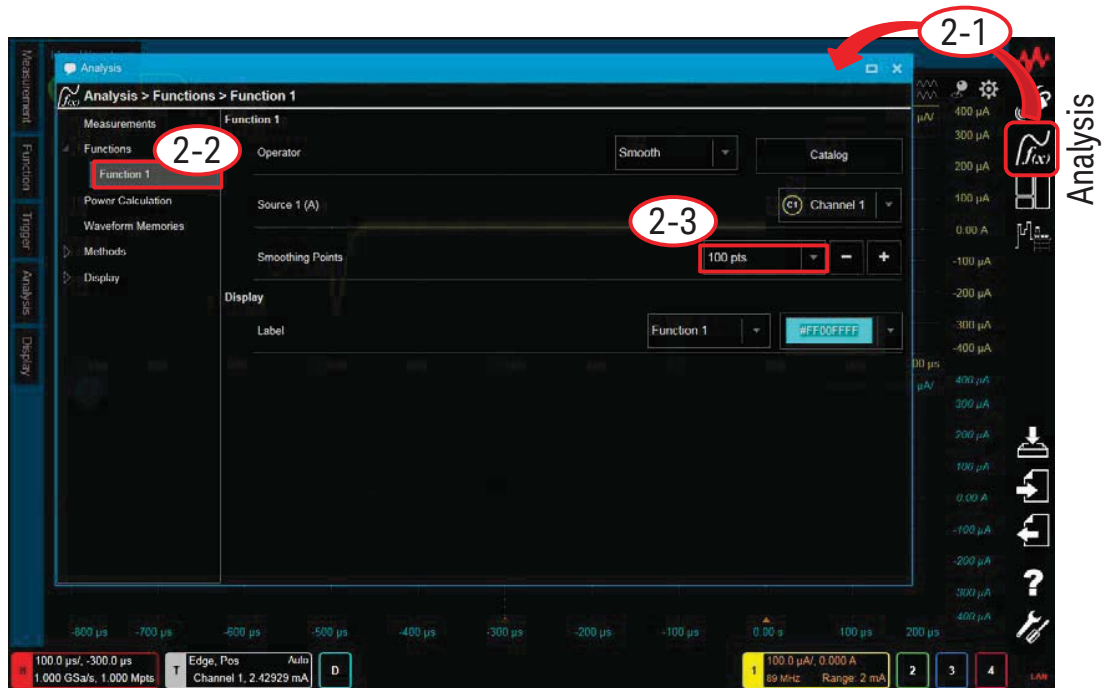


Keysight Technologies CX3300 Current Waveform Analyzer

Quick Operation Guide

How to perform the FFT (Fast Fourier Transform)


1. How to apply the FFT
 - 1-1. Click the “Analysis” icon on the Tool Palette.
 - 1-2. Drag and drop the “Spectrum” icon on the Analysis tool palette onto the chosen area of interest on the waveform.

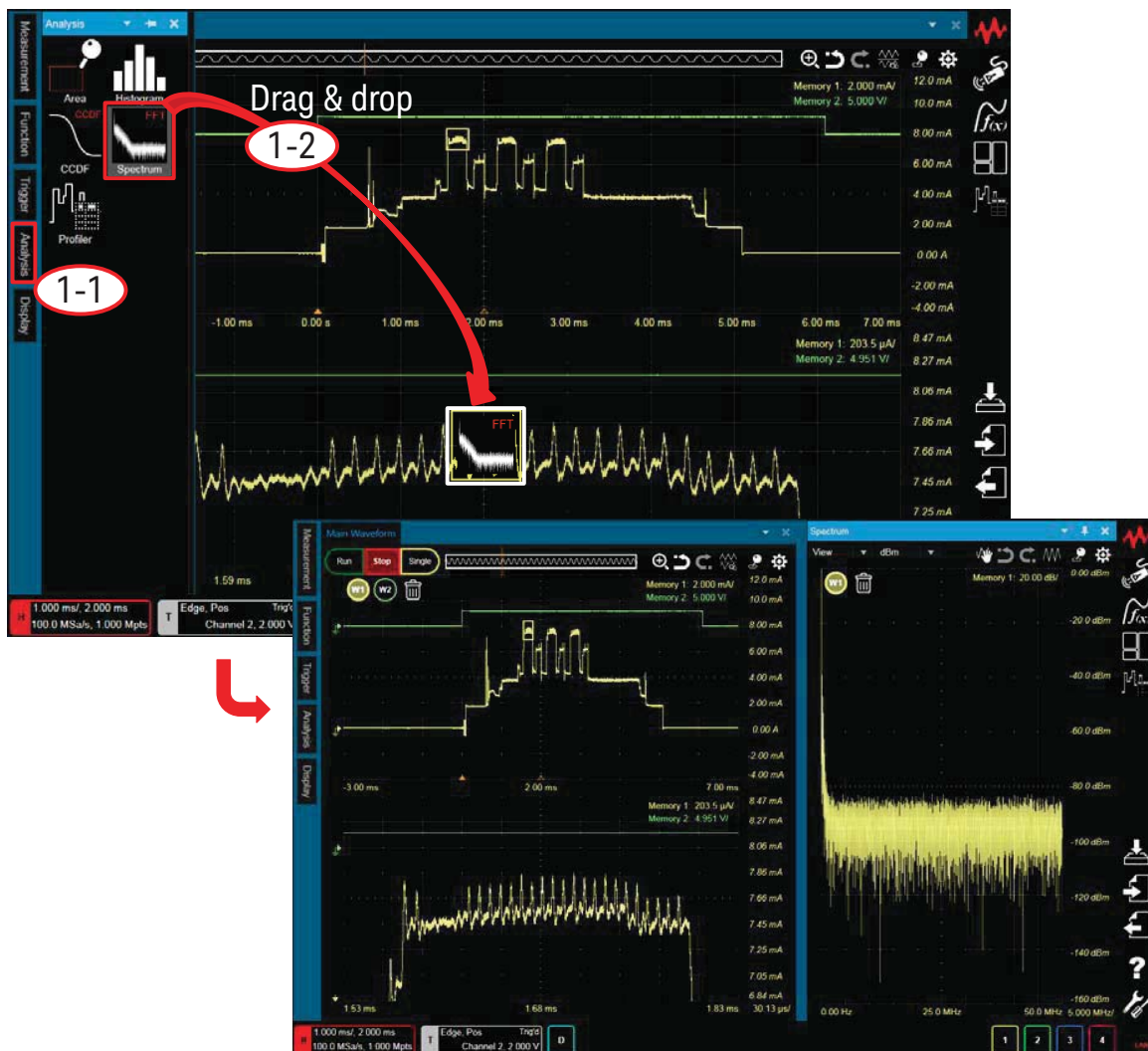

Note:

- If the Zoom window is displayed in the lower half, the dragged icon is always effective for the zoomed waveform.

2. How to change the FFT settings ?
 - 2-1. Double-click anywhere on the horizontal axes in the Spectrum graph view.
 - 2-2. Un-check the "Auto" box.
 - 2-3. Click this button until the Scale is 100.0KHz.

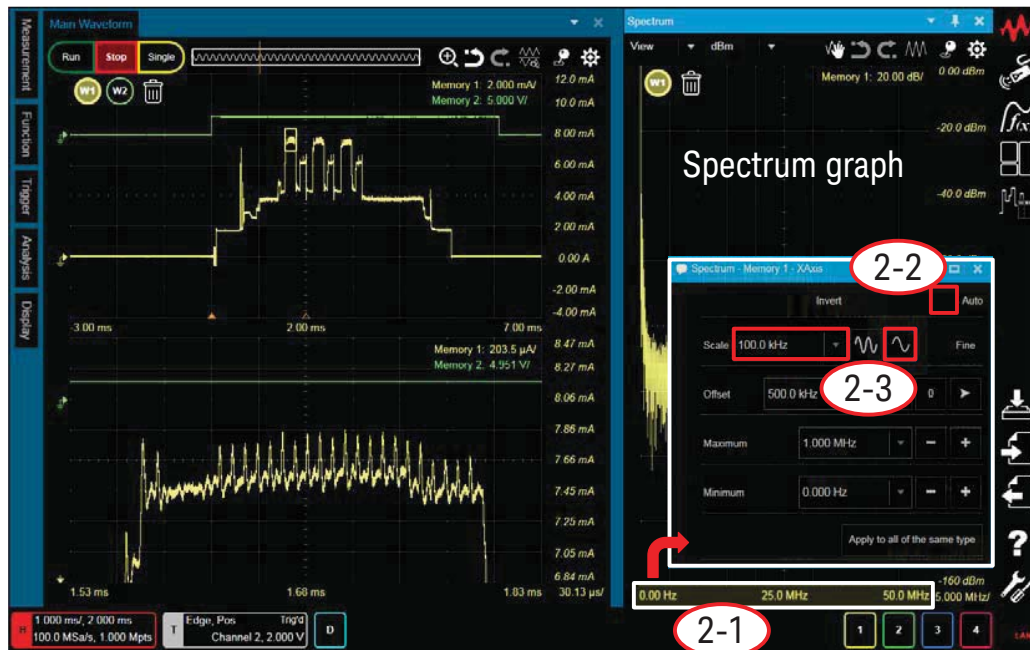


3. How to adjust the chosen area of interest on the waveform ?

The CX3300 has two ways in which to adjust the chosen area of interest.

- a) Gating Function > FFT
- b) Waveform Zoom > FFT

Refer to the next page for further details.

a) Gating Function > FFT

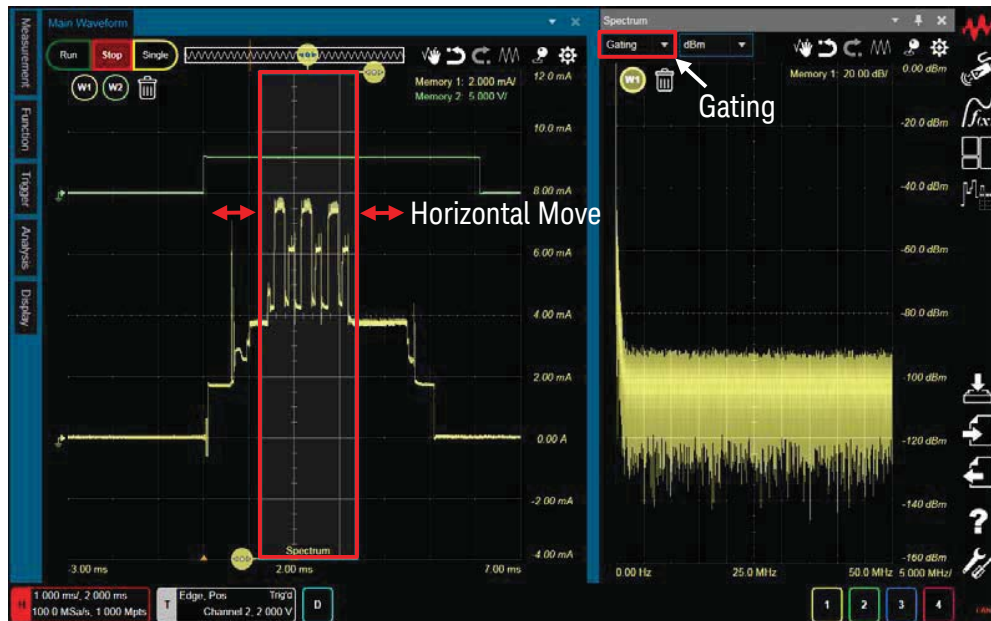

b) Waveform Zoom > FFT



Note:

- The result of the FFT is changed in real time by adjusting the chosen area of interest on the waveform.

Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to analyze the waveform in the specified area with an Area marker


There are four types of Area markers as follows; Area, Area (+width), Area (-width) and Area (period).

Using the mouse drag and drop option, the Area marker allows you to freely change the area range.


Using Area, Area (+width), Area (-width) and Area (period) marker, the area is automatically determined according to where you drag and drop the marker.

[Apply the “Area” marker]

1. Turn the Marker ON.
 - 1-1. Click the “Display” icon on the Tool Palette.
 - 1-2. Select “Markers”.
 - 1-3. Drag and drop the “Area” icon on the Display tool palette onto the specified area on the waveform.

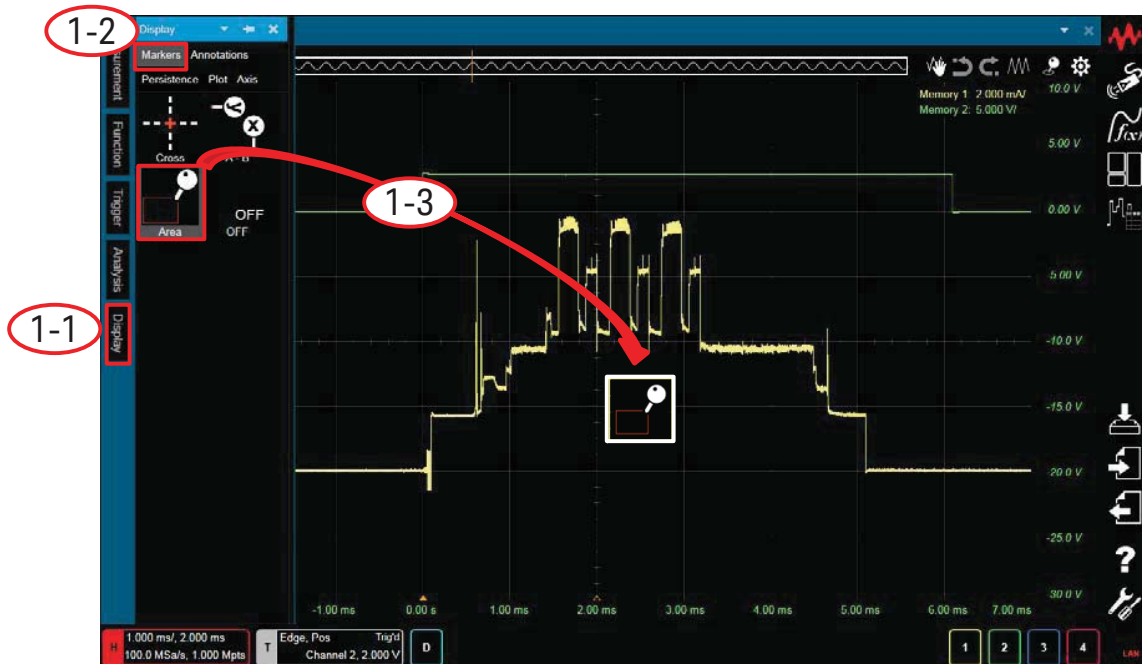


- Click and hold the mouse button, and horizontally drag the Start or Stop line over the specified area of interest.

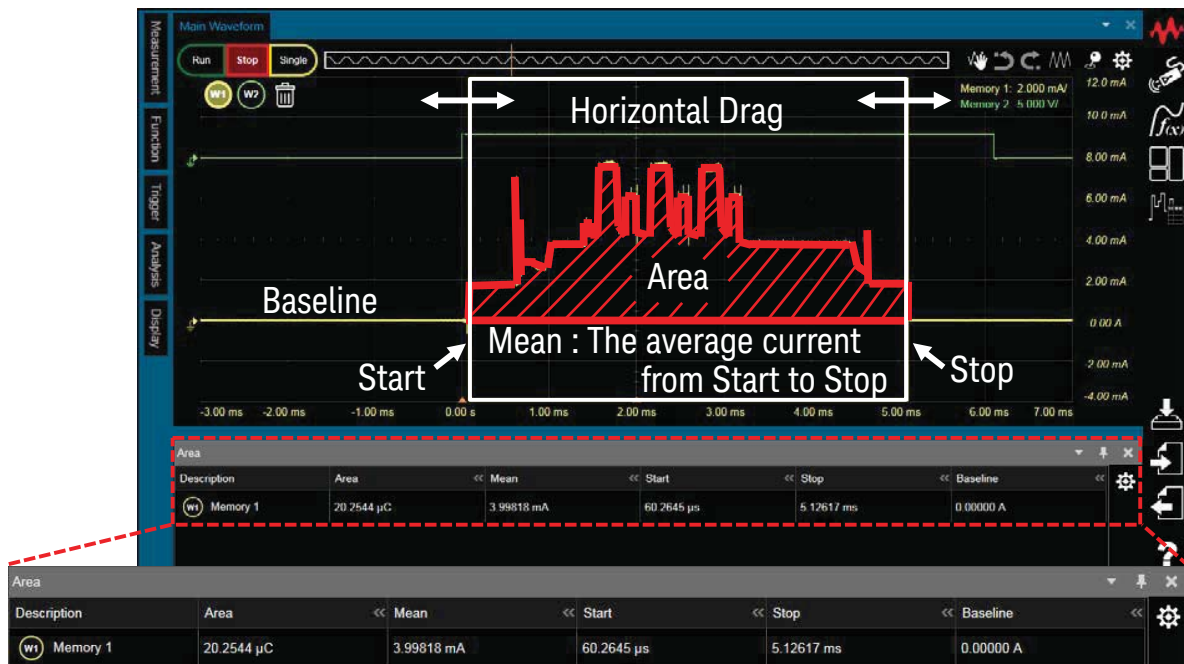

[Apply the “Area (+width), Area (-width), Area (period) marker”]

- Apply the “Area (+width) marker”
For example, as with the Area marker, drag and drop the “Area (+width)” icon on the Display tool palette onto the area of interest.

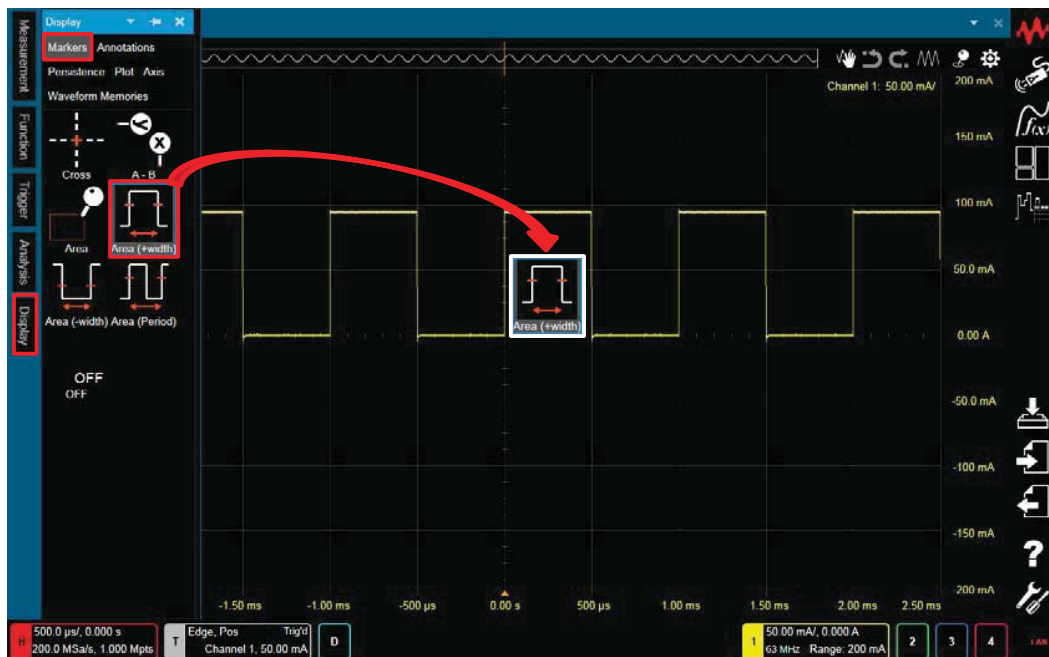

- The result of applying the "Area (+ width), Area (-width), Area (period) marker"

The area is automatically determined according to where you drop the marker.




Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to apply the Power and Current Profiler

1. Drag and drop the Display icon onto the target waveform.
 - 1-1. Click the “Analysis” icon on the Tool Palette.
 - 1-2. Drag and drop the “Profiler” icon on the Analysis tool palette onto the chosen area on the waveform.

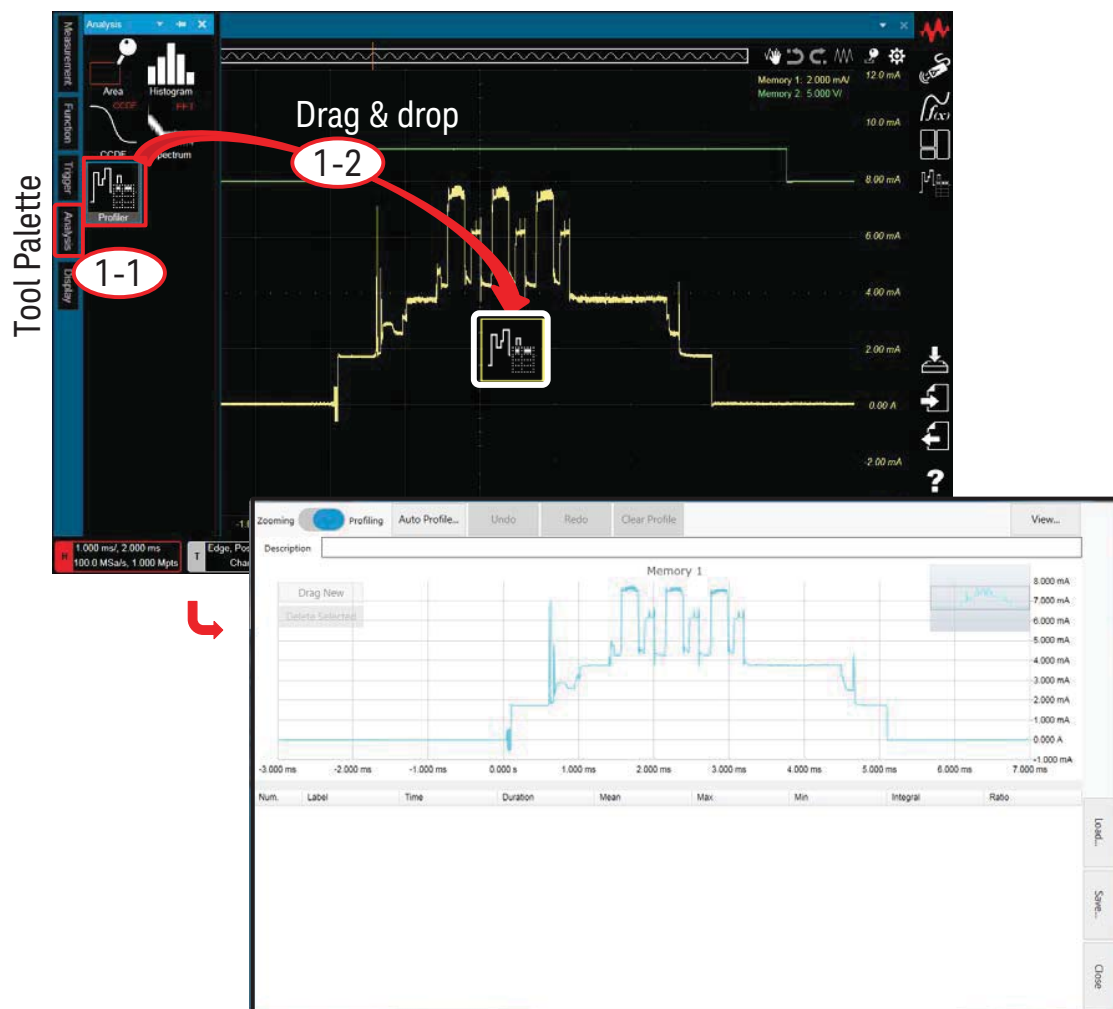


2. Perform the Auto Profile (Optional)

- The Auto Profile function can be used to help divide the segments into levels to view the vertical level difference. You can also manually adjust the levels of the segments.

2-1. Click the “Auto Profile” button to open the Auto Profile window.

2-2. Check the Condition

- Set the number of levels that you want to divide the segments into for the Number of Y Levels.
- Set a resolution that is lower than the minimum resolution that you want to divide into levels for the Y Step Resolution.
- Set the time which you want not to divide the different segments for X Step Resolution.


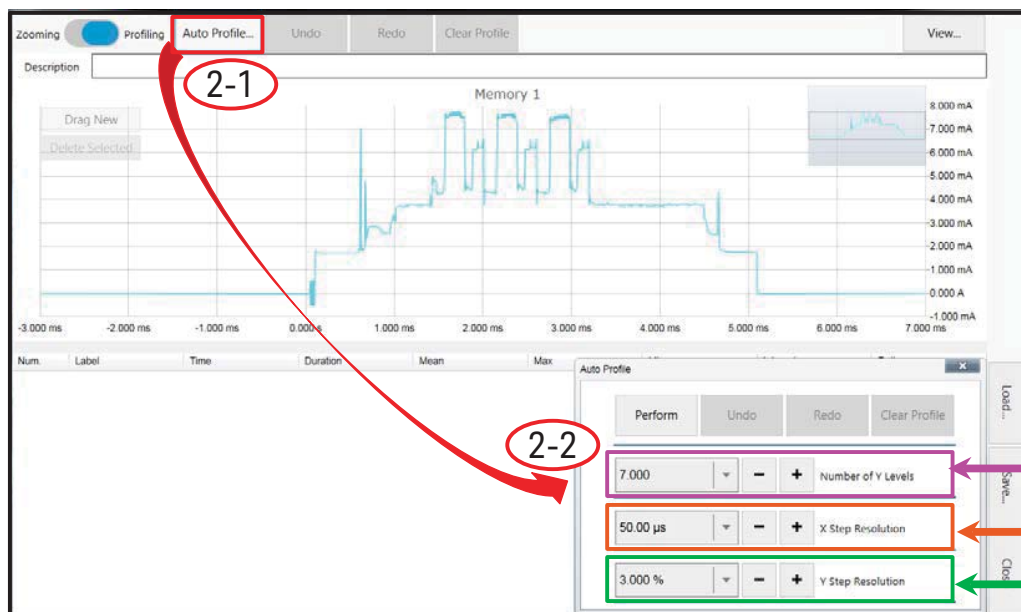


Example




2. Perform Auto Profile (Continued)

2-3. Click the “Perform” button.

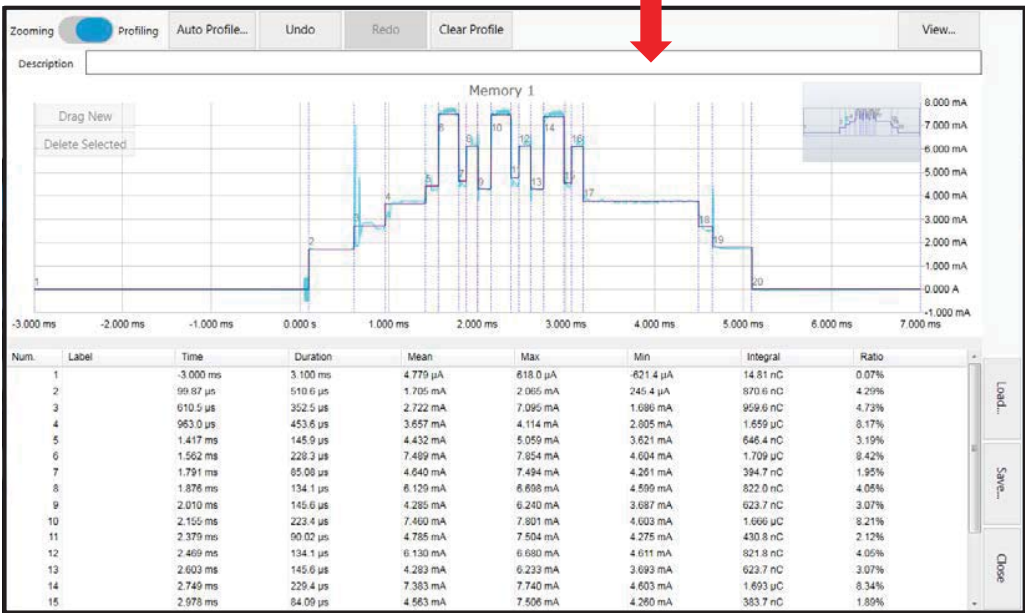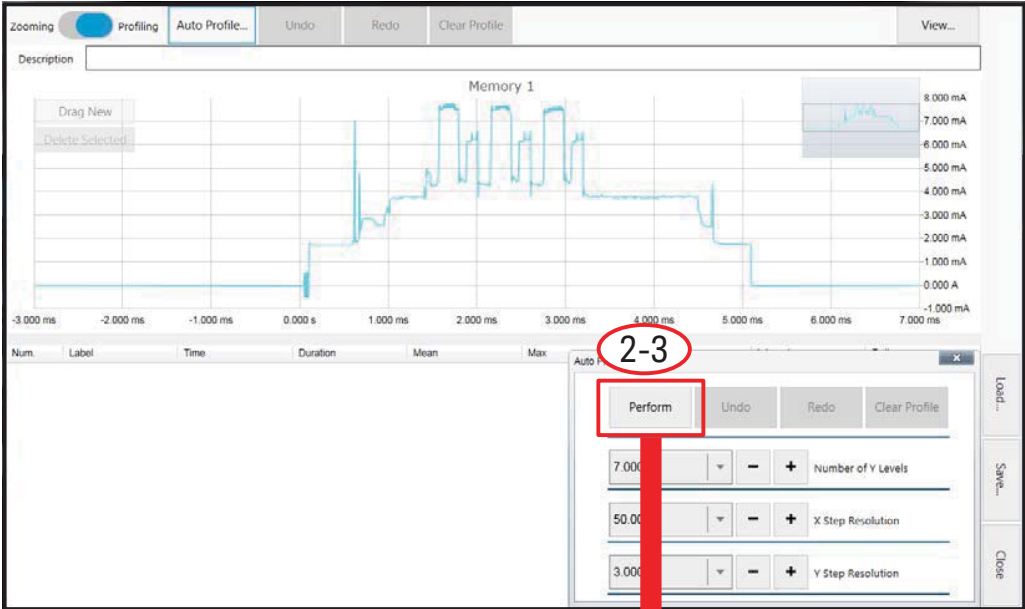

3. Adjust the position of the line of the segments manually.

3-1 Select "Zoom".

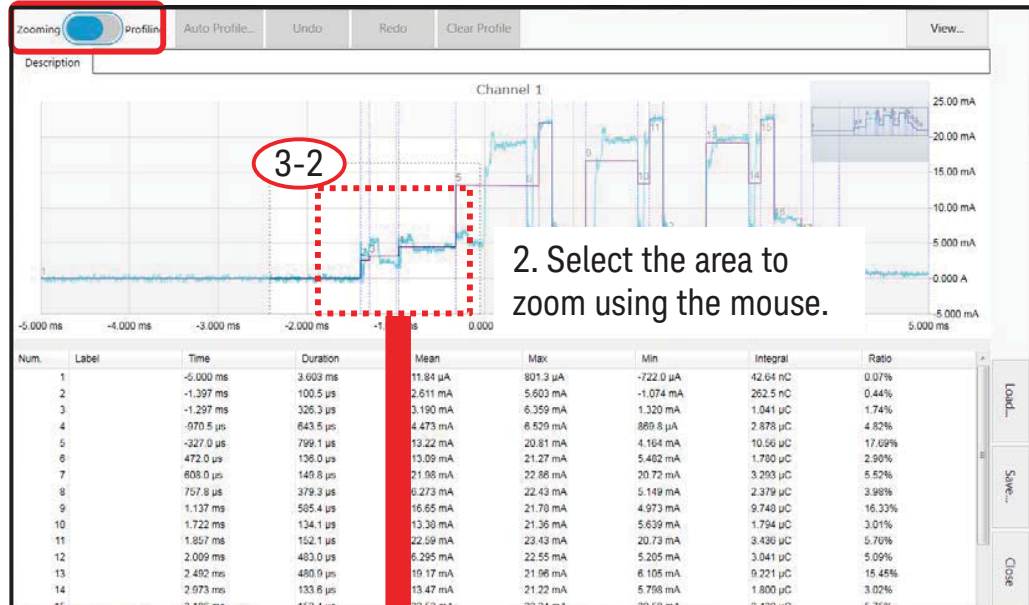

3-3 Select "Profile".



3-5 Reset zoom

Keysight Technologies

CX3300 Current Waveform Analyzer


Quick Operation Guide

How to apply the XY analysis


1. How to apply the XY analysis

1-1. Click the “Analysis” icon on the Tool Palette.


1-2. Drag and drop the “XY” icon on the Analysis tool palette on the channel set as Source1.

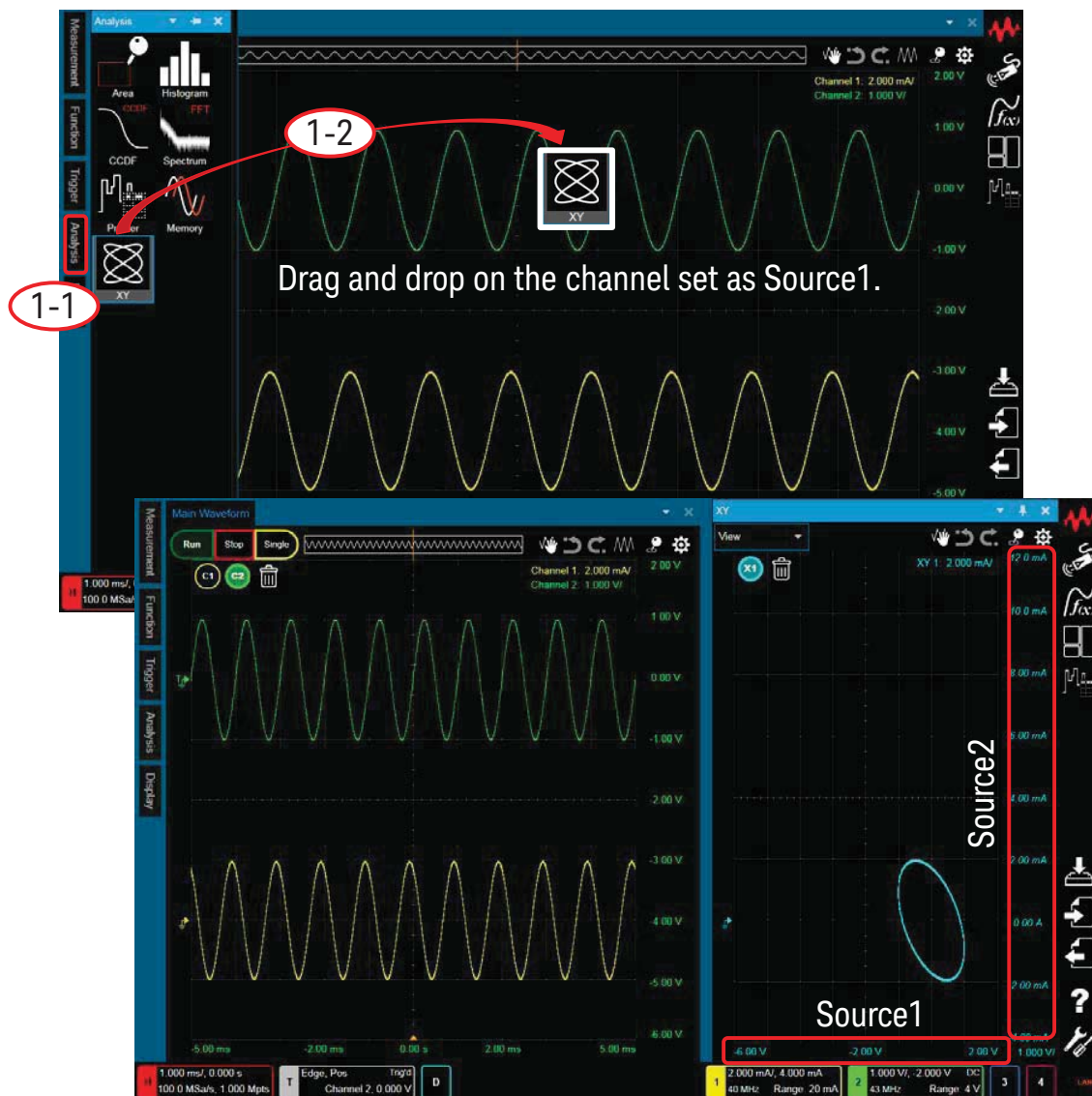


2. How to change the target channels
 - 2-1. Double-click "X1" icon in the XY graph view.
 - 2-2. Set the target channel as Source1 and Source2 on the XY1 window.

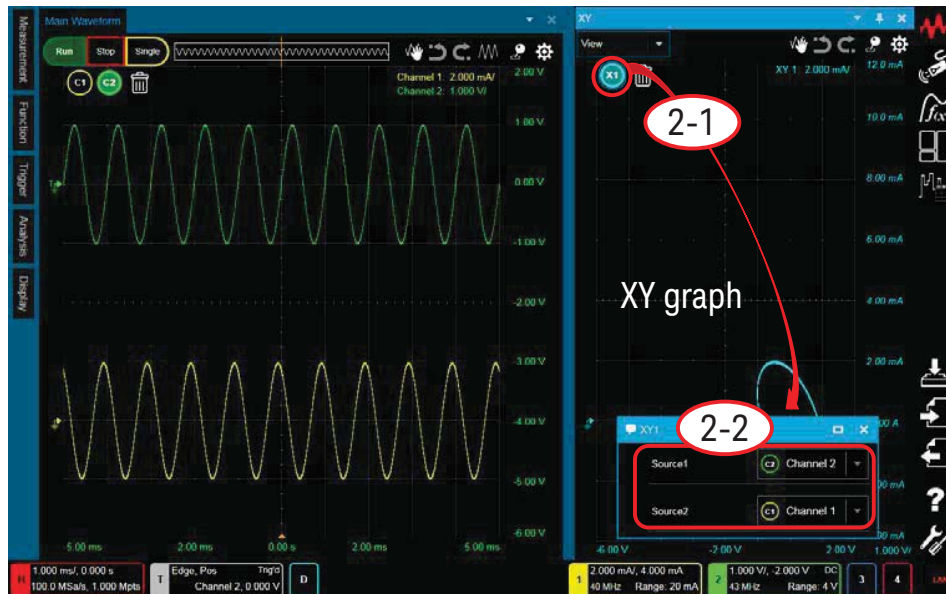

3. How to adjust the scale of the graph
 - 3-1. Double-click anywhere on the horizontal or vertical axes in the XY graph view.
 - 3-2. Un-check the "Auto" box.
 - 3-3. Adjust the scale on this area.

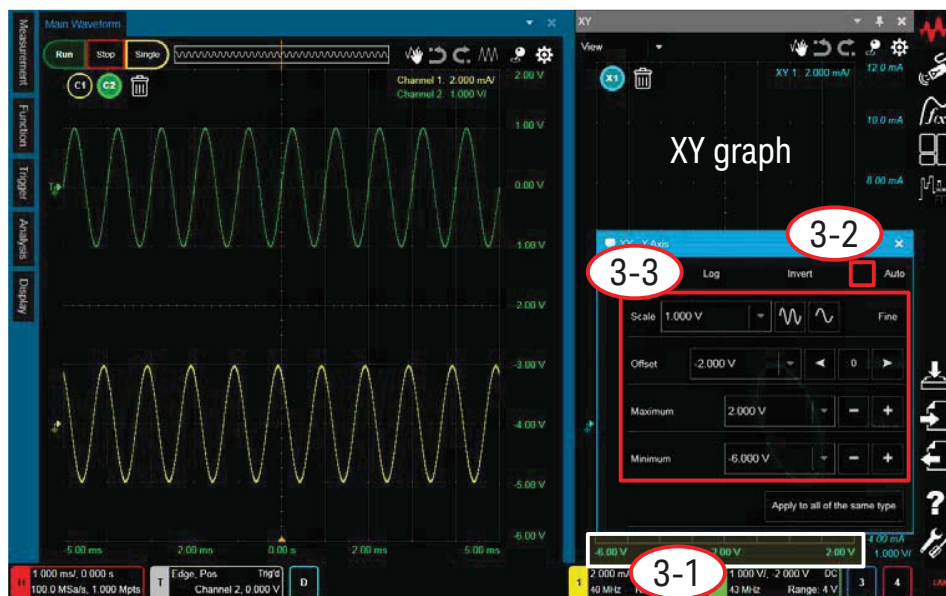

4. How to adjust a specific area on the waveform.
- 4-1. Select “Gating” from “All/View/Gating” at the top left on the XY graph.
 - 4-2. The result of the XY graph is changed in real time by adjusting the specified area on the waveform.

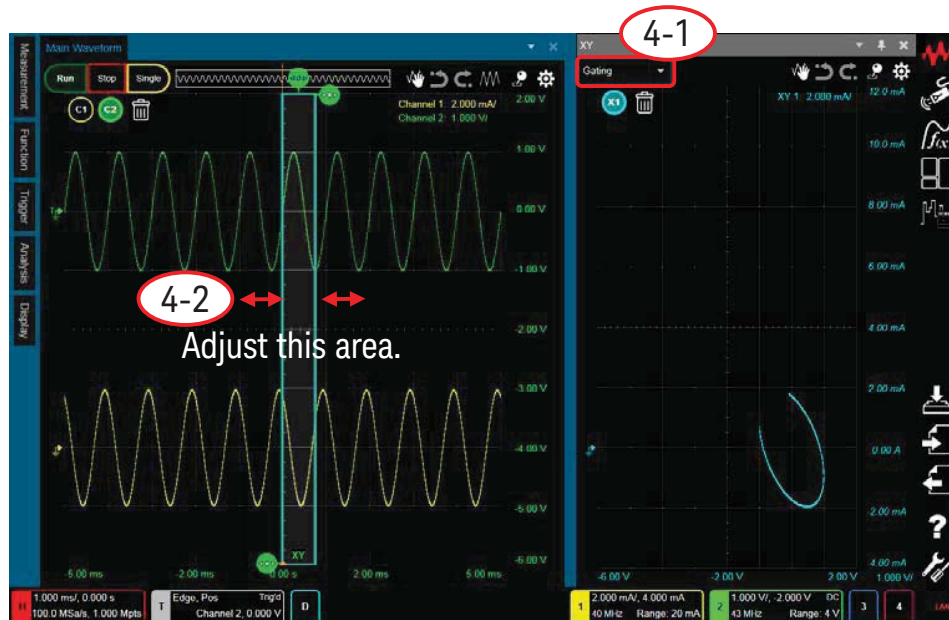


Keysight Technologies CX3300 Current Waveform Analyzer


Quick Operation Guide

How to change the measurement units (Clone to Ampere Hour) in Area analysis

1. Select “Markers” in the Display Tool Pallet, then drag and drop the “Area” icon on the interest area on the waveform.

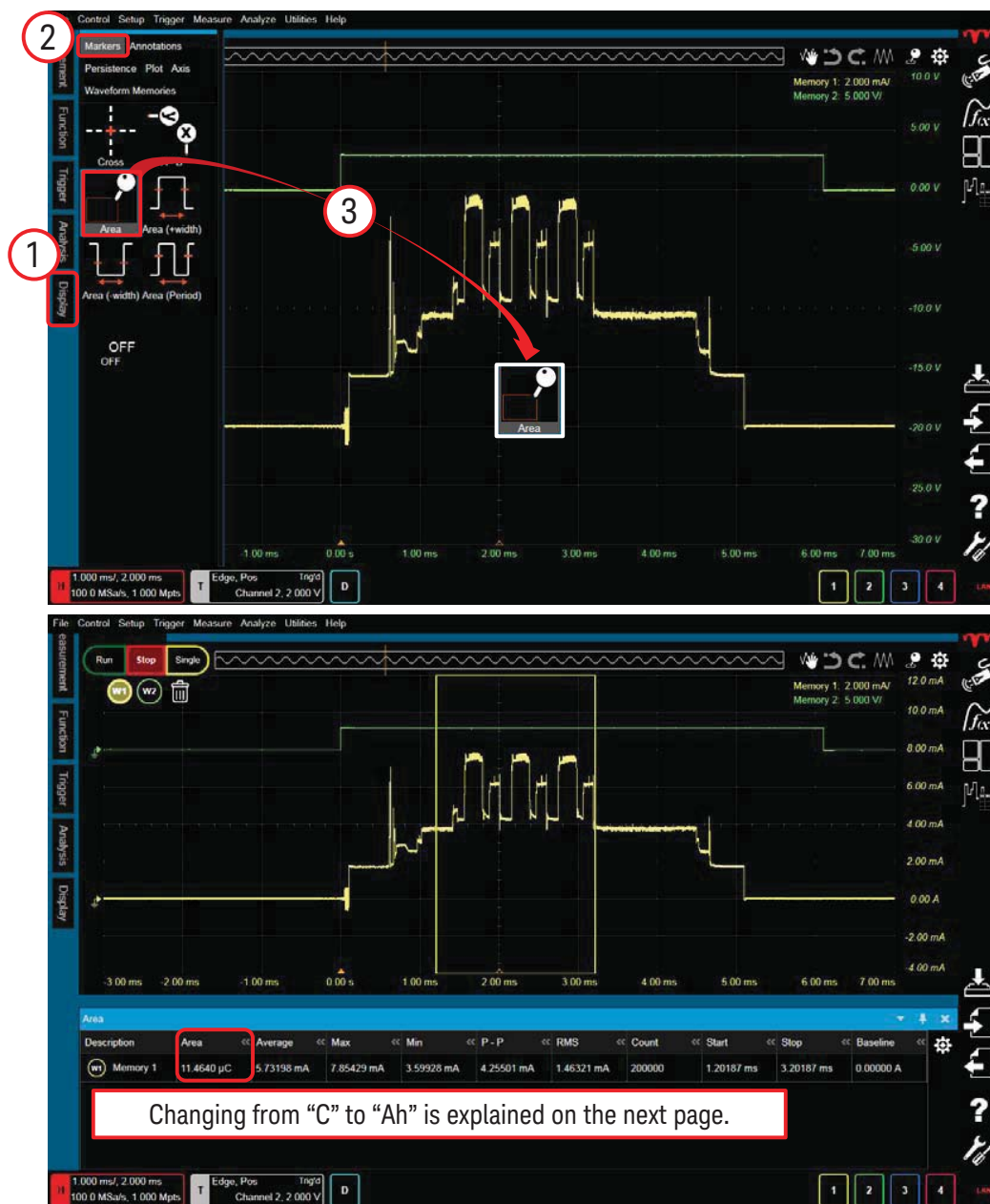

2. Select “Utilities > Configuration > User Preference” window.

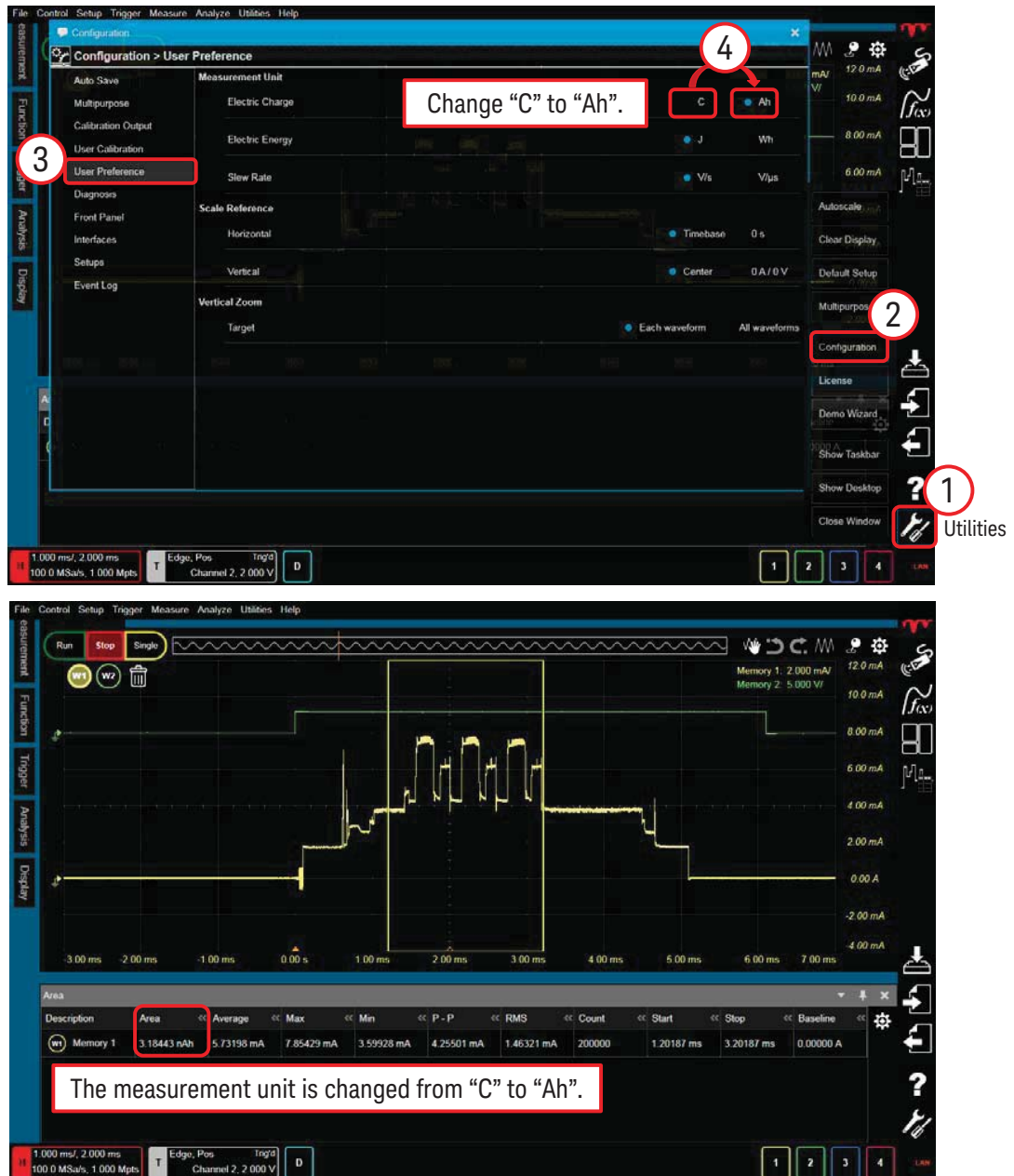

Note:

- The setting is also retained in CX3300 power cycle and [Default Setup].
- It returns to the factory condition with [Factory Default]¹.

1. [Factory Default] can be selected on the “Utilities > Configuration > Setups” window.

Keysight Technologies


CX3300 Current Waveform Analyzer

Quick Operation Guide

How to perform User calibration

User calibration is effective for making more accurate measurements by allowing you to perform error corrections of measurement data.

1. Before starting the user calibration, connect the CX1100 to the CX3300 analog input channel.

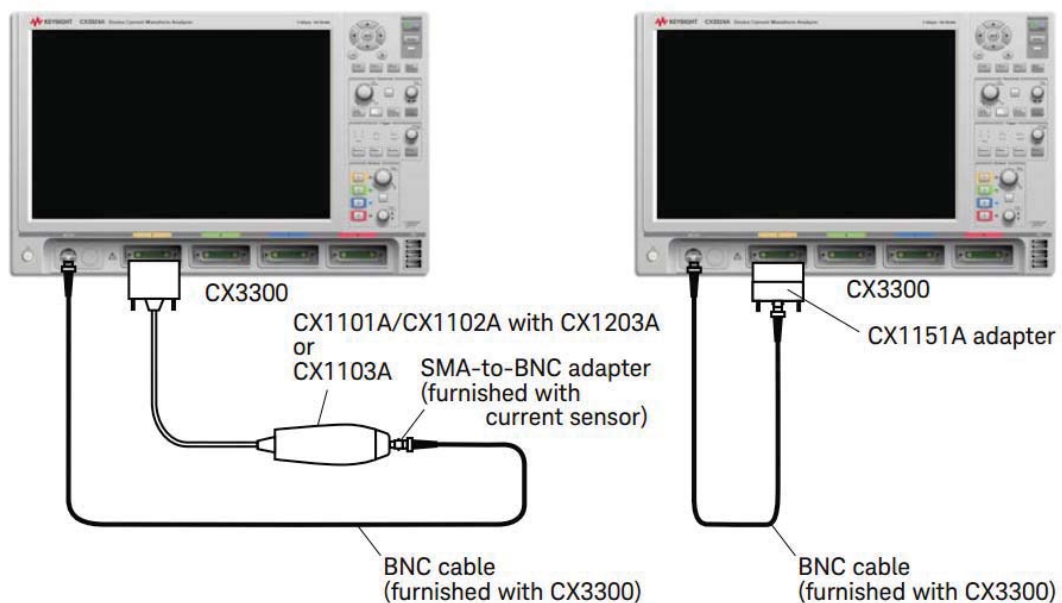

Note:


- If you connect the CX1101A/CX1102A, attach the CX1203A sensor head to it and set the CX1203A slide switch to 0 Ω .

2. How to perform User calibration ?

- 2-1. Click on the “Utilities” icon to open the Utilities menu.
- 2-2. Click “Configuration” to open the “Configuration” dialog box.
- 2-3. Click “User Calibration” to display the “Configuration>User Calibration” screen.
- 2-4. Select the “Gain/Offset” or “Offset” radio button for the channel under user calibration.
- 2-5. Click the “Execute” button of the channel.



- 2-6. The dialog box which asks you to connect the CX1100 input to the Aux Out terminal is shown. If you are ready, click the “OK” button. The user calibration will start.

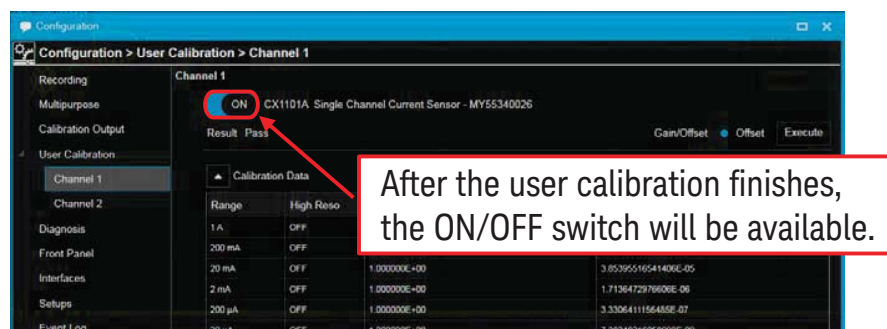

Note:

- The data is not cleared by turning the CX3300 off.
- Once you disconnect the CX1100 from the CX3300, the user calibration data for this channel will be cleared.


Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to check the firmware revisions

1. Click on the "About Device Current Waveform Analyzer" icon.
2. The "About CX3300 Device Current Waveform Analyzer" window will open.
3. The Firmware revisions are displayed.

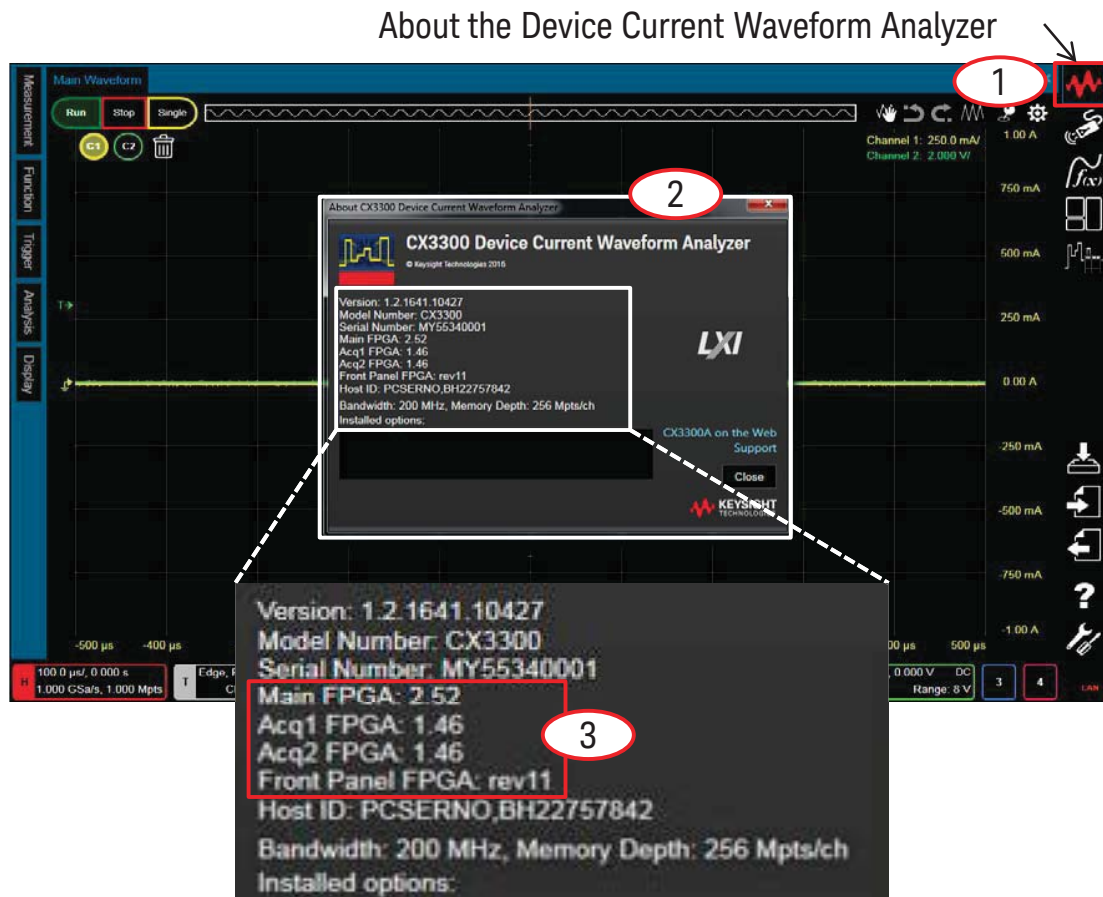

Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

How to simulate the 50 MHz or 100 MHz bandwidth option on the mainframe with the 200 MHz option.

The mainframe has three maximum bandwidth options, 50 MHz, 100 MHz and 200 MHz. In order to simulate the 50 MHz or 100 MHz option on the mainframe with the 200 MHz option, follow the procedure below.

1. Click on the "Setting" icon.
2. The "Setting" dialog box opens.
3. Click on "Setting>Channels>Channel 1" to open the channel 1 setup panel.
4. Set the Bandwidth Limit setting to "ON".
5. Set the "Bandwidth Limit" to 50 MHz for checking the performance of the 50 MHz option and to 100 MHz for checking the performance of the 100 MHz option.




Keysight Technologies

CX3300 Current Waveform Analyzer

Quick Operation Guide

Where to find the Online help

1. Click the “Help” icon on the Side Bar.
2. Click the “On-line Help (.chm)” button to open the online help.

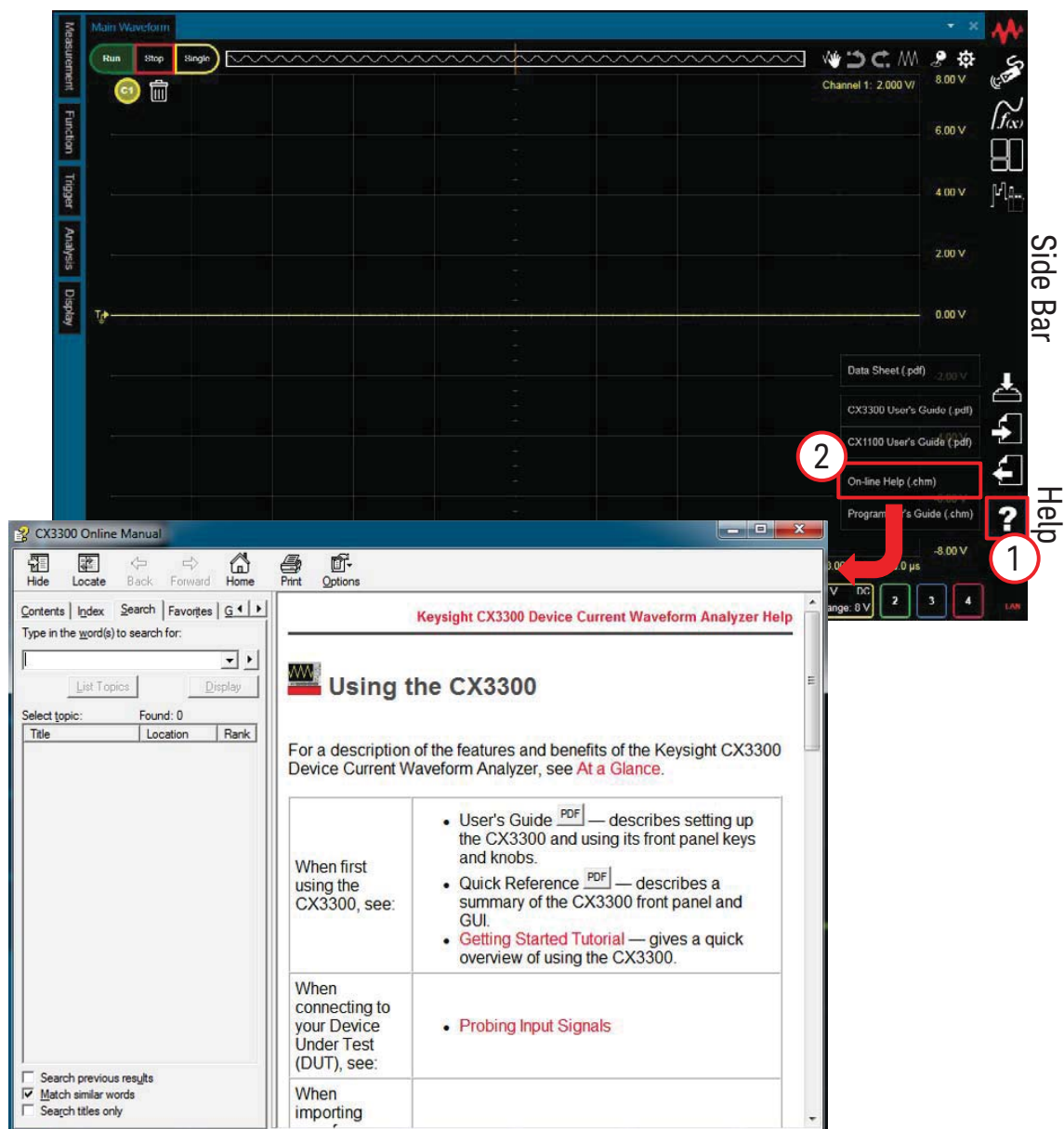

Keysight Technologies


CX3300 Current Waveform Analyzer

Quick Operation Guide

How to display the menu bar

There are two ways to display the menu bar.

1. Hover the mouse arrow over the upper edge of the monitor.
2. Click the blue bar at the top of the monitor, or click any other area.

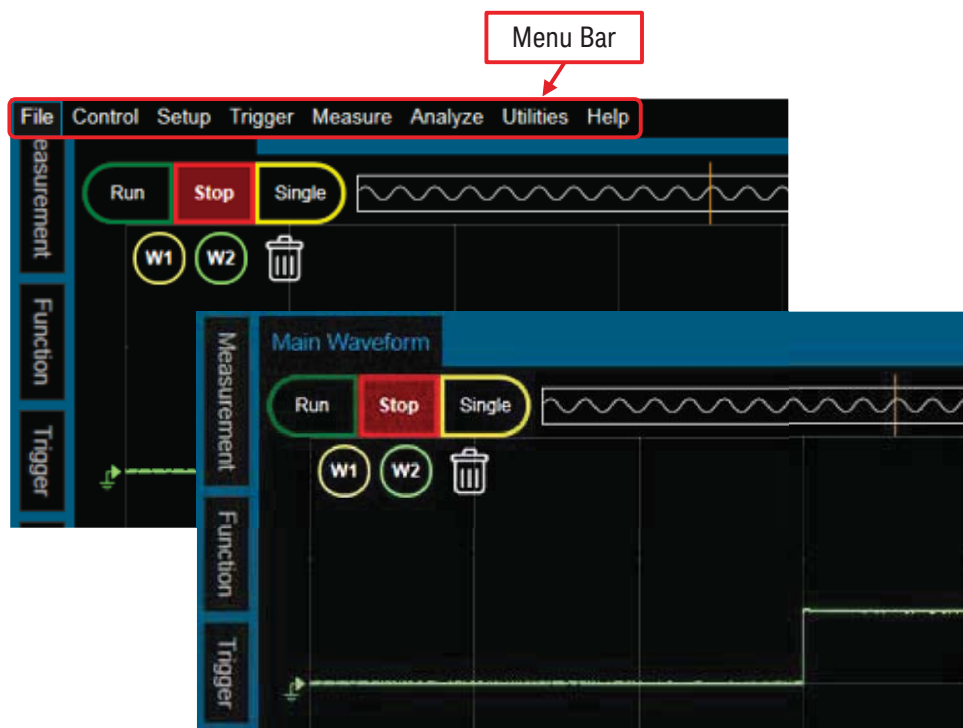

Note:

- By selecting "Utilities > Menu Bar" from the Menu Bar, it is displayed at all times. By doing the same operation again, you can hide the Menu Bar.
- The setting is also retained in CX3300 power cycle and [Default Setup].
- It returns to the factory condition with [Factory Default]¹.

1. [Factory Default] can be selected on the "Utilities > Configuration > Setups" window.

